

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [REVITIVE](#) /

› [New Revitive Medic Circulation Booster | Fight Leg Aches & Pains and Swelling | Drug-Free Relief from Persistent Foot/Leg Problems & Body Aches User Manual](#)

REVITIVE 5573AQ

REVITIVE Medic Circulation Booster User Manual

Model: 5573AQ

INTRODUCTION

The REVITIVE Medic Circulation Booster is a medical device designed to alleviate leg and foot discomfort. It utilizes advanced Electromagnetic Muscle Stimulation (EMS) technology to improve circulation, reduce swelling, and relieve muscle stiffness and cramping. This manual provides comprehensive instructions for safe and effective use of your device.

Many individuals experience leg aches, pain, cramps, or swollen feet and ankles due to poor circulation. Conditions such as diabetes, osteoarthritis, high cholesterol, high blood pressure, age, or inactivity can impact blood flow in the legs and feet, leading to these issues. The REVITIVE Medic is engineered to address these concerns.



Image: The REVITIVE Medic Circulation Booster, a circular device with footpads and a central display, shown from a top-down perspective.



Image: A diagram illustrating common areas of leg aches, pain, and swelling, highlighting how poor circulation can affect the lower limbs.

SAFETY INFORMATION

Important: Please read all safety warnings before using the REVITIVE Medic Circulation Booster.

- Do not use if you are fitted with an electronic implant, such as a pacemaker or AICD.
- Do not use if you are being treated for, or suspect you have, a deep vein thrombosis (DVT).
- Do not use if you are pregnant.
- Caution is recommended when using with metallic implants, as direct contact may cause discomfort.
- Do not apply foot-pads or body pads directly on open wounds, rashes, swollen, red, infected, or inflamed areas, or skin eruptions (such as phlebitis, thrombophlebitis, varicose veins, cellulitis).
- Do not apply on or close to malignant tumors.
- Do not apply on areas treated with radiotherapy within the past 6 months.
- Apply stimulation only to normal, intact, healthy skin.

- If you do not feel any stimulation while the device is on, consult your doctor as you may have an undiagnosed medical condition.

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

WHAT'S IN THE BOX

Your REVITIVE Medic Circulation Booster package includes the following components:

- The rechargeable REVITIVE Medic Circulation Booster device.
- One charging cable for internal battery charging.
- One remote control.
- One set of body pads (for use on other body areas like shoulders or back).



Image: The REVITIVE Medic Circulation Booster device, its remote control, and a set of body pads, laid out on a surface.

SETUP

1. **Charge the Device:** Before first use, ensure the REVITIVE Medic is fully charged using the provided charging cable. The display will indicate the battery level.
2. **Placement:** Place the device on a flat surface in front of a comfortable chair or sofa.
3. **Remote Control:** Ensure the remote control has working batteries. It allows for convenient adjustment of settings during use.



Image: The REVITIVE Medic Circulation Booster positioned on a wooden floor, ready for use, with a small plant in the background.

OPERATING INSTRUCTIONS

The REVITIVE Medic is designed for ease of use and can be operated cordlessly. It is recommended to use the device for 30-60 minutes per day for 6-8 weeks for optimal benefits.

1. **Positioning:** Sit comfortably and place your bare feet on the footpads of the device. Ensure full contact with both footpads for effective stimulation.

2. **Power On:** Press the power button on the device or the remote control to turn it on.
3. **Adjust Intensity:** Use the '+' and '-' buttons on the device or remote to adjust the intensity level (0-99). Start at a comfortable level, such as 30, and gradually increase as your feet become accustomed to the sensation. The display shows the current intensity.
4. **Set Time:** The device typically defaults to a 30-minute session. You can adjust the session duration if needed, though 30 minutes is generally recommended.
5. **Body Pad Use (Optional):** If using body pads for other areas, connect them to the device and apply them to the desired area. Note that foot relief and body pad use cannot be performed simultaneously.
6. **Session Completion:** The device will automatically turn off after the set time.



Image: A woman sitting on a sofa with her bare feet placed on the REVITIVE Medic Circulation Booster, demonstrating typical use.

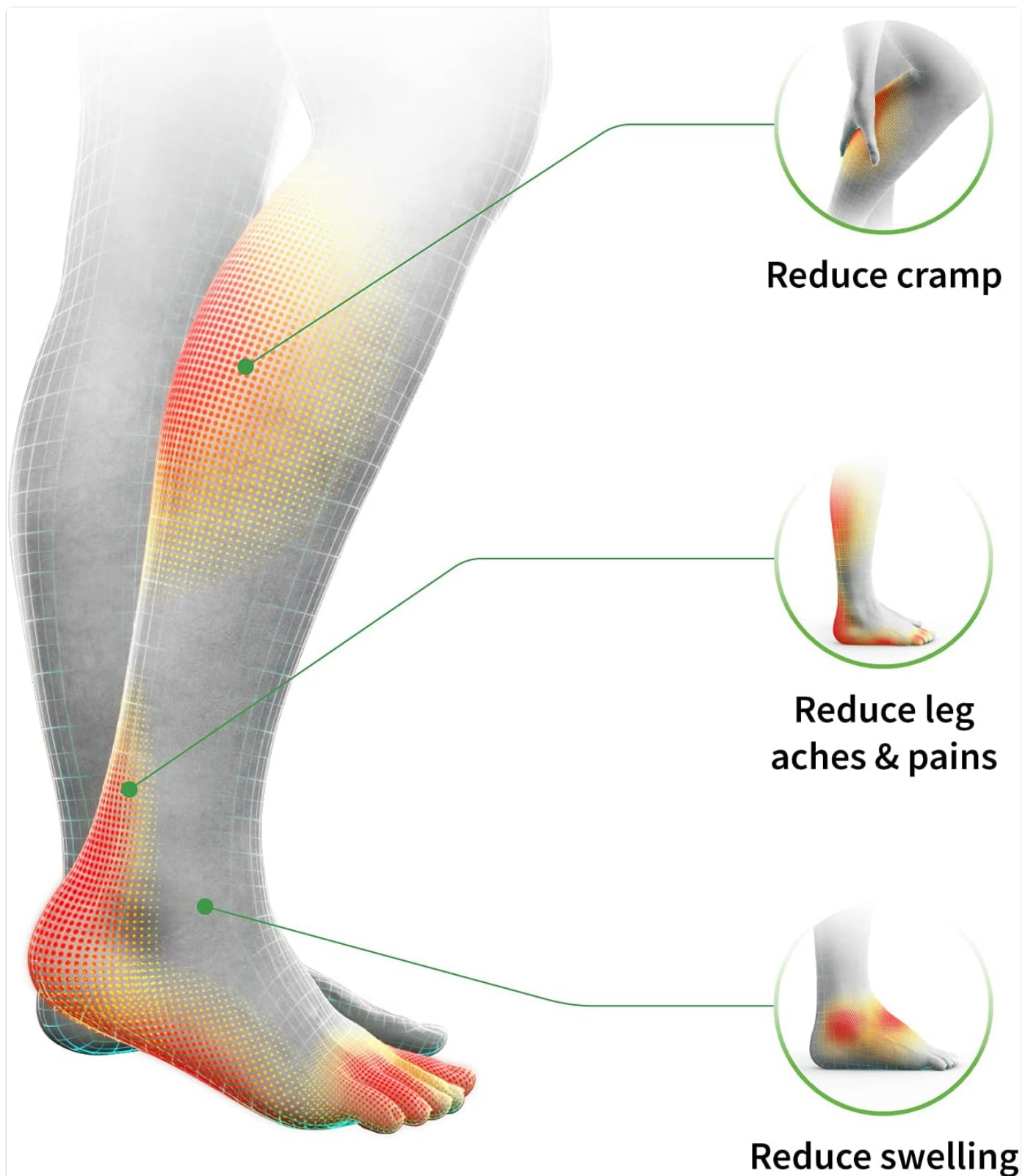


Image: A woman with body pads attached to her upper back, indicating the use of the device for targeted pain relief beyond the feet.

Understanding the Technology

The REVITIVE Medic uses OxyWave Technology, which is clinically proven to increase circulation. It sends gentle electrical impulses to stimulate calf muscles, pumping fresh, oxygen-rich blood throughout your legs and feet. The patented IsoRocker system allows your feet and ankles to rock, replicating natural heel-toe exercises to maximize the benefit of each therapy session.

Your browser does not support the video tag.

Video: An official product video from Actegy Health Inc. demonstrating how the REVITIVE Medic helps relieve foot and leg pain by stimulating circulation.

MAINTENANCE

- **Cleaning:** Wipe the device with a soft, damp cloth after each use. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device in a cool, dry place away from direct sunlight when not in use.
- **Battery Care:** Recharge the device regularly to maintain battery health, even if not used daily.
- **Pad Replacement:** Replace body pads when they lose their adhesiveness or become worn.

TROUBLESHOOTING

Problem	Possible Cause	Solution
No stimulation felt.	Feet not making full contact; intensity too low; dry skin.	Ensure both feet are fully on the pads. Increase intensity. Apply a small amount of lotion to feet if skin is dry.
Device not turning on.	Low battery; power button not pressed correctly.	Charge the device. Press and hold the power button firmly.
Remote control not working.	Dead batteries; remote out of range.	Replace remote batteries. Ensure remote is pointed at the device and within close range.

SPECIFICATIONS

- Product Dimensions:** 14.8 x 14.8 x 3.7 inches
- Item Weight:** 4.59 Pounds
- Model Number:** 5573AQ
- Batteries:** 1 Lithium Ion battery required (included)
- Controls Type:** Remote
- Display Type:** Built-in Digital Display
- Color:** White and Black
- Manufacturer:** Actegy Health, Inc.

WARRANTY & SUPPORT

The REVITIVE Medic Circulation Booster comes with a 2-Year Limited Warranty upon registration. For detailed warranty information, product support, or to purchase replacement parts, please refer to the official user manual PDF or visit the REVITIVE website.

You can access the full User Manual [here](#).

For additional support, you may visit the [REVITIVE Store on Amazon](#).