

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

> [comfheat](#) /

> [Comfheat USB Heating Pad User Manual](#)

comfheat CFH003

Comfheat USB Heating Pad User Manual

Brand: comfheat

Model: CFH003

INTRODUCTION

Thank you for choosing the Comfheat USB Heating Pad. This portable heating pad is designed to provide soothing warmth and pain relief for various body parts, ideal for travel, outdoor activities, or simply relaxing at home. This manual provides essential information for safe and effective use of your heating pad.



Image: The Comfheat USB Heating Pad in grey, showing its compact size and attached USB power cable.

SAFETY INFORMATION

Please read all safety warnings and instructions before using this product to avoid injury or damage.

- Avoid prolonged use directly on one area of the skin to prevent burns.
- Always ensure proper electrical connections.
- Do not leave the heating pad unattended while in use.
- This USB heating pad cannot be charged. It requires continuous connection to a 5V/2A USB power source to operate.
- Do not use if the pad or cord is damaged.
- Keep out of reach of children and pets.
- Do not fold or crush the heating pad during use.
- Do not use while sleeping.
- Consult a physician before use if you have any medical conditions, especially those affecting circulation or sensation.

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

PRODUCT FEATURES

- **Convenient USB-Powered Heating:** Measures 12"x16" with a 5.9ft long USB cord. Plugs into any 5V/2A USB power source (power bank, laptop, etc.).
- **3 Heat Settings:** High (130 - 140°F), Medium (112 - 122°F), and Low (100 - 104°F).
- **Auto Shut-Off:** Built-in intelligent control chip ensures stable heating and prevents overheating; automatically turns off after 45 minutes for safety and energy saving.
- **Versatile Pain Relief:** Provides soothing warmth for sciatica, muscle soreness, and chronic pain. Suitable for various body parts and situations like travel, camping, or post-workout recovery.
- **Compact & Lightweight:** Easy to carry and store, weighing only 0.33 lbs.

3-Level Temperature

comfheat®

Warms you up, improves circulation and relieves stiffness and pain



High
130 - 140°F



Medium
112 - 122°F



Low
94 - 104°F

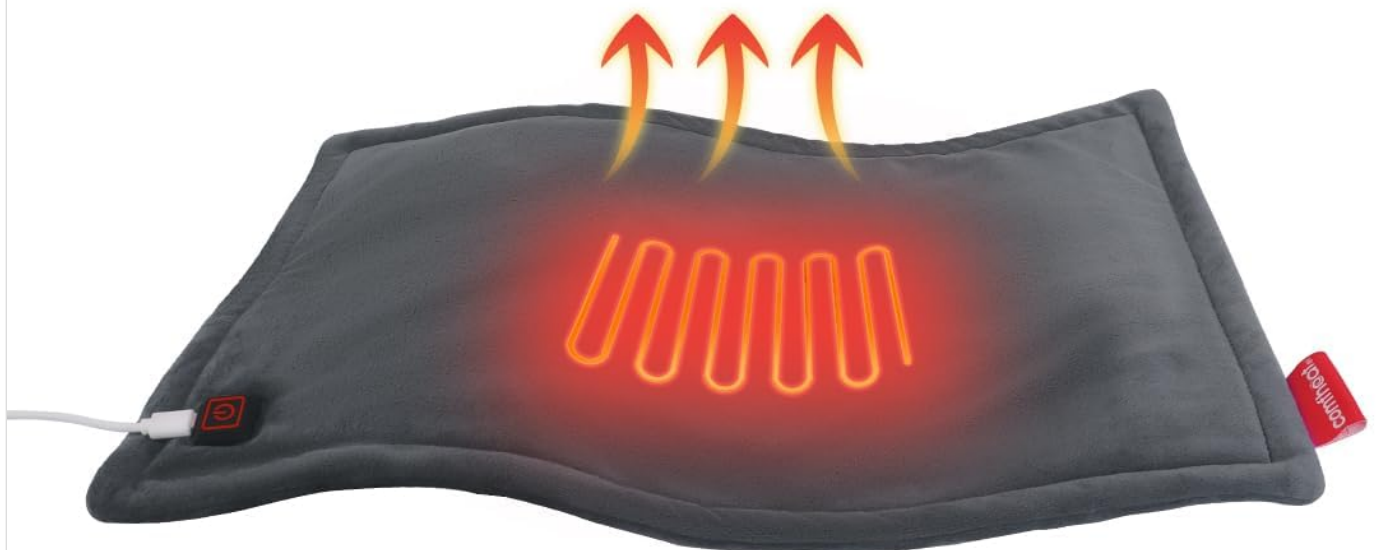


Image: Visual representation of the three heat settings and their corresponding temperature ranges.

PACKAGE CONTENTS

The Comfheat USB Heating Pad package includes:

- 1 x Comfheat USB Heating Pad (Grey, 16" x 12")
- 1 x Integrated USB Power Cord (5.9 ft)
- 1 x Storage Pouch

Note: Power bank and USB adapter are not included and must be purchased separately.



Image: The heating pad's compact size and the included storage pouch for portability.

SETUP

To set up your Comfheat USB Heating Pad for use:

1. **Connect to Power Source:** Insert the USB plug of the heating pad's cord into a compatible 5V/2A USB power

source. This can be a power bank, a laptop's USB port, a car's USB port, or a USB wall adapter.

2. **Position the Pad:** Place the heating pad on the desired area of your body for heat therapy. The pad is versatile and can be used on the back, abdomen, shoulders, or other areas requiring warmth.

comfheat

Convenient USB Powered

Note: Need to plug **2A** port to work



Optional Power Source (Not included)



PC USB Port



Power Bank



Cigarette Lighter



USB Adapter

Image: Demonstrates various ways to power the USB heating pad, including laptops, power banks, and wall adapters.

Your browser does not support the video tag.

Video: This video demonstrates how to connect the USB heating pad to a power bank or a wall adapter, highlighting its portable design. It also shows how to turn on the device and cycle through temperature settings, and illustrates its use for shoulder, abdomen, and back pain relief.

OPERATING INSTRUCTIONS

Follow these steps to operate your Comfheat USB Heating Pad:

1. **Turn On:** Once connected to a power source, long press the power button located on the pad's control unit. The red light will illuminate, indicating high temperature mode.
2. **Adjust Temperature:** Short press the power button to cycle through the temperature settings:
 - Red Light: High Temperature (130 - 140°F)
 - Blue Light: Medium Temperature (112 - 122°F)
 - Green Light: Low Temperature (100 - 104°F)
3. **Automatic Shut-Off:** The heating pad is designed to automatically turn off after approximately 45 minutes of continuous use for safety and energy conservation. If you wish to continue use, simply long press the power button again to restart.
4. **Turn Off:** To manually turn off the heating pad, long press the power button until the light indicator turns off. Alternatively, unplug the USB power cord from the source.



Image: Detailed view of the control button and its light indicators for temperature selection.



comfheat®

For Camping Use



In the Car



Warms up
Your Pet

Image: Examples of the heating pad's versatile use for personal comfort and pain relief in different environments.

MAINTENANCE

Proper care will extend the life of your Comfheat USB Heating Pad:

- **Cleaning:** Unplug the heating pad before cleaning. Wipe the surface with a damp cloth and mild detergent. Do not submerge the pad in water or machine wash. Allow to air dry completely before next use.
- **Storage:** When not in use, store the heating pad in its provided storage pouch in a cool, dry place, away from direct sunlight and excessive heat.
- **Handling:** Avoid sharp objects that could puncture or damage the heating elements inside the pad. Do not twist or excessively bend the cord.

TROUBLESHOOTING

If you encounter issues with your Comfheat USB Heating Pad, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Heating pad does not turn on.	Not properly connected to power, insufficient power supply (less than 5V/2A), or power button not long-pressed.	Ensure the USB plug is fully inserted into a 5V/2A power source. Long press the power button for 3 seconds to turn on. Try a different USB port or power source.
Heating pad turns off automatically before 45 minutes.	This is normal operation for safety and energy saving.	Long press the power button to restart the heating cycle.
Pad is not warm enough.	Low temperature setting selected, or power source is not providing sufficient current.	Short press the power button to cycle to a higher temperature setting (red light for high). Ensure your USB power source provides at least 5V/2A.
Pad feels too hot.	High temperature setting selected.	Short press the power button to cycle to a lower temperature setting (blue for medium, green for low). Discontinue use if discomfort persists.

SPECIFICATIONS

Feature	Detail
Model Number	CFH003
Dimensions	16" x 12" (40.64 cm x 30.48 cm)
Power Input	5V/2A USB (DC)
Wattage	10 watts
Temperature Settings	High (130-140°F), Medium (112-122°F), Low (100-104°F)
Auto Shut-Off	45 minutes
Material	Grey fabric (soft plush)
Item Weight	0.18 Kilograms (approx. 6.38 ounces)
Cord Length	5.9 ft

WARRANTY AND SUPPORT

The Comfheat USB Heating Pad comes with a **1-year warranty** and hassle-free replacement. If you experience any issues with your product, please reach out to the manufacturer, Comfheat, for assistance. They are committed to ensuring your complete satisfaction and will assist you within 12 hours.

For support, please refer to the contact information provided with your product packaging or visit the official Comfheat store on Amazon: [Comfheat Amazon Store](#).

