

SINGCALL B09D7RCNW9

SINGCALL Caregiver Pager Wireless Call Button System Instruction Manual

Model: B09D7RCNW9 | Brand: SINGCALL

INTRODUCTION

This manual provides detailed instructions for the setup, operation, and maintenance of your SINGCALL Caregiver Pager Wireless Call Button System. This system is designed to provide a reliable and convenient way for individuals needing assistance to alert their caregivers.

PRODUCT OVERVIEW

Components

- **1 Mobile Receiver:** A portable device that displays call alerts.
- **2 Call Buttons:** Wireless buttons used to send call alerts.

Key Features

- **Wrist Pager Receiver:** Can be worn around the neck with a lanyard.
- **Multi-functional Receiver:** Includes calendar, clock, and alarm clock functions.
- **Rechargeable Battery:** Receiver offers extended standby time (2-3 days).
- **Fast Delete Function:** Allows quick deletion of individual or all registered numbers from the receiver.
- **Encryption Function:** Receiver protects against accidental data loss.
- **Waterproof Call Button:** Emergency call buttons are designed to be waterproof.
- **Wide Compatibility:** SOS Transmitters are compatible with all SINGCALL receivers.
- **Long Transmission Distance:** Approximately 500 feet wireless range.
- **Portable Call Button:** Easy to carry, can be worn, taped, or placed conveniently.


Product Images



The complete SINGCALL Caregiver Pager system, featuring one mobile receiver and two call buttons.




An elderly person and a caregiver demonstrating the use of the SINGCALL pager system in a home setting.

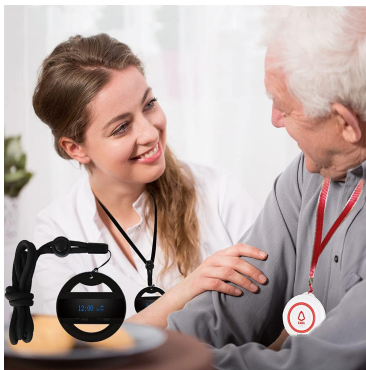

A close-up view of the SINGCALL call button and receiver being worn by an elderly person and caregiver.




An overhead diagram showing the effective wireless range of the system across multiple rooms in a building.

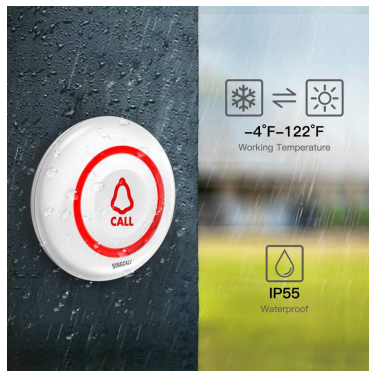

The call button's durability, highlighting its IP55 waterproof rating and operational temperature range.

Compact Size

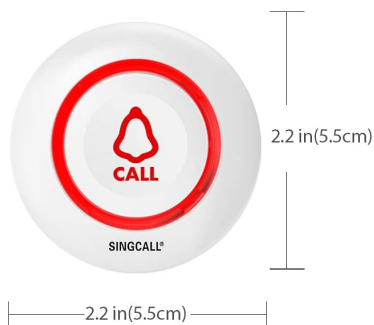

A visual representation of the compact dimensions of the SINGCALL call button.



The SINGCALL call button shown with its red lanyard, ready for neck wear.



The SINGCALL mobile receiver displayed with its black lanyard for convenient carrying.

SETUP INSTRUCTIONS

1. Charging the Receiver

Before first use, fully charge the mobile receiver using the provided charging cable. The battery indicator on the OLED display will show charging status.

2. Pairing Call Buttons with Receiver

To ensure the call buttons communicate with the receiver, they must be paired. Follow these steps:

1. Turn on the mobile receiver.
2. Enter the pairing mode on the receiver (refer to the receiver's on-screen menu for 'Pairing' or 'Registration' option).
3. Press the call button you wish to pair. The receiver should confirm successful pairing with a display message or sound.
4. Repeat for all additional call buttons.

For a visual guide on the pairing process, please refer to the video below:

Your browser does not support the video tag.

Video demonstrating the pairing process between the SINGCALL mobile receiver and call buttons.

3. Wearing and Placement

- **Receiver:** Use the provided lanyard to wear the receiver around the neck for easy access by the caregiver.
- **Call Buttons:** Call buttons can be worn around the neck, placed in a pocket, or mounted on a wall near the user's bed, bathroom, or common areas. Ensure they are within easy reach.

OPERATING INSTRUCTIONS

1. Making a Call

To send an alert, simply press the 'CALL' button on the wireless call button. The button will illuminate briefly, indicating the signal has been sent.

2. Receiving a Call Alert

When a call button is pressed, the mobile receiver will:

- Display the calling number on its OLED screen.
- Emit an audible alert (if enabled).
- Vibrate (if enabled).

Acknowledge the call by attending to the individual who initiated the alert.

3. Using Receiver Functions

- **Clock/Calendar:** The receiver displays the current time and date. Adjust these settings via the 'FUN' (Function) button.
- **Alarm Clock:** Set alarms through the receiver's menu for reminders.
- **Delete Registered Numbers:** Use the 'ENT' (Enter) button and 'FUN' button to navigate the menu

and delete specific call button registrations or clear all.

MAINTENANCE

1. Battery Replacement (Call Buttons)

The call buttons use CR2 batteries. When the button's indicator light fails to illuminate or the range decreases significantly, replace the battery. Carefully open the button casing, replace the battery with a new CR2 battery, and ensure the casing is sealed properly to maintain its waterproof integrity.

2. Cleaning

Wipe the receiver and call buttons with a soft, damp cloth. Do not use harsh chemicals or abrasive cleaners. Ensure the devices are dry before use or charging.

TROUBLESHOOTING

- **No Signal/Receiver Not Alerting:**
 - Ensure the call button battery is not depleted.
 - Verify the call button is paired with the receiver. Re-pair if necessary.
 - Check if the receiver's volume or vibration settings are enabled.
 - Ensure the call button and receiver are within the specified operating range (approx. 500ft in open air, less indoors due to obstacles).
- **Receiver Not Charging:**
 - Check the charging cable and adapter for damage.
 - Ensure the charging port on the receiver is clean and free of debris.
- **Incorrect Time/Date on Receiver:**
 - Access the receiver's menu using the 'FUN' button to adjust time and date settings.

PRODUCT SPECIFICATIONS



Feature	Detail
Control Method	Touch
Noise Level	9 Decibels
Size (Call Button)	2.2 inches (5.5 cm) diameter
Item Weight	7.2 ounces (total package)
Package Dimensions	4.3 x 3.6 x 2.5 inches
Batteries (Call Buttons)	2 CR2 batteries required (included)
Display Style (Receiver)	OLED
Battery Cell Type	Lithium Manganese Dioxide

Power Source	Battery Powered
Voltage	3.7 Volts
Working Temperature	-4°F to 122°F
Waterproof Rating (Call Button)	IP55

WARRANTY AND SUPPORT

For warranty information and technical support, please refer to the product packaging or contact SINGCALL customer service directly. Keep your purchase receipt for warranty claims.

Related Documents - B09D7RCNW9

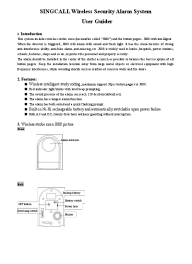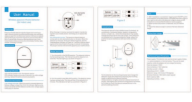	<p>SINGCALL J008 Wireless Security Alarm System User Guide</p> <p>User guide for the SINGCALL J008 Wireless Security Alarm System, detailing its features, operation instructions for pairing pagers, volume and alarm time settings, and technical specifications.</p>
	<p>SINGCALL Wireless Passive Infrared Detector User Manual</p> <p>Comprehensive user manual for the SINGCALL wireless passive infrared detector (indoor use), covering features, appearance, anti-dismantle function, mode settings, installation, testing, detection range, working indicators, functions, and technical specifications. Learn how to set up and use your PIR motion sensor.</p>