

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [REVITIVE](#) /

› [Revitive ProHealth Circulation Booster User Manual](#)

REVITIVE 6302-5572AQ-US

Revitive ProHealth Circulation Booster User Manual

Model: 6302-5572AQ-US | Brand: REVITIVE

INTRODUCTION

The Revitive ProHealth Circulation Booster is an FDA-cleared medical device designed to alleviate symptoms of tired, aching, and heavy-feeling legs and feet, with or without occasional foot and ankle swelling. It utilizes advanced Electromagnetic Muscle Stimulation (EMS) technology to improve blood circulation in the lower extremities.



Image: The Revitive ProHealth Circulation Booster, a compact, oval-shaped device with a white central panel and black textured footpads.

IMPORTANT SAFETY INFORMATION

Please read all instructions carefully before using the Revitive ProHealth Circulation Booster. Failure to do so may result in injury or damage to the device.

Contraindications:

Revitive is not suitable if you are:

- Fitted with an electronic implanted device such as a heart pacemaker or AICD.
- Being treated for, or have the symptoms of, an existing Deep Vein Thrombosis (DVT).
- Pregnant.



Image: A warning symbol overlaid on the Revitive device, emphasizing the importance of reviewing contraindications before use.

WHAT COMES IN THE BOX

Upon opening your Revitive ProHealth package, you should find the following items:

- One (1) Revitive ProHealth Circulation Booster unit.
- One (1) Power cable.
- One (1) Remote control.
- One (1) User manual (this document).
- Information regarding a 2-Year Limited Warranty (with Registration).

SETUP AND OPERATION

The Revitive ProHealth is designed for easy setup and use. Follow these steps to begin your therapy:

1. Placement and Power Connection:

Place the Revitive ProHealth unit on a flat, stable surface in front of a comfortable chair. Ensure it is within reach of a wall outlet. Connect the power cable to the unit and then plug it into a wall outlet. The device is corded and must be used while plugged in.



Image: The Revitive ProHealth device positioned on the floor, ready for use while seated.

2. Seated Position and Foot Placement:

Sit comfortably in a chair. Place both bare feet on the textured footpads of the device. Ensure full contact for optimal performance. **Do not stand on the device.**



Images: (Left) Close-up of feet being placed on the device. (Right) A user seated, operating the device with the remote control.

3. Adjusting Intensity:


Use the remote control to turn on the device and adjust the intensity level. The device offers 99 intensity levels. It is recommended to start at intensity level 30 and gradually increase it over time to maximize benefits, based on your foot sensitivity. You should feel a comfortable muscle contraction.



Image: A close-up of the Revitive ProHealth display showing intensity level '40' and the remote control used for adjustments.

4. Recommended Usage:

For optimal benefits, it is recommended to use the Revitive ProHealth for 30 to 60 minutes per day for at least 6 weeks. The device's OxyWave Technology sends gentle electrical waves to stimulate calf muscles, pumping oxygen-rich blood through your legs and feet. The patented IsoRocker system also allows for a 15-degree rocking motion, replicating heel-toe exercises to further enhance circulation.

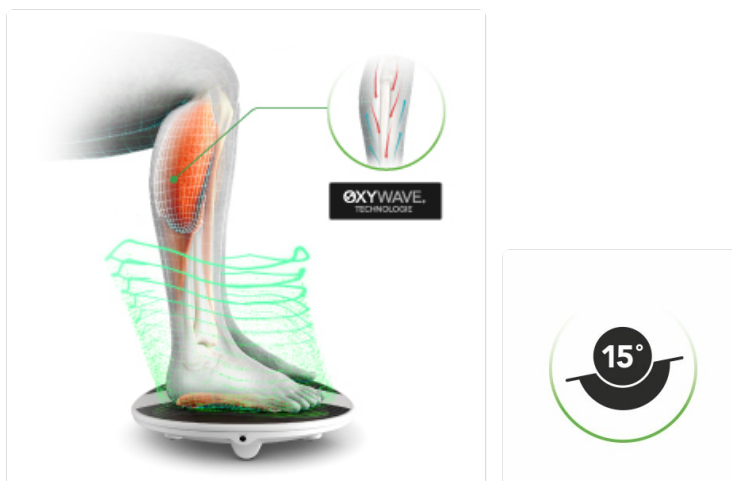

Images: (Left) A visual representation of OxyWave Technology stimulating leg muscles to improve circulation. (Right) An icon depicting the 15-degree IsoRocker movement for ankle exercise.

MAINTENANCE

To ensure the longevity and optimal performance of your Revitive ProHealth Circulation Booster, follow these simple maintenance guidelines:

- **Cleaning:** Wipe the device with a soft, damp cloth after each use. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device in a cool, dry place away from direct sunlight and extreme temperatures.
- **Power Cable:** Inspect the power cable regularly for any signs of damage. If damaged, discontinue use and contact customer support.

TROUBLESHOOTING

If you encounter any issues with your Revitive ProHealth Circulation Booster, please refer to the following common troubleshooting tips:

Problem	Possible Cause	Solution
Device does not turn on.	Not properly plugged in; power outlet issue.	Ensure power cable is securely connected to the device and a working wall outlet. Try a different outlet.
No sensation or weak sensation.	Feet not making full contact; intensity too low; dry skin.	Ensure both feet are fully on the footpads. Increase intensity gradually. Lightly dampen your feet with a damp towel for better conductivity.
Remote control not working.	Dead batteries; obstruction.	Replace remote control batteries. Ensure no obstructions between remote and device.

For issues not listed here, or if troubleshooting steps do not resolve the problem, please contact customer support.

SPECIFICATIONS

Feature	Detail
Brand	REVITIVE

Feature	Detail
Model Number	6302-5572AQ-US
Power Source	Corded Electric
Material	Plastic, Rubber
Item Weight	1.11 Kilograms
Use For	Feet, Legs
FDA Clearance	Yes (Medical Device)

WARRANTY AND SUPPORT

The Revitive ProHealth Circulation Booster comes with a 2-Year Limited Warranty upon registration. For warranty claims, technical support, or further assistance, please refer to the contact information provided in your product packaging or visit the official REVITIVE website.

A comprehensive User Manual in PDF format is also available for download:[Download User Manual \(PDF\)](#)