

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

Manuals.plus /

› [ocallcoo](#) /

› ocallcoo Electric Tooth Cleaner and Tartar Remover User Manual - Model B09FS68W6P

ocallcoo B09FS68W6P

ocallcoo Electric Tooth Cleaner and Tartar Remover User Manual

Model: B09FS68W6P

1. PRODUCT OVERVIEW

The ocallcoo Electric Tooth Cleaner and Tartar Remover is designed to effectively remove dental plaque, calculus, tartar, and stains using advanced technology. This device features multiple cleaning modes, an integrated LED light, and comes with a comprehensive set of dental tools to support your oral hygiene routine.



Image 1.1: ocallcoo Electric Tooth Cleaner and Tartar Remover with included accessories.

2. SAFETY INFORMATION

- This device is not intended for use by children without adult supervision.
- If you experience any discomfort, pain, or adverse reactions while using the device, discontinue use immediately and consult a healthcare professional.
- Do not press the cleaning head hard against teeth or gums. Gentle pressure is sufficient.
- Avoid using the device on gums or soft skin to prevent injury. The device is designed to vibrate only when touching hard objects like teeth.
- Do not disassemble the product. This will void the warranty and may pose a safety risk.
- Keep out of reach of children and store in a cool, dry place away from direct sunlight.

3. PACKAGE CONTENTS

Your ocallcoo Electric Tooth Cleaner kit includes:

- 1 x Plaque Remover Device
- 2 x Replaceable Cleaning Heads (1 pointed, 1 flat)

- 2 x Stainless Steel Dental Tools
- 1 x Mouth Mirror
- 1 x Wrench (for attaching/detaching heads)
- 1 x USB Charging Cable
- 1 x Protective Cover
- 1 x User Manual

Advanced All-in-one Tooth Cleaner Set



Image 3.1: All-in-one dental cleaning kit components.

4. KEY FEATURES

- **High-Frequency Vibration:** The device utilizes ultrasonic technology with a vibration frequency of up to 40 kHz, reaching 2,000,000 vibrations per minute, for efficient plaque and tartar removal.
- **Safe & Painless Operation:** The cleaning head only vibrates when it touches hard objects (teeth). It automatically stops vibrating upon contact with gums or soft skin, ensuring safety and preventing discomfort.
- **4 Cleaning Modes:** Choose from Soft, Normal, Middle, and Strong modes to suit various sensitivity levels and cleaning needs.
- **Built-in LED Light:** An integrated LED light illuminates the oral cavity, providing clear visibility for

precise cleaning.

- **Portable & Rechargeable:** The device is compact and comes with a USB charging cable, offering up to 20 days of standby time on a full charge, making it ideal for travel.
- **IPX6 Waterproof:** The device has an IPX6 waterproof rating, allowing for direct washing with water for easy cleaning and maintenance.

FOUR Safe

Cleaning Modes



Soft

Gentle cleaning for beginners & sensitive teeth



Normal

For regular cleaning, oral hygiene



Middle

Deep cleaning for massive tartar and plaque



Strong

For stubborn calculus and hard tartar stains

Note: For beginners, we would recommend using the Soft Mode and light touch until you are used to the feel of using it



Image 4.1: Illustration of the four safe cleaning modes.

Electric Teeth Cleaning Kit

100% Safe for Gums

* Only work on teeth or hard objects.

Stop vibrating when touching gums or soft skin


2,000,000
Vibrations/Minute


40kHz
Vibration Frequency

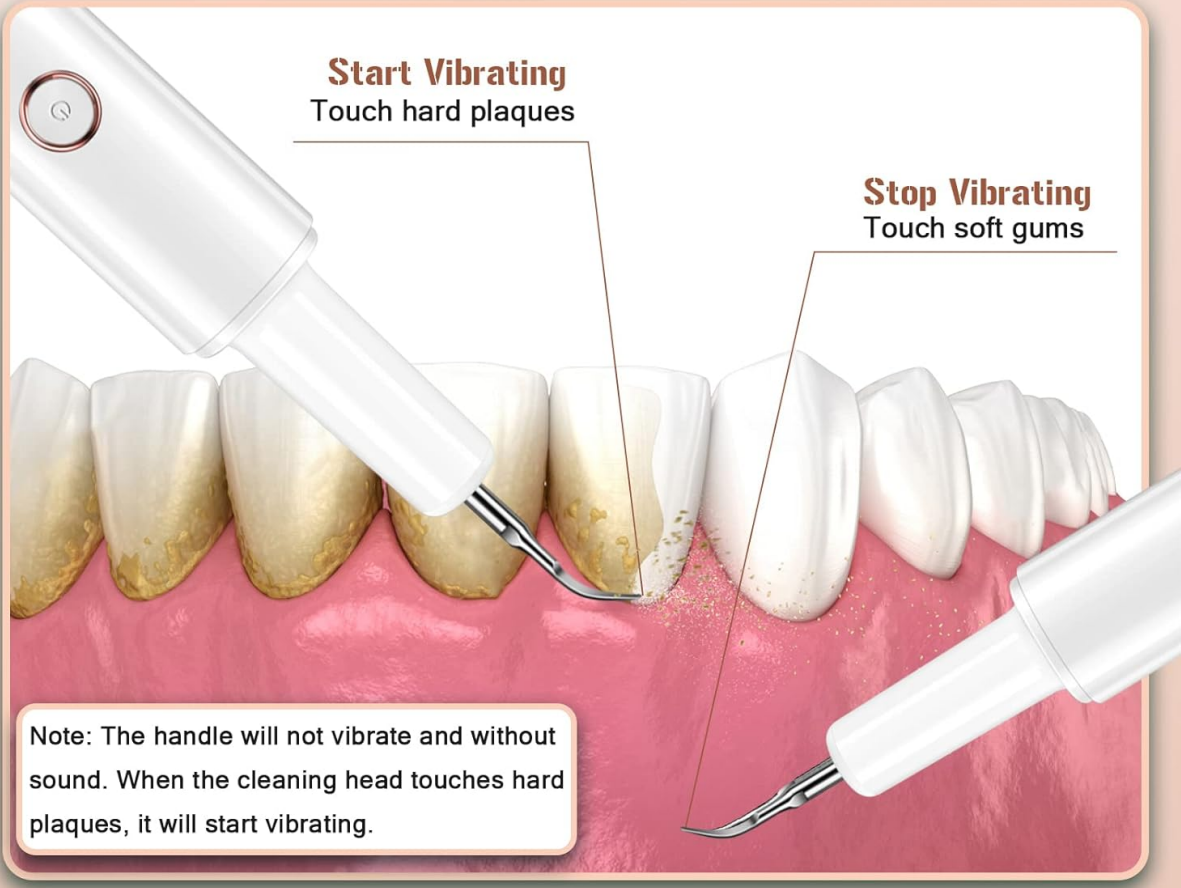


Image 4.2: The device's safety feature: stops vibrating on soft surfaces.

5. SETUP

1. **Charge the Device:** Before first use, fully charge the device using the provided USB charging cable. Connect the cable to the device and a compatible USB power source (e.g., power bank, adapter, laptop).
2. **Attach Cleaning Head:** Select the appropriate cleaning head (pointed for small calculus, flat for large calculus). Use the provided wrench to securely tighten the head onto the device by turning clockwise.

TWO Replacement Cleaning Heads



Tighten the head with the wrench



Pointed Head

Clean small calculus in the gums



Flat Head

Clean large calculus on the teeth surface

Image 5.1: Attaching the pointed or flat cleaning head using the wrench.

A Full Charge = 20 Days of Standby Time



IPX6 Waterproof

Image 5.2: USB charging and waterproof feature.

6. OPERATING INSTRUCTIONS

1. **Power On:** Press and hold the power button to turn on the device.
2. **Select Mode:** Press the power button repeatedly to cycle through the 4 modes: Soft, Normal, Middle, and Strong. For beginners or those with sensitive teeth, start with the Soft Mode.
3. **Prepare for Use:** Wet the tip of the device. You may apply a small amount of toothpaste if desired.
4. **Cleaning:** Place the tip against your teeth at a 90-degree angle. Move it gently along the gum line and over areas with plaque or stains. The device will vibrate when in contact with hard surfaces. Use the mouth mirror for better visibility of hard-to-reach areas.
5. **Auto Shutdown:** The device has an auto-shutdown function after 10 minutes of continuous use.
6. **Power Off:** Press and hold the power button to turn off the device after use.

Built-in LED Light

Achieve precise results



Image 6.1: Using the built-in LED light for precise cleaning.

Official Product Video: Ecardy Calculus Remover Model 2

Your browser does not support the video tag.

Video 6.1: Demonstration of the Ecardy Calculus Remover, showcasing its operation and effectiveness.

Official Product Video: Pearl Glide Water Flow Demo with Narration

Your browser does not support the video tag.

Video 6.2: A narrated demonstration of the Pearl Glide Water Flow feature, highlighting its precision and benefits.

7. MAINTENANCE

- **Cleaning the Device:** The device is IPX6 waterproof. After each use, remove the cleaning head and rinse both the device body and the cleaning head thoroughly under running water.
- **Storage:** Ensure the device and cleaning heads are completely dry before storing them in the protective cover. Store in a cool, dry place.
- **Replacing Heads:** Replace cleaning heads as needed, or if they show signs of wear.

8. TROUBLESHOOTING

- **Device not vibrating:** The device is designed to only vibrate when the cleaning head touches a hard object (like a tooth). It will not vibrate or make sound when touching soft tissue or when held in the air. Ensure the device is powered on and the cleaning head is in contact with a tooth.
- **Device not charging:** Ensure the USB cable is securely connected to both the device and a working power source. Check if the power source is active.
- **Ineffective cleaning:** Ensure you are using the appropriate cleaning mode for your needs. For stubborn calculus, use the Strong mode. Ensure the cleaning head is properly attached and not worn out.

9. SPECIFICATIONS

Brand	ocallcoo
Model Number	B09FS68W6P
Power Source	Battery Powered (1 Lithium Ion battery required)
Special Feature	Rechargeable
Waterproof Rating	IPX6
Vibration Frequency	Up to 40 kHz (2,000,000 Vibrations/Minute)
Modes	Soft, Normal, Middle, Strong
Package Dimensions	9.61 x 3.82 x 1.54 inches
Item Weight	11.68 ounces
UPC	755419323867

10. WARRANTY AND SUPPORT

For information regarding product warranty, returns, or technical support, please refer to the contact details provided with your purchase or visit the ocallcoo official website. The product is eligible for a refund/replacement within 2026 days of purchase, as per Amazon's return policy.