



Manuals.plus /

› KENZZI /

› KENZZI IPL Hair Removal Handset User Manual

KENZZI B09G166B9M

KENZZI IPL Hair Removal Handset User Manual

Model: B09G166B9M

1. INTRODUCTION

The KENZZI IPL Hair Removal Handset is designed for permanent hair reduction using Intense Pulsed Light (IPL) technology. This device offers a convenient and effective solution for hair removal in the comfort of your home. It is suitable for various body parts, including the upper lip, chin, underarms, legs, arms, and bikini area.

This manual provides essential information on how to safely and effectively use your KENZZI IPL Handset, including setup, operation, maintenance, and troubleshooting. Please read all instructions carefully before first use.



Image 1: The KENZZI IPL Hair Removal Handset, a compact white device with a gold-colored flash window.

2. How IPL TECHNOLOGY WORKS

Intense Pulsed Light (IPL) is a non-invasive procedure that uses high-intensity pulses of visible light. The light energy is absorbed by the melanin (pigment) in the hair follicle. This absorption converts light energy into heat, which damages the hair follicle and inhibits future hair growth. The process stimulates the follicles to go into a resting phase, leading to natural shedding of hair over time.



Image 2: Diagram showing IPL light penetrating the skin to hair follicles, stimulating them to enter a resting phase and inhibit growth.

3. SAFETY INFORMATION AND IMPORTANT WARNINGS

Your safety is paramount. Adhere to the following guidelines to ensure safe operation of the KENZZI IPL Handset.

3.1 General Safety Precautions

- Always ensure the device is in full contact with your skin before flashing. The handset has a built-in safety feature that prevents flashing in the air.
- Use the device in a well-lit room to minimize the perceived brightness of the flash.
- Avoid looking directly at the flash during use. While the device is designed with safety in mind, prolonged direct exposure should be avoided.
- Keep the device dry. Do not use near water or in wet environments.
- Keep out of reach of children.

3.2 Contraindications (Do Not Use If)

- You have tattoos in the treatment area.
- You have had significant sun exposure less than 48 hours prior to treatment.
- You have used retinol-based skin treatments less than 2 weeks prior to treatment.
- You have applied fake tans before use.
- You have medical equipment installed (e.g., pacemaker).
- You are pregnant or nursing.
- You have light sensitivity conditions.
- Your skin tone or hair color is incompatible with IPL technology (refer to Section 3.3).

3.3 Skin Tone and Hair Color Compatibility


IPL technology is most effective on fair to warm tan skin tones with dark hair. It works by targeting the melanin in the hair follicle.

For IPL to work effectively and safely, there needs to be a difference between your hair color and skin tone. Using it on dark skin tones with dark hair, or light hair with lighter skin tones, can render the device unsafe and ineffective.

Skin Tone Chart

IPL is suitable for all **fair** to **warm tan** skin tones. Refer to the skin tone chart below to see whether or not IPL is suitable for you.

FAIR TO LIGHT TO FAIRLY MEDIUM

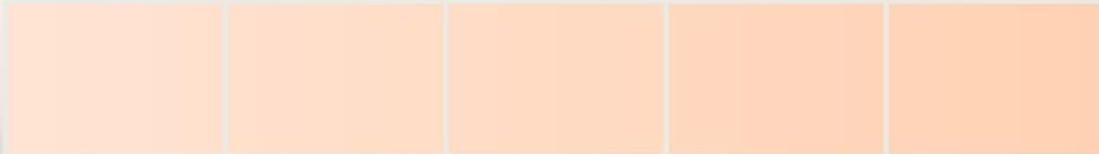

FAIRLY MEDIUM TO MEDIUM



MEDIUM TAN TO WARM TAN



For IPL to work effectively, there needs to be a difference between your hair color and skin tone. Using it on **dark hair** and **darker skin tones** and **light hair** and **lighter skin tones** can render this device unsafe and ineffective.

Image 3: Skin Tone Chart illustrating compatibility for IPL treatment. IPL is suitable for fair to warm tan skin tones. Darker skin tones or very light hair may not be suitable.

4. PRODUCT COMPONENTS

The KENZZI IPL Hair Removal Handset typically includes:

- KENZZI IPL Handset
- Power Adapter
- User Manual (this document)

Designed for Safety, Efficiency & Luxury

KENZZI is salon expertise at the palm of your hand, in the comfort of your home.



Image 4: The KENZZI IPL Handset connected to its power adapter.

5. SETUP

1. **Unpack:** Carefully remove the KENZZI IPL Handset and all accessories from the packaging.
2. **Inspect:** Check the device and power adapter for any visible damage. Do not use if damaged.
3. **Connect Power:** Plug the power adapter into the handset, then plug the adapter into a suitable electrical outlet. The device will power on.

6. OPERATING INSTRUCTIONS

6.1 Preparation Before Treatment

1. **Clean Skin:** Ensure the treatment area is clean, dry, and free from any lotions, oils, or deodorants.
2. **Shave Hair:** Shave the hair in the treatment area. Do not wax, pluck, or epilate, as IPL requires the hair follicle to be present beneath the skin surface.

3. **Perform Skin Test:** Before your first full treatment, perform a skin test on a small area of skin you intend to treat. Start with the lowest intensity level (Level 1) and gradually increase if comfortable, up to the highest setting that causes no adverse reactions. Wait 24 hours to check for any skin reactions.

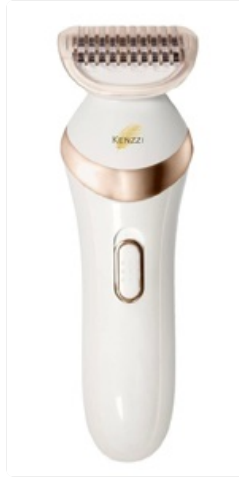


Image 5: Exfoliate skin before treatment.



Image 6: Shave the treatment area.

6.2 Performing the Treatment

1. **Select Intensity Level:** The KENZZI IPL Handset features 5 intensity levels. Adjust the intensity based on your comfort level and skin test results. Start with Level 1 and increase as tolerated.
2. **Position Device:** Place the flash window firmly against your skin in the desired treatment area. Ensure full contact for the device to flash.
3. **Initiate Flash:** Press the flash button. You will see a bright flash of light.
4. **Move to Next Area:** Lift the device and move it to the next adjacent area, ensuring slight overlap to cover all hair.
5. **Glide Mode:** For larger areas, the device may have a glide mode for automatic flashing when continuously moved across the skin. Refer to the device's specific button functions for activation.
6. **Complete Treatment:** Continue until the entire desired area has been treated.

Designed for a Lasting Luxurious Comfort

- ✓ Professional Salon Technology
- ✓ 500,000 Flashes
- ✓ Can Use Up to 10 Years
- ✓ Use from Head to Toe
- ✓ No Replacements or Refills Needed
- ✓ No Goggles or Sunglasses Required

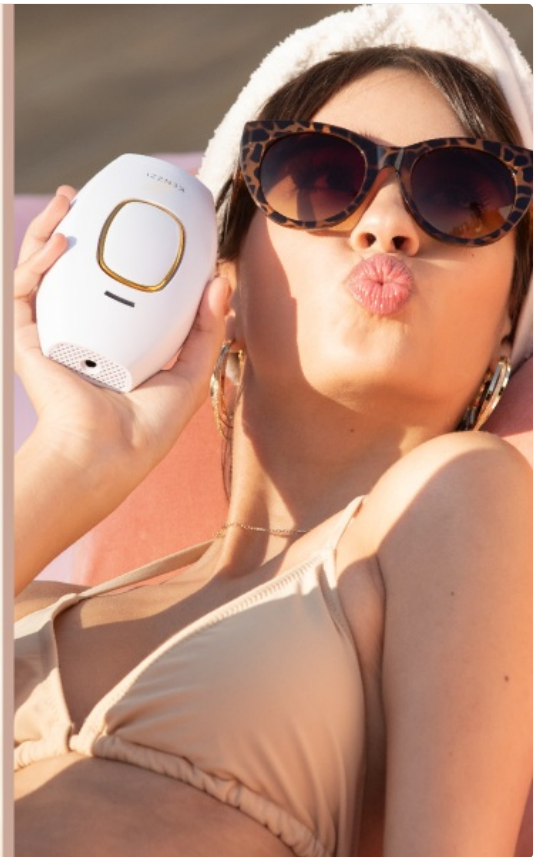


Image 7: Explanation of the 5-level intensity feature and safety precautions for the KENZZI IPL Handset.



Image 8: Applying the KENZZI IPL Handset to the skin.

6.3 Applicable Treatment Areas

The KENZZI IPL Handset is safe to use on various body parts, including:

- Legs
- Chest
- Back
- Upper lip
- Chin
- Underarms
- Arms
- Bikini line



FDA Cleared IPL Technology

KENZZI is safe to use on:



Legs



Chest



Back



Upper lip



Chin



Underarms



Arms



Bikini line

Image 9: Visual guide to body areas where the KENZZI IPL Handset can be safely used.

6.4 Post-Treatment Care

- Apply a soothing lotion or moisturizer to the treated area if desired.
- Avoid direct sun exposure on treated areas for at least 48 hours. Use sunscreen if sun exposure is unavoidable.



Image 10: Apply soothing lotion after treatment.

6.5 Treatment Schedule and Expected Results

For optimal results, it is recommended to use the KENZZI IPL Handset once a week for 12 weeks. Visible hair reduction can typically be observed within 3-4 weeks. Full benefits are usually achieved between 10-12 weeks. After the initial 12-week period, touch-up treatments can be performed every other month as needed to maintain results.

Before & After

Results You Can See & Feel

We recommend using KENZZI once a week for 12 weeks to achieve significant hair reduction. Results may vary based on individual hair regrowth and across different body areas.



3-4 weeks to see visible results



6-8 weeks



10-12 weeks for full benefits

Image 11: Visual representation of hair reduction progress over 3-4 weeks, 6-8 weeks, and 10-12 weeks of consistent use.

7. MAINTENANCE

7.1 Cleaning the Device

- Ensure the device is unplugged and cooled down before cleaning.
- Wipe the flash window and the exterior of the device with a soft, dry, lint-free cloth.
- Do not use abrasive cleaners, alcohol, or harsh chemicals, as these can damage the device.

7.2 Storage

- Store the device in a cool, dry place, away from direct sunlight and moisture.
- Keep the device in its original packaging or a protective case to prevent damage.

8. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Not properly plugged in; power outlet issue.	Ensure the power adapter is securely connected to both the device and a working electrical outlet. Try a different outlet.
Device does not flash.	Flash window not in full contact with skin; device not ready.	Ensure the flash window is pressed firmly and completely against the skin. Wait for the "Ready" indicator light (if applicable, typically bluish-pink) to appear before attempting to flash. If the light remains green, reposition the device.
Skin irritation or discomfort.	Intensity level too high; sensitive skin.	Reduce the intensity level. Ensure skin is properly prepared (shaved, clean, dry). If irritation persists, discontinue use and consult a healthcare professional.
No visible hair reduction after several weeks.	Inconsistent use; incompatible skin/hair type; incorrect technique.	Ensure consistent weekly treatments. Verify your skin tone and hair color are compatible (refer to Section 3.3). Ensure the entire area is treated with slight overlap. Results vary by individual.

9. SPECIFICATIONS

- **Model:** B09G166B9M
- **Brand:** KENZZI
- **Technology:** Intense Pulsed Light (IPL)
- **Intensity Levels:** 5
- **Flashes:** Up to 500,000 flashes (estimated 10+ years of treatments)
- **Product Dimensions:** 9.06 x 5.83 x 3.86 inches
- **Weight:** 1.54 Pounds
- **Manufacturer:** KENZZI

10. WARRANTY AND SUPPORT

For warranty information and customer support, please refer to the documentation included with your purchase or visit the official KENZZI website.

If you experience any issues or have questions regarding your KENZZI IPL Hair Removal Handset, please contact KENZZI customer service for assistance.

