

## THRIVE MD-450(IV)

# THRIVE MD-450 Kneading Massager User Manual

Model: MD-450 (Ivory)

## 1. INTRODUCTION

Thank you for purchasing the THRIVE MD-450 Kneading Massager. This device is a controlled medical device designed to provide a firm kneading massage to various parts of your body, including the neck, shoulders, waist, and thighs. Please read this manual thoroughly before use to ensure safe and effective operation.



Figure 1.1: THRIVE MD-450 Kneading Massager (Ivory)

## 2. SAFETY PRECAUTIONS

**WARNING: To reduce the risk of electric shock, burns, fire, or injury to persons:**

- Always unplug the massager from the electrical outlet immediately after using and before cleaning.
- Do not use while bathing or in a shower.
- Do not place or store appliance where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Never use pins or other metallic fasteners with this appliance.
- Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been

dropped or damaged, or dropped into water.

- Keep the cord away from heated surfaces.
- Never use while sleeping or drowsy.
- Never drop or insert any object into any opening.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not use on sensitive skin areas or where circulation is poor.
- This product is not a toy. Keep out of reach of children.

### 3. PRODUCT FEATURES

The THRIVE MD-450 Kneading Massager is equipped with advanced features to provide a comprehensive massage experience:

- **Six Kneading Balls:** Features six massage balls (three upper, three lower) that operate simultaneously to firmly grasp and relieve muscle stiffness.
- **Targeted Massage Areas:** Effectively massages the neck, shoulders, waist, and thighs.
- **Adjustable Grip and Pressure:** The direction of the grips and the applied pressure can be adjusted to find the optimal massage position.
- **Multiple Modes:** Includes two automatic modes and a continuous mode with two speed settings (Fast and Slow).
- **Auto-Timer:** Equipped with a 10-minute auto-timer for safety and to prevent overuse.

# THRIVE

つかみもみマッサージャー MD-450

6つのもみボールが  
首・肩・腰・太もものコリを  
捕らえてマッサージ

メーカー  
**1年保証**

医療機器認証

Figure 3.1: Targeted massage for neck and shoulders.

# つかみもみマッサージャー

上もみボールと下もみボールが同時に動作する  
つかみもみマッサージ。  
グリップの向きや力加減で最適な位置でほぐせます。

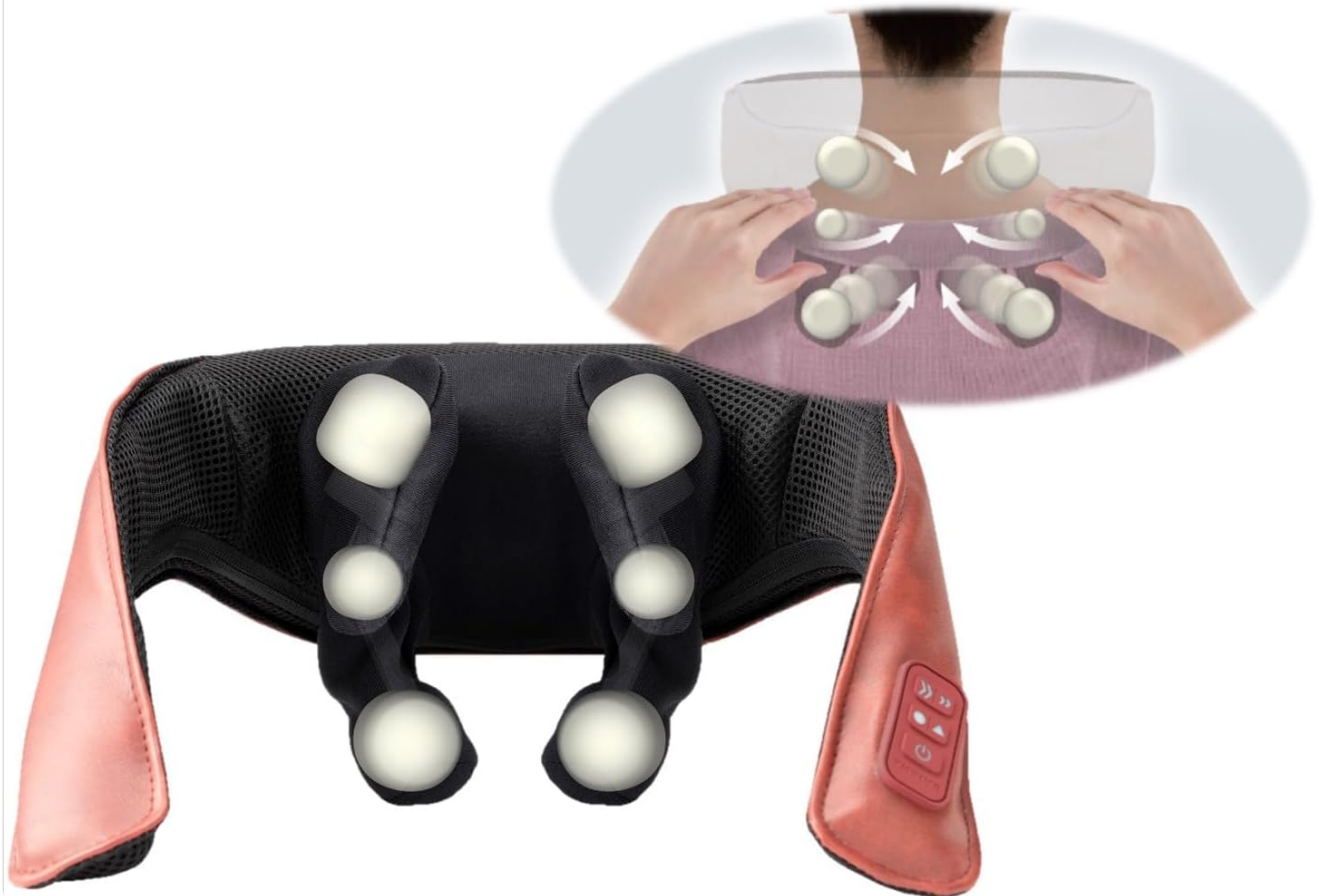


Figure 3.2: Illustration of the six kneading balls in action.

## 4. SETUP

The THRIVE MD-450 is designed for easy setup:

1. Unpack the massager and all its components from the box.
2. Ensure the power cord is free from kinks or damage.
3. Plug the power cord into a standard AC100V 50/60Hz electrical outlet.
4. The device is now ready for use.

## 5. OPERATING INSTRUCTIONS

Follow these steps to operate your THRIVE MD-450 Kneading Massager:

### 5.1 Control Panel

The control panel is located on one of the arm straps for easy access during use.



## 自動モード 2種類

速度の切り替えと停止を、それぞれ異なる動きでマッサージします。

### 連続モード 【速】【遅】

2段階で速度が切り替わり、一定の速度で動作します。

使いすぎや切り忘れを防止する10分オートタイマー付き。



●自動Ⅰ ▲自動Ⅱ >>速 >>遅

Figure 5.1: Control Panel Overview

- **Power Button** ( ): Press to turn the massager ON or OFF.
- **Mode Button** (▲): Cycles through the two automatic modes (Automatic I and Automatic II).
- **Speed Buttons** (»): Adjusts the speed in continuous mode between Fast (» Fast) and Slow (» Slow).

### 5.2 Usage for Different Body Parts

The massager can be used on various parts of the body. Adjust the position and hold the grips as shown:



首・肩、腰、太ももをギュッとつかみもみ。



Figure 5.2: Usage examples for neck, waist, and thighs.

- **Neck/Shoulders:** Place the massager around your neck and shoulders, holding the arm straps to adjust pressure and position.
- **Waist:** Position the massager on your lower back or waist area, using the arm straps to secure it.
- **Thighs:** Place the massager on your thighs. You can hold the main unit or the straps to keep it in place.

The massager will automatically turn off after approximately 10 minutes of continuous operation to prevent overuse.

## 6. MAINTENANCE

Proper maintenance ensures the longevity and hygiene of your massager:

- **Cleaning:** Always unplug the massager before cleaning. Wipe the surface with a soft, slightly damp cloth. Do not use abrasive cleaners or immerse the unit in water.
- **Storage:** Store the massager in a cool, dry place away from direct sunlight and extreme temperatures. Avoid placing heavy objects on top of the massager.
- **Cord Care:** Do not wrap the power cord tightly around the unit. Avoid twisting or bending the cord excessively.

## 7. TROUBLESHOOTING

If you encounter any issues with your massager, please refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Massager does not turn on.	Not plugged in correctly. Power outlet issue. Damaged cord/plug.	Ensure the power cord is securely plugged into a working outlet. Try a different outlet. Do not use if cord or plug is damaged; contact customer support.
Massage feels weak or inconsistent.	Incorrect positioning. Overuse protection activated.	Adjust the massager's position and the grip pressure. The unit may have automatically shut off after 10 minutes. Wait a few minutes before restarting.
Unusual noise during operation.	Internal malfunction.	Turn off the unit immediately and unplug it. Do not attempt to repair. Contact customer support.

If the problem persists after trying these solutions, please contact our customer support center.

## 8. SPECIFICATIONS

Feature	Detail
Model Number	MD-450(IV)
Product Dimensions	37 x 35.5 x 15.5 cm (Length x Width x Height)
Product Weight	1.6 kg
Color	Ivory
Rated Power	AC100V 50/60Hz 20W
Auto-Timer	20 minutes (internal timer, though feature bullets mention 10 min auto-off for safety)
Kneading Speed (Approx.)	Fast: 15 times/minute, Slow: 7 times/minute
Materials	Synthetic leather, Polyester, Nylon
Country of Origin	China
Water Resistance Level	Non-waterproof

## 9. WARRANTY AND SUPPORT

Daito Electric Industry Co., Ltd. is committed to ensuring your peace of mind when using THRIVE products.





Figure 9.1: THRIVE MD-450 in different colors.

## 9.1 Warranty Information

This product comes with a **1-year manufacturer's warranty** from the date of purchase. Please retain your proof of purchase for warranty claims. The warranty covers defects in materials and workmanship under normal use.

## 9.2 Customer Support

For any inquiries, repairs, or support regarding your THRIVE product, please contact our Customer Consultation Center:

**Phone Number:** 0570-064-762

**Operating Hours:** 9:00 AM - 5:00 PM JST (Monday - Friday, excluding holidays)


Calls from general phones are charged at local call rates.








Figure 9.2: Warranty and Support Contact Information.

© 2024 Daito Electric Industry Co., Ltd. All rights reserved.

## Related Documents - MD-450(IV)

	<p><a href="#">Medcursor MD-LM01 Air Compression Leg Massager with Heat - User Manual</a></p> <p>Comprehensive user manual for the Medcursor MD-LM01 Air Compression Leg Massager with Heat. Learn about product features, safety instructions, operation, troubleshooting, and care.</p>
---	---

 <p>The image shows the cover of the Medcursor MD-KP02 Knee Massager with Heat User Manual. It features the Medcursor logo at the top, followed by the title 'USER MANUAL' and 'Knee Massager with Heat'. Below the title is a line drawing of the knee massager. At the bottom, there is a green banner with the text 'READ THIS MANUAL Thoroughly before using and save it for future reference.'</p>	<p><a href="#">Medcursor MD-KP02 Knee Massager with Heat User Manual - Operation, Care, Troubleshooting</a></p> <p>Comprehensive user manual for the Medcursor MD-KP02 Knee Massager with Heat. Learn about operation, safety precautions, care, maintenance, troubleshooting, and warranty information for this cordless massager.</p>
 <p>The image shows the cover of the Medcursor Foot Massager User Manual. It features the Medcursor logo at the top, followed by the title 'MEDCURSOR FOOT MASSAGER' and 'USER MANUAL'. Below the title is a photograph of the foot massager. At the bottom, there is a green banner with the text 'READ THIS MANUAL Thoroughly before using and save it for future reference.'</p>	<p><a href="#">Medcursor Foot Massager User Manual - Model MD-52230</a></p> <p>Comprehensive user manual for the Medcursor Foot Massager (Model MD-52230). Learn how to operate, maintain, and troubleshoot your device. Includes safety warnings and warranty information.</p>
 <p>The image shows the cover of the NURSAL MD-82103A-3 Massage Cushion: Operating Guide and Features. It features a photograph of the massage cushion and a diagram showing its various massage modes.</p>	<p><a href="#">NURSAL MD-82103A-3 Massage Cushion: Operating Guide and Features</a></p> <p>Comprehensive operating guide for the NURSAL MD-82103A-3 Massage Cushion, detailing features, safety precautions, operating procedures, and technical specifications for effective use.</p>
 <p>The image shows the cover of the Medcursor MD-HM02 Hand Massager User Manual. It features the Medcursor logo at the top, followed by the title 'USER MANUAL' and 'HAND MASSAGER'. Below the title is a photograph of the hand massager. At the bottom, there is a green banner with the text 'READ THIS MANUAL Thoroughly before using and save it for future reference.'</p>	<p><a href="#">Medcursor MD-HM02 Hand Massager User Manual</a></p> <p>Comprehensive user manual for the Medcursor MD-HM02 Hand Massager, covering important safeguards, operation instructions, care, maintenance, troubleshooting, and warranty information.</p>
 <p>The image shows the cover of the DR-HO'S Neck Pain Pro User Manual and Pain Relief Guide. It features the 'neck pain pro' logo at the top, followed by a photograph of a man wearing the neck pain pro device. Below the photograph is a green banner with the text 'READ THIS MANUAL Thoroughly before using and save it for future reference.'</p>	<p><a href="#">DR-HO'S Neck Pain Pro User Manual and Pain Relief Guide</a></p> <p>Comprehensive guide to the DR-HO'S Neck Pain Pro system, detailing its features, setup, usage instructions, pad placement for various body areas, and important safety information for effective pain relief.</p>