

Oclean FlowBlue

Oclean Flow Electric Sonic Toothbrush User Manual

Model: FlowBlue

INTRODUCTION

Thank you for choosing the Oclean Flow Electric Sonic Toothbrush. This manual provides detailed instructions for the safe and effective use of your new toothbrush. Please read it thoroughly before first use and keep it for future reference.





The Oclean Flow Electric Sonic Toothbrush, shown with its sleek blue packaging.

WHAT'S IN THE BOX

Upon unboxing your Oclean Flow Electric Sonic Toothbrush, please ensure all components are present:

- 1 x Oclean Flow Electric Sonic Toothbrush Handle
- 1 x Brush Head
- 1 x Charging Cable
- 1 x User Manual

ES GIBT KEINEN PLANETEN B

So halten wir das Oclean Flow Paket kompakt, wiederverwendbar und wiederverwertbar

Die Oclean Flow Box enthält:



All items included in the Oclean Flow Electric Sonic Toothbrush package.

SETUP AND FIRST USE

Follow these steps to set up your Oclean Flow toothbrush for the first time:

1. Remove the toothbrush and all accessories from the packaging.
2. Before initial use, fully charge the toothbrush. Connect the charging cable to the toothbrush and a suitable power source.
3. Moisten the brush head and apply a small amount of toothpaste.
4. Place the brush head in your mouth before turning on the device to prevent splashing.



The Oclean Flow toothbrush connected to its charging cable.

OPERATING INSTRUCTIONS

The Oclean Flow toothbrush is designed for effective and convenient oral care.

Brushing Procedure:

1. Press the power button to start brushing.
2. Brush for two minutes. The smart timer will remind you every 30 seconds to change the brushing area, ensuring all quadrants of your mouth are cleaned thoroughly.
3. After brushing, rinse the toothbrush head under running water and allow it to air dry.



Visual representation of the 2-minute smart timer and 30-second quadrant reminders.

Brushing Modes:

The Oclean Flow offers 5 distinct cleaning modes to suit various oral care needs. You can switch between modes by pressing the mode button on the handle.

- **Morning Mode:** Gentle intensity for a fresh start.
- **Night Mode:** Higher intensity for thorough cleaning before bed.
- **Standard Mode:** Balanced cleaning for daily use.
- **Whitening Mode:** Designed to help remove surface stains.
- **Gentle Mode:** Very soft cleaning, ideal for sensitive teeth and gums.

5 BÜRSTENMODI



Morgens

Intensität des Bürstens: ★★★



Nacht

Intensität des Bürstens: ★★★★★



Standard Reinigung

Intensität des Bürstens: ★★★★★



Aufhellen

Intensität des Bürstens: ★★★★★



Sanft

Intensität des Bürstens: ★★



Detailed overview of the 5 available brushing modes and their intensity levels.

Pressure Sensor:


The Oclean Flow is equipped with a pressure sensor that warns you if you apply too much pressure while brushing. This helps protect your gums and tooth enamel from damage.

Suitable for Braces:

The Oclean Flow Sonic Toothbrush is safe for use with braces, fillings, crowns, and veneers. Its gentle yet effective cleaning action helps maintain oral hygiene without causing damage.

GEEIGNET FÜR MENSCHEN MIT ZAHNSPANGE

Machen Sie Ihr Zähneputzen zu einer Reinigungsmassage, ohne Zahnschmelz und Zahnfleisch zu verletzen

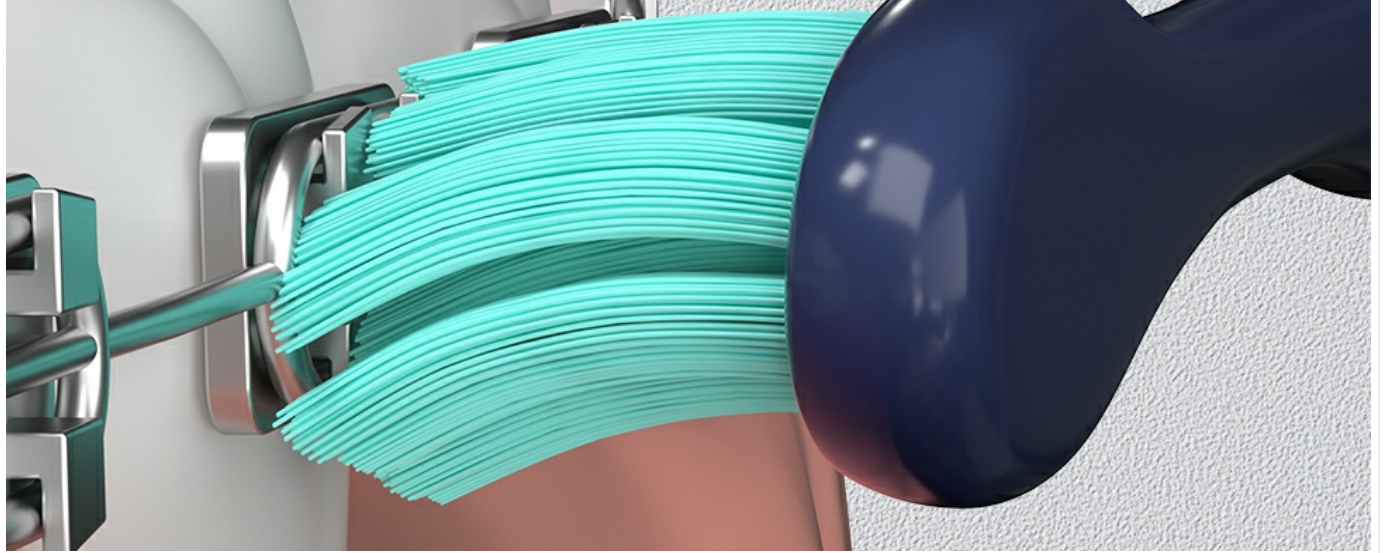

The Oclean Flow toothbrush is designed to effectively clean around braces without causing discomfort or damage.


BATTERY LIFE AND CHARGING

The Oclean Flow boasts an impressive battery life, making it ideal for travel and daily use without frequent recharging.

- **180-Day Battery Life:** A full charge provides up to 180 days of continuous use, requiring only two charges per year.
- **Charging Time:** Approximately 5 hours for a full charge.

SUPER LANGE AKKULAUFZEIT

180 Tage Akkulaufzeit

1  Jahr = **2**  Aufladungen



The Oclean Flow's exceptional 180-day battery life, requiring minimal charging.



oclean

1 = **180**
volle Ladung = Tage Nutzung



A visual summary of the Oclean Flow's long-lasting battery performance.

**180 Tage lange Batterielebensdauer
Ideal für Reisen**



The Oclean Flow's extended battery life makes it an ideal companion for travel.

BRUSH HEAD AND CLEANING PERFORMANCE

The Oclean Flow utilizes advanced technology for superior cleaning.

- **DuPont Brush Heads:** Equipped with precision bristles that penetrate deeper into interdental spaces and cover more surface area for a cleaner mouth than manual toothbrushes.
- **Magnetic Levitation Motor:** Delivers up to 76,000 movements per minute (38,000 vibrations in both left and right directions) for fast and effective plaque removal.
- **Gentle Care:** Despite its powerful cleaning, the toothbrush is designed to gently care for gums, making it ideal for sensitive teeth.

PATENTIERTER MAGNETSCHWEBEBAHNMOTOR


bis zu

76,000

Bewegungen pro Minute



The patented magnetic levitation motor ensures powerful and efficient cleaning.

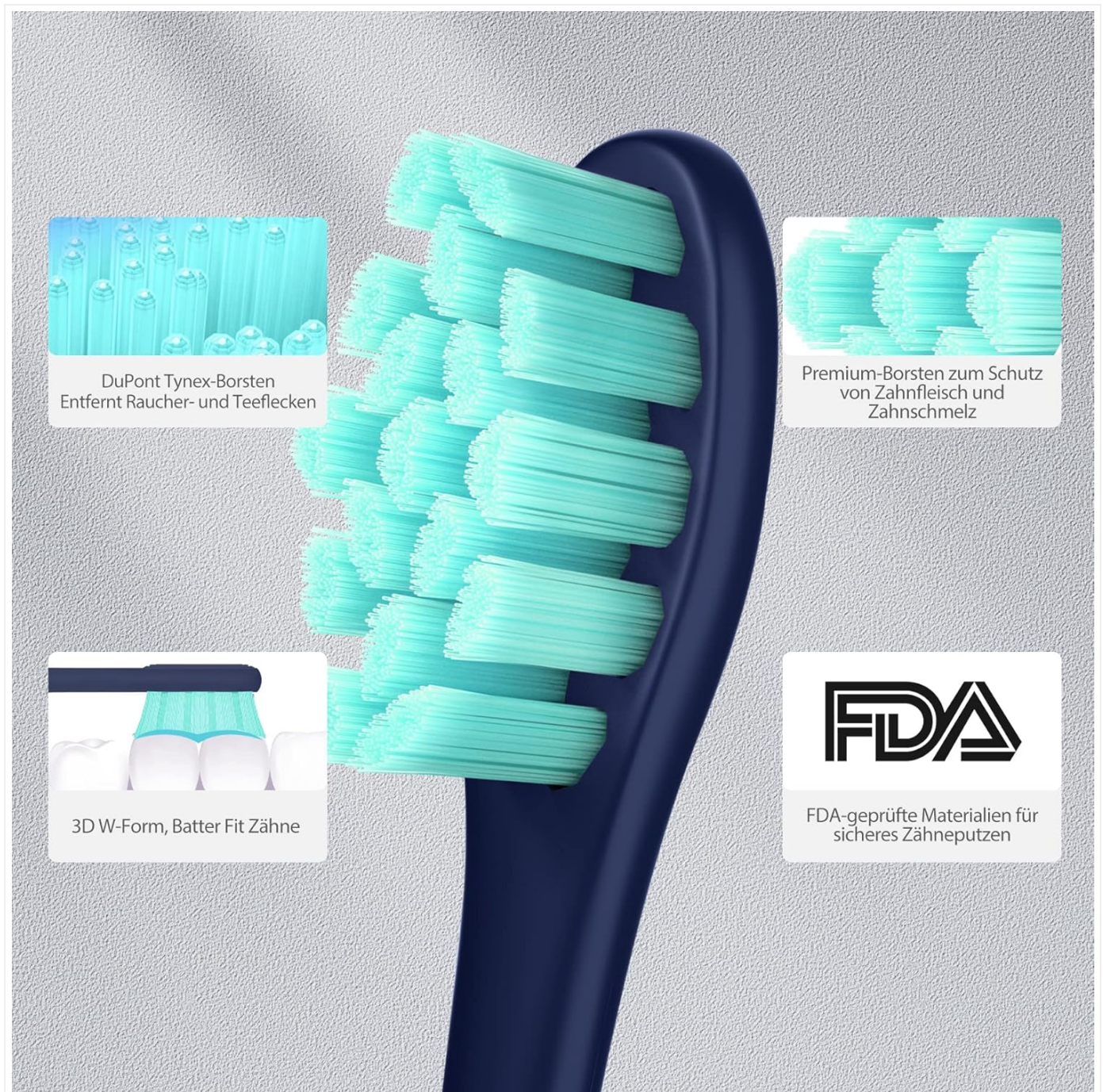

Key features of the DuPont brush head, designed for effective and safe cleaning.

4 WOCHEN ANWENDUNG



Vor

Nachher



Visual demonstration of the whitening effect after 4 weeks of consistent use.

MAINTENANCE

Proper maintenance ensures the longevity and hygiene of your Oclean Flow toothbrush.

- **Daily Cleaning:** After each use, rinse the brush head and handle thoroughly under running water.
- **Deep Cleaning:** Periodically remove the brush head and clean the area where it connects to the handle to prevent toothpaste residue buildup.
- **Brush Head Replacement:** For optimal hygiene and cleaning performance, replace the brush head every 3 months, or sooner if the bristles become frayed.
- **Storage:** Store the toothbrush in a dry place, away from direct sunlight.

Waterproof Rating:

The Oclean Flow is IPX7 waterproof, meaning it can be submerged in up to 1 meter of water for 30 minutes without damage. This allows for easy cleaning under running water and safe use in the shower.

IPX7 WATERPROOF



The Oclean Flow's IPX7 waterproof design allows for easy cleaning and safe use in wet environments.

TROUBLESHOOTING

If you encounter any issues with your Oclean Flow toothbrush, please refer to the following general advice:

- **Toothbrush not turning on:** Ensure the toothbrush is fully charged.
- **Weak vibrations:** Recharge the toothbrush. If the issue persists, the battery may be nearing the end of its life cycle.
- **Unusual noise:** Ensure the brush head is securely attached. If the noise continues, contact customer support.

For more complex issues or if troubleshooting steps do not resolve the problem, please contact Oclean customer support.

SPECIFICATIONS

| Feature | Specification |
|---------|---------------|
|---------|---------------|

| | |
|-----------------------|---|
| Product Name | Sonic Electric Toothbrush |
| Model Number | FlowBlue |
| Brand | Oclean |
| Battery Life | Up to 180 days |
| Charging Time | Approx. 5 hours |
| Brushing Modes | 5 (Morning, Night, Standard, Whitening, Gentle) |
| Vibrations per Minute | Up to 76,000 |
| Waterproof Rating | IPX7 |
| Dimensions | 2.5 x 2.5 x 24.5 cm |
| Weight | 135 grams |
| Power Source | Battery-powered, Electric |
| Age Range | Adults; Teenagers |

| VERGLEICH MIT ANDEREN | | |
|-----------------------|----------------------|---|
| | Oclean Flow | Others |
| Batterielebensdauer | 180 Tage | Etwa 14 Tage |
| Reinigungsleistung | 76000 Bewegungen/min | 17.600~62.000 Bewegungen/min |
| Modus | 5 | 1~3 |
| Motor | Oclean Maglev Motor | Rotationsmotor Elektromotor Hohlbechermotor |
| Batterie | 2500mAh | 400~600mAh |

A comparison of Oclean Flow's specifications against other electric toothbrushes.

Safety Warnings

Please read and adhere to the following safety warnings to prevent injury or damage to the product:

- Use the product strictly according to the instructions provided in this manual.
- Do not disassemble, repair, or modify the device yourself. Contact authorized service personnel for any repairs.

- Keep the toothbrush out of reach of children and individuals with reduced physical, sensory, or mental capabilities unless supervised by a responsible adult.
- Do not use the toothbrush if it is damaged or malfunctioning.
- Avoid exposing the charging port to water when not in use or during charging. Ensure the charging port cover is securely closed.
- If you experience any discomfort or pain while using the toothbrush, discontinue use and consult a dental professional.
- Only use Oclean original or recommended replacement brush heads.

WARRANTY AND SUPPORT

Oclean is committed to providing high-quality products and services. Your Oclean Flow Electric Sonic Toothbrush comes with a warranty.

- **Warranty:** Oclean will provide a free replacement if there is a quality problem within 2 years from the date of purchase.
- **Customer Support:** If you have any questions, concerns, or require assistance with your product, please do not hesitate to contact Oclean customer support.

For contact information, please refer to the official Oclean website or the contact details provided in your product packaging.