

COMFIER CF-5001-US

COMFIER CF-5001-US 2-in-1 Foot and Calf Massager Instruction Manual

MODEL: CF-5001-US

1. Product Overview

The COMFIER CF-5001-US is a versatile 2-in-1 foot and calf massager designed to provide a comprehensive massage experience. It features Shiatsu, compression, vibration, and heat functions. The unique design allows it to function as both a foot and calf massager and a heated ottoman, integrating seamlessly into your living space.

Key Features:

- **Dual Functionality:** Transforms from a foot and calf massager into a heated ottoman (19.6" x 15.3").
- **Massage Modes:** Offers Shiatsu, compression, vibration, and gentle heat for feet, ankles, and calves.
- **Adjustable Settings:** Three compression intensity options and two leg massage modes.
- **Ergonomic Design:** Features an adjustable support bar for a 20-degree inclined massage angle.
- **Hygiene:** Includes removable and washable foot covers for easy cleaning.
- **Foot Size Compatibility:** Accommodates foot sizes up to men's size 13.
- **Automatic Shut-off:** Automatically powers off after 15 minutes of use.

Multi Foot Massage Machine

Footrest for chair or high bed
(heated Ottoman design)



Footstool, Sofa Tea Stool
(Provide warm seat)



Foot Massager
(Shiatsu, compression,
Rolling, vibration)

Figure 1: The COMFIER massager shown in its compact ottoman form and expanded for foot and calf massage.

2. Safety Information

Please read all instructions and safety warnings before using this appliance to prevent injury or damage. Keep this manual for future reference.

- Always unplug the appliance from the electrical outlet immediately after use and before cleaning.
- Do not use while bathing or in a shower.
- Do not place or store appliance where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Never use pins or other metallic fasteners with this appliance.
- Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.

- Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
- Keep the cord away from heated surfaces.
- Do not use outdoors.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- To disconnect, turn all controls to the “OFF” position, then remove plug from outlet.
- This appliance is designed for household use only.
- Consult a doctor before use if you have certain medical conditions, including but not limited to pregnancy, diabetes, heart disease, or if you have a pacemaker.

3. Product Components

The COMFIER CF-5001-US massager includes the main unit, which functions as both an ottoman and a massager, and a power adapter. The massager features integrated controls for ease of use.

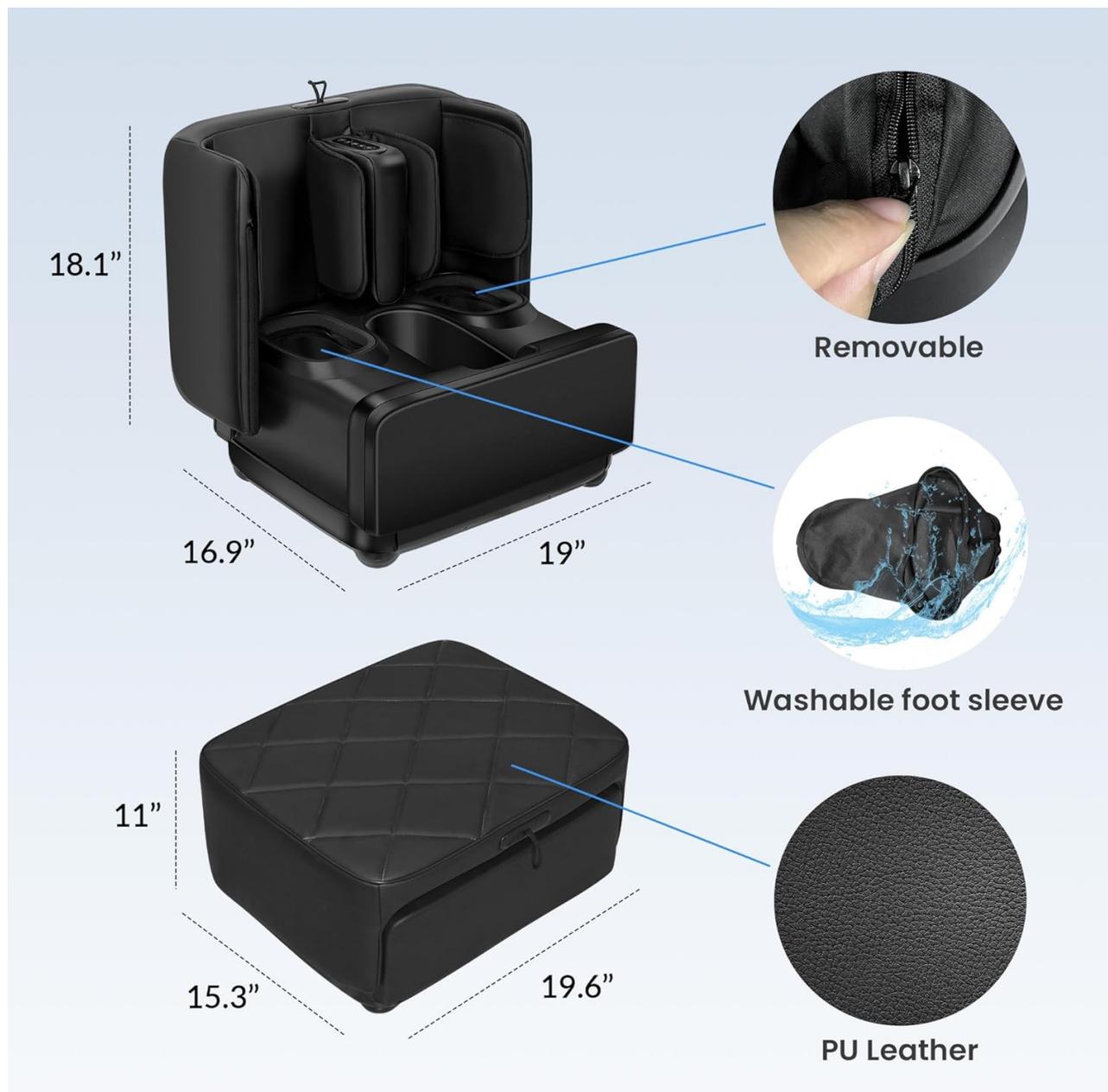


Figure 2: Key components and dimensions of the massager, highlighting the removable foot sleeve and PU leather construction.

4. Setup Instructions

1. **Unpack:** Carefully remove the massager from its packaging.
2. **Positioning:** Place the massager on a flat, stable surface near an electrical outlet.
3. **Power Connection:** Connect the power adapter to the massager's power port, typically located on the side or rear of the unit. Plug the adapter into a suitable electrical outlet.
4. **Convert to Massager Mode:** If the unit is in ottoman mode, lift the top cover to reveal the foot and calf massage compartments. Ensure the cover is securely positioned to support the calf massage section.
5. **Adjust Angle (Optional):** Pull up the support bar at the bottom of the massager to achieve a 20-degree inclined angle for enhanced comfort during use.



Figure 3: The adjustable support bar allows for a comfortable 20-degree inclined massage angle.

5. Operating Instructions

The control panel is conveniently located on the top of the massager for easy access.

1. **Power On/Off:** Press the **On/Off Button** to start or stop the massager.
2. **Foot Shiatsu:** Press the **Foot Shiatsu Button** to activate or change the foot massage mode (M1/M2/M3).
3. **Leg Mode:** Press the **Leg Mode Button** to activate or change the leg massage mode (M1/M2/Turn off).
4. **Heat Function:** Press the **Heat Button** to turn on the gentle heating function. The heat function provides soothing warmth to your feet and calves.
5. **Pressure Intensity:** Press the **Pressure Intensity Button** to adjust the compression intensity (L1/L2/L3).
6. **Placement:** Sit comfortably and place your feet into the foot compartments and your calves into the leg compartments.
7. **Automatic Timer:** The massager will automatically shut off after 15 minutes.

CUSTOMIZABLE CALF COMPRESSION

3 Intensity Levels for Perfect Fit



Best for 32-42 cm calves



Figure 4: The control panel allows selection of various massage modes and intensity levels.

Official Product Video: COMFIER Shiatsu Foot Massage Machine

Video 1: This video demonstrates the operation and features of the COMFIER Shiatsu Foot Massage Machine, including its transformation and various massage functions.

6. Maintenance and Care

- **Cleaning:** The foot covers are removable and washable. Unzip the covers from the foot compartments and hand wash them with mild detergent. Allow them to air dry completely before reattaching.
- **Exterior Cleaning:** Wipe the exterior of the massager with a soft, damp cloth. Do not use abrasive cleaners or solvents.
- **Storage:** When not in use, fold the massager back into its ottoman form for compact storage. Store in a cool, dry place away from direct sunlight.
- **Avoid Liquids:** Do not immerse the massager in water or any other liquid.

7. Troubleshooting

Problem	Possible Cause	Solution
Massager does not turn on.	Power cord not securely plugged in.	Ensure the power cord is firmly connected to both the massager and the electrical outlet.
No massage function.	Specific massage mode not activated.	Press the corresponding button (Foot Shiatsu, Leg Mode) to activate the desired massage function.
Heat function not working.	Heat function not activated.	Press the Heat Button to turn on the heat. Note that the heat is designed to be gentle and soothing, not high-temperature.

Problem	Possible Cause	Solution
Massage feels too intense/weak.	Intensity level needs adjustment.	Adjust the Pressure Intensity using the dedicated button.

8. Specifications

- **Product Dimensions:** 19 x 16.9 x 18.1 inches
- **Item Weight:** 20 Pounds
- **Item Model Number:** CF-5001-US
- **Manufacturer:** COMFIER
- **Power Source:** Corded Electric
- **Target Use Body Part:** Ankle, Feet, Leg
- **Massager Form:** Foot Massager Machine
- **Massage Therapy Technique:** Shiatsu
- **Material Type:** Leather
- **UPC:** 652805972510

9. Warranty and Support

Manufacturer Warranty: This product comes with a 1-year manufacturer warranty.

For warranty claims, technical support, or any inquiries regarding your COMFIER CF-5001-US massager, please refer to the contact information provided in the product packaging or visit the official COMFIER website. You may also visit the [COMFIER Store on Amazon](#) for additional resources.