Manuals+

Q & A | Deep Search | Upload

manuals.plus /

- > Philips Sonicare /
- > Philips Sonicare 2100 Series Electric Toothbrush Instruction Manual (Model HX3661/04)

Philips Sonicare HX3661/04

Philips Sonicare 2100 Series Electric Toothbrush Instruction Manual

Model: HX3661/04 Brand: Philips Sonicare

1. Introduction

The Philips Sonicare 2100 Series Electric Toothbrush is designed to provide a powerful yet gentle cleaning experience. Utilizing advanced sonic technology, this toothbrush helps remove plaque more effectively than a manual toothbrush. It features two intensity settings, a Smartimer, and a QuadPacer to guide you through your brushing routine, ensuring comprehensive oral care.

2. WHAT'S IN THE BOX

Your Philips Sonicare 2100 Series package includes the following components:

- 1x Philips Sonicare 2100 Series electric toothbrush handle
- 1x C1 Simply Clean brush head
- 1x USB charging stand (power adapter not included)



Figure 1: The Philips Sonicare 2100 Series package contents, showing the toothbrush handle, a C1 Simply Clean brush head, and the USB charging stand.

3. SETUP

3.1 Charging the Toothbrush

Before first use, fully charge your Philips Sonicare toothbrush. The toothbrush provides up to 14 days of regular brushing on a single full charge (based on two two-minute brushings per day on standard mode).

- 1. Connect the USB charging stand to a compatible USB power adapter (not included) and plug it into a wall outlet.
- 2. Place the toothbrush handle onto the charging stand. The battery indicator light on the handle will illuminate to show charging status.
- 3. Allow the toothbrush to charge for at least 24 hours for a full charge.



Figure 2: The Philips Sonicare 2100 Series toothbrush handle placed on its USB charging stand.

3.2 Attaching the Brush Head

Align the brush head with the metal shaft on the toothbrush handle and push it down firmly until it clicks into place. Ensure the brush head is securely attached before use.



Figure 3: A close-up view of the brush head in action, demonstrating its gentle yet effective plaque removal.

4. OPERATING INSTRUCTIONS

4.1 Powering On and Off

Press the power button once to turn the toothbrush on. Press it again to turn it off.

4.2 Advanced Sonic Technology

The Philips Sonicare 2100 Series utilizes advanced sonic technology that pulses fluid between teeth and along the gumline, providing a thorough and effective clean.

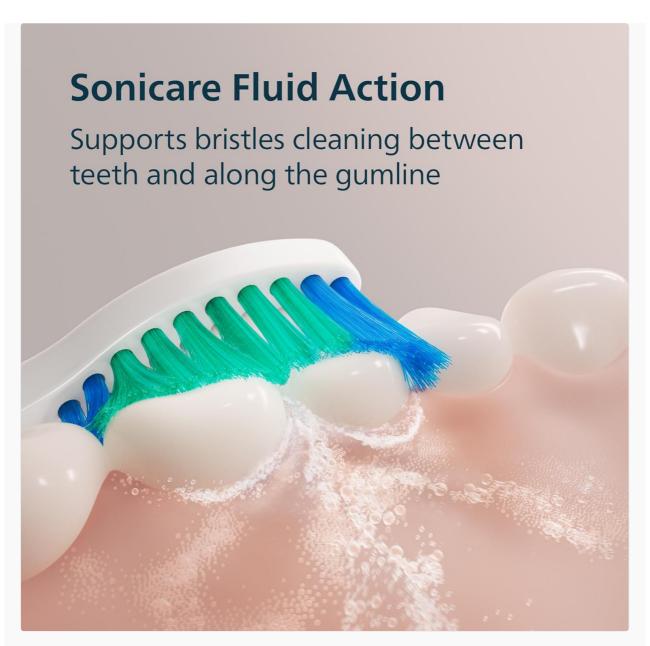


Figure 4: Visual representation of Sonicare Fluid Action, illustrating how bristles and fluid work together for effective cleaning.

4.3 Intensity Settings

The toothbrush offers two intensity settings: low and high. You can select your preferred intensity by pressing the mode button (if available, or by cycling through power button presses depending on the model variation). The EasyStart feature gradually increases the brushing power over the first 14 uses to help you adapt to the Sonicare experience.



Figure 5: The toothbrush handle displaying the two intensity settings and the 14-day battery life icon.

4.4 Smartimer and QuadPacer

The built-in Smartimer ensures you brush for the dentist-recommended 2 minutes. The QuadPacer feature prompts you to move to a new section of your mouth every 30 seconds, ensuring even and thorough cleaning across all four quadrants of your mouth.



Figure 6: Visual explanation of SmarTimer and QuadPacer, guiding users to complete a full 2-minute brushing cycle.

4.5 Official Product Video: Philips Sonicare 2000 Series Overview

Video 1: This video provides an overview of the Philips Sonicare 2000 Series, demonstrating features such as sonic technology, intensity settings, QuadPacer, Easy-Start, and USB charging, which are also applicable to the 2100 Series.

5. MAINTENANCE

5.1 Replacing Brush Heads

For optimal performance and hygiene, dental professionals recommend replacing your brush head every 3 months. To replace, simply pull the old brush head off the handle and snap a new one into place.

PHILIPS SOTTICATE 2100 Poorer To Confidential Confidenti

packaging may vary

Figure 7: The C1 Simply Clean brush head, which should be replaced every 3 months for effective cleaning.

5.2 Cleaning the Toothbrush

After each use, rinse the brush head and bristles thoroughly under running water. Wipe the handle with a damp cloth. Do not submerge the handle in water for extended periods.

6. TROUBLESHOOTING

- **Toothbrush not turning on:** Ensure the toothbrush is fully charged. If the issue persists, try cleaning the charging contacts.
- **Reduced brushing power:** This may indicate a low battery. Place the toothbrush on the charger. If the battery is fully charged and power is still low, the brush head might be worn out and needs replacement.
- **Brush head not staying on:** Ensure the brush head is correctly aligned and pushed down firmly onto the metal shaft. If the brush head is old, it might be worn and require replacement.

7. SPECIFICATIONS

Feature Detail

Model Number	HX3661/04
Dimensions	2.48 x 4.29 x 8.66 inches
Weight	9.07 ounces
Power Source	Battery Powered (1 Lithium Ion battery included)
Special Features	Brushing Timer, Ergonomic Design, Rechargeable
Item Firmness Description	Extra Soft

8. SAFETY INFORMATION

Please refer to the safety information provided in the packaging for detailed warnings and precautions regarding the use and maintenance of your Philips Sonicare 2100 Series Electric Toothbrush.

9. WARRANTY AND SUPPORT

For warranty information, product registration, or technical support, please visit the official Philips Sonicare website or contact their customer service. Details can typically be found in the product packaging or on the manufacturer's website.

Related Documents - HX3661/04



Philips Sonicare Brilliant Clean HX3792/21: Smart Electric Toothbrush Overview

Explore the Philips Sonicare Brilliant Clean HX3792/21 sonic electric toothbrush. This document details its smart screen features, advanced sonic technology, C3 Premium Plaque Control brush head, multiple brushing modes (Clean, White), and technical specifications for a personalized and gentle oral care experience.



Philips Sonicare HX3792/11 Advanced Clean Electric Toothbrush - Intelligent Brushing and Gum Care

Discover the Philips Sonicare HX3792/11 Advanced Clean electric toothbrush. Featuring smart display, pressure sensor, advanced sonic technology, and 5 brushing modes for a gentle, personalized, and superior clean. Protect your gums and optimize your oral hygiene.



Philips Sonicare 4100 Series HX3689/44: Electric Toothbrush for Whiter Teeth and Gum Protection

Discover the Philips Sonicare 4100 Series HX3689/44 electric toothbrush. Features advanced Sonicare technology for whiter teeth in 1 week, gum protection with Pressure Sensor, SmartTimer, and a 14-day battery life. Learn about its specifications and benefits.



Philips Sonicare Series 6000 Rechargeable Sonic Toothbrush HX7400/03 - Advanced Gum and Teeth Care

Discover the Philips Sonicare Series 6000 HX7400/03 rechargeable sonic toothbrush. Experience 1000% more plaque removal, gentle care for sensitive teeth and gums, and advanced features like a pressure sensor and SmarTimer.



Philips Sonicare Series 6100 Rechargeable Sonic Toothbrush - Product Overview

Discover the Philips Sonicare Series 6100 rechargeable sonic toothbrush. Featuring advanced Sonic technology, personalized cleaning modes, a smart pressure sensor, and long-lasting battery life for a gentle yet effective clean.



Philips Sonicare Series 6100 HX7406/01: Innovative Sonic Toothbrush for Sensitive Teeth

Discover the Philips Sonicare Series 6100 HX7406/01 sonic toothbrush. Featuring advanced Sonicare technology, 2 cleaning modes, 3 intensity levels, and a pressure sensor for gentle yet effective cleaning. Ideal for sensitive teeth and gums.