

SNOW 075-0035

SNOW Electric Toothbrush with Sonic Technology - LED Whitening - User Manual

Model: 075-0035

1. INTRODUCTION

The SNOW Electric Toothbrush combines sonic technology with LED light brush heads for an enhanced oral care experience. This manual provides essential information for the safe and effective use, maintenance, and troubleshooting of your device.

2. SAFETY INFORMATION

Warning:

- Do not allow the charging base or cord to sit in water.
- For use with 5V only.

Always follow basic safety precautions when using electrical products, especially when children are present.

3. PACKAGE CONTENTS

Your SNOW Electric Toothbrush package includes the following items:

- SNOW Electric Toothbrush Handle
- 2 LED Brush Heads
- Charging Base
- USB Charging Adapter
- User Manual

BUILDS HEALTHIER BRUSHING HABITS

The SNOW Teeth Whitening Electric Toothbrush has a built-in 2-minute timer, and it even tells you (with a little buzz) exactly when to switch sides.



Image: The SNOW Electric Toothbrush package contents, showing the toothbrush handle, two brush heads, charging base, and USB charging adapter neatly arranged.

4. SETUP

To prepare your toothbrush for first use:

1. Carefully remove all components from the packaging.
2. Align the brush head with the metal shaft on the toothbrush handle.
3. Push the brush head down firmly until it clicks into place.



Image: The SNOW Electric Toothbrush handle with a brush head securely attached, ready for use.

5. CHARGING THE TOOTHBRUSH

Ensure your toothbrush is fully charged before initial use to maximize battery life.

- Place the toothbrush handle onto the charging base.
- Connect the USB charging adapter to the charging base and plug it into a 5V power source.
- The indicator light on the handle will illuminate red while charging and turn green when fully charged.
- A full charge typically takes approximately 4 hours and provides up to 14 days of battery life.
- It is normal for the toothbrush to be warm when charging.



**SONIC
TECHNOLOGY
TO CLEAN TEETH**



**2MIN TIMER
FOR BETTER
BRUSHING**



**EASY
TO USE**



**EASY TO USE
LED TOOTHBRUSHES**

INTELLIGENT BRUSHING:

The 30 second quad pacer and smart timer helps ensure you reach the dental professional recommended 2 minute brushing time and encourages brushing in each quadrant of the mouth.

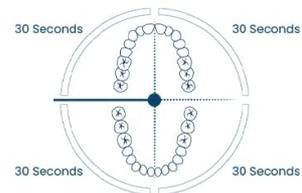


Image: The SNOW Electric Toothbrush standing upright on its charging base, indicating it is connected and charging.

6. OPERATING INSTRUCTIONS

Brushing Modes

The SNOW Electric Toothbrush offers four unique brushing modes to cater to different dental needs:

- **Clean:** Standard mode for daily plaque removal.
- **White:** Designed to remove surface stains and enhance brightness.
- **Polish:** For a brighter, more polished tooth surface.
- **Sensitive:** Gentle cleaning for sensitive teeth and gums.

To select a mode, press the power button (snowflake icon) to turn on the toothbrush. Continue pressing the power button to cycle through the modes. The corresponding LED indicator will light up for the selected mode.



WHITENING ELECTRIC LED TOOTHBRUSH



Sonic technology for better cleaning
Vibration Frequency: 39600 rpm



Blue LED light trusted by millions for whitening support



2-minute timer with 30-second intervals for easy, intelligent brushing



4 unique brushing modes for all your smile's needs



Long-lasting battery life with charging base



Have the cleanest, brightest, and healthiest smile possible.

INTELLIGENT BRUSHING:

The 30 second quad pacer and smart timer helps ensure you reach the dental professional recommended 2 minute brushing time and encourages brushing in each quadrant of the mouth.

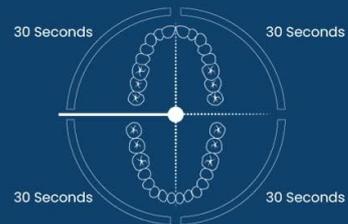


Image: A detailed view of the SNOW Electric Toothbrush handle, highlighting the four brushing modes: Clean, White, Polish, and Sensitive, each with an indicator light.

Intelligent Brushing Timer

The toothbrush features a built-in 2-minute timer with 30-second interval pauses. This helps ensure you brush each quadrant of your mouth for the dentist-recommended duration.

- Divide your mouth into four quadrants: upper right, upper left, lower right, and lower left.
- Brush one quadrant until the toothbrush briefly pauses, signaling you to move to the next quadrant.
- The toothbrush will automatically turn off after 2 minutes.



Image: A visual guide demonstrating the 30-second interval timer, dividing the mouth into four sections for balanced brushing.

7. MAINTENANCE AND CLEANING

Proper maintenance ensures the longevity and hygiene of your toothbrush.

- After each use, rinse the brush head thoroughly under running water for a few seconds to remove any toothpaste residue or debris.
- Occasionally, remove the brush head and rinse the metal shaft and the inside of the brush head to prevent toothpaste buildup, which can cause bacteria growth.
- Store the toothbrush in an upright position to allow proper air drying.
- If not using for an extended period, unplug the unit from the outlet and store it in a cool, dry area, away from direct sunlight.
- Replace brush heads every 3 months, or sooner if bristles become frayed.

WELCOME TO THE FUTURE OF BRUSHING

Meet the 2022 NEW BEAUTY® Award Winner! Our **LED electric toothbrush** is called the **"Apple of Oral Care."** It not only cleans; it also whitens! Equipped with **blue LED** light for whitening support and sonic technology for better cleaning, this toothbrush leaves your smile whiter and brighter. And, with **4 unique modes**, your new toothbrush customizes each brushing session to exactly what you need: Clean, Whiten, Polish, and Sensitive.



Image: The SNOW Electric Toothbrush with its blue LED light activated, demonstrating its operational state during a brushing cycle.

8. TROUBLESHOOTING

If you encounter any issues with your SNOW Electric Toothbrush, please refer to the following common solutions:

Problem	Possible Cause	Solution
Toothbrush does not turn on.	Low battery or not charged.	Place on charging base and ensure indicator light is red. Charge for 4 hours.
Indicator light not working.	Loose connection or power issue.	Check power outlet and USB connection. Ensure toothbrush is correctly seated on the base.
Reduced brushing power.	Low battery or worn brush head.	Recharge the toothbrush. Replace the brush head if bristles are splayed or worn.
Unusual noise during operation.	Loose brush head or debris.	Remove and reattach brush head. Clean any debris from the shaft and brush head.

If the problem persists after trying these solutions, please contact SNOW customer support.

9. SPECIFICATIONS

- **Product Dimensions:** 2.63 x 4 x 10.88 inches
- **Item Weight:** 14.56 ounces
- **Model Number:** 075-0035
- **Batteries:** 1 Lithium Ion battery (included)
- **Power Source:** Battery Powered
- **Item Firmness Description:** Medium
- **Special Feature:** LED Whitening, Sonic Technology, Built-in Timer

10. WARRANTY AND SUPPORT

Your SNOW Electric Toothbrush comes with a **2-year warranty** from the date of purchase, covering manufacturing defects. Additionally, SNOW offers a **60-day money-back guarantee** for your satisfaction. For warranty claims, technical support, or any questions regarding your product, please visit the official SNOW website or contact their customer service department.