

OasisSpace CTS001

OasisSpace Cold Therapy Machine User Manual

Model: CTS001

Brand: OasisSpace

1. INTRODUCTION

The OasisSpace Cold Therapy Machine is engineered to provide effective relief from localized pain and swelling following surgery or injury. By circulating cold water, it helps to reduce inflammation, alleviate joint pain, and accelerate recovery. This system is designed for ease of use and quiet operation, making it an ideal solution for continuous cold therapy.

Key Features:

- **Ice Therapy Relief:** Provides soothing, constant circulation of cold water to reduce pain and swelling.
- **Customizable Therapy:** Offers five different flow rates and five time adjustment modes (30-45-60-75-90 minutes).
- **Durable & Leak Proof:** Constructed with durable composite material and reinforced seams to prevent leaks.
- **Low Noise:** Equipped with a high-quality diaphragm pump for quiet operation.
- **Transparent Bucket:** Allows easy monitoring of ice levels without opening the lid.

2. WHAT'S IN THE BOX

Your OasisSpace Cold Therapy Machine package includes the following components:

- OasisSpace Cold Therapy Machine Main Unit (CTS001)
- 2 Flexible Therapy Pads (one for knee/shoulder, one for ankle/other areas)
- Insulated Pipe/Hose
- Adjustable Elastic Straps (various sizes for different body parts)
- Power Adapter
- User Manual



Image: All components included in the OasisSpace Cold Therapy Machine package.

3. SETUP INSTRUCTIONS

Follow these steps to set up your cold therapy machine:

1. **Prepare the Unit:** Press the red button on the main unit to open the top cover.
2. **Add Ice and Water:** Fill the transparent bucket with ice cubes up to the indicated "Ice Mark" line. Then, add cold water until it covers the sensor tube on the lid (the unit has a 6-liter capacity).



Image: Proper filling of the unit with ice and water.

3. **Close the Cover:** Insert the buckle of the top cover into the groove and press down firmly to close it tightly.
4. **Connect the Hose:** Connect one end of the insulated hose to the connectors on the main unit's cover. Ensure a secure connection by pushing the silver buttons until they click into place.
5. **Connect the Therapy Pad:** Connect the other end of the hose to the desired flexible therapy pad. Again, ensure a secure connection by pushing the silver buttons until they click.
6. **Plug in Power:** Plug the power adapter into the main unit and then into a wall outlet.

Video: A step-by-step guide on how to set up the OasisSpace Cold Therapy Machine, including filling with ice and water, and connecting the hoses and pads.

4. OPERATING INSTRUCTIONS

Once the unit is set up, follow these steps to operate it:

1. **Power On:** Flip the power switch on the main unit to the "On" position. The LED display will illuminate.
2. **Select Flow Rate:** The display will show "LPM" (liters per minute) and the current flow rate. Press and hold the "M" (Mode) button for a couple of seconds until "LPM" flashes. Use the up/down arrow buttons to adjust the flow rate to your desired intensity.
3. **Set Timer:** After selecting the flow rate, press the "M" button again. The "Mins" (Minutes) display will flash. Use the up/down arrow buttons to set the desired therapy duration (30, 45, 60, 75, or 90 minutes).
4. **Start Therapy:** Once the settings are confirmed, press the Play/Pause button to start the cold therapy. You will hear the pump circulating water, and bubbles may be visible, indicating proper operation.
5. **Adjust Application:** Secure the therapy pad to the desired body part using the adjustable elastic straps. Ensure the pad is comfortably but securely placed for optimal contact.

Video: An overview of the OasisSpace Cold Therapy Machine in operation, demonstrating its quietness and effective cold circulation.

5. USAGE & APPLICATION

The flexible therapy pads and adjustable straps allow for versatile application on various body parts:

- **Knee:** Position the pad around the knee, ensuring the central opening aligns with the kneecap. Secure with straps.
- **Shoulder:** Place the pad over the shoulder joint and secure it with the appropriate straps.
- **Ankle:** Wrap the pad around the ankle, ensuring full coverage of the affected area, and fasten with straps.
- **Other Areas:** The flexible design allows for application on the waist, back, or elbow as needed.

SMOOTH & QUIET

Give you a relaxed recovery

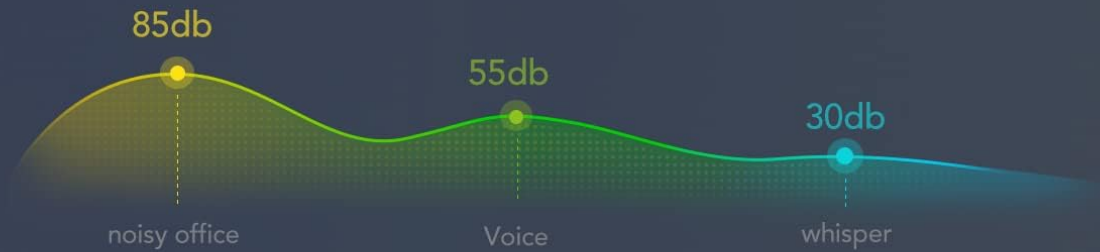


Image: Versatile application of the therapy pads on different body areas.

6. MAINTENANCE

To ensure the longevity and proper functioning of your OasisSpace Cold Therapy Machine:

- **Drain Water:** After each use, carefully drain all water from the therapy pad and the insulated pipe to prevent leakage or damage.
- **Cleaning:** Hand wash the pads in warm water using mild soap rinse thoroughly. Air dry completely before storage.
- **Storage:** Store the unit and its components in a cool, dry place away from direct sunlight.

7. TROUBLESHOOTING

If you encounter any issues with your cold therapy machine, refer to the following common solutions:

- **No Water Circulation / Pump Not Working:** Ensure the unit is filled with enough water to cover the sensor tube on the lid (up to 6 liters). If bubbles are produced in the water during normal operation or water flows back out when opening the lid, it indicates the pump is running. Check all connections for tightness.
- **Leaking:** Ensure all connections (hose to unit, hose to pad) are securely fastened. When turning off the machine, make sure to drain water from the pad and pipe to prevent residual water from causing issues.
- **Pad Not Getting Cold Enough:** Verify that sufficient ice is used and that the water level is adequate. Ensure the pad is properly secured to the body part for efficient heat transfer.

8. SPECIFICATIONS

Product Dimensions	11.4 x 10.6 x 11.2 inches
Net Weight of Base Unit	6.8 lbs
Item Model Number	CTS001
Manufacturer	OasisSpace
Water Capacity	Up to 6 Liters

9. WARRANTY & SUPPORT

Your OasisSpace Cold Therapy Machine comes with a **One Year Warranty**. Your satisfaction is our greatest aim. Please contact us if you encounter any difficulties during use or if the machine malfunctions within the first year of purchase.

For additional support, you may refer to the official User Guide (PDF) available online or contact OasisSpace customer service.