

[Manuals.plus](#) /

› [CURECURE](#) /

› CURECURE Electric Heating Pad User Manual

CURECURE B09NVSZM97

CURECURE Electric Heating Pad User Manual

Model: B09NVSZM97

INTRODUCTION

Thank you for choosing the CURECURE Electric Heating Pad. This manual provides important information regarding the safe operation, maintenance, and care of your heating pad. Please read all instructions carefully before use and retain this manual for future reference.

This heating pad is designed to provide soothing warmth for temporary relief of minor aches and pains associated with muscle soreness, stiffness, and tension in areas such as the back, shoulders, neck, abdomen, and legs.

IMPORTANT SAFETY INSTRUCTIONS

WARNING: To reduce the risk of burns, electric shock, fire, and personal injury, this product must be used in accordance with the following instructions:

- Do not use while sleeping.
- Do not use on an infant.
- Do not use on an invalid, a sleeping or unconscious person, or a person with poor blood circulation unless carefully attended.
- Do not use on areas of insensitive skin.
- Do not use in an oxygen-enriched environment or near equipment that stores or emits oxygen.
- Never pull this pad by the supply cord.
- Do not use pins or other metallic means to fasten this pad.
- Do not sit on, or crush pad. Avoid sharp folds.
- Never leave the appliance unattended, especially if children are present.
- Do not use this pad with liniment, salve, or ointment preparations that contain heat-producing ingredients, as skin burns could result.
- Do not use this product if the cord or plug is damaged.

- Unplug from outlet when not in use.
- This pad is equipped with an automatic shut-off feature for safety.
- Do not use for purposes other than those for which it is intended.

PRODUCT OVERVIEW

The CURECURE Electric Heating Pad is designed for targeted heat therapy. It features a soft flannel material, multiple heat settings, and an automatic shut-off for safety and convenience.

Key Features:

- **Large Coverage Area:** 24" x 33" size for comprehensive coverage of shoulders, neck, back, and other body parts.
- **4 Heat Settings:** Adjustable temperature from approximately 108°F to 140°F (42°C to 60°C).
- **Automatic Shut-Off:** 2-hour auto shut-off for safety and energy saving.
- **Fast Heating:** Provides rapid warmth for quick relief.
- **Secure Fit:** Designed with snap fasteners, weighted edges, and adjustable straps for stability during use.
- **Soft & Skin-Friendly Material:** Made from microplush flannel for comfort.
- **Machine Washable:** Easy to clean (controller must be removed).
- **Safety Certified:** ETL certified with overheating protection.

Components:

- Heating Pad
- Controller with power cord
- Adjustable straps

4 Heating Setting

±108°F ~ ±140°F

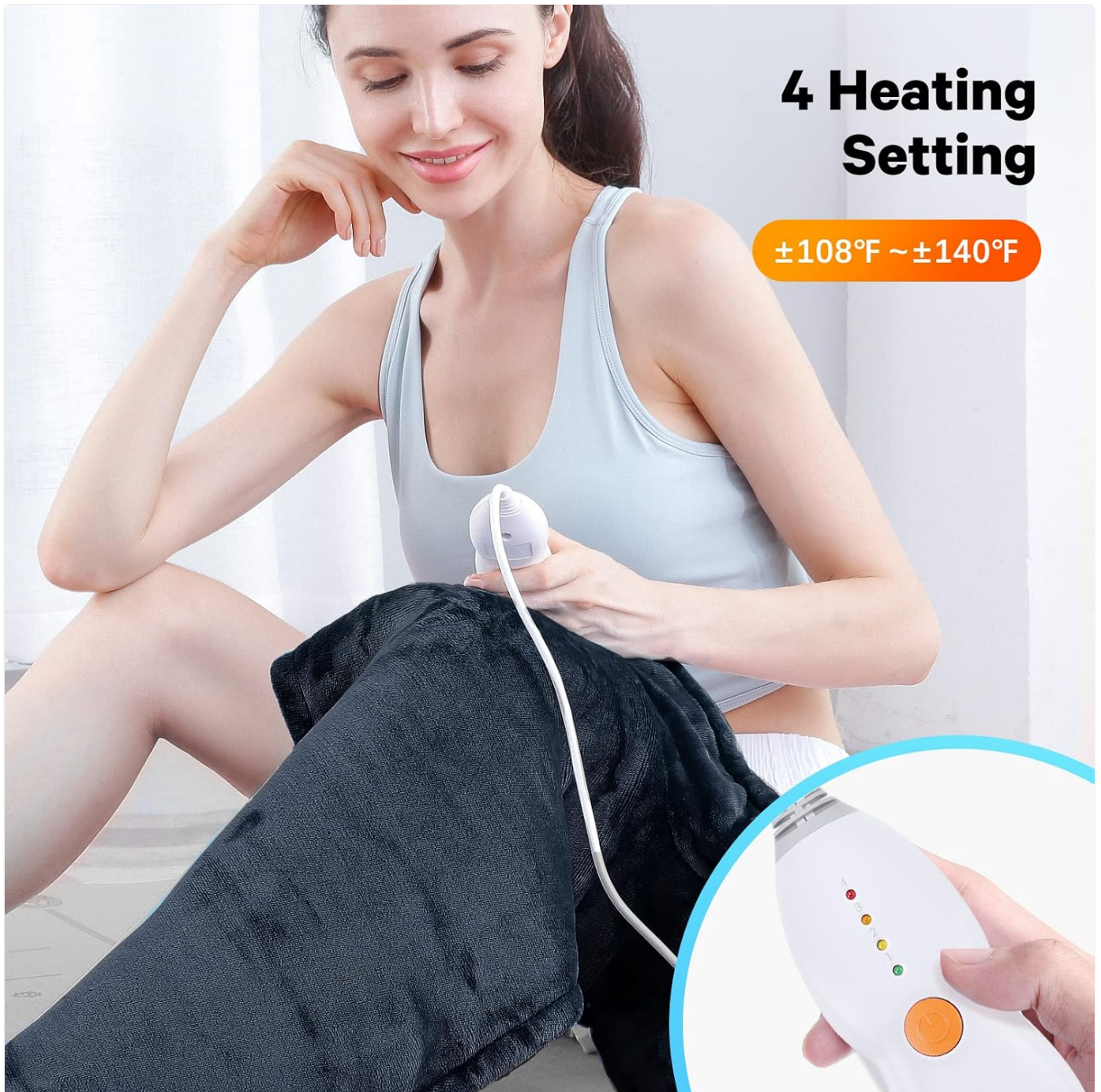


Image: The CURECURE Electric Heating Pad with its attached controller. The controller displays indicators for 4 heat settings, ranging from approximately 108°F to 140°F.



Image: A detailed diagram highlighting the heating pad's design features: neck contour for warmth retention, hands-free magnetic clasps, full back heating capability, and flexible wearing with an adjustable 65 cm strap.

SETUP

1. **Unpack:** Carefully remove the heating pad and controller from its packaging.
2. **Inspect:** Check the heating pad and power cord for any signs of damage. Do not use if damaged.
3. **Connect:** Ensure the controller is securely plugged into the heating pad.
4. **Position:** Place the heating pad on the desired body area (shoulders, back, neck, etc.). Use the snap fasteners and adjustable straps to secure it comfortably in place.
5. **Plug In:** Insert the power plug into a standard 120V AC electrical outlet.



Image: A woman demonstrating how to secure the heating pad around her shoulders using the integrated button design for a snug fit.

OPERATING INSTRUCTIONS

1. **Power On:** Press the power button on the controller to turn on the heating pad. The indicator light will illuminate.
2. **Select Heat Setting:** Press the heat setting button to cycle through the 4 available heat levels. The corresponding indicator light will show the selected setting.
 - Setting 1: Low heat
 - Setting 2: Medium-low heat
 - Setting 3: Medium-high heat
 - Setting 4: High heat
3. **Automatic Shut-Off:** The heating pad is designed to automatically shut off after 2 hours of continuous use. This is a safety feature. To restart, simply press the power button again.
4. **Power Off:** To turn off the heating pad manually, press the power button until all indicator lights are off, or unplug

the unit from the electrical outlet.

5. **Wet Heat Therapy (Optional):** Lightly dampen the heating pad with water using a spray bottle. Do not soak the pad. Ensure the controller remains dry.

Soft Heating Pad

At home or in the office,
anytime, anywhere
Warm your back, shoulder and neck



Image: A woman comfortably using the soft heating pad on her back and shoulder while relaxing on a couch, illustrating typical use.



Image: An illustration demonstrating the versatile application of the heating pad on different body parts, including the shoulder, back, abdomen, and knee, for targeted relief.

MAINTENANCE AND CLEANING

Proper care ensures the longevity and safe operation of your heating pad.

1. **Disconnect Power:** Always unplug the heating pad from the electrical outlet before cleaning.
2. **Remove Controller:** Detach the controller from the heating pad by gently pulling the connector. The controller is NOT washable.
3. **Machine Wash:** The fabric heating pad is machine washable. Wash on a gentle cycle with cold water. Do not bleach.
4. **Tumble Dry:** Tumble dry on low heat or line dry. Ensure the heating pad is completely dry before reattaching the controller and storing or using. Do not iron.
5. **Storage:** When not in use, store the heating pad in a cool, dry place. Avoid sharp folds or placing heavy objects on the pad.



Image: A heating pad being placed into a washing machine, illustrating its machine washable feature after the controller has been removed.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Heating pad does not heat.	◦ Not plugged in. ◦ Controller not securely connected to pad. ◦ Automatic shut-off activated. ◦ Unit is damaged.	◦ Ensure power cord is firmly plugged into a working outlet. ◦ Check that the controller connector is fully inserted into the pad. ◦ Press the power button to restart the unit. ◦ If still not working, discontinue use and contact customer support.
Heating pad heats unevenly.	◦ Pad is folded or creased during use.	◦ Ensure the pad is laid flat or contoured smoothly to the body without sharp folds.
Indicator lights are off, but unit is plugged in.	◦ Automatic shut-off activated. ◦ Power outage or faulty outlet.	◦ Press the power button to restart. ◦ Check the outlet with another appliance.

SPECIFICATIONS

- **Model:** B09NVSZM97
- **Brand:** CURECURE
- **Dimensions:** 24" x 33" (approx. 61 cm x 84 cm)
- **Material:** Flannel (Microplush)
- **Heat Settings:** 4 levels (approx. 108°F - 140°F / 42°C - 60°C)
- **Auto Shut-Off:** 2 hours
- **Power Input:** 120V AC (Standard US outlet)

- **Certifications:** ETL
- **Color:** Dark Blue
- **UPC:** 799513718778

WARRANTY AND SUPPORT

No specific warranty information was provided in the product details. For warranty inquiries or customer support, please refer to the contact information provided with your purchase or visit the official CURECURE website.

Manufacturer: CURECURE

