

YUNCHI K3

YUNCHI K3 Sonic Toddler Electric Toothbrush Instruction Manual

Model: K3 Sonic (LMYC-K3-004HL)

INTRODUCTION

This manual provides essential instructions for the safe and effective use, maintenance, and troubleshooting of your YUNCHI K3 Sonic Toddler Electric Toothbrush. Designed for children aged 0-12 years, this toothbrush features sonic cleaning technology, an LED timer, and soft bristles to promote healthy oral hygiene habits.

PRODUCT FEATURES

- **Sonic Cleaning Technology:** Equipped with 16,000 vibrations per minute for effective cleaning of baby and toddler teeth.
- **Food-Safe Materials:** Made from safe ABS and TPE materials, BPA-free.
- **Soft Bristles:** Selected soft bristles to protect children's tooth enamel and tender gums.
- **2-Minute Smart Timer with 30-Second Pulse Reminder:** Encourages dentist-recommended brushing time by signaling to move to a different quadrant of the mouth every 30 seconds.
- **LED Light:** Soft and glare-free LED indicator helps parents observe children's dental health and ensures thorough cleaning.
- **IPX7 Waterproof Design:** Allows for easy cleaning under running water.
- **Easy-to-Grip Silicone Handle:** Non-slip design for small hands.
- **Suction Cup Base:** For easy and stable storage on smooth surfaces.
- **Battery Powered:** Uses 1 AAA battery (not included), providing up to 150 days of use.
- **Interchangeable Brush Heads:** Includes 4 replacement heads (2 small, 2 large) suitable for different age ranges (Mini: 0-3 years, Medium: 3-6 years, Big: 4-12 years).

PACKAGE CONTENTS

- 1 x YUNCHI K3 Electric Toothbrush Handle
- 4 x Replacement Brush Heads (2 small, 2 large)

- 1 x Protective Bristle Cover
- 1 x Instruction Manual



Image: The YUNCHI K3 Electric Toothbrush handle shown with four different sized replacement brush heads, illustrating the complete package contents.

SETUP

1. Install Battery:

Unscrew the bottom cap of the toothbrush handle. Insert one AAA battery according to the polarity markings (+/-). Securely screw the bottom cap back on. A small screwdriver may be required to open the battery compartment, which is typically included with the product.



LED Light
Convenient to see your baby, kids teeth



High-Quality Nylon Bristles
Long time use



Food Grade Silicone Material
Safe & Durable



Easy to Install Battery
1*AAA battery needed

Battery not included in package



Sonic Technology
Deep Cleaning



YUNCHI®

Baby & Kids Electric Toothbrush
Your Kids Oral Health Guard



Image: A close-up view of the YUNCHI K3 toothbrush handle, highlighting the battery compartment at the bottom, indicating where to insert the AAA battery.

2. Attach Brush Head:

Select the appropriate brush head size for your child's age (Mini: 0-3 years, Medium: 3-6 years, Big: 4-12 years). Gently push the chosen brush head onto the metal shaft at the top of the toothbrush handle until it clicks into place.

Brush with baby's first tooth

Yunchi your best oral care friends,
accompany your life



Image: The YUNCHI K3 toothbrush handle displayed alongside three different sizes of replacement brush heads, labeled for age ranges 0-3 months, 3-6 years, and 4-12 years, demonstrating the variety available.

OPERATING INSTRUCTIONS

1. Apply Toothpaste:

Apply a pea-sized amount of children's toothpaste to the brush head.

2. Power On:

Press the power button on the handle to turn on the toothbrush. The LED light will illuminate, and the sonic vibrations will begin.



Image: A toddler smiling while holding and using the YUNCHI K3 electric toothbrush, demonstrating its ease of use for young children.

3. **Brushing Technique:**

Guide your child to gently move the brush head across all surfaces of their teeth. The 2-minute timer will ensure adequate brushing time, with a pulse reminder every 30 seconds to encourage moving to a different quadrant of the mouth (upper right, upper left, lower right, lower left).


SCIENTIFIC BRUSHING



X2 power
modes



2 minute timer

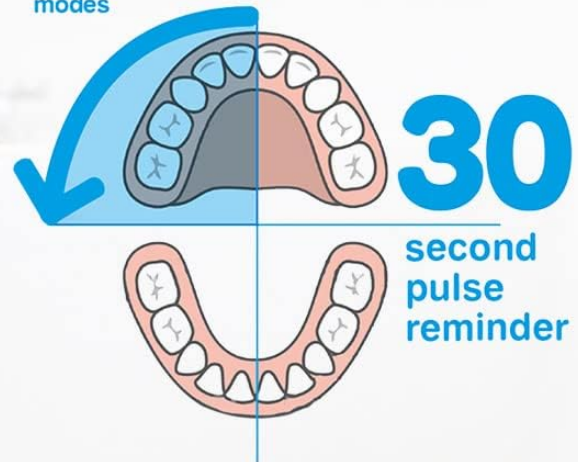

HOW TO BRUSH YOUR TEETH

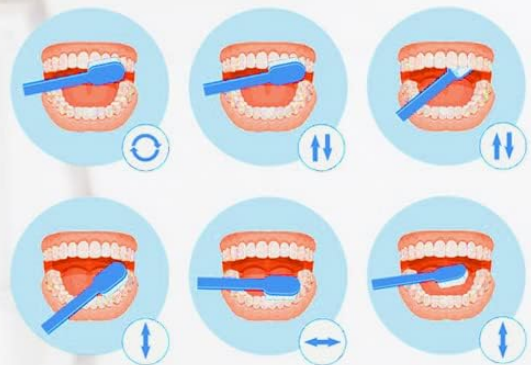


Image: A visual guide illustrating the scientific brushing method, featuring a 2-minute timer and 30-second pulse reminders for thorough cleaning of all tooth quadrants.

4. Power Off:

After 2 minutes, the toothbrush will automatically turn off. You can also manually turn it off by pressing the power button again.

MAINTENANCE

1. Cleaning the Toothbrush:

After each use, rinse the brush head thoroughly under running water. The IPX7 waterproof design allows the entire toothbrush to be rinsed. Wipe the handle with a damp cloth if necessary.

IPX7

Waterproof design



16000

Powerful Sonic
Vibration Movements



Image: The YUNCHI K3 toothbrush being held under running water, illustrating its IPX7 waterproof capability for easy cleaning.

2. Replacing Brush Heads:

Dentists recommend replacing brush heads every 3 months, or sooner if bristles become frayed. To replace, gently pull the old brush head off the shaft and push a new one into place.

3. Battery Replacement:

When the toothbrush performance decreases, it's time to replace the AAA battery. Follow the battery installation steps in the "Setup" section.

4. Storage:

Store the toothbrush upright using its suction cup base on a clean, smooth surface. Ensure the protective bristle cover is in place when not in use.



Making Brushing Teeth Fun & Easy

Dental Advice: Healthy teeth should start from baby



2 Minutes Smart Timer



Easy to Grip

Food Grade Silicone Handle – Non-Slip



Suction Cup

Easy to stand or suck on smooth surface
(mirror, table)

Image: The YUNCHI K3 toothbrush showcasing its ergonomic, easy-to-grip handle and the practical suction cup base for stable upright storage.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Toothbrush does not turn on or vibrate.	Low or dead battery; incorrectly installed battery.	Replace the AAA battery with a new one, ensuring correct polarity. Check if the battery compartment is securely closed.
Vibration is weak or inconsistent.	Low battery power.	Replace the AAA battery.
Brush head feels loose or falls off.	Brush head not properly attached.	Ensure the brush head is pushed firmly onto the metal shaft until it clicks.
Bristles are coming out.	Normal wear and tear; excessive biting on bristles.	Replace the brush head. Monitor child's brushing habits to prevent excessive biting.

SPECIFICATIONS

Feature	Detail
Model Name	K3 Sonic
Item Model Number	LMYC-K3-004HL
Brand	YUNCHI
Age Range Description	Baby (0-12 years with appropriate brush heads)
Power Source	1 AAA Battery (not included)
Operation Mode	Automatic (Sonic Vibration)
Vibration Frequency	16,000 strokes/minute
Timer	2-minute timer with 30-second interval pulse
Waterproof Rating	IPX7
Bristle Material	Nylon (Soft)
Handle Material	Acrylonitrile Butadiene Styrene (ABS) & TPE (BPA Free)
Included Components	Toothbrush handle, 4 replacement heads, protective bristle cover, instruction manual
Dimensions	7.32 x 3.82 x 1.46 inches
Weight	3.21 ounces (0.09 Kilograms)
UPC	643375703291

WARRANTY AND SUPPORT

For any questions, concerns, or support regarding your YUNCHI K3 Sonic Toddler Electric Toothbrush, please contact YUNCHI customer service. We are committed to resolving any issues and ensuring your satisfaction. Contact information may be found on the product packaging or through the retailer where the product was purchased.

© 2023 YUNCHI. All rights reserved.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page’s manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question