

YUNCHI AK-CS-02

YUNCHI AK-CS-02 Ultrasonic Physiotherapy Massager

User Manual

1. INTRODUCTION

Thank you for choosing the YUNCHI AK-CS-02 Ultrasonic Physiotherapy Massager. This device is designed to provide targeted relief for muscle and joint pain through deep-penetrating sound waves. It features adjustable intensity levels and a timer for controlled usage. Please read this manual carefully before operating the device to ensure safe and effective use.

2. SAFETY INSTRUCTIONS

- Always apply an appropriate amount of gel to the probe or skin before turning on the device. Do not use the device on dry skin or without gel.
- Avoid keeping the probe in one place for longer than one minute to prevent discomfort or potential skin irritation.
- This device is not a vibrating or heated massager. Its operation involves ultrasonic waves, which may not be felt directly.
- Keep the device out of reach of children.
- Consult a healthcare professional before use if you have any medical conditions, are pregnant, or have implanted electronic devices.

3. PACKAGE CONTENTS

Upon unpacking, please ensure all items are present:

- YUNCHI AK-CS-02 Physiotherapy Device
- Ultrasonic Probe
- Power Adapter
- Carrying Case (may vary by package)
- User Manual



Image: The YUNCHI AK-CS-02 device, its power adapter, and the ultrasonic probe.



Image: A grey carrying case, often included for convenient storage and transport of the device and its accessories.

4. DEVICE COMPONENTS

Familiarize yourself with the main components and controls of your AK-CS-02 device:

Bequeme Touch-Steuerung

3 Intensitätsstufen

Niedrig, Mittel, Hoch

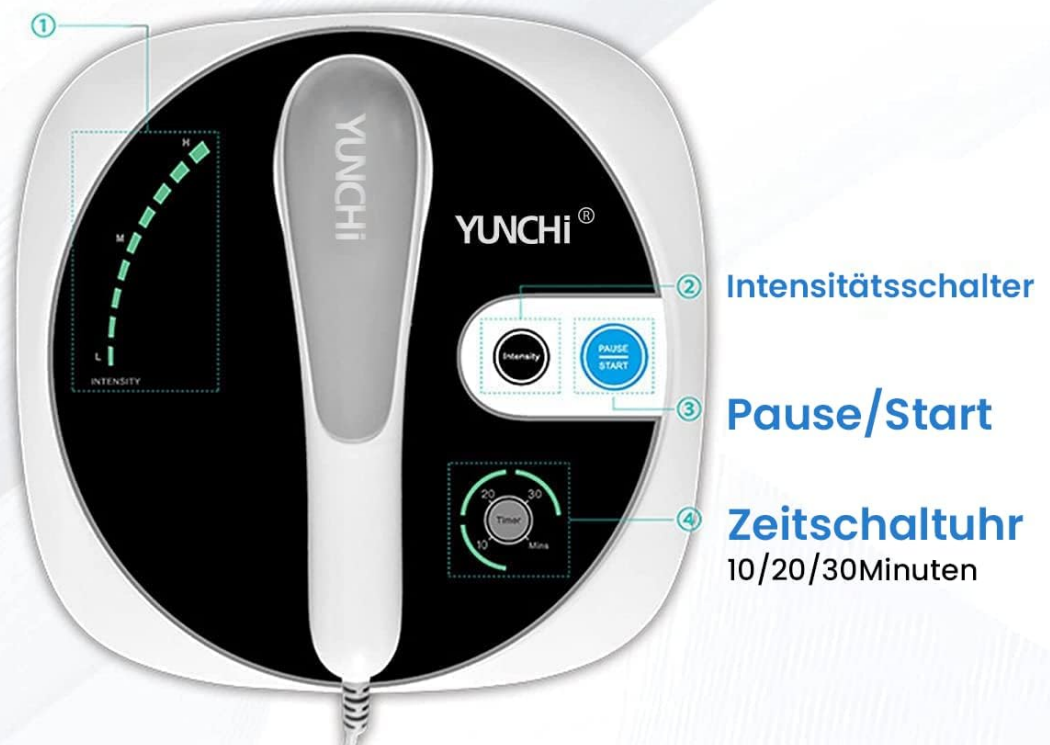


Image: Control panel diagram illustrating the key functions:

1. **Intensity Levels:** Displays the selected intensity (Low, Medium, High).
2. **Intensity Switch:** Button to cycle through intensity levels.
3. **Pause/Start:** Button to begin or pause the treatment.
4. **Timer:** Settings for 10, 20, or 30 minutes of operation.



Image: A close-up view of the metal probe head, which emits the ultrasonic waves.

5. SETUP

1. Connect the power adapter to the device's power input port and then plug it into a suitable electrical outlet.
2. Ensure the ultrasonic probe is securely connected to the device's output port.
3. Prepare the treatment area by ensuring it is clean and dry.

6. OPERATING INSTRUCTIONS

1. **Apply Gel:** Apply an appropriate amount of ultrasonic gel evenly to the metal probe head or directly to the skin area you intend to treat. This is crucial for effective ultrasonic transmission and to prevent skin irritation.
2. **Power On:** Press the 'Start' button to turn on the device.
3. **Select Intensity:** Use the 'Intensity Switch' button to cycle through the three intensity levels: Low, Medium, or High. Choose the level that feels most comfortable and effective for your needs.
4. **Set Timer:** Use the 'Timer' button to select your desired treatment duration: 10, 20, or 30 minutes. The device will automatically shut off after the selected time.

5. **Begin Treatment:** Gently place the probe head onto the skin area with the applied gel.
6. **Move the Probe:** Slowly move the probe in small, circular motions over the treatment area. Do not hold the probe stationary in one spot for more than one minute.
7. **End Treatment:** The device will automatically stop after the set timer. You can also press the 'Pause/Start' button to stop the treatment manually at any time.

7. CHECKING DEVICE FUNCTIONALITY

The YUNCHI AK-CS-02 is an ultrasonic device, not a vibrating or heated massager. Therefore, you may not feel any direct sensation when it is operating. To verify that the device is working correctly:

- Apply a small drop of water to the metal surface of the probe head.
- Turn on the device and select an intensity level.
- Observe the water drop. If the device is functioning, the ultrasonic waves will cause a visible fogging effect or slight ripples on the water surface.

8. MAINTENANCE

- After each use, wipe the probe head clean with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.
- Store the device in its carrying case in a cool, dry place, away from direct sunlight and extreme temperatures.
- Ensure the device is unplugged before cleaning or storing.

9. TROUBLESHOOTING

- **Device does not turn on:** Ensure the power adapter is securely connected to both the device and a working power outlet.
- **No sensation during use:** This is normal for an ultrasonic device. Refer to Section 7 for how to check functionality.
- **Skin irritation:** Ensure sufficient ultrasonic gel is applied. Do not hold the probe in one spot for too long. If irritation persists, discontinue use and consult a healthcare professional.

10. SPECIFICATIONS

Feature	Detail
Model Name	AK-CS-02
Brand	YUNCHI
Power Source	Battery Powered (via adapter)
Material	Plastic
Intensity Levels	3 (Soft, Medium, Strong)
Timer Settings	10, 20, 30 minutes
Massager Form	Handheld

Product Dimensions (L x W x H)	24 x 24 x 11 cm
Target Body Part	Whole Body
Specific Product Uses	General muscle and joint pain relief
Wireless	No
UPC	643375704076
ASIN	B09QKPGX46

11. WARRANTY AND SUPPORT

For warranty information, technical support, or customer service inquiries, please refer to the contact details provided with your purchase or visit the official YUNCHI website. Keep your proof of purchase for warranty claims.

© 2023 YUNCHI. All rights reserved.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page’s manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question