

COMFIER CF-5250

COMFIER Foot Massager CF-5250 User Manual

Model: CF-5250

1. INTRODUCTION

Thank you for choosing the COMFIER Foot Massager CF-5250. This manual provides essential information for the safe and effective use of your new foot massager. Please read it thoroughly before operation and retain it for future reference. The COMFIER Foot Massager CF-5250 is designed to provide comprehensive foot massage, incorporating deep kneading Shiatsu, rolling, air compression, and soothing heat functions to help relieve tired muscles, plantar fasciitis, and neuropathy pain.

2. IMPORTANT SAFETY INFORMATION

To reduce the risk of electric shock, burns, fire, or injury, always follow these basic safety precautions:

- Do not use while sleeping.
- Keep away from water.
- Not suitable for use if you have certain medical conditions. Consult your physician before use if you have any pre-existing medical conditions, are pregnant, or have a pacemaker.
- Always follow the instructions provided with the product.
- Unplug the appliance from the outlet when not in use and before cleaning.
- Do not operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- Do not use pins or other metallic fasteners with this appliance.
- Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.

3. PRODUCT OVERVIEW

The COMFIER Foot Massager CF-5250 is designed for optimal foot relaxation and therapy. Key features include:

- Deep Kneading Shiatsu Massage
- Rolling Massage
- Air Compression Massage (3 adjustable levels)
- Soothing Heat Function (2 adjustable levels, up to 131°F)

- Adjustable Timer (15/30 minutes)
- Easy-to-use Touch Panel Controls
- Detachable and Washable Foot Sleeves
- Accommodates foot sizes up to US Men's 13



Figure 3.1: Top-down view of the COMFIER Foot Massager CF-5250, showing the foot inserts and control panel.



MAX.: US Size 13



Fits Whole Family

Figure 3.2: The foot massager is designed to fit various foot sizes, up to US Men's size 13, making it suitable for most users.

Comprehensive & Easy to use



Touch panel



Figure 3.3: Close-up of the intuitive touch panel controls, allowing easy adjustment of settings.

4. SETUP

1. **Unpacking:** Carefully remove the foot massager from its packaging. Ensure all components, including the power adapter and user manual, are present.
2. **Placement:** Place the foot massager on a flat, stable surface. Ensure there is enough space around the unit for proper ventilation.
3. **Power Connection:** Connect the power adapter to the DC input port on the massager. Plug the other end of the adapter into a standard electrical outlet (AC 100-240V, 50/60Hz).
4. **Foot Sleeves:** Ensure the detachable foot sleeves are properly installed inside the foot chambers. These can be removed for cleaning.

5. OPERATING INSTRUCTIONS

The COMFIER Foot Massager CF-5250 features a user-friendly control panel. Insert your feet into the foot chambers before starting the massage.

5.1 Power On/Off

- Press the **Power** button on the control panel to turn the unit on. The default massage mode will activate.
- Press the **Power** button again to turn the unit off.

5.2 Massage Modes and Intensity

- **Kneading & Rolling:** Press the **Kneading** button to cycle through 3 intensity levels (Low, Medium, High) for the deep kneading and rolling massage. This targets the soles of your feet.
- **Air Compression:** Press the **Air Intensity** button to adjust the air compression level (Low, Medium, High). This applies pressure to the instep, heel, and ankles.



Figure 5.1: Illustration of the deep kneading and rolling massage mechanisms targeting the soles of the feet.

360° Compression Massage



Figure 5.2: Visual representation of the 360-degree air compression massage, showing how air bags inflate to apply pressure around the foot.

5.3 Heat Function

- Press the **Heat** button to activate the soothing heat function. Press again to cycle through 2 heat levels (approx. 113°F to 131°F) or turn it off.
- The heat function can be used independently or in combination with massage modes.

Heat Up in 1 min

Max. temperature of 131 °F
(Warm enough)



2 Levels

113 °F ~ 131 °F



Figure 5.3: The heat function quickly warms up, providing soothing warmth to cold or tired feet.

5.4 Timer Setting

- Press the **Timer** button to select a massage duration of 15 or 30 minutes. The unit will automatically shut off after the selected time.

6. MAINTENANCE

6.1 Cleaning

- Always unplug the massager before cleaning.
- Wipe the exterior of the unit with a soft, damp cloth. Do not use abrasive cleaners or solvents.
- The detachable foot sleeves can be unzipped and removed for washing. Hand wash with mild detergent and air dry completely before reattaching.

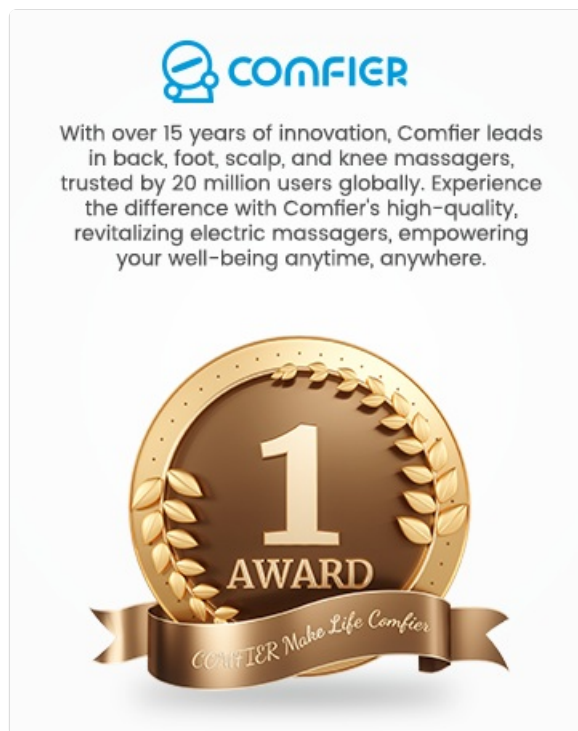


Figure 6.1: The foot massager features detachable and washable foot sleeves for easy cleaning and hygiene.

6.2 Storage

- Store the massager in a cool, dry place when not in use.
- Avoid placing heavy objects on the unit.
- Do not wrap the power cord around the appliance.

7. TROUBLESHOOTING

Problem	Possible Cause	Solution
Unit does not turn on.	Power cord not properly connected. No power from outlet.	Ensure power cord is securely plugged into the unit and a working electrical outlet. Try a different outlet.
Massage intensity feels too weak or too strong.	Incorrect intensity setting.	Adjust the Kneading or Air Intensity settings using the respective buttons on the control panel.
Heat function is not working or not warm enough.	Heat function not activated or set to lower level.	Press the Heat button to activate or cycle through the heat levels. Allow a few minutes for the heat to build up.
Unit automatically shuts off too soon.	Timer setting. Overheating protection activated.	Check the Timer setting. If the unit feels hot, unplug it and allow it to cool down for 30 minutes before re-using.

8. SPECIFICATIONS

Feature	Detail
Model Number	CF-5250
Product Dimensions	17.5 x 13.6 x 7.5 inches

Feature	Detail
Item Weight	7.42 Pounds
Power Source	Corded Electric
Material	Plastic
Max Foot Size	US Men's 13
Heat Temperature	Up to 131°F (55°C)

9. WARRANTY AND SUPPORT

The COMFIER Foot Massager CF-5250 is covered by a 1-year warranty, which is extendable to 3 years. For warranty claims, technical support, or any inquiries, please contact COMFIER customer service through their official website or the contact information provided with your purchase receipt.

For more information, visit the official COMFIER Store: [COMFIER Store on Amazon](#).





10. LEGAL DISCLAIMER





Statements regarding this product have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is not a substitute for professional medical advice, diagnosis, or treatment. Always seek the advice of your physician or other qualified health provider with any questions you may have regarding a medical condition.



© 2025 COMFIER. All rights reserved.

Related Documents - CF-5250

<div><p>COMFIER FUß-MASSAGEGERÄT EE FOOT MASSAGER GB MASSEUR DE PIEDS FR MASSAGGIATORE PIEDI IT MASSAJADOR DE PIES ES</p></div> <div><p>CE User Manual CF-5913</p></div>	<p>Comfier CF-5913 Shiatsu Foot Massager User Manual</p> <p>User manual for the Comfier CF-5913 Shiatsu Foot Massager Machine. Learn how to use this kneading foot massager with heat for pain relief, plantar fasciitis, and general relaxation. Includes operation, safety, maintenance, and warranty details.</p>
<div><p>COMFIER Shiatsu Neck and Back Massager with Heat</p></div> <div><p>User Manual CF-2303-KR</p></div>	<p>COMFIER CF-2303-KR Shiatsu Neck and Back Massager with Heat User Manual</p> <p>Comprehensive user manual for the COMFIER CF-2303-KR Shiatsu Neck and Back Massager with Heat, covering safety instructions, features, technical data, setup, operation, controller functions, care, cleaning, troubleshooting, and warranty information.</p>

<div><div><div>COMFIER</div><div>Shiatsu massage pillow with heat</div></div><div></div><div><div>CE</div><div>User Manual</div><div>CF-1228</div></div></div>	<p>Comfier CF-1228 Shiatsu Massage Pillow with Heat User Manual</p> <p>User manual for the Comfier CF-1228 Shiatsu massage pillow with heat, detailing safety instructions, operating procedures, features, maintenance, disposal, technical specifications, and warranty information.</p>
<div><div><div>COMFIER</div><div>Leg and Foot Massager</div></div><div></div><div><div>CE</div><div>User Manual</div><div>CF-F25H</div></div></div>	<p>Comfier CF-F25H Leg and Foot Massager User Manual</p> <p>Comprehensive user manual for the Comfier CF-F25H Leg and Foot Massager, detailing safety instructions, product features, setup and operation, maintenance, troubleshooting, warranty information, and FCC compliance.</p>
<div><div><div>COMFIER</div><div>Shiatsu massagecushion with heat Shiatsu massage cushion with heat Coussin de massage Shiatsu avec chaleur Cuscino per massaggio Shiatsu riscaldato Cojín masajeador Shiatsu con calefacción</div></div><div></div><div><div>CE</div><div>User Manual</div><div>CF-2606MC</div></div></div>	<p>COMFIER CF-2606MC Shiatsu Massage Cushion with Heat User Manual</p> <p>Comprehensive user manual for the COMFIER CF-2606MC Shiatsu massage cushion with heat. Includes safety instructions, operating guide, features, maintenance, and warranty information.</p>
<div><div><div>COMFIER</div><div>Shiatsu Neck and Back Massager</div></div><div></div><div><div>CE</div><div>User Manual</div><div>CF-2309A</div></div></div>	<p>Comfier Shiatsu Neck and Back Massager CF-2309A User Manual</p> <p>User manual for the Comfier Shiatsu Neck and Back Massager, model CF-2309A. Includes safety instructions, features, technical data, setup and operation, controller instructions, care and cleaning, troubleshooting, and warranty information.</p>