

## MASDERM Professional 2-in-1 Radiofrequency Device and Firming Cream Kit (ASIN: B09RBFK98H)

# MASDERM Professional 2-in-1 Radiofrequency Device and Firming Cream User Manual

Model: Professional 2-in-1 Radiofrequency Device and Firming Cream Kit (ASIN: B09RBFK98H)

Brand: MASDERM

## 1. INTRODUCTION

This manual provides detailed instructions for the safe and effective use of your MASDERM Professional 2-in-1 Radiofrequency Device and Firming Body Cream. This system is designed to assist in body contouring, skin firming, and improving skin elasticity and hydration through advanced radiofrequency technology combined with a specialized firming cream.



Image 1.1: The MASDERM Professional 2-in-1 Radiofrequency Device and Firming Body Cream.

## Safety Information

- Read all instructions before use.
- Do not use on broken, irritated, or infected skin.
- Consult a healthcare professional if you have any pre-existing medical conditions or concerns.
- Keep out of reach of children.
- Ensure the device is clean before and after each use.
- Always use the device with sufficient cream to ensure proper gliding and prevent skin irritation.

## 2. PACKAGE CONTENTS

Your MASDERM Professional 2-in-1 Radiofrequency Pack includes:

- 1 x MASDERM Professional Radiofrequency Device
- 1 x MASDERM Firming Body Cream (1 Liter)
- 1 x Power Adapter (if applicable, not specified in input but common for devices)
- 1 x User Manual (this document)

## 3. GETTING STARTED

### Initial Preparation

1. **Charge the Device:** If the device is battery-operated, ensure it is fully charged before first use. Connect the power adapter to the device and a suitable power outlet.
2. **Cleanse Skin:** Ensure the area of skin you intend to treat is clean and dry.
3. **Apply Cream:** Apply a generous amount of the MASDERM Firming Body Cream to the target area. The cream acts as a conductor for the radiofrequency and ensures smooth gliding of the device.

## 4. OPERATING INSTRUCTIONS

### Application of Cream and Device

For optimal results, apply the firming cream daily, both morning and night. When using the radiofrequency device, apply the device over the cream. It is recommended to use the device 3 nights per week, for approximately 10 minutes per zone.

## FUNCIONES

### Luz verde

**EMS (Electroestimulación Muscular)** ejercita el movimiento de la musculatura para favorecer una mayor tensión en la piel. Piel más firme y menos flácida.

### Luz azul

**Slim (Anticelulitis)** promueve la absorción de la crema en la piel y quema grasa corporal. Piel del cuerpo más lisa, con menos volumen y más tersa.

### Luz roja

**RF (Radiofrecuencia)** reafirma la piel y suaviza la zona promoviendo la circulación sanguínea del cuerpo, mejorando así su firmeza y elasticidad. Piel del cuerpo más firme.

### Luz lila

**SLIM+RF** quema grasa corporal y, al mismo tiempo, promueve la circulación de la piel. Notas la piel más firme, suave y con menos volumen.



Image 4.1: Demonstrating the application of the firming cream followed by the radiofrequency device.

## Device Operation and Modes

Turn on the device and select the desired mode. It is recommended to start with the lowest power setting for any function and gradually increase as your skin adapts. Always ensure there is sufficient cream on the skin for the device to glide properly.

## FUNCIONES

### Luz verde

**EMS (Electroestimulación Muscular)** ejercita el movimiento de la musculatura para favorecer una mayor tensión en la piel. Piel más firme y menos flácida



### Luz azul

**SLIM (Anticelulitis)** promueve la absorción de la crema en la piel y quema grasa corporal. Piel del cuerpo más lisa, con menos volumen y más tersa



### Luz roja

**RF (Radiofrecuencia)** reafirma la piel y suaviza la zona promoviendo la circulación sanguínea del cuerpo, mejorando así su firmeza y elasticidad. Piel del cuerpo más firme



### Luz lila

**EMS (Electroestimulación Muscular)** ejercita el movimiento de la musculatura para favorecer una mayor tensión en la piel. Piel más firme y menos flácida



Image 4.2: The device indicating different operational modes via colored LED lights.

**Green Light (EMS; Fitness Mode):** This mode utilizes Electrical Muscle Stimulation to exercise muscle movement, promoting greater tension in the skin. It helps build muscle tone, improves oxygen utilization, and has soothing effects on the skin, leading to firmer and less flaccid skin.

**Blue Light (SLIM; Body Sculpting Mode):** This mode promotes skin absorption and activation, firming the skin by aiding in collagen synthesis. It significantly improves oily, sensitive, and/or acne-prone skin, contributing to smoother skin with reduced volume and increased firmness.

**Red Light (RF; Radiofrequency Mode):** This mode enhances skin tension and elasticity, restoring its radiance and vitality. It maintains optimal skin conditions long-term, improves luminosity, and reduces wrinkles by penetrating to the deepest layers of the skin.

**Lilac Light (Combined Mode: RF + Body Sculpting):** This mode combines radiofrequency with body sculpting. It restores skin firmness, eliminates toxins, inhibits bacteria, and improves acne marks, resulting in firmer, smoother skin with reduced volume.

## Treatment Areas and Technique

Move the device in slow, circular motions over the target area. Ensure continuous contact with the skin. The device is suitable for various body areas including thighs, abdomen, and arms.



Image 4.3: Application on the thigh.



Image 4.4: Application on the abdomen.



Image 4.5: Application on the arm.



Image 4.6: Example of visible improvements after consistent use.

## 5. CARE AND MAINTENANCE

### Cleaning the Device

- After each use, wipe the device head with a soft, damp cloth.
- Do not use abrasive cleaners or immerse the device in water.
- Ensure the device is completely dry before storing.

### Storage

- Store the device and cream in a cool, dry place, away from direct sunlight.
- Keep the cream container tightly closed when not in use.

## 6. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device not turning on	Low battery; Not properly connected to power (if wired)	Charge the device fully; Ensure power adapter is securely connected.
Device not gliding smoothly	Insufficient cream applied	Apply more MASDERM Firming Body Cream to the treatment area.

Problem	Possible Cause	Solution
Skin irritation after use	Power setting too high; Sensitive skin; Allergic reaction to cream	Reduce power setting; Discontinue use if irritation persists and consult a dermatologist. Review cream ingredients for known allergens.
No visible results	Inconsistent use; Incorrect technique; Insufficient duration	Ensure consistent use as per instructions (3 times/week, 10 min/zone); Review application technique; Allow sufficient time for results (typically visible within 1 month).

## 7. TECHNICAL SPECIFICATIONS

**Brand:** MASDERM

**Product Dimensions (L x W x H):** 15 x 15 x 15 cm (5.9 x 5.9 x 5.9 inches)

**Item Weight:** 2.09 Kilograms (4.61 pounds)

**Cream Volume:** 1000 Milliliters

**Cream Form:** Cream

**Skin Type:** All

**Product Benefits:** Anti-cellulite, Firming

**Specific Uses for Product:** Cellulite, Wrinkles (Flaccidity)

**ASIN:** B09RBFK98H

## 8. WARRANTY AND SUPPORT

### Warranty Information



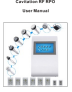

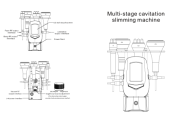

This product is covered by a standard manufacturer's warranty. Please refer to your purchase documentation or contact MASDERM customer service for specific warranty terms and conditions.

### Customer Support

For any questions, technical assistance, or support regarding your MASDERM Professional 2-in-1 Radiofrequency Device and Firming Cream, please contact MASDERM customer support through their official website or the retailer from whom you purchased the product.

*Thank you for choosing MASDERM.*

**Related Documents - Professional 2-in-1 Radiofrequency Device and Firming Cream Kit (ASIN: B09RBFK98H)**

<p>Masderm</p>  <p>RF Body Device DISPOSITIVO RADIOFRECUENCIA CORPÓREA</p>	<p><a href="#">Masderm RF Body Device User Manual</a></p> <p>Comprehensive user manual for the Masderm RF Body Device, detailing its functions, usage, and maintenance. Learn how to effectively use the RF, EMS, and SLIM technologies for skin firming, cellulite reduction, and body contouring.</p>
 <p>3-in-1 40K Ultrasonic Cavitation Machine USER MANUAL V1.1 NuLifeBeauty.com</p>	<p><a href="#">NuLifeBeauty 3-in-1 40K Ultrasonic Cavitation Machine User Manual and Guide</a></p> <p>Comprehensive user manual for the NuLifeBeauty 3-in-1 40K Ultrasonic Cavitation Slimming Machine. Learn how to use, safety precautions, benefits, treatment modes, and troubleshooting for effective body contouring and skin rejuvenation.</p>
 <p>Cavitation RF RPO User Manual</p>	<p><a href="#">Cavitation RF RPO User Manual: Comprehensive Guide to Beauty Device Operation</a></p> <p>This user manual provides detailed instructions for operating the Cavitation RF RPO beauty device, covering its various functions including 40K cavitation, RF treatments (body, face, eye), EMS, Laser Lipolysis, and Cold therapy. It outlines contraindications, notices, package contents, and step-by-step usage guides for achieving optimal skin and body sculpting results.</p>
 <p>Eris Lipo Cavitation Machine 6 in 1 Professional USER MANUAL</p>	<p><a href="#">Eris Lipo Cavitation Machine 6 in 1 Professional User Manual</a></p> <p>Comprehensive user manual for the Rejuva Fresh Eris Lipo Cavitation Machine 6 in 1 Professional, detailing product information, specifications, warnings, treatment selection, and recommended positions for effective use.</p>
 <p>Multi stage cavitation slimming machine</p>	<p><a href="#">Beemyi CS-119 Multi-stage Cavitation Slimming Machine User Manual</a></p> <p>User manual for the Beemyi CS-119 multi-stage cavitation slimming machine. Provides detailed instructions on operation, therapy handle modes, treatment plans, application areas, safety precautions, and contraindications for effective body shaping and skin tightening.</p>
 <p>40K SLIMMING MACHINE "IntuitionSlimmer"</p> <p>In order to use this product correctly, please read the manual carefully and keep it properly.</p>	<p><a href="#">FREYARA 40K Slimming Machine Instruction Manual - Body Contouring and Skin Tightening Guide</a></p> <p>Comprehensive instruction manual for the FREYARA 40K Slimming Machine (LW-202). Learn about its various heads, functions for body contouring, skin tightening, wrinkle reduction, and safe usage guidelines for effective beauty treatments.</p>