

## Escali XLSCALE

# Escali High Capacity Digital Bathroom Scale (Model: XLSCALE) Instruction Manual

Brand: Escali | Model: XLSCALE

[Overview](#)

[Setup](#)

[Operation](#)

[Maintenance  
Support](#)

[Troubleshooting  
Media](#)

[Specifications](#)

[Warranty &](#)

## 1. PRODUCT OVERVIEW

The Escali High Capacity Digital Bathroom Scale (Model: XLSCALE) provides reliable and accurate body weight measurement. It features an extra-high capacity of 562 pounds, a slip-resistant glass surface, and a large, backlit display for clear reading. The scale incorporates instant-on technology and an audible alert to confirm weight readings.

### Key Features:

- **High Capacity:** Supports up to 562 lb (255 kg).
- **Anti-Slip Surface:** Etched glass platform provides secure footing.
- **Extra-Large Backlit Display:** Ensures clear readability in various lighting conditions.
- **Instant-On Technology:** Activates automatically upon stepping on the scale.
- **Audible Measurement Alert:** Beeps to confirm a stable weight reading.
- **Automatic Shut-Off:** Conserves battery life by turning off after inactivity.

## 2. SETUP

### 2.1 Unpacking and Initial Preparation

1. Carefully remove the scale from its packaging.
2. **Remove Protective Wedges:** Locate and remove any small cardboard or plastic wedges inserted between the scale's components for shipping protection. Failure to remove these may prevent the scale from functioning correctly.
3. Ensure the scale is placed on a hard, flat surface. Avoid carpets or uneven flooring, as this can affect accuracy.

## 2.2 Battery Installation

The scale requires 4 AAA batteries, which are included. The battery compartment is located on the underside of the scale.

1. Open the battery compartment cover on the bottom of the scale.
2. Insert the 4 AAA batteries, ensuring correct polarity (+/-).
3. Close the battery compartment cover securely.

## 2.3 Unit Selection

Most scales have a switch on the underside to select between pounds (lb), kilograms (kg), or stones (st). Adjust this switch to your preferred unit of measurement.

## 3. OPERATING INSTRUCTIONS

---

### 3.1 First Use and Recalibration

Before the first use, or if the scale has been moved, it is recommended to recalibrate it for accurate readings.

1. Gently step onto the scale until the display activates and shows a reading.
2. Step off the scale. The display should return to "0.0" or turn off. This completes the recalibration process.

### 3.2 Taking a Weight Measurement

1. Place the scale on a hard, flat surface.
2. Step onto the center of the scale platform with bare feet. Stand still and distribute your weight evenly.
3. The display will show your weight. Once the reading stabilizes, it may blink three times and/or emit an audible beep to confirm the measurement.
4. You can then step off the scale. The weight will remain displayed for approximately 10 seconds before the scale automatically shuts off.

*Note: For consistent readings, always weigh yourself at the same time of day, under similar conditions.*

## 4. MAINTENANCE

---

### 4.1 Cleaning

- Wipe the scale surface with a soft, damp cloth.
- Do not use abrasive cleaners or immerse the scale in water.

### 4.2 Battery Replacement

- If the display shows a "Lo" indicator or becomes dim, replace all 4 AAA batteries.
- Ensure proper disposal of old batteries according to local regulations.

### 4.3 Storage

- Store the scale in a cool, dry place.
- Avoid extreme temperatures and direct sunlight.

## 5. TROUBLESHOOTING

---

Problem	Possible Cause	Solution
Scale does not turn on.	Batteries are dead or incorrectly installed. Protective shipping wedges not removed.	Check battery orientation and replace if necessary. Ensure all shipping wedges are removed.
Inaccurate or inconsistent readings.	Scale is on an uneven or soft surface (e.g., carpet). Scale needs recalibration.	Place the scale on a hard, flat surface. Step on and off the scale to recalibrate.
"Lo" or "Err" displayed.	Low battery or overload.	Replace batteries. Ensure weight does not exceed 562 lbs (255 kg).

## 6. SPECIFICATIONS

Feature	Detail
Brand	Escali
Model Number	XLSCALE
Weight Limit	562 Pounds (255 kg)
Readout Accuracy	+/- 0.1 lb or +/- 0.1 kg
Display Type	LCD with Backlight
Material	Tempered Glass
Product Dimensions	13 x 1.5 x 13 inches
Power Source	4 AAA batteries (included)
Special Features	Auto Shut Off, Instant-On, Audible Alert


## 7. WARRANTY & SUPPORT

This Escali product typically comes with a manufacturer's warranty. For specific warranty details, including duration and coverage, please refer to the documentation included with your purchase or visit the official Escali website. For customer support, product inquiries, or troubleshooting assistance, please contact Escali customer service directly.

**Note:** Some retailers may offer extended protection plans. Check your purchase details for additional coverage options.

## 8. PRODUCT MEDIA

### 8.1 Product Images

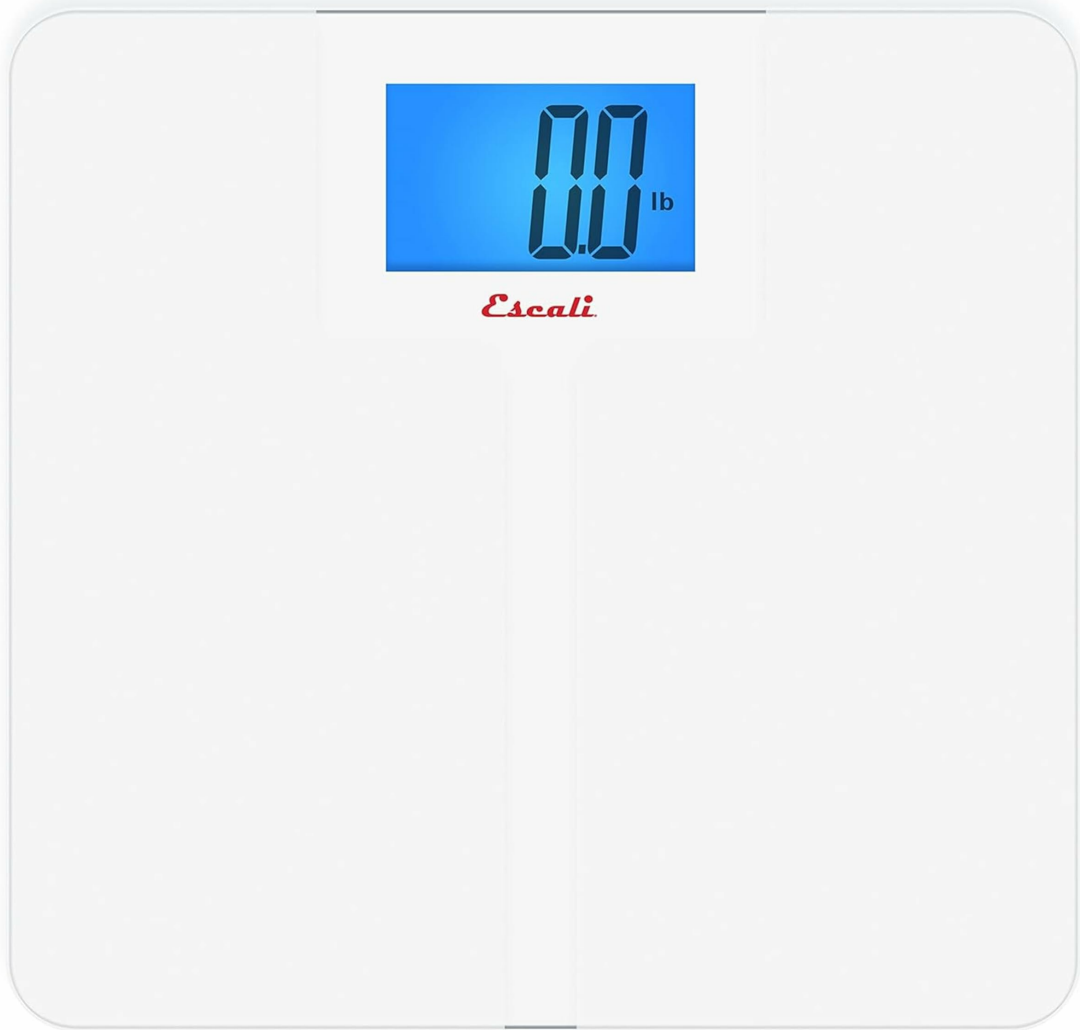

Top-down view of the Escali High Capacity Digital Bathroom Scale, showing the large display and clean white surface.



Angled view of the Escali High Capacity Digital Bathroom Scale, highlighting its sleek design and glass surface.



Close-up of the scale's digital display, showing the "0.0 lb" reading and the Escali logo.



The Escali scale positioned on a wooden floor in a bathroom setting, demonstrating its aesthetic integration.



The Escali scale shown alongside fitness accessories like a dumbbell, emphasizing its role in health tracking.



An illustrative diagram highlighting the scale's key features: secure etched surface, backlit display, 562lb capacity, and smart power display.



This high capacity bathroom scale is

## perfect for anyone

who wants a reliable and sturdy digital scale. The 562-pound capacity is among the highest available. A slip resistant glass surface ensures safe and accurate readings, while the durable construction and larger size make the scale stable and less prone to tipping.



A graphic detailing the scale's features (560lb capacity, large stable platform, extra-large backlit display, slip-resistant etched glass surface, audible measurement alert, brilliant white platform finish) and specifications (562lb capacity, 0.1lb increments, multiple units, dimensions, instant-on, display hold, auto-off, 4-AAA batteries included).



An image illustrating the dimensions of the Escali High Capacity Digital Bathroom Scale, showing its length and width.

## 8.2 Official Product Videos

Your browser does not support the video tag.

This video, titled "Escali: Measuring What's Next", provides an overview of Escali's commitment to quality and innovation in measuring devices, including kitchen scales, food thermometers, and body scales. It also highlights their social responsibility initiatives.

