

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [SYNLETT](#) /

› [SYNLETT Caregiver Pager Wireless Call Buttons User Manual](#)

## SYNLETT D091G1

# SYNLETT Caregiver Pager Wireless Call Buttons User Manual

Model: D091G1

## INTRODUCTION

---

The SYNLETT Caregiver Pager Wireless Call Button system is designed to provide reliable and portable assistance for caregivers monitoring elderly individuals, patients, or anyone requiring immediate attention. This system enhances safety and peace of mind by allowing the user to alert their caregiver with a simple press of a button.

Key features include:

- Portable caregiver pager for freedom of movement.
- Easy and simple operation with a single SOS call button.
- Flexible and versatile mounting options for call buttons (belt clip, wall mount, lanyard).
- Two distinct alert tones (DingDong and Alarm) to ensure calls are not missed.
- Effective wireless signal range up to 500 feet in open areas.
- Long-lasting battery life for both pagers and call buttons.

## PACKAGE CONTENTS

---

Your SYNLETT Caregiver Pager system includes the following components:

- 1 x Portable Caregiver Pager Receiver
- 2 x Wireless Call Buttons
- 2 x Wall Mount Brackets for Call Buttons
- 2 x Lanyards for Call Buttons
- All necessary batteries (4 AA batteries included for the pager, button batteries for call buttons).
- Mounting hardware (screws, adhesive patches).
- User Manual

# COMES WITH ALL YOU NEED



Figure 1: All components included in the SYNLETT Caregiver Pager system, including the pager, two call buttons, wall mounts, lanyards, batteries, and user manual.

## SETUP GUIDE

### 1. Battery Installation

The pager receiver requires 4 AA batteries, which are included. The call buttons come with pre-installed button cell batteries.

1. **For the Pager Receiver:** Open the battery compartment on the back of the pager. Insert 4 AA batteries, ensuring correct polarity (+/-). Close the battery compartment securely.
2. **For Call Buttons:** The call buttons are ready to use with pre-installed batteries. If replacement is needed in the future, refer to the Maintenance section.

### 2. Pager Placement

The pager receiver is portable and can be carried by the caregiver using its built-in belt clip or placed on a flat surface. For

optimal performance, ensure the pager is within the wireless signal range of the call buttons.

**IT IS NOT A DOORBELL**

# EASY TO CARRY AS YOU GO

Take it with you - a true caregiver pager

## 2 Alert Tones

Simple and Discreet  
pocket size



-  DingDong - ring once  
*less disturbing*
-  Alarm - ring for 20s  
*emergency mode*

## Portable



Belt Clip

## 3-Color Indicator

Differentiate calls by reading the light color



Adjustable Volume  
low, medium and high



**Vs**



doorbell

The receiver needs to stay plugged in.  
Inflexible...

It is just a doorbell, NOT a pager.  
It rings once only. Calls are easily missed.

Too many tone options. Most are useless or irrelevant.  
Go through the endless selection of ringtones.

Figure 2: The portable pager can be easily carried with its belt clip, offering convenience for caregivers.

### 3. Call Button Mounting Options

The call buttons offer multiple flexible mounting options to suit various needs and locations:

- **Screw Mounting:** Use the provided wall mount brackets and screws to fix the call button securely to a wall or bedside.
- **Adhesive Patch Mounting:** For a less permanent solution, use the included adhesive patches to attach the call button bracket to a clean, dry surface.
- **Lanyard Wear:** The call button can be worn around the neck as a pendant using the provided lanyard, ensuring it is always within reach. The lanyard features a breakaway safety clasp.

# FLEXIBLE AND LIGHTWEIGHT

Reliable call button with safe feature



## Screw Mounting



## Adhesive Patch Mounting



## Hang with Lanyard



Figure 3: Various methods for mounting or wearing the call button for optimal accessibility.

## OPERATING INSTRUCTIONS

### 1. Initiating a Call

To alert the caregiver, simply press the large red "Push for Help" button on the call button. The button will send a wireless signal to the pager receiver.

# PORTABLE CAREGIVER PAGER

Push the call button when in need of help or assistance.

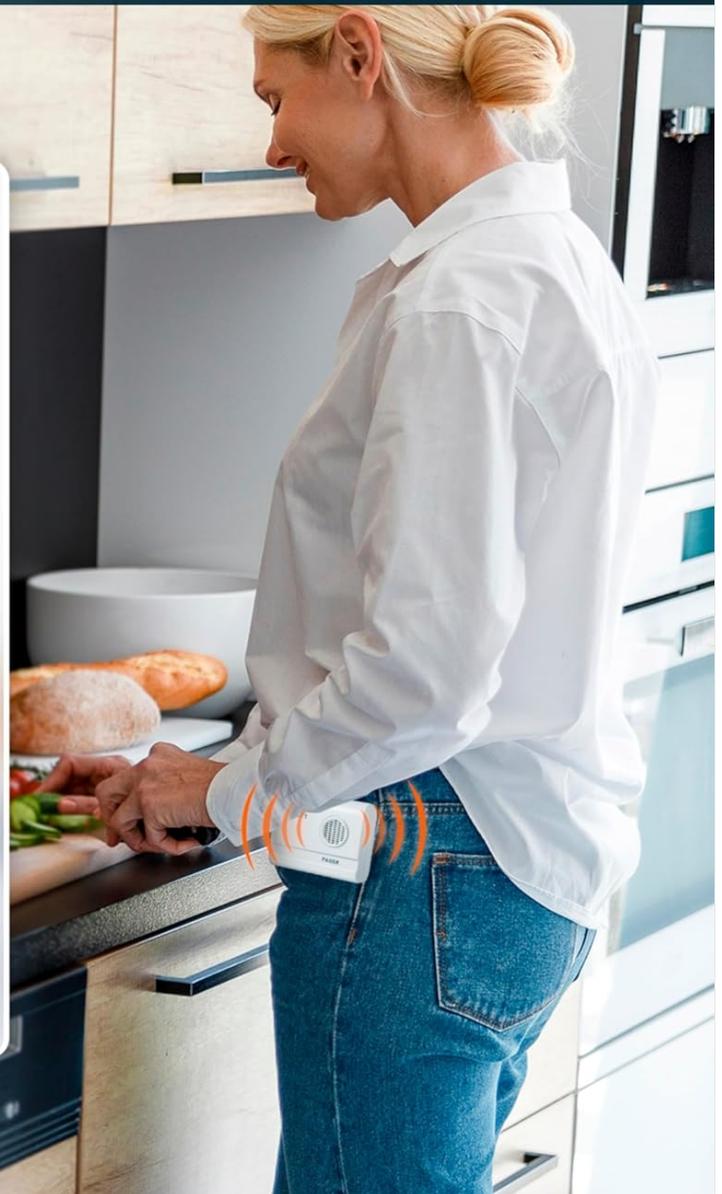


Figure 4: A user activating the call button to send an alert to the caregiver's pager.

## 2. Pager Alert Tones and Volume

The pager receiver offers two alert modes and adjustable volume settings:

- **DingDong Mode:** For a less disturbing alert, switch the pager to "DINGDONG". It will chime once when a call button is pressed.
- **Alarm Mode:** For continuous alerting in emergency situations, switch the pager to "ALARM". It will emit a continuous siren until the alert is acknowledged or the button is pressed again.
- **Volume Control:** Adjust the volume using the dial on the side of the pager (low, medium, high).

Each call button can be differentiated by a 3-color indicator on the pager, allowing the caregiver to identify which button was pressed (e.g., red, green, blue). This feature helps in quickly locating the person needing assistance if multiple buttons are in use.

## 3. Wireless Range

The wireless signal range between the pager and call buttons is up to 500 feet in open areas. Environmental factors such as concrete or brick walls and heavy electronic equipment may affect the signal range. It is recommended to test the range in your specific environment after setup.

# MAKE CAREGIVING EASY

The illustration is divided into four quadrants. The top-left quadrant shows a caregiver sitting on a hospital bed, with yellow concentric arcs representing a wireless signal emanating from a pager on their chest. The top-right quadrant shows a caregiver in green scrubs holding a white pager device, with red concentric arcs representing a wireless signal. The bottom-left quadrant shows a bathroom with a toilet and sink, with yellow concentric arcs representing a wireless signal emanating from a call button on the wall. The bottom-right quadrant is a grey background with text.

Caregiver Pager

Wireless Range Up to 500ft\*

\*Environmental factors such as concrete or brick walls and heavy electronic equipment may affect the signal range.

Figure 5: Illustration of the effective wireless range of the SYNLETT Caregiver Pager system.

## MAINTENANCE

### 1. Battery Replacement

- **Pager Receiver:** The pager uses 4 AA alkaline batteries. When the sound becomes softer or the pager stops responding, it's time to replace the batteries. Open the battery compartment on the back and replace all 4 batteries with new ones, ensuring correct polarity.
- **Call Buttons:** A blinking LED light on the call button indicates a low battery. To replace the button cell battery (CR2032, typically), carefully open the call button casing using a small screwdriver. Replace the old battery with a new one, observing polarity, and securely close the casing.

Note: Always dispose of used batteries responsibly according to local regulations.

## 2. Cleaning

To clean the pager and call buttons, use a soft, dry cloth. For stubborn dirt, a slightly damp cloth with mild soap can be used, followed by immediate drying. Do not use abrasive cleaners or immerse the devices in water.

## TROUBLESHOOTING

| Problem                                      | Possible Cause   | Solution   |
|--|--|--|
| Pager does not sound when button is pressed. | <ul style="list-style-type: none"><li>• Pager batteries are low or dead.</li><li>• Call button battery is low or dead.</li><li>• Pager is out of range.</li><li>• Pager is switched off.</li></ul> | <ul style="list-style-type: none"><li>• Replace pager batteries.</li><li>• Check for blinking LED on call button; replace battery if needed.</li><li>• Move pager closer to the call button.</li><li>• Ensure pager is switched to "DINGDONG" or "ALARM" mode.</li></ul> |
| Pager alert sound is too quiet or too loud.  | Volume setting is incorrect.   | Adjust the volume dial on the side of the pager.   |
| Call button does not stay in its holder.     | Holder may be slightly oversized or button not seated correctly.   | Ensure the button is firmly clicked into the holder. If issues persist, consider using screw mounting or the lanyard option for secure placement.  |
| Alarm sound is too aggressive.               | The "ALARM" mode is designed for emergency.  | Switch the pager to "DINGDONG" mode for a single chime alert. Adjust volume to the lowest setting if "ALARM" mode is necessary but less intrusive sound is desired.  |

## SPECIFICATIONS

- **Model Number:** D091G1
- **Product Dimensions:** 6.14 x 3.94 x 2.44 inches (Pager)
- **Item Weight:** 8.8 ounces
- **Power Source:** Battery Powered
- **Pager Batteries:** 4 AA batteries (included)
- **Call Button Batteries:** Button cell (pre-installed, typically CR2032)
- **Control Method:** Touch
- **Mounting Type:** Wall Mount, Portable (belt clip), Lanyard
- **Maximum Wireless Range:** Up to 500 Feet (in open area)
- **Manufacturer:** SYNLETT
- **UPC:** 768497814479

## IMPORTANT SAFETY INFORMATION AND DISCLAIMERS

- **CHOKING HAZARD:** Small parts. Not for children under 3 years.
- **Supervision:** This product is not a substitute for full-time supervision and monitoring of your loved one or patient. It

is designed to augment personal supervision by a caregiver.

- **Read Instructions:** Always read the instructions and legal disclaimer prior to using the product.
- **Voltage:** The system operates on 12 volts.

## WARRANTY AND SUPPORT

---

**Warranty:** This SYNLETT product comes with a 1-year warranty from the date of sale, covering manufacturing defects.

**Customer Support:** For any questions, technical assistance, or warranty claims, please contact SYNLETT customer support. SYNLETT is committed to providing cost-effective and easy-to-use call solutions to improve care for your loved ones.

For more information or to contact support, please visit the official SYNLETT brand store: [SYNLETT Store on Amazon](#).