

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [Philips](#) /
- › [Philips Lumea IPL Hair Removal Device BRI945/00 User Manual](#)

Philips BRI945/00

Philips Lumea IPL Hair Removal Device BRI945/00 User Manual

Model: BRI945/00

1. INTRODUCTION

This manual provides essential information for the safe and effective use of your Philips Lumea IPL Hair Removal Device, model BRI945/00. Please read all instructions carefully before first use and retain this manual for future reference.

1.1 Product Overview

The Philips Lumea IPL device utilizes Intense Pulsed Light (IPL) technology for long-term hair reduction. It features SenseIQ technology, a SmartSkin Sensor, and specialized attachments for different body areas.



Image 1.1: Philips Lumea IPL Hair Removal Device BRI945/00, showing the device, its packaging, and the included body and face attachments.

1.2 What's in the Box

- 1x Philips Lumea IPL hair removal device 8000 Series (BRI945/00)
- 2x Attachments (Body and Face)
- 1x Storage Pouch
- 1x Power Adapter
- 1x Cleaning Cloth
- 1x User Manual (this document)

2. SAFETY INFORMATION

Before using the device, it is crucial to read and understand all safety warnings and instructions. Failure to do so may result in injury or damage to the device.

- **Skin Contact Sensors:** The device is equipped with skin contact sensors to prevent unintentional flashing. The device will only emit a flash when it is in full contact with the skin.
- **Eye Protection:** Goggles are not required during treatment due to the integrated skin contact sensors. However, avoid looking directly into the light flash.
- **Skin Tone Suitability:** The device is suitable for skin tones I-V (dark blonde, brown, and black hair). It is not suitable for very dark skin tones or very light hair colors (red, light blonde, white/grey). Refer to the skin tone chart in the Operating section.
- **Medical Conditions:** Consult a healthcare professional before use if you have any underlying medical conditions, skin conditions, or are taking medications that may increase skin sensitivity.
- **Pregnancy and Breastfeeding:** Do not use the device if you are pregnant or breastfeeding.
- **Children:** Keep the device out of reach of children.
- **Cleaning:** Always unplug the device before cleaning. Use the provided cleaning cloth.



Built for safety
Skin contact sensors
prevent unintentional
flashing



Image 2.1: The device features skin contact sensors to ensure safe operation by preventing flashes unless in full contact with the skin.

3. SETUP

1. **Prepare Skin:** Ensure the treatment area is clean, dry, and shaved. Do not use waxing, epilating, or plucking methods immediately before treatment, as the hair root needs to be present for IPL to be effective.
2. **Connect Power:** Plug the power adapter into the device and then into a wall outlet. The device is corded and requires continuous power during use.
3. **Select Attachment:** Choose the appropriate attachment for the body area you intend to treat. The device comes with specific attachments for the body and face.
4. **Power On:** Press the power button to turn on the device.

4. OPERATING INSTRUCTIONS

4.1 SmartSkin Sensor and Light Settings

The SmartSkin Sensor automatically detects your skin tone and recommends the most comfortable light setting. There are five adjustable light settings.

1. After powering on, place the device on your skin. The SmartSkin Sensor will analyze your skin tone.
2. The device will indicate the recommended light setting. You can manually adjust this setting if desired, but always start with the recommended level.
3. Ensure the light setting is comfortable. If you experience discomfort, reduce the setting.



Personalized light setting

Integrated sensors measure your skin tone for a comfortable treatment

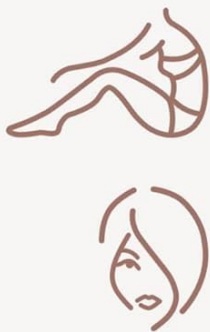


Image 4.1: The SmartSkin Sensor automatically detects skin tone to suggest an optimal light setting for comfortable and effective treatment.

4.2 Using Attachments

The device includes uniquely curved attachments designed to fit different body contours and optimize treatment for specific areas.

- **Body Attachment:** Use for larger areas like legs, arms, and abdomen.
- **Face Attachment:** Use for sensitive areas like the upper lip, chin, and jawline.



Intelligent curved attachments

Designed for your body's curves

Image 4.2: The device comes with specialized attachments, each designed to conform to different body curves for targeted treatment.



Image 4.3: Example of device application on the armpit area.

4.3 Treatment Procedure

1. Place the device window firmly against your skin. The "ready to flash" indicator light will illuminate.
2. Press the flash button. A flash of light will be emitted.
3. Move the device to the next adjacent area, ensuring slight overlap to cover all areas. Do not treat the same area multiple times in one session.
4. Continue until the entire desired area has been treated.

4.4 Treatment Schedule

For optimal results, follow the recommended treatment schedule:

- **Initial Phase:** Perform 4 treatments, one every two weeks.
- **Maintenance Phase:** After the initial phase, perform touch-up treatments monthly to maintain results.

Suitable on skin tones I–V

Suitable for dark blonde, brown and black hair and for skin tones I–V.

Suitable skin tones



Suitable body hair colors



Treatment schedule

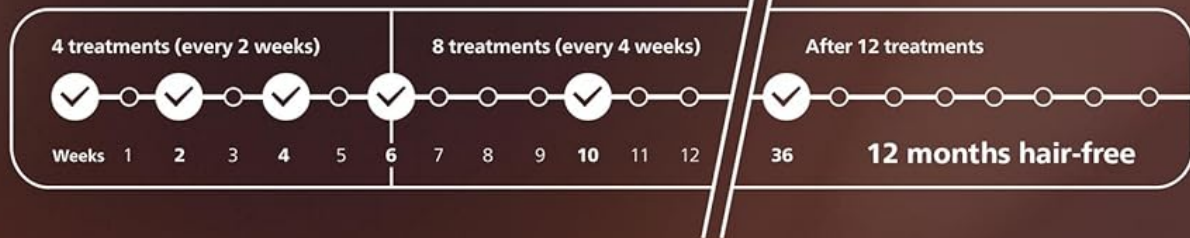


Image 4.4: This diagram illustrates suitable skin tones (I-V) and hair colors (dark blonde to black) for effective IPL treatment, along with the recommended treatment schedule for initial and maintenance phases.



Image 4.5: Visual representation of hair reduction progress: before first treatment, after three treatments, and after full treatment.

The Philips Lumea IPL app can assist you in creating a personalized treatment plan, tracking sessions, and receiving reminders.

5. MAINTENANCE

- **Cleaning the Device:** After each use, unplug the device and allow it to cool down. Use the provided cleaning cloth to gently wipe the treatment window and the exterior of the device. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device and its attachments in the provided storage pouch in a cool, dry place, away from direct sunlight and moisture.

6. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Not properly plugged in.	Ensure the power adapter is securely connected to both the device and the wall outlet.
Device does not flash.	Skin contact sensors not fully engaged.	Ensure the treatment window is in full and firm contact with the skin. The "ready to flash" indicator should be lit.
Discomfort during treatment.	Light setting is too high for your skin tone or sensitivity.	Reduce the light intensity setting. If discomfort persists, discontinue use and consult a healthcare professional.
No hair reduction after several treatments.	Incorrect treatment schedule or unsuitable hair/skin type.	Review the treatment schedule and ensure consistent application. Verify your skin tone and hair color are suitable for IPL treatment (refer to Image 4.4).

7. SPECIFICATIONS

Model Number: BRI945/00

Brand: Philips

Technology: Intense Pulsed Light (IPL) with SenseIQ

Power Source: Corded Electric

Attachments: Body and Face

Skin Tone Sensor: SmartSkin Sensor

Light Settings: 5 adjustable settings

Product Dimensions: 16.5 x 10.5 x 22 cm

Item Weight: 1.1 kg

Country of Origin: Slovenia

Batteries: 1 Lithium Ion battery required (included, likely for internal sensor/memory, not primary power)

8. WARRANTY AND SUPPORT

The Philips Lumea IPL Hair Removal Device BRI945/00 typically comes with a **2-year warranty**. Please retain your original receipt as proof of purchase for warranty claims.

For further assistance, troubleshooting, or warranty inquiries, please contact Philips customer support or visit the official Philips website. Refer to the packaging or the Philips website for specific contact details in your region.

You can also utilize the Philips Lumea App for personalized guidance and support related to your device.