

## Ejoyous Ejoyous3pkgid7e65135-02

# Ejoyous Microdermabrasion Filter Instruction Manual

Model: 15mm Replacement Cotton Filters

## INTRODUCTION

---

This instruction manual provides guidance for the proper use and maintenance of Ejoyous Microdermabrasion Replacement Cotton Filters. These filters are designed to be used with various beauty machines, such as blackhead removal machines and dermabrasion machines, to ensure hygienic operation and protect the machine's internal components.

## SETUP AND INSTALLATION

---

The Ejoyous Microdermabrasion Filters are designed for easy installation into compatible beauty machines. Follow these steps to properly install a new filter:

1. Ensure your beauty machine is turned off and disconnected from its power source before handling any components.
2. Locate the filter housing or chamber on your microdermabrasion wand or device. This is typically a removable part at the tip or base of the wand.
3. Carefully unscrew or detach the component that holds the filter. Remove any used or old filter pads.
4. Take one new Ejoyous cotton filter (15mm diameter for this model) from the package.
5. Place the new filter pad securely into the designated slot or chamber. Ensure it lies flat and covers the opening completely to prevent debris from entering the machine.
6. Reattach the filter housing or component, ensuring it is tightened properly to prevent air leaks during operation.



Image: Disassembled microdermabrasion wand with a filter pad ready for insertion. This illustrates where the filter is placed within the device.



Image: A hand demonstrating the insertion of a cotton filter into the tip of a microdermabrasion wand. This shows the practical step of placing the filter.

## OPERATION

Once the filter is correctly installed, your beauty machine is ready for operation. The filter's primary function is to prevent impurities, dead skin cells, and other debris generated during treatments from being drawn into the machine's main motor. This protection is crucial for maintaining the machine's performance and extending its lifespan.

- Always ensure a clean, new filter is in place before each treatment session.
- Operate your beauty machine according to its specific manufacturer's instructions.
- The filter acts as a barrier, collecting particles and ensuring only clean air passes through the machine's internal mechanisms.

# Filter Cotton



Specially designed for microdermabrasion machines' filtration

Image: A magnified view of the filter cotton material, highlighting its texture designed for effective filtration in microdermabrasion machines.

## MAINTENANCE AND REPLACEMENT

Regular replacement of the cotton filters is essential for optimal performance and hygiene. These filters are disposable and designed for single-use or replacement after each treatment.

- **Frequency:** It is strongly recommended to replace the filter after each treatment session. This prevents cross-contamination and ensures maximum filtration efficiency.
- **Disposal:** Used filters should be disposed of properly according to local waste disposal guidelines. The cotton material is designed to be easily decomposed, reducing waste.
- **Inspection:** Before each use, visually inspect the filter housing and ensure no debris has bypassed the filter. If any damage or excessive wear is observed in the housing, consult your beauty machine's manual or manufacturer.



Image: A quantity of new, unused cotton filters, illustrating the supply provided for regular replacement.

## TROUBLESHOOTING

These cotton filters are passive components and typically do not require troubleshooting. However, if you experience issues related to filtration or machine performance, consider the following:

- **Reduced Suction:** If your beauty machine experiences reduced suction, first check if the filter is clogged. Replace it with a new one. Also, ensure the filter is correctly seated and the housing is sealed.
- **Debris in Machine:** If you observe debris inside the machine's internal components, it may indicate an improperly installed filter or a filter that has been used beyond its recommended lifespan. Ensure proper installation and timely replacement.
- **Incorrect Size:** Ensure you are using the correct filter size (15mm for this product) for your specific beauty machine model. Using an incorrect size can lead to ineffective filtration or damage.

## SPECIFICATIONS

Brand	Ejoyous
-------	---------

<b>Model Number</b>	Ejoyous3pkgid7e65135-02
<b>Material</b>	Cotton
<b>Color</b>	White
<b>Diameter (Current Product)</b>	15mm (0.59 inches)
<b>Available Diameters</b>	10mm, 15mm, 20mm, 25mm, 30mm, 40mm, 50mm
<b>Quantity per Pack</b>	100 pieces
<b>Reusability</b>	Disposable
<b>Package Dimensions</b>	3.54 x 3.54 x 0.39 inches
<b>Package Weight</b>	0.11 ounces (approx. 2-33g for various sizes)

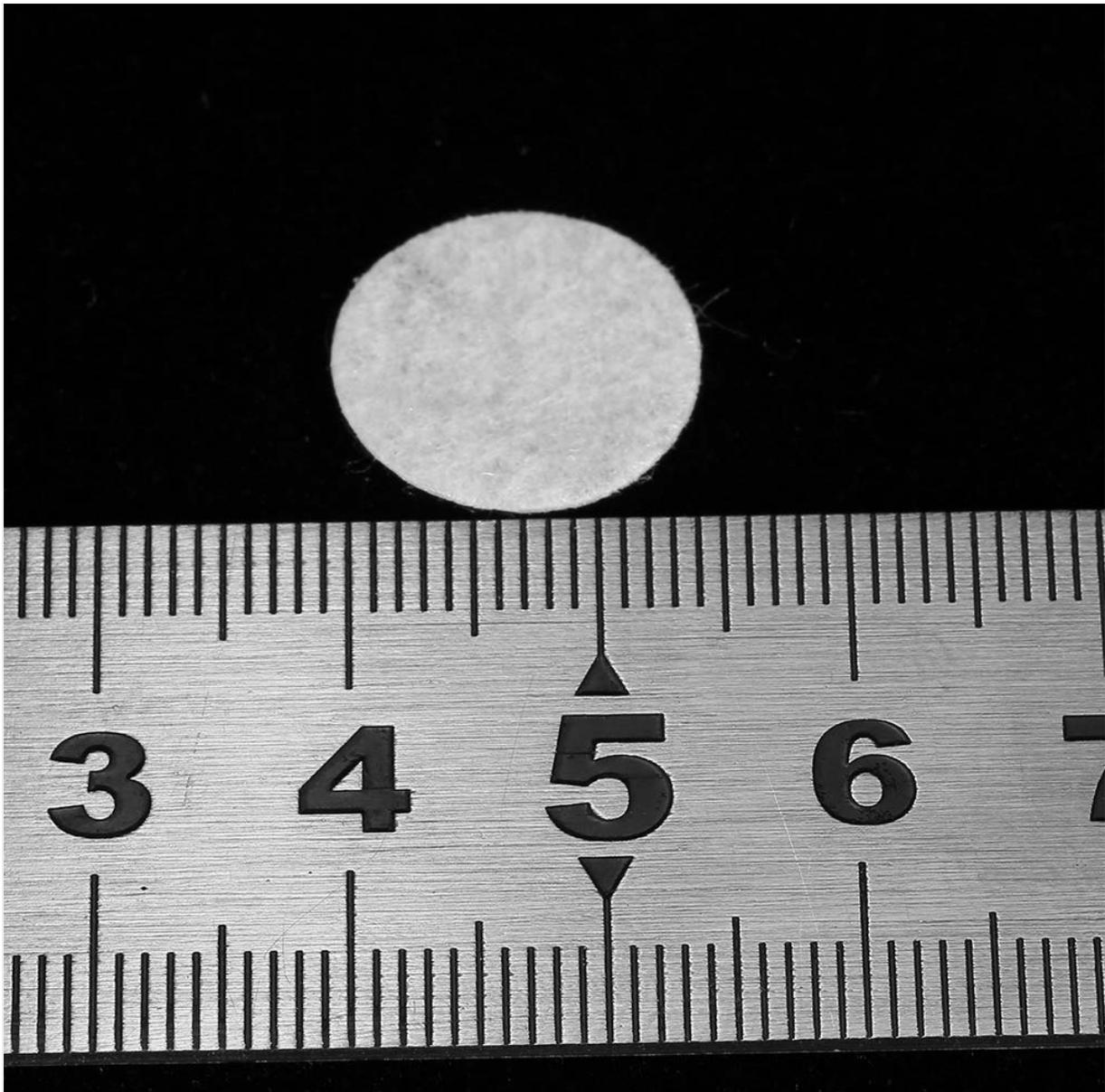


Image: A single 15mm filter pad displayed next to a ruler, indicating its precise diameter.

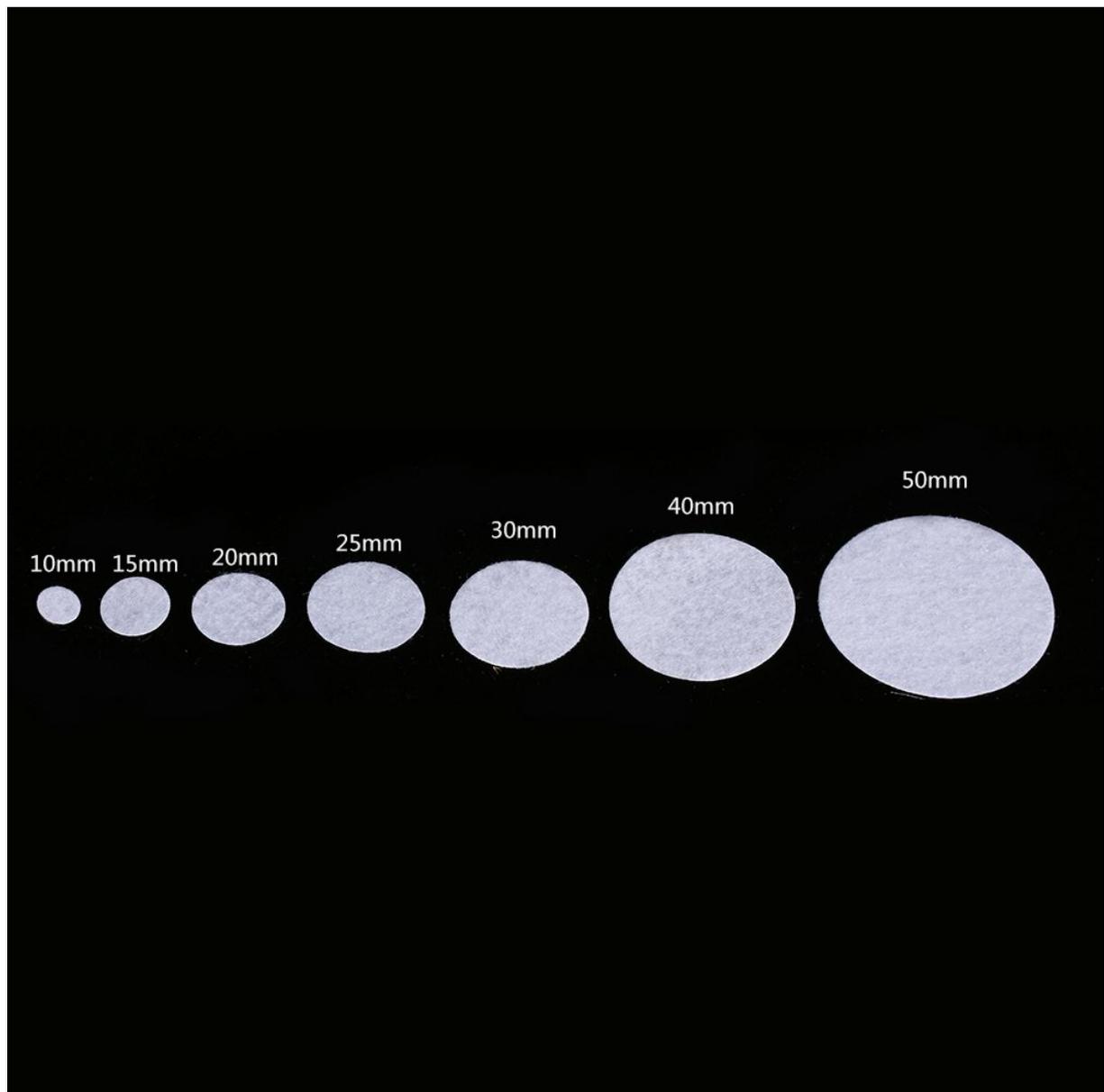


Image: An array of Ejoyous filter pads in different diameters (10mm, 15mm, 20mm, 25mm, 30mm, 40mm, 50mm), demonstrating the range of available sizes.

## WARRANTY INFORMATION

---

Warranty information for Ejoyous Microdermabrasion Filters is not specified in the product details. For any warranty-related inquiries, please contact the manufacturer directly.

## CUSTOMER SUPPORT

---

For further assistance, questions, or concerns regarding Ejoyous Microdermabrasion Filters, please refer to the contact information provided by your retailer or the Ejoyous brand directly. You may also visit the [Ejoyous Store on Amazon](#) for more product information.