



Manuals.plus /

› **Magic Wand /**

› Magic Wand Mini HV-135 Cordless Massager User Manual

Magic Wand HV-135

Magic Wand Mini HV-135 Cordless Massager User Manual

Model: HV-135 | Brand: Magic Wand

INTRODUCTION

The Magic Wand Mini is a compact, cordless personal massager designed for various applications including neck, hand, back, muscle, and head massage. This manual provides essential information for the safe and effective use, maintenance, and care of your device.

It features a powerful 6,000 RPM brushless motor, multiple speed settings, and a soft silicone head with a flexible neck for comfortable use. The rechargeable battery offers up to 2.5 hours of cordless operation.



Image: The Magic Wand Mini HV-135 massager, shown in white with blue control buttons and a soft white silicone head. This image illustrates the overall design of the product.

SAFETY INFORMATION

- This device is registered with the FDA as a medical device.
- Do not use in or near water.
- Keep out of reach of children.
- Always follow safe practices when using personal care devices.
- Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

SETUP AND CHARGING

1. **Initial Charge:** Charge the device fully before its first use.
2. **Recharging:** Recharge the device after every 2.5 hours of cordless use.
3. **Charging Indicator:** Refer to the device's indicator lights for charging status (specifics not provided, general instruction).

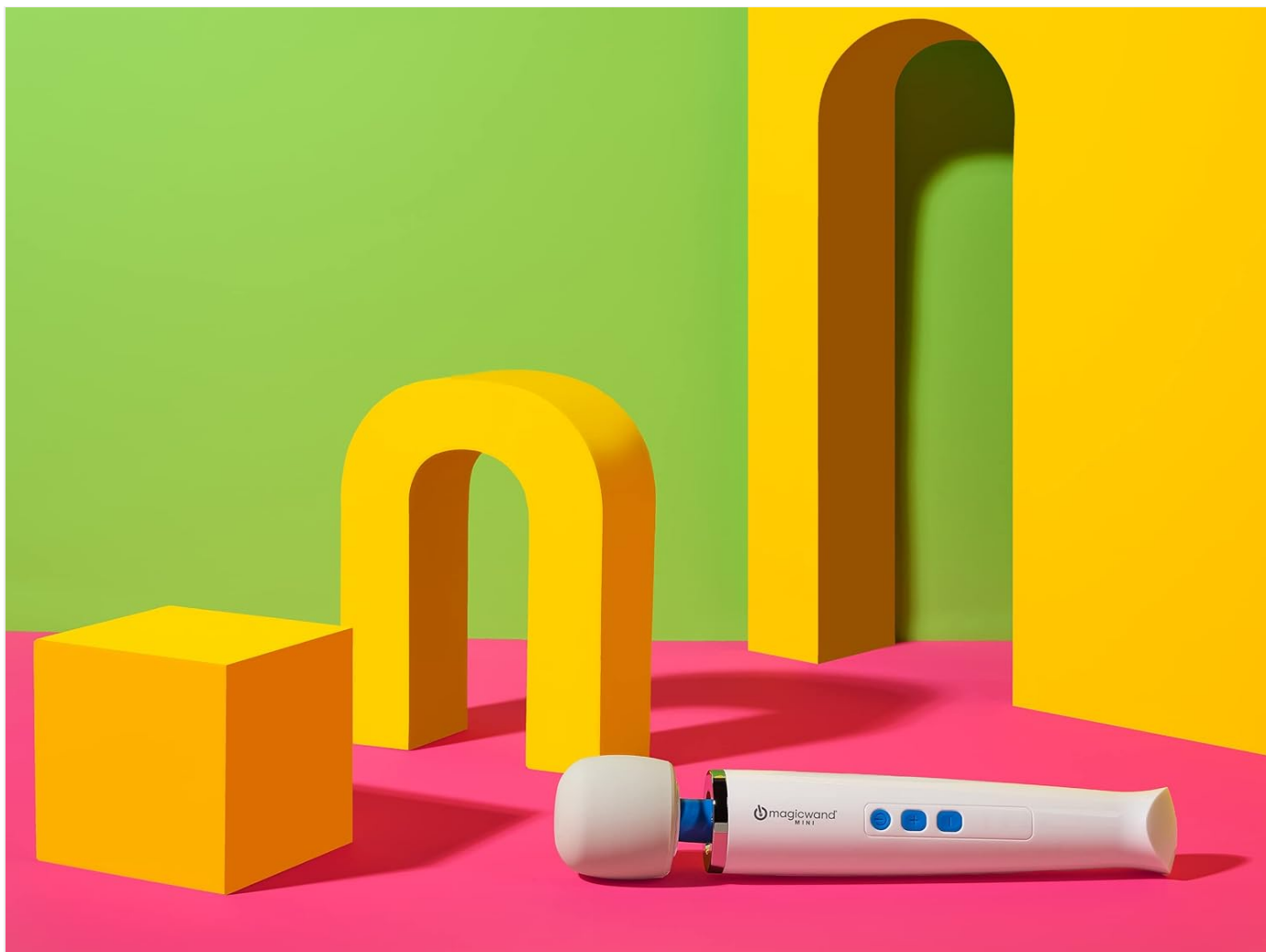


Image: The power adapter and charging cord for the Magic Wand Mini. The cord is detachable from the body massager for easy storage.

OPERATING INSTRUCTIONS

1. **Preparation:** Apply a small amount of personal lubricant to the silicone head if desired.
2. **Power On:** Press the power button to turn on the device.
3. **Adjusting Speed:** Use the "+" and "-" buttons to cycle through the 3 available vibration intensity levels.
4. **Application:** Place the massager head against the desired area and apply gentle pressure.
5. **Power Off:** To turn off the device, press and hold the power button for 3 seconds.



High Capacity Motor



Up to 2.5H of Full Charged Battery



3-Speeds (3500-6000RPM)



Soft Silicone Head



 magicwand[®]
MINI

Image: The Magic Wand Mini, highlighting its key features including a high-capacity motor, 2.5 hours battery life, 3 speeds (3500-6000 RPM), and a soft silicone head.

MAINTENANCE AND STORAGE

- **Cleaning:** Clean the device after each use with a toy cleaner spray and a soft cloth.
- **Storage:** Store the massager in a cool, dry place away from direct sunlight.

TROUBLESHOOTING

If you experience any issues with your Magic Wand Mini, please ensure the device is fully charged and cleaned according to the instructions. For persistent problems or technical support, please refer to the warranty and support section for contact information.

SPECIFICATIONS

Feature	Detail
Product Dimensions	12.2 x 2.36 x 1.97 inches
Weight	1.07 Pounds
Item Model Number	HIT135
Batteries	1 Nonstandard Battery required
Manufacturer	Magic Wand
Color	Mini (White with blue accents)
Material	Silicone (head)
Model Name	HP-HIT135
Speeds	3 intensity levels (3500-6000 RPM)
Run Time	Up to 2.5 hours (cordless)

WARRANTY AND SUPPORT

The Magic Wand Mini HV-135 includes a **1-year warranty**. For warranty claims, technical assistance, or customer support, please visit the official Magic Wand website or contact their customer service department. You can find more information and support resources by visiting the Magic Wand Store on Amazon.

