

[Manuals.plus](#) /

> [beatXP](#) /

> beatXP SmartPlus Prime Weight Machine Instruction Manual

beatXP SmartPlus Prime

beatXP SmartPlus Prime Weight Machine Instruction Manual

Model: SmartPlus Prime

1. INTRODUCTION

This manual provides instructions for the proper use and maintenance of your beatXP SmartPlus Prime Weight Machine. This device is designed to measure body weight and analyze various body composition parameters through bioelectrical impedance analysis (BIA) technology. Please read this manual thoroughly before use to ensure accurate readings and safe operation.

2. SAFETY INFORMATION

- Do not use this device if you have a pacemaker or other implanted medical device.
- Do not use this device if you are pregnant.
- Place the scale on a hard, flat surface for accurate measurements. Avoid carpets or soft surfaces.
- The scale surface can be slippery when wet. Ensure the surface is dry before stepping on it.
- Do not jump on the scale.
- Keep the scale away from water and extreme temperatures.
- This scale is for household use only and not intended for medical diagnosis or treatment.
- The first reading after moving the scale may be inaccurate. Please ignore it and take a second reading.

3. PACKAGE CONTENTS

- beatXP SmartPlus Prime Weight Machine
- Batteries (pre-installed or included separately)
- Instruction Manual

4. PRODUCT OVERVIEW

The beatXP SmartPlus Prime features a sleek design with a durable tempered glass top and high-precision sensors for accurate body composition analysis. It includes a VA LED screen for clear display of key parameters.



Image: Components of the beatXP SmartPlus Prime Weight Machine. This image highlights the durable tempered glass, enlarged electrode pads, anti-slip pads, and high-precision sensors that contribute to the scale's thoughtful design and durable construction.

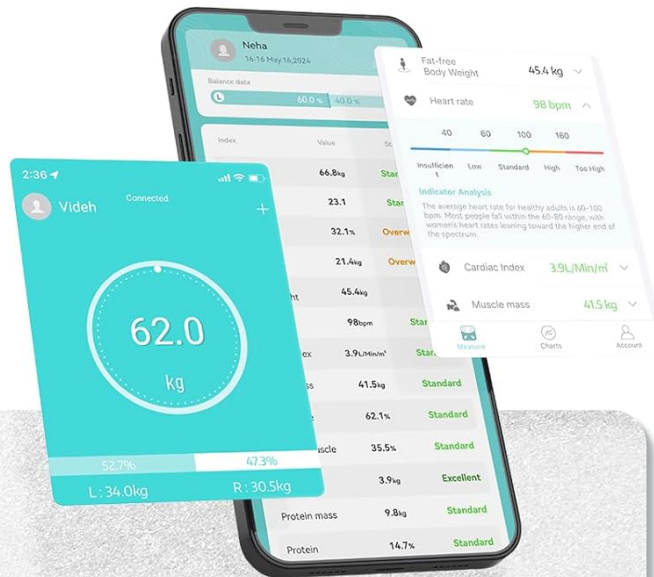
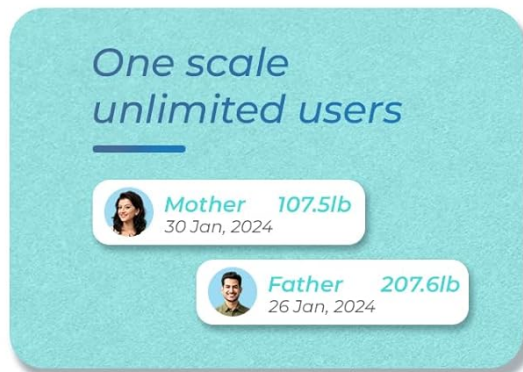


Image: Underside of the scale with anti-skid padding and battery compartment. This image also shows the scale's connectivity with the Fitdays app on a smartphone, emphasizing high accuracy and durability.

5. SETUP

5.1 Battery Installation

If batteries are not pre-installed, open the battery compartment on the underside of the scale and insert the provided batteries, ensuring correct polarity.

5.2 Placing the Scale


Place the scale on a hard, flat, and stable surface. Avoid placing it on carpets, rugs, or uneven flooring, as this can affect measurement accuracy.

5.3 App Download and Pairing

1. Download the **Fitdays** app from your smartphone's app store (Google Play Store for Android or Apple App Store for iOS). You can also scan the QR code provided in the product packaging or on the product images.
2. Enable Bluetooth on your smartphone.
3. Open the Fitdays app and create an account or log in.
4. Step onto the scale briefly to activate it. The scale will automatically attempt to pair with the app via Bluetooth.
5. Follow the on-screen prompts in the app to complete the pairing process.

Seamlessly Sync with Your Preferred Health Apps


Scan the QR code to download **Fitdays** app.





GET IT ON
Google Play


Download on the
App Store


Sync with


Google Fit


Apple Health


Fitbit


Samsung Health

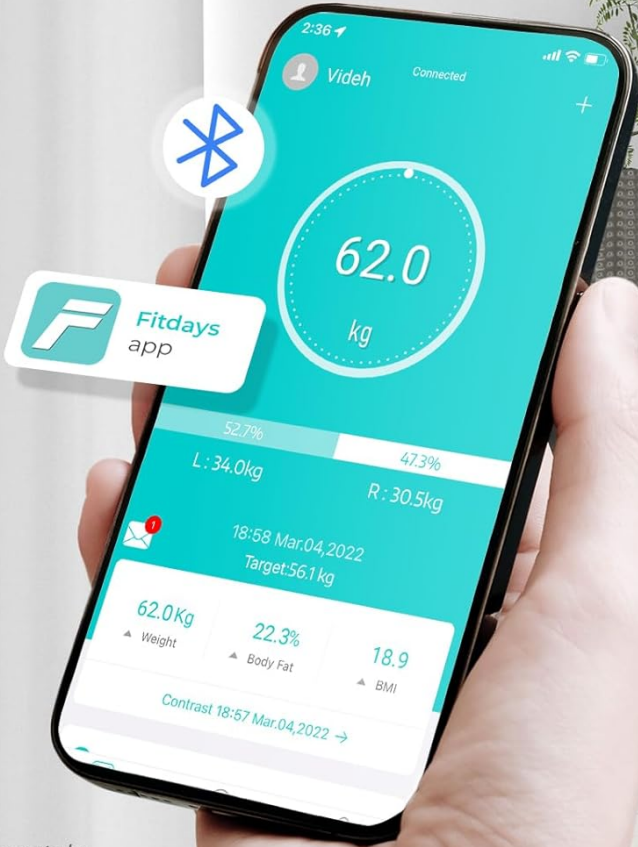

Note: Mobile phone need to be purchased separately.

Image: Seamless synchronization with health apps. This image illustrates how to download the Fitdays app via QR code and its compatibility with various health platforms like Apple Health, Google Fit, Fitbit, and Samsung Health.

For more information on the Fitdays app, visit: <http://www.fitfeel.cn/fitdays/>

6. OPERATION

6.1 Taking a Measurement

1. Ensure the scale is on a hard, flat surface.
2. Step onto the scale barefoot with your feet evenly placed on the electrode pads. Stand still.
3. The scale's VA LED display will show your weight and other key parameters.
4. If connected to the Fitdays app, your data will automatically sync to your smartphone for detailed analysis.



Image: Enhanced visibility of the scale's display. This image shows a user taking a measurement, with the large, dark LCD screen clearly displaying weight and body composition percentages. A crucial note reminds users to place the scale on a hard, flat surface for accuracy.

6.2 Understanding the Display

The VA LED screen displays your weight and up to six other health parameters. It also indicates low battery and overload conditions.

6.3 Using the Fitdays App

The Fitdays app provides a comprehensive analysis of 15 essential body parameters, including:

- Body Weight
- Body Mass Index (BMI)
- Body Fat Percentage
- Fat-Free Body Weight
- Subcutaneous Fat Percentage
- Visceral Fat
- Body Water Percentage
- Skeletal Muscle Percentage
- Muscle Mass
- Bone Mass
- Protein Percentage
- Basal Metabolic Rate (BMR)
- Body Age
- Heart Rate
- Heart Index

The app allows you to track your progress with detailed charts and graphs, set fitness goals, and manage data for up to 24 different users, making it suitable for families.

Video: [beatXP SmartPlus Prime BMI Weighing Scale Overview](#). This video demonstrates the scale's features, including its ability to measure 15 key parameters, high precision sensors, connectivity to fitness apps like Fitbit, Google Fit, Apple Health, and Samsung Health, and its sleek design. It also shows how to download the Fitdays app and track progress.

7. SPECIFICATIONS

Parameter	Value
Brand	beatXP
Model Number	BXPBMI1002_TEL
Weight Limit	180 Kilograms
Accuracy	0.5 kg
Measurement Type	kg, lbs
Display Type	VA LED
Material Type	Acrylonitrile Butadiene Styrene (ABS)
Item Weight	1.25 Kilograms
Connectivity	Bluetooth

Body Parameters Measured	15 (Body Weight, BMI, Body Fat Percentage, Fat-Free Body Weight, Subcutaneous Fat Percentage, Visceral Fat, Body Water Percentage, Skeletal Muscle Percentage, Muscle Mass, Bone Mass, Protein Percentage, Basal Metabolic Rate, Body Age, Heart Rate, Heart Index)
--------------------------	---

8. TROUBLESHOOTING

8.1 Inaccurate Readings

- Ensure the scale is placed on a hard, flat surface. Avoid carpets.
- Make sure your feet are clean, dry, and making full contact with the electrode pads.
- The first reading after moving the scale may be inaccurate. Take a second reading.
- Check battery levels and replace if necessary.

8.2 App Connectivity Issues

- Ensure Bluetooth is enabled on your smartphone.
- Make sure the Fitdays app is open and running in the foreground.
- Restart your smartphone and the scale.
- Ensure the scale is within the Bluetooth range of your smartphone.
- If issues persist, try unpairing and re-pairing the device in the Fitdays app settings.

8.3 Display Not Working

- Check if the batteries are correctly installed and have sufficient charge. Replace if needed.
- Ensure there is no obstruction on the scale surface.

9. CARE AND MAINTENANCE

- Clean the scale surface with a soft, damp cloth. Do not use abrasive cleaners or immerse the scale in water.
- Store the scale in a cool, dry place.
- Remove batteries if the scale will not be used for an extended period.

10. WARRANTY AND SUPPORT

Your beatXP SmartPlus Prime Weight Machine comes with an **18-month replacement warranty** on manufacturing defects. To activate your warranty or for any support inquiries, please contact beatXP customer care.

Manufacturer Contact Information: GHV Medical Anchor Private Limited, Plot N0 1453. ST Tower Old Delhi Gurgaon Road -122008, Contact - 91-7303762766, email support@beatxp.com

