



Manuals.plus /

› YUNCHI /

› YUNCHI Sonic Electric Toothbrush Y2 User Manual

## YUNCHI Y2

# YUNCHI Sonic Electric Toothbrush Y2 User Manual

Model: YC-Y2-08B

## INTRODUCTION

---

This manual provides detailed instructions for the safe and effective use of your YUNCHI Sonic Electric Toothbrush, Model Y2. Please read this manual thoroughly before first use and retain it for future reference.

## SAFETY INFORMATION

---

- Not suitable for children under 3 years old. Adult supervision is required for children.
- Do not use this device for purposes other than oral hygiene.
- Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.
- Ensure the charging port cover is securely closed when not charging to maintain waterproof integrity.

## PACKAGE CONTENTS

---

The YUNCHI Sonic Electric Toothbrush Y2 package includes:

- 1 x YUNCHI Electric Toothbrush Handle
- 8 x Charcoal Replacement Toothbrush Heads
- 1 x USB Charging Cable
- 1 x User Manual



Image: YUNCHI Sonic Electric Toothbrush Y2, showing the handle, eight charcoal brush heads, and the product packaging.

## SETUP

### 1. Attaching the Brush Head

1. Align the brush head with the shaft of the toothbrush handle.
2. Gently push the brush head onto the shaft until it clicks into place. Ensure it is securely attached.

# YUNCHI

## Sonic Electric Toothbrush Idea for adult & Kids

DEEP CLEANING BETWEEN TEETH  
FAST & EFFECTIVE

30 Days long lasting  
Travel case



**3D BRISTLE  
CUTTING PROCESS  
40,000 VPM**

Image: Close-up view of the YUNCHI Sonic Electric Toothbrush Y2 brush head, highlighting the 3D bristle cutting process and charcoal bristles.

## 2. Charging the Toothbrush

Before first use, fully charge the toothbrush. A full charge takes approximately 4 hours and provides up to 30 days of normal usage (2 minutes, twice per day).

1. Open the waterproof cover at the bottom of the toothbrush handle.
2. Insert the provided USB charging cable into the Type-C charging port.
3. Connect the other end of the USB cable to a compatible power source (e.g., power bank, phone charger, PC).
4. The indicator lights will show charging progress. Once fully charged, disconnect the cable and securely close the waterproof cover.

**IP X7**  
**WATER PROOF**  
Can be used while shower

YUNCHI

White  
Clean  
Sensitive  
Polish  
Massage

Power bank  
Phone charger  
PC

**POWERFUL BATTERY**  
**4 HOURS FAST CHARGE LAST UP TO 30 DAYS**  
**ONCE CHARGE FOR 60 DAYS**

**GREAT FOR TRAVEL**

Image: YUNCHI Sonic Electric Toothbrush Y2 demonstrating its IPX7 waterproof feature and various USB charging options like power bank, phone charger, and PC.

## OPERATING INSTRUCTIONS

### 1. Power On/Off and Mode Selection

The toothbrush features 5 brushing modes: White, Clean, Sensitive, Polish, and Massage.

1. Press the power button once to turn on the toothbrush. It will start in the default 'White' mode.
2. Press the power button repeatedly to cycle through the different modes. The indicator light next to the selected mode will illuminate.
3. To turn off the toothbrush, press and hold the power button for 2 seconds, or cycle through all modes until it turns off automatically.



Image: YUNCHI Sonic Electric Toothbrush Y2 handle, highlighting the power button and the five distinct brushing modes: White, Clean, Sensitive, Polish, and Massage.

## 2. Brushing Technique and Smart Timer

The Y2 toothbrush incorporates a 2-minute smart timer and a 30-second interval reminder to ensure thorough cleaning.

1. Wet the brush head and apply a small amount of toothpaste.
2. Place the brush head against your teeth at a slight angle towards the gum line.
3. Press the power button to start. Gently guide the brush head over each tooth surface.
4. The toothbrush will pause briefly every 30 seconds to signal you to move to the next quadrant of your mouth.
5. After 2 minutes, the toothbrush will automatically turn off, indicating the completion of the recommended brushing time.

# GREAT FOR ADULTS & KIDS



2 MIN SMART TIMER



30S INTERVAL PAUSE



Never miss a zone  
Auto stop after 2 mins

Image: A visual representation of the YUNCHI Sonic Electric Toothbrush Y2's 2-minute smart timer and 30-second interval pause feature, guiding users to brush each quadrant of their mouth evenly.

### 3. Benefits of Charcoal Bristles and Sonic Technology

- The charcoal-infused DuPont nylon bristles are designed to help remove surface stains, contributing to whiter teeth and fresher breath.
- The sonic technology delivers up to 40,000 vibrations per minute, providing a deep clean that effectively removes plaque.
- End-rounded bristles ensure safe and gentle cleaning, minimizing irritation to gums.

# ACTIVATED CHARCOAL BRISTLES

For whiter teeth, healthier gums  
fresher breath




SONIC TECHNOLOGY  
**40000**  
Micro-brushes per minute



## Charcoal infused bristles

helps remove surface stains  
whiten your teeth

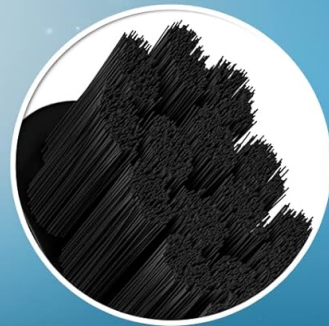

## End-Rounded Bristles

Safe & Extra Gentle on Gums  
Avoid Bleeding Gums

Image: Close-up of the YUNCHI Sonic Electric Toothbrush Y2's activated charcoal bristles and a graphic indicating 40,000 micro-brushes per minute sonic technology.

## MAINTENANCE

### 1. Cleaning the Toothbrush

After each use, rinse the brush head thoroughly under running water. The toothbrush handle is IPX7 waterproof, allowing it to be rinsed directly under the faucet or used in the shower.

### 2. Replacing Brush Heads

Dentists recommend replacing brush heads every 3 months. The YUNCHI Y2 comes with 8 replacement heads, providing up to 2 years of use.

1. To remove a brush head, gently twist it counter-clockwise and pull it off the shaft.
2. To attach a new brush head, align it with the shaft and twist clockwise until it locks into place.



Image: YUNCHI Sonic Electric Toothbrush Y2 with a display of eight activated charcoal brush heads, indicating their long-term usability.

## TROUBLESHOOTING

- **Toothbrush not turning on:** Ensure the toothbrush is fully charged. If the issue persists, try charging with a different USB port or cable.
- **Weak vibrations:** The battery may be low. Charge the toothbrush fully. If the issue continues, try a different brushing mode.
- **Brush head not staying on:** Ensure the brush head is correctly aligned and twisted into place.

## SPECIFICATIONS

Feature	Detail
Model Number	YC-Y2-08B
Vibrations Per Minute (VPM)	Up to 40,000
Brushing Modes	5 (White, Clean, Sensitive, Polish, Massage)
Timer	2-minute smart timer with 30-second interval reminders
Waterproof Rating	IPX7
Charging Type	USB Type-C
Battery Life	Up to 30 days on a 4-hour charge
Bristle Material	DuPont Nylon with Activated Charcoal
Power Source	Battery Powered (1 Lithium Polymer battery included)
Item Weight	7.2 ounces
Age Range	Adults and Kids

## WARRANTY AND SUPPORT

---

The YUNCHI Sonic Electric Toothbrush Y2 comes with a 1-year warranty.

For any questions or support, please contact the seller through your Amazon account:

1. Log in to your Amazon account.
2. Go to 'Your Orders'.
3. Find the order ID for your YUNCHI Electric Toothbrush.
4. Click 'Contact seller'.

## PRODUCT VIDEOS

---

### Electric Toothbrush for Children & Adults

Your browser does not support the video tag.

Video: A demonstration of the YUNCHI Sonic Electric Toothbrush Y2, showcasing its features and ease of use for both children and adults. The video covers unboxing, brush head attachment, mode selection, and brushing technique.

### Kids Electric Toothbrush Sonic Brush Charcoal Toothbrush

Your browser does not support the video tag.

Video: This video focuses on the YUNCHI Sonic Electric Toothbrush Y2's suitability for children, highlighting the gentle nature of the charcoal bristles and the various modes available for different sensitivities. It demonstrates the product's design and functionality.