

[Manuals.plus](#) /

› [Curaprox](#) /

› Curaprox Interdental Brush Set CPS Prime Start Instruction Manual


Curaprox CPS Prime Start

Curaprox Interdental Brush Set CPS Prime Start Instruction Manual

Model: CPS Prime Start (20-12015-500)

1. INTRODUCTION

The Curaprox Interdental Brush Set CPS Prime Start is designed for comprehensive interdental cleaning, helping to prevent gum inflammation and bad breath. This set includes five different brush sizes (CPS 06, CPS 07, CPS 08, CPS 09, CPS 11) and two ergonomic holders (UHS 409 and UHS 470) to ensure effective cleaning of various interdental spaces.

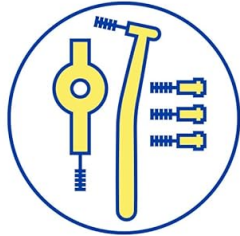

Starter Set



Efficient cleaning
thanks to
umbrella effect



Swiss made -
premium quality



CURAPROX

Figure 1: Curaprox Interdental Brush Set CPS Prime Start packaging.

2. SAFETY INFORMATION

- **Choking hazard:** Keep out of reach of children.
- Not intended for use by children who may not understand or follow safety instructions.
- Do not use for purposes other than oral hygiene.
- Replace brushes regularly to maintain effectiveness and safety.

3. PRODUCT COMPONENTS AND SETUP

The set includes:

- 5 Interdental Brushes: 1 x CPS 06, 1 x CPS 07, 1 x CPS 08, 1 x CPS 09, 1 x CPS 11
- 1 UHS 409 Brush Holder
- 1 UHS 470 Brush Holder

3.1 Attaching a Brush to the Holder

The Curaprox holders feature a convenient click system for easy brush attachment. Simply insert the metal wire end of the interdental brush into the opening of the holder until it clicks securely into place.



Figure 2: Interdental brushes and holders with click system.



Figure 3: All components of the Curaprox CPS Prime Start set.

4. OPERATING INSTRUCTIONS

Follow these steps for effective interdental cleaning:

1. Remove the interdental brush set from its packaging.
2. Identify the appropriate brush size for your interdental spaces. The mixed set allows you to find the perfect fit for each space.
3. Gently insert the brush between your teeth, starting from the gum line. Do not force the brush into tight spaces.
4. Move the brush back and forth a few times, then remove it. A single cleaning movement (once in and out) is often sufficient due to the ultra-thin and extra-strong surgical wire (CURAL).
5. Repeat this process for all interdental spaces, using the appropriate brush size for each.
6. Rinse the brushes and holder as needed during and after use.
7. Store the brushes in the holder provided when not in use.

For optimal oral hygiene, use the interdental brushes once a day.



Figure 4: Demonstrating the use of a Curaprox interdental brush.

5. MAINTENANCE

To ensure the longevity and hygiene of your Curaprox Interdental Brush Set:

- After each use, rinse the interdental brush and holder thoroughly under running water.
- Allow the brushes to air dry before storing them in the provided holder.
- Replace interdental brushes regularly, typically every 1-2 weeks, or when the bristles show signs of wear, bending, or damage. Regular replacement ensures optimal cleaning effectiveness and prevents potential injury to gums.

6. TROUBLESHOOTING

- **Brush does not fit between teeth:** If a brush feels too tight or requires force to insert, do not proceed. Select a smaller brush size from the mixed set. Forcing a brush can cause gum injury.
- **Gum bleeding:** Slight bleeding may occur initially if gums are inflamed. Continue gentle daily use. If bleeding persists for more than a few days, consult your dental professional.
- **White scars on teeth:** Inadequate removal of bacteria can lead to demineralization (white scars) on tooth surfaces, especially around braces. Consistent and thorough interdental cleaning is crucial to

prevent this.

7. SPECIFICATIONS

Product Dimensions	2.36 x 0.55 x 5.98 inches
Item Weight	0.01 ounces
Item Model Number	20-12015-500
Manufacturer	Curaden
Age Range (Description)	Adult
Special Feature	Interdental Brush
Power Source	Manual
Item Firmness Description	Soft

8. WARRANTY AND SUPPORT

This product is subject to a 30-day return policy for refund or replacement, in accordance with standard retail terms.

For further assistance, product inquiries, or support, please contact Curaprox customer service or visit their official website. Refer to the product packaging for specific contact details.