

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [Dioche](#) /

› [Dioche Paraffin Wax Bath Machine 2.3L Kit Instruction Manual](#)

Dioche 2.3L Kit

Dioche Paraffin Wax Bath Machine 2.3L Kit Instruction Manual

Model: 2.3L Kit

1. IMPORTANT SAFETY INFORMATION

Please read all instructions carefully before using the Dioche Paraffin Wax Bath Machine. Failure to follow these safety guidelines may result in electric shock, fire, or personal injury.

- **Electrical Safety:** Ensure the power outlet matches the voltage requirements of the machine. Do not operate with a damaged cord or plug. Keep the machine away from water or wet surfaces. Unplug when not in use and before cleaning.
- **Heat Safety:** The wax and internal components become very hot during operation. Always use the provided heat insulation board. Do not touch hot surfaces directly. Allow wax to cool sufficiently before use to avoid burns. Test wax temperature on a small area of skin before full application.
- **Wax Handling:** Use only paraffin wax specifically designed for wax bath machines. Do not add water or other liquids to the wax. Do not overheat the wax.
- **Supervision:** Keep out of reach of children and pets. Do not leave the machine unattended while in operation.
- **Medical Conditions:** Consult a physician before use if you have diabetes, circulatory problems, skin conditions, or open wounds.

2. PRODUCT COMPONENTS

Familiarize yourself with the parts of your Dioche Paraffin Wax Bath Machine:

- **Main Unit:** The primary body of the machine containing the wax reservoir and heating element.
- **Lid:** Covers the wax reservoir during heating and storage.
- **Heat Insulation Board:** Placed at the bottom of the wax reservoir to prevent direct contact with the heating element and ensure even melting.
- **Switch Button:** Controls power and temperature settings.
- **Indicator Light:** Illuminates to show power status or heating activity.
- **Power Cord:** Connects the machine to an electrical outlet.

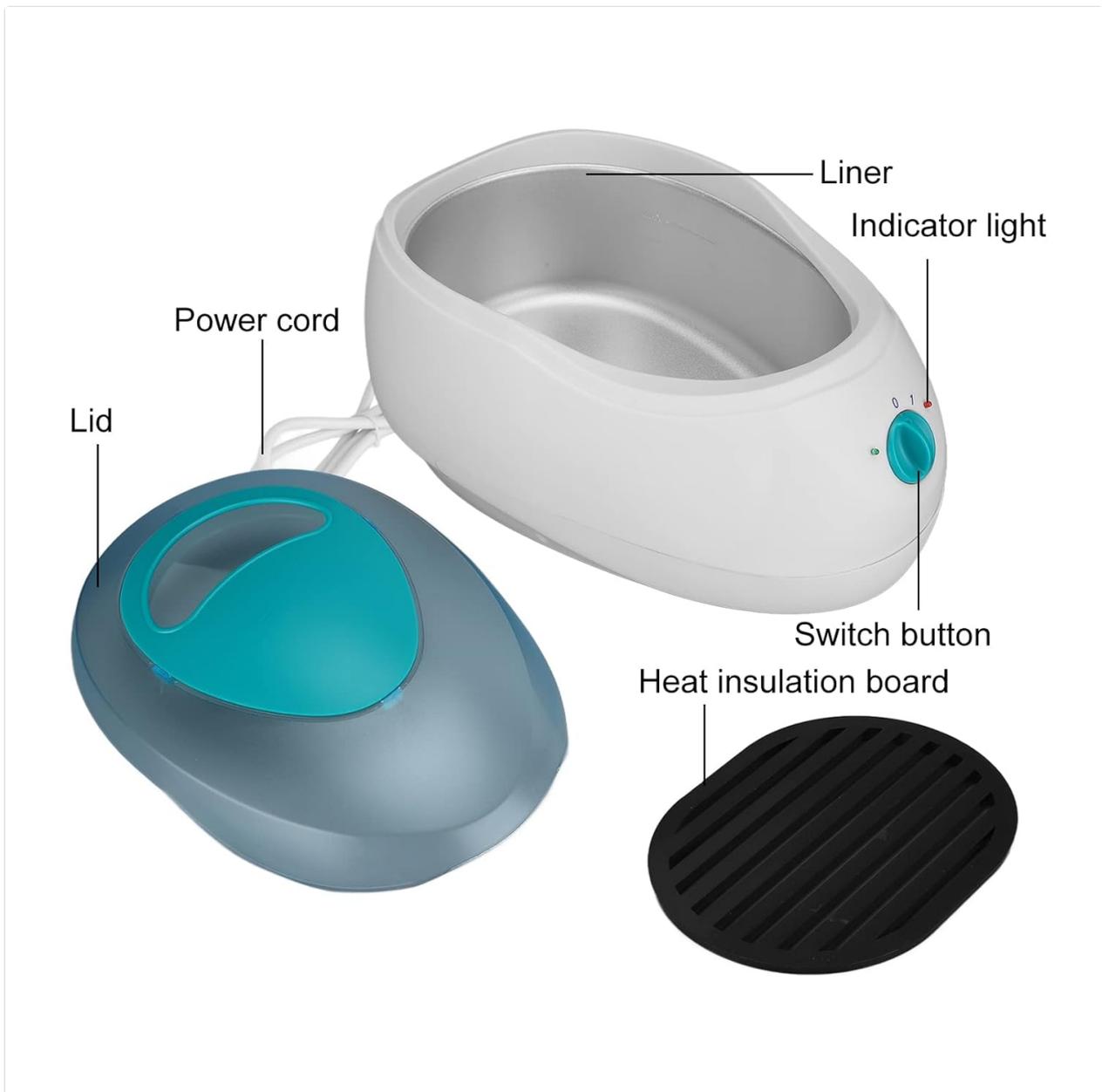


Image: Overview of the Dioche Paraffin Wax Bath Machine with its key components labeled for easy identification.

3. SETUP INSTRUCTIONS

1. **Unpack:** Carefully remove all components from the packaging.
2. **Placement:** Place the machine on a stable, heat-resistant, and level surface. Ensure it is away from water sources and flammable materials.
3. **Insert Heat Insulation Board:** Place the heat insulation board at the bottom of the wax reservoir. This prevents the wax from directly touching the heating element.
4. **Add Paraffin Wax:** Place appropriate paraffin wax blocks into the reservoir. The machine has a 2.3L capacity. Do not overfill.
5. **Connect Power:** Plug the power cord into a suitable electrical outlet.

With heat insulation board
To prevent burns



Image: The interior of the wax bath machine, highlighting the heat insulation board designed to prevent burns and ensure safe operation.

4. OPERATING INSTRUCTIONS

4.1 Melting the Wax

1. **Turn On:** Rotate the switch button to the '1' position to turn on the machine. The indicator light will illuminate.
2. **Set Temperature:** Adjust the temperature dial to the desired setting. For initial melting, a higher setting can be used.
3. **Melt Wax:** The machine features quick melt technology, dissolving paraffin blocks evenly within approximately 20 minutes. Keep the lid closed during melting to retain heat.
4. **Check Consistency:** Once melted, rotate the switch to a lower setting to maintain a comfortable temperature. The wax should be fully liquid and smooth.



Image: The paraffin wax bath machine demonstrating its quick melt capability, showing that wax blocks can melt evenly within 20 minutes.

Switch

Simple operation

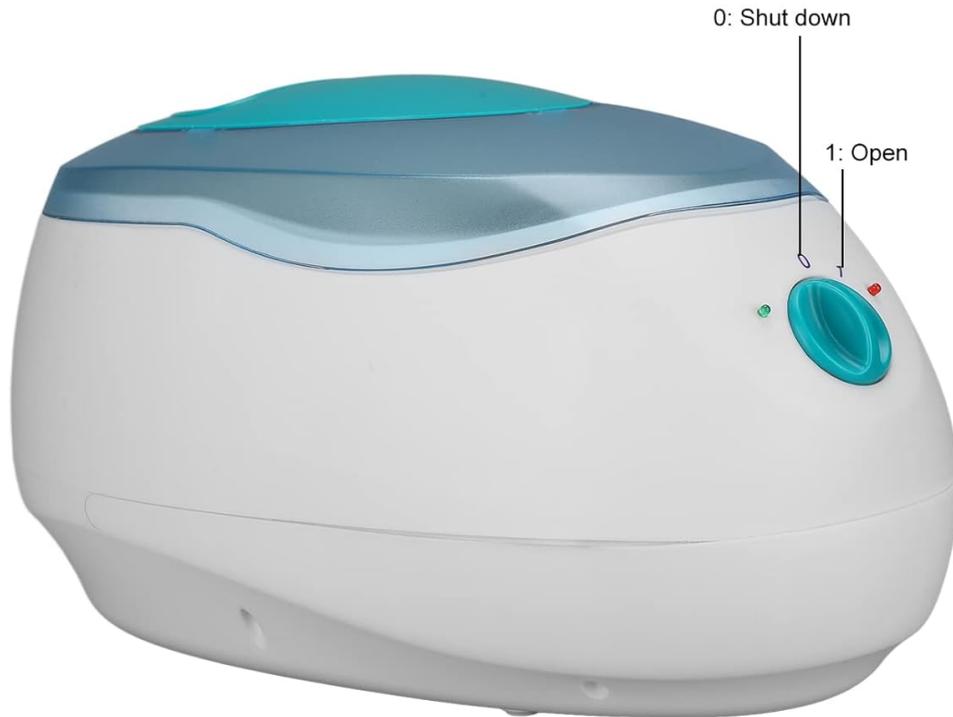


Image: Close-up of the simple operation switch on the machine, indicating positions for '0: Shut down' and '1: Open' (power on).

4.2 Applying the Wax

Before applying, ensure the wax has cooled to a safe and comfortable temperature. Test a small amount on your wrist or inner arm.

1. **Prepare Skin:** Clean and dry the area to be treated (hands, feet, elbows).
2. **Dip or Brush:** Slowly dip the hand or foot into the melted wax, or use a brush to apply layers. Allow each layer to dry slightly before re-dipping/re-applying. Repeat 3-5 times to build a thick layer.
3. **Cover:** For enhanced heat retention, cover the waxed area with plastic liners or a towel.
4. **Relax:** Leave the wax on for 10-20 minutes.
5. **Remove Wax:** Gently peel off the hardened wax. Dispose of used wax.

Your browser does not support the video tag.

Video: Demonstrates the process of melting paraffin wax, applying it to the hand, and removing it after treatment. This video shows the full application process.

Your browser does not support the video tag.

Video: Shows the initial setup of the paraffin wax machine, including placing the heat insulation board and turning on the unit.



Image: A hand coated in paraffin wax, illustrating the application process for moisturizing and soothing skin.

5. MAINTENANCE

5.1 Cleaning the Machine

1. **Unplug and Cool:** Always unplug the machine and allow the wax to cool and solidify completely before cleaning.
2. **Remove Solidified Wax:** Once the wax is solid, it can usually be lifted out in one piece. If not, gently scrape it out with a non-abrasive tool.
3. **Wipe Interior:** Use a soft, damp cloth to wipe the interior of the wax reservoir. Do not use abrasive cleaners or immerse the unit in water.
4. **Clean Exterior:** Wipe the exterior with a soft, dry cloth.

5.2 Wax Storage

If you plan to reuse the wax, ensure it is free from debris before allowing it to solidify. Store the machine with the lid closed in a cool, dry place.

6. TROUBLESHOOTING

- **Machine Not Turning On:**

Check: Ensure the power cord is securely plugged into both the machine and a working electrical outlet. Verify the switch button is in the '1' (on) position.

- **Wax Not Melting:**

Check: Confirm the machine is powered on and the temperature setting is appropriate for melting. Ensure the lid is closed to retain heat. Allow sufficient time (approx. 20 minutes for initial melt).

- **Wax Too Hot/Cold:**

Adjust: Use the temperature dial to fine-tune the wax temperature. Always test the wax temperature before application.

- **Uneven Wax Melting:**

Ensure: The heat insulation board is correctly placed at the bottom of the reservoir. Break larger wax blocks into smaller pieces for more even melting.

7. SPECIFICATIONS

- **Model:** Dioche Paraffin Wax Bath Machine 2.3L Kit
- **Capacity:** 2.3 Liters
- **Dimensions:** Approximately 14.9 x 10.4 x 8.3 inches
- **Weight:** Approximately 2.65 Pounds
- **Power:** Standard electrical outlet (voltage details typically found on product label)
- **Features:** Quick melt technology, heat insulation board, transparent protective cover.

8. WARRANTY AND SUPPORT

For warranty information or technical support, please refer to the contact details provided with your purchase or visit the official Dioche website. Keep your purchase receipt as proof of purchase.