

NuFACE TRINITY+ Starter Kit

NuFACE TRINITY+ Microcurrent Facial Device Kit User Manual

Model: TRINITY+ Starter Kit

INTRODUCTION

The NuFACE TRINITY+ Microcurrent Facial Device Kit is designed to provide at-home microcurrent treatments for facial toning and skin tightening. This device utilizes low-level electrical currents to stimulate facial muscles, aiming to improve contour and reduce the appearance of fine lines and wrinkles. This manual provides essential instructions for the safe and effective use of your TRINITY+ device and accompanying activators.



Image: The NuFACE TRINITY+ device, shown with the Hydrating Aqua Gel and Firming + Brightening Silk Creme activators.

Safety Information

The NuFACE TRINITY+ device is intended for external use only. Please read and understand all safety warnings and instructions before operating the device. Consult a healthcare professional if you have any medical conditions or concerns regarding microcurrent treatment.

- Do not use on broken or irritated skin.
- Avoid direct contact with eyes.
- Keep out of reach of children.
- Do not use if pregnant, have a pacemaker, or any electronic implanted device.
- If you experience any discomfort or adverse reactions, discontinue use immediately and consult a physician.

WHAT'S INCLUDED

Your NuFACE TRINITY+ Microcurrent Facial Device Kit typically includes the following items:

- NuFACE TRINITY+ Microcurrent Facial Device
- NuFACE Hydrating Aqua Gel Activator
- NuFACE Firming + Brightening Silk Creme Activator
- Charging Cradle and Power Adapter
- User Manual & Quick Start Guide (this document serves as the detailed manual)

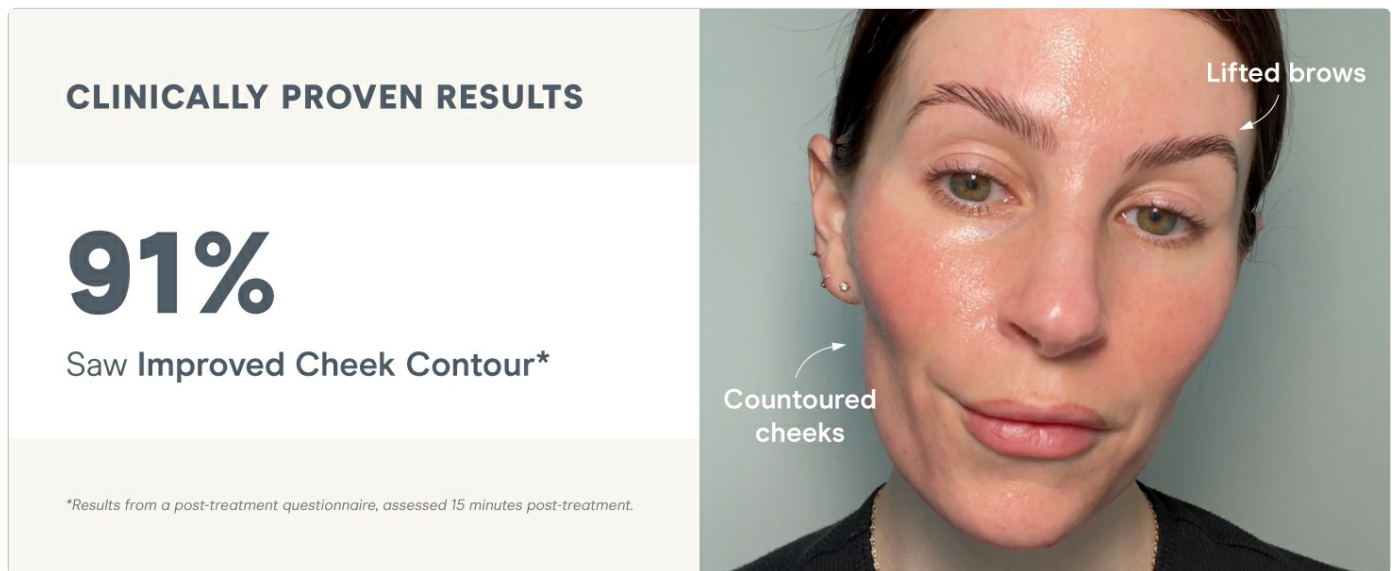


Image: An overview of the NuFACE TRINITY+ Starter Kit components.

SETUP

1. **Charge the Device:** Before first use, fully charge your NuFACE TRINITY+ device for 12 hours. Place the device in its charging cradle and connect the power adapter to an electrical outlet. The indicator lights on the device will illuminate to show charging status.
2. **Download the NuFACE Smart App:** For guided treatments and personalized routines, download the NuFACE Smart App from your device's app store. The app provides step-by-step tutorials and allows you to track your progress. You can find a QR code for the app on the product packaging or by visiting qrco.de/bbWhvO.
3. **Prepare Your Skin:** Ensure your skin is clean and free of oil-based cleansers or makeup before applying the microcurrent activator.



Image: The NuFACE Smart App, offering guided tutorials and progress tracking features.


OPERATING INSTRUCTIONS

Follow these steps for a typical microcurrent treatment session:

1. **Cleanse Skin:** Begin with clean, dry skin. Use an oil-free cleanser to remove any impurities.
2. **Apply Activator:** Apply a mask-like layer of your chosen NuFACE Microcurrent Activator (Aqua Gel or Silk Creme) to the treatment area. The activator ensures proper microcurrent conduction. Apply in sections as you treat, ensuring the skin remains hydrated throughout the session.
3. **Power On:** Press the power button on your TRINITY+ device to turn it on.
4. **Perform Glides and Holds:** Using the device, perform slow, upward glides and holds across the targeted zones of your face and neck. Refer to the NuFACE Smart App for detailed, step-by-step guided treatments for specific areas like cheeks, brows, and neck.
5. **Treatment Duration:** A typical session lasts 5 minutes. For optimal results, use the device 5 minutes a day, 5 days a week.
6. **Power Off:** Press the power button to turn off the device after your treatment.
7. **Finish:** Gently massage any remaining activator into your skin or wipe off excess. Continue with your regular skincare routine.



Image: Visual guide on applying activator and performing glides with the TRINITY+ device.

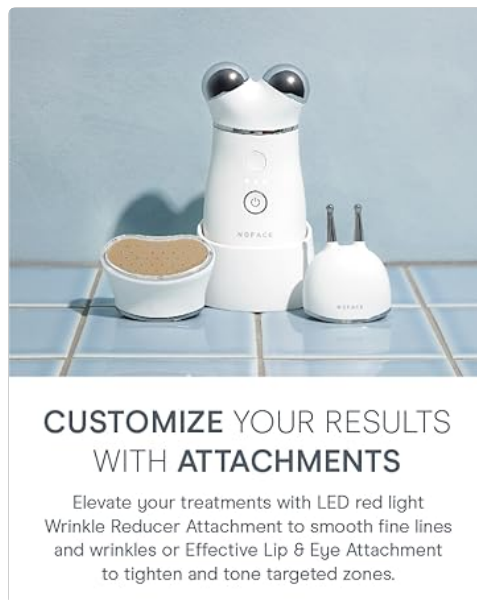

CUSTOMIZE YOUR RESULTS WITH ATTACHMENTS


Elevate your treatments with LED red light
Wrinkle Reducer Attachment to smooth fine lines
and wrinkles or Effective Lip & Eye Attachment
to tighten and tone targeted zones.

Image: A user demonstrating the application of the NuFACE TRINITY+ device on the neck area.

DEVICE MODES AND FEATURES

The NuFACE TRINITY+ device features 3-Depth Technology, offering various modes for targeted treatments:

- **Skin-Tightening Mode:** Designed to smooth lines and firm skin.
- **Instant-Lift Mode:** Aims to sculpt and contour for immediate, visible results.
- **Pro-Toning Mode:** Provides deeper toning to the muscle for long-lasting effects.

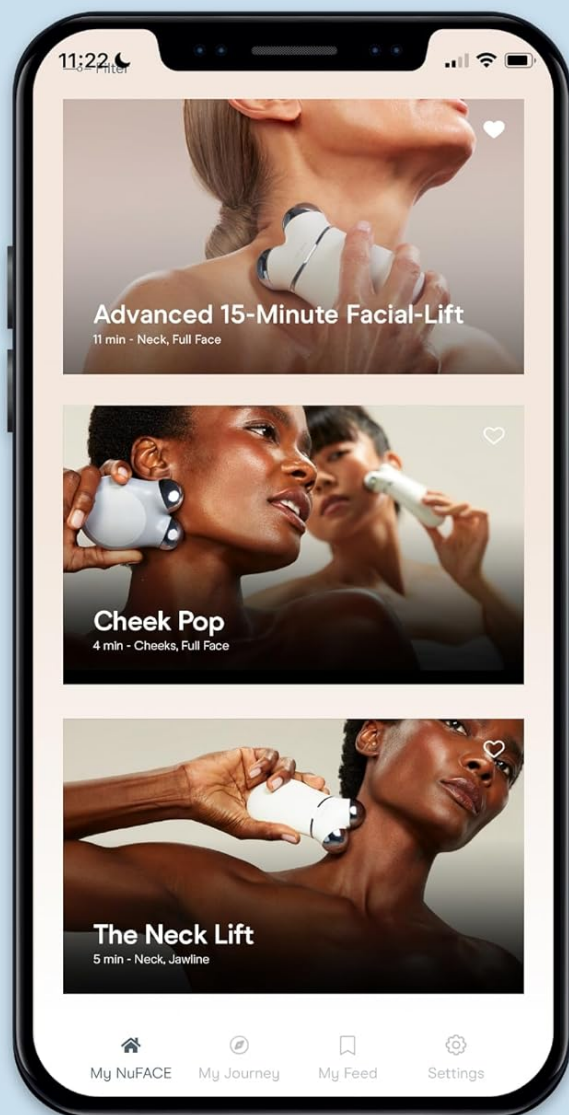

NuFACE[®] SMART APP

**FOLLOW
STEP-BY-STEP TUTORIALS**

TRACK PROGRESS

**CUSTOMIZE RESULTS WITH
APP-EXCLUSIVE TREATMENTS**

SET TREATMENT REMINDERS

Image: The NuFACE TRINITY+ device highlighting its 3-Depth Technology modes.

The TRINITY+ device is also compatible with various attachments (sold separately) to customize your treatments, such as the Wrinkle Reducer Attachment for fine lines and the Effective Lip & Eye Attachment for targeted zones.

NuFACE SMART APP

The NuFACE Smart App enhances your treatment experience by providing:

- Step-by-step guided treatment tutorials.
- App-exclusive treatments.
- Progress tracking with a Selfie Tracker.
- Treatment reminders.

Download the app for free to optimize your NuFACE routine. The app can be accessed via this link: [NuFACE Smart App](#).

MAINTENANCE

- **Cleaning:** After each use, wipe the device's spheres and body with a damp, soft cloth. Do not submerge the device in water. Avoid using alcohol or abrasive cleaners.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and extreme temperatures. Keep it in its charging cradle to ensure it's ready for your next treatment.
- **Battery Care:** The device contains a Lithium Ion battery. For optimal battery life, avoid fully discharging the device frequently. Charge as needed.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device not turning on	Low battery or not charged.	Ensure the device is fully charged. Connect to the charging cradle and power adapter.
No sensation during treatment	Insufficient activator applied or dry skin.	Apply a generous, mask-like layer of NuFACE Microcurrent Activator to the skin. Reapply if skin becomes dry during treatment.
Device lights not illuminating	Device is off or battery is completely drained.	Press the power button to turn on. If still unresponsive, place on charger for several hours.

If you encounter issues not listed here, please refer to the NuFACE Smart App for additional support or contact customer service.

SPECIFICATIONS

- **Model Number:** TRINITY+ Starter Kit
- **Power Source:** 1 Lithium Ion battery (included)
- **Product Dimensions:** 5.41 x 5.16 x 9.02 inches
- **Weight:** 2.2 Pounds
- **Manufacturer:** NuFACE
- **FDA Status:** FDA Cleared

INGREDIENTS (ACTIVATORS)

The microcurrent activators are formulated to conduct the microcurrent from your NuFACE device to your facial muscles. Key ingredients include:

- **Aqua Gel Activator:** Water, Propanediol, Sea Silt Extract, Phenoxyethanol, Carbomer, Sodium Citrate, Hydroxyethylcellulose, Glycerin, Hydrolyzed Hyaluronic Acid, Caprylyl Glycol, Ethylhexylglycerin, Sodium Hydroxide, Betaine, Adenosine, Tremella Fuciformis Sporocarp Extract, Olivine Extract, Benzyl Alcohol, Potassium Sorbate, Glyceryl Caprylate, Glyceryl Undecylenate, Monosodium Phosphate, Disodium Phosphate, Sodium Metabisulfite, Tocopherol, Sodium Nitrate.
- **Silk Creme Activator:** (Specific ingredients not provided in the input, but generally includes hydrating and firming agents.)

Note: Always check the product packaging for the most current and complete ingredient list.

WARRANTY AND SUPPORT

For information regarding product warranty, returns, or technical support, please visit the official NuFACE website or contact their customer service department. Keep your proof of purchase for warranty claims.





Manufacturer: NuFACE


Website: www.nuface.com (example link, actual link may vary)



© 2024 NuFACE. All rights reserved. This manual is for informational purposes only.

Related Documents

	<p>NuFACE 5-Minute Facial-Lift: At-Home How-To Guide</p> <p>Learn how to tone, lift, and contour your skin in just 5 minutes with the NuFACE 5-Minute Facial-Lift at-home treatment. This guide covers supplies needed, steps to follow, and microcurrent contraindications.</p>
	<p>NuFACE 5-Minute Facial-Lift: At-Home How-To Guide</p> <p>Learn how to perform the NuFACE 5-Minute Facial-Lift at home for toned, lifted, and contoured skin. This guide covers essential supplies, microcurrent contraindications, and step-by-step application instructions for the NuFACE Trinity and Mini devices.</p>
	<p>NuFACE 5-Minute Facial-Lift: At-Home How-To Guide</p> <p>Learn how to perform a 5-minute facial lift at home with NuFACE Trinity, Mini, or Mini+ devices. This guide covers supplies, steps, and microcurrent contraindications for a toned, lifted, and contoured complexion.</p>
	<p>NuFACE Trinity Lip & Eye Attachment: At-Home How-To Guide</p> <p>Learn how to use the NuFACE Trinity Device with the Effective Lip & Eye Attachment for targeted microcurrent treatments around the eyes, lips, and nasolabial folds. Includes supplies, contraindications, and step-by-step application instructions.</p>

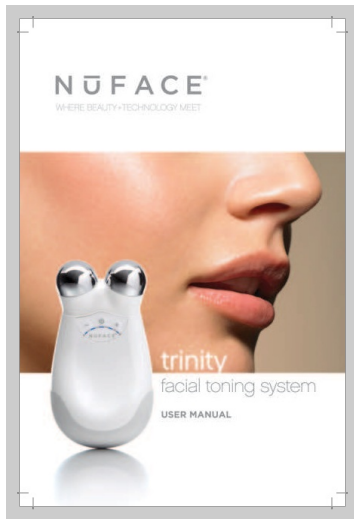

[\[pdf\]](#) User Manual Accessories

User Manual Amazon NuFACE Trinity Starter Kit Facial Toning Device with Hydrating Leave On Gel Primer Wrinkles Contour Professional 2 Fl Oz Beauty Personal Care A1smLpRJzmL m media amazon images | |||

WELCOME TO NuFACE Congratulations on receiving your NuFACE Trinity Facial Toning Device. We look forward to serving your entire skin health needs. To assure your immediate success, and prior to using your NuFACE Trinity Device, please read this entire user manual and view the online instructional v...

lang:en score:33 filesize: 1.44 M page_count: 12 document date: 2022-09-20