

Voraiya Light Therapy Lamps

Voraiya Light Therapy Lamp 10000 Lux User Manual

Model: Light Therapy Lamps

1. INTRODUCTION AND OVERVIEW

The Voraiya Light Therapy Lamp is designed to simulate natural sunlight, providing 10,000 Lux of UV-free light. This device is intended to support well-being by enhancing energy levels and creating a comfortable atmosphere, particularly in environments with limited natural light. It features adjustable brightness, multiple color temperatures, and a built-in timer for personalized use.

2. SAFETY INFORMATION

Please read all instructions carefully before using the Voraiya Light Therapy Lamp. Keep this manual for future reference.

- This lamp emits 10,000 Lux of UV-free light, making it safe for eyes and skin.
- Avoid direct prolonged staring into the light source.
- Keep the device away from water or damp environments.
- Do not attempt to disassemble or repair the device yourself. Contact customer service for assistance.
- Ensure the power adapter is correctly connected to a suitable power source.

3. PACKAGE CONTENTS

Your package includes the following items:

- 1 x Voraiya Light Therapy Lamp
- 1 x Power Adapter
- 1 x User Manual

4. PRODUCT FEATURES

- **10,000 Lux Brightness:** Emits a powerful, UV-free light to simulate natural sunlight.
- **3 Color Temperatures:** Offers three distinct temperature settings (2700K-6000K) for customizable lighting.
- **5 Brightness Levels:** Adjustable brightness from 10% to 100% (10%-25%-50%-70%-100%) to suit individual preferences and environments.
- **4 Timer Settings:** Built-in timer with options for 15, 30, 45, or 60 minutes.
- **Memory Function:** Automatically restores previous settings (color, brightness, timer) for convenience.
- **Enhanced Design:** Modern and attractive design suitable for various home and office settings.

5. SETUP

1. Unpack the Voraiya Light Therapy Lamp and all included accessories.
2. Place the lamp on a stable, flat surface where it can be positioned at arm's length from you.
3. Connect the power adapter to the lamp's power input port.
4. Plug the power adapter into a standard electrical outlet.

6. OPERATING INSTRUCTIONS

Power On/Off

Press the power button (usually indicated by a circle with a vertical line) on the lamp to turn it on or off.

Adjusting Brightness

The lamp offers 5 levels of brightness: 10%, 25%, 50%, 70%, and 100%. Use the '+' and '-' buttons to cycle through these brightness levels. For light therapy, the brightest setting is often recommended.

4 LEVEL BRIGHTNESS ADJUSTMENTS



Image: The Voraiya Light Therapy Lamp demonstrating its four adjustable brightness levels, from 25% to 100% intensity.

Changing Color Temperature

The lamp provides three color temperature settings: Daylight (5500-6000K), Natural Light (4000-4500K), and Warm Light (2700-3000K). Long-press the color temperature button (often indicated by a sun icon or 'O' key) for 2 seconds to switch between these modes. The memory function will recall your last setting.



Image: The Voraiya Light Therapy Lamp showcasing its three distinct color temperature options, allowing users to choose between different light moods.

Setting the Timer

Select your desired session duration using the timer buttons: 15, 30, 45, or 60 minutes. The lamp will automatically turn off after the selected time. The memory function will also remember your last timer setting.

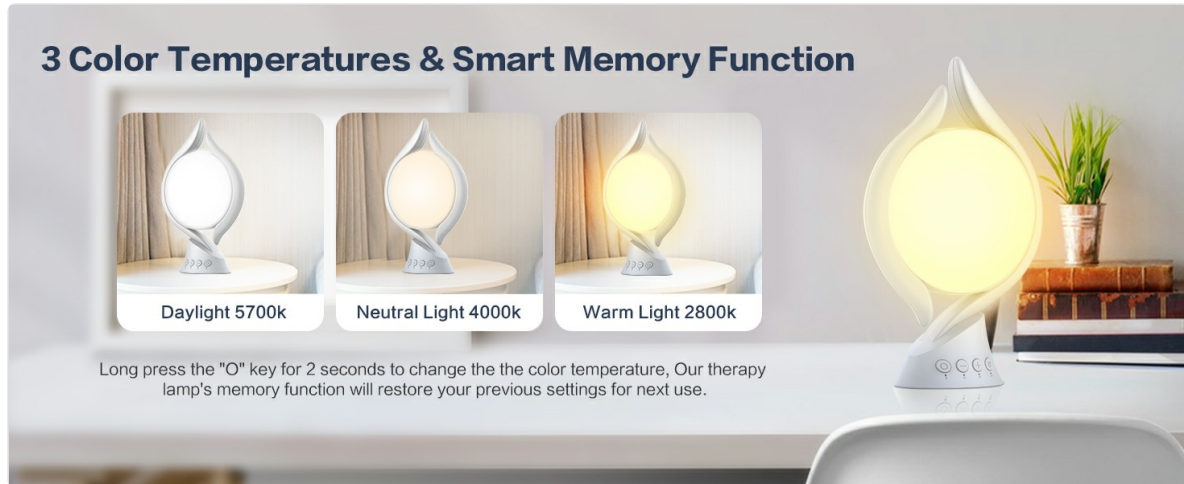


Image: The Voraiya Light Therapy Lamp illustrating its timer functionality, with options for 15, 30, 45, and 60-minute sessions.

Memory Function

The lamp is equipped with a memory function that automatically saves your last used brightness, color temperature, and timer settings. This eliminates the need to reset your preferences each time you use the lamp.

7. USAGE RECOMMENDATIONS

For optimal results, it is recommended to use the light therapy lamp for 15-30 minutes daily, preferably in the morning within the first hour of waking up. Position the lamp at arm's length from you, ensuring the light is directed towards your eyes but not directly staring into it. Adjust brightness and color temperature to your comfort level.

Your browser does not support the video tag.

Video: A demonstration of the Voraiya Light Therapy Lamp, highlighting its features and ease of use in a home setting.

Your browser does not support the video tag.

Video: An overview of the convenient Voraiya Therapy Light, showcasing its 5 brightness settings and 3 color temperatures.

8. MAINTENANCE

To ensure the longevity and optimal performance of your Voraiya Light Therapy Lamp, follow these maintenance guidelines:

- **Cleaning:** Disconnect the power adapter before cleaning. Use a soft, dry cloth to wipe the lamp's

surface. Do not use abrasive cleaners or solvents.

- **Storage:** When not in use for extended periods, store the lamp in a cool, dry place away from direct sunlight and extreme temperatures.
- **Handling:** Handle the lamp with care to avoid dropping or damaging it.

9. TROUBLESHOOTING

If you encounter any issues with your Voraiya Light Therapy Lamp, please refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Lamp does not turn on	Power adapter not connected or faulty outlet.	Ensure the power adapter is securely plugged into the lamp and a working electrical outlet. Try a different outlet.
Brightness/Color Temperature cannot be adjusted	Buttons not pressed correctly or temporary glitch.	Ensure you are using the correct buttons for adjustment (short press for brightness, long press for color temperature). Unplug and replug the lamp to reset.
Timer function not working	Timer button not activated.	Press the timer button to select a duration. The lamp should indicate the active timer setting.

10. SPECIFICATIONS

Feature	Detail
Light Source	LED
Luminous Flux	10000 Lux
Color Temperature	2700K-6000K (3 settings)
Brightness Levels	5 (10%-25%-50%-70%-100%)
Timer Settings	15/30/45/60 minutes
Power Source	Corded Electric
Wattage	12 watts
Product Dimensions	11.8 x 8.55 x 1.29 inches
Item Weight	1.63 pounds
Material	Plastics
Switch Style	Touch

11. WARRANTY AND SUPPORT

The Voraiya Light Therapy Lamp comes with a worry-free 12-month warranty. For any product inquiries, technical support, or warranty claims, please contact our responsive customer service team. Refer to the contact information provided in your purchase documentation or on the official Voraiya website.