



Manuals.plus /

› Nulift /

› Nulift NL04 Radio Frequency Skin Tightening Device User Manual

Nulift NL04

Nulift NL04 Radio Frequency Skin Tightening Device User Manual

Model: NL04 (BM12-0W)

1. INTRODUCTION

Thank you for choosing the Nulift NL04 Radio Frequency Skin Tightening Device. This device is designed to provide a multi-functional skincare experience at home, targeting various skin concerns such as fine lines, wrinkles, sagging skin, and puffiness. Please read this manual thoroughly before first use to ensure safe and effective operation.

2. IMPORTANT SAFETY INFORMATION

To ensure safe operation and prevent injury or damage to the device, please observe the following precautions:

- Do not use the device on broken, irritated, or infected skin.
- Consult a physician before use if you have any underlying medical conditions, are pregnant, or have electronic implants.
- Keep the device away from water and do not immerse it in any liquid.
- Do not use the device near flammable materials.
- Keep out of reach of children.
- Use only the provided charging cable and adapter.
- Discontinue use immediately if you experience discomfort or adverse reactions.

3. PRODUCT OVERVIEW

3.1 Package Contents

The Nulift NL04 device package typically includes:

- Nulift NL04 Radio Frequency Skin Tightening Device
- RF Base Gel (100ml)
- USB Charging Cable
- User Manual



Figure 1: Nulift NL04 device, RF base gel, and USB charging cable as typically packaged.

3.2 Device Components

The Nulift NL04 device features a multi-functional treatment head and intuitive controls:

- **Treatment Head:** Features multiple contact points for various modes (RF, TS-Lifting, Cold, Warm, Vibration).
- **Power Button:** Long press to turn on/off.
- **Mode Selection Button:** Short press to cycle through skincare modes.
- **Energy Level Button:** Adjusts intensity levels (1, 2, 3).
- **Vibration Button:** Activates or deactivates vibration mode.

4. SETUP

4.1 Charging the Device

1. Connect the USB charging cable to the device's charging port.
2. Plug the other end of the USB cable into a suitable power adapter (not included) or a USB port.
3. The device will indicate charging status (e.g., battery icon on display). Charge fully before first use.

4.2 Applying Conductive Gel

Before each treatment, it is essential to apply a conductive gel to the skin. This ensures proper contact and effective delivery of the device's functions.

1. Cleanse your face thoroughly and pat dry.
2. Apply a sufficient amount of the Nulift RF Base Gel (or another suitable conductive gel) to the treatment area. Ensure an even layer for optimal conductivity.

5. OPERATING INSTRUCTIONS

5.1 Basic Operation

1. **Power On:** Long press the power button (first button) to turn on the Nulift NL04 device.
2. **Select Mode:** Short press the mode selection button (second button) to choose your desired skincare mode.
3. **Set Energy Level:** Use the energy level button (third button) to set the desired intensity level (1, 2, or 3).
4. **Activate Vibration (Optional):** The fourth button controls the vibration mode.
5. **Apply to Skin:** Place the treatment head on the intended area, ensuring full contact with the skin. Move the treatment head in slow circular or back-and-forth motions.


Your browser does not support the video tag.

Video 1: Demonstration of turning on the Nulift NL04 device and selecting modes.

5.2 Understanding Modes


The Nulift NL04 device combines several advanced technologies for comprehensive skincare:

HOW IT WORKS

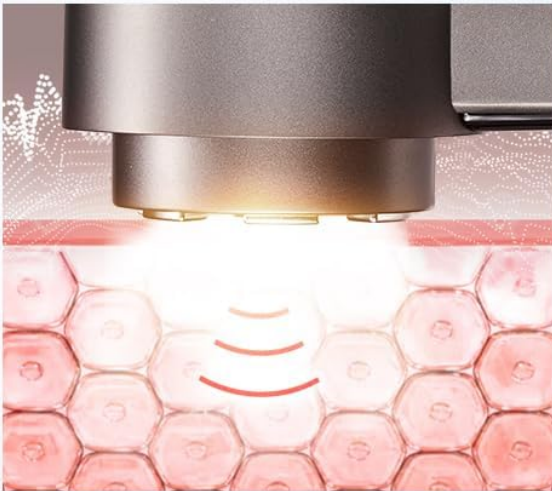

IAP mode

Radio frequency energy increases collagen production, plump skin, along with Microcurrent, lift firm and tone skin, improves skin laxity, wrinkles, facial contour...

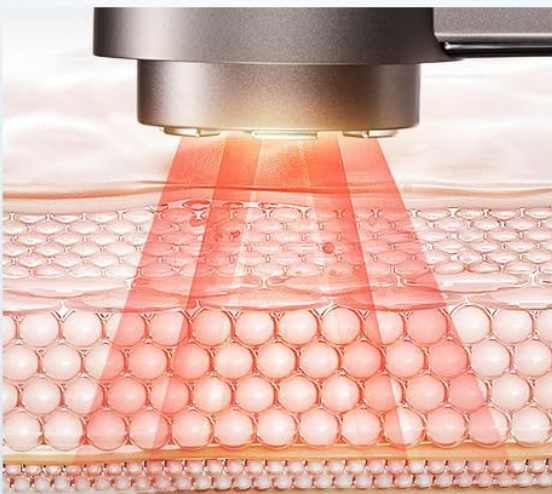

TS-Lifting mode

Microcurrent stimulates the muscles, lift & tone skin, smooth skin, reduce fine lines and wrinkles, prevent sagging...

Figure 2: Illustration of IAP (Radio Frequency) and TS-Lifting (Microcurrent) modes.

- **IAP Mode (Radio Frequency):** Increases collagen production, plumps skin, improves skin laxity, reduces wrinkles, and refines facial contour.
- **TS-Lifting Mode (Microcurrent):** Stimulates facial muscles, lifts and tones skin, smooths fine lines and wrinkles, and helps prevent sagging.

Cold Massage

Cool and calm skin, shrinks pores, locks in moisture, minimizes puffiness, penetrates the skin, refines the skin



Warm Massage

Open pores & increased skin absorption, wake up skin vitality water molecules



Vibration mode

Relax tension points, smooth lines & wrinkles, stimulate microcirculation and enhance product absorption

Figure 3: Visual representation of Cold, Warm, and Vibration modes.

- **Cold Massage Mode:** Cools and calms skin, shrinks pores, locks in moisture, minimizes puffiness, and refines skin texture.
- **Warm Massage Mode:** Opens pores, increases skin absorption of skincare products, and awakens skin vitality.
- **Vibration Mode:** Relaxes tension points, smooths lines and wrinkles, stimulates microcirculation, and enhances product absorption.

5.3 Usage Frequency and Recommendations

For optimal results, use the Nulift NL04 device regularly 3-5 times a week. You can combine modes for a comprehensive treatment.



Figure 4: Example of daily skincare routine using Ts-Lifting and Cold modes to depuff and tone skin. For enhanced skincare product absorption, use the Warm mode. For shrinking pores and retaining moisture, use the Cooling mode. The Vibration mode is beneficial for relaxing tension points and stimulating microcirculation.

ADVANTAGES OF OVER OTHER BRANDS

- ✓ The most high quality rechargeable RF device on the market
- ✓ Deep hydrating and massage skin, as do salon's facial
- ✓ Combines the most advanced tech: Radio frequency, TS-Lifting, heating, cooling and vibration
- ✓ De-puff, shrink pores, fine line and wrinkle removal
- ✓ Prevent sagging and anti wrinkle
- ✓ 1.57 inches medical grade Large treatment head
- ✓ Visible results in 8 - 12 weeks



Figure 5: Key benefits of using the Nulift NL04 device, including deep hydration, wrinkle reduction, and anti-sagging effects.

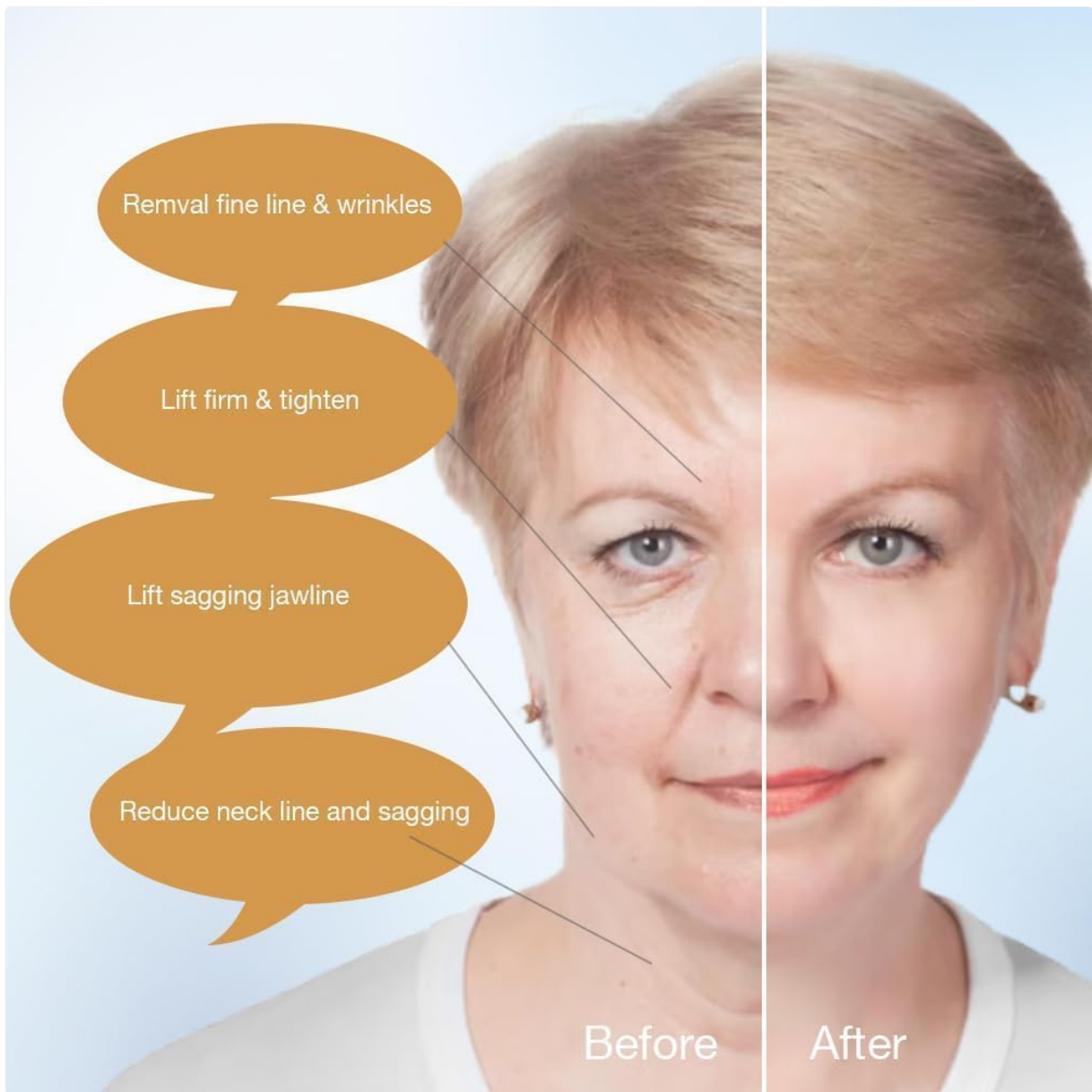


Figure 6: Illustrative results of using the device, showing improvements in fine lines, skin firmness, and jawline definition.



Figure 7: The Nulift NL04 device is suitable for use in various personal care routines.



Figure 8: The portable design allows for comfortable use in various settings.

6. MAINTENANCE AND CARE

6.1 Cleaning the Device

After each use, clean the treatment head to maintain hygiene and device longevity:

- Wipe the treatment head with a soft, damp cloth or a cotton pad.
- Do not use abrasive cleaners, alcohol, or harsh chemicals.
- Ensure the device is completely dry before storing.

6.2 Storage

Store the device in a cool, dry place, away from direct sunlight and extreme temperatures. Keep it in its original packaging or a protective case when not in use.

7. TROUBLESHOOTING

If you encounter any issues with your Nulift NL04 device, please refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Device does not turn on	Low battery; Device not charged	Charge the device fully. Ensure the charging cable is properly connected.
No sensation or weak effect	Insufficient conductive gel; Low energy level; Poor skin contact	Apply more conductive gel. Increase the energy level. Ensure the treatment head has full contact with the skin.
Device suddenly stops during use	Low battery; Automatic shutdown feature	Recharge the device. The device may have an automatic shutdown after a certain period of use; restart if necessary.

8. SPECIFICATIONS

Feature	Detail
Brand	Nulift
Model Name	NL04
Item Model Number	BM12-0W
Material	Plastic
Item Weight	0.64 Kilograms (1.41 Pounds)
Color	Black
Package Dimensions	7.72 x 6.18 x 2.32 inches
Manufacturer	MLAY
Date First Available	June 6, 2022

9. WARRANTY AND CUSTOMER SUPPORT

9.1 Warranty Information

The Nulift NL04 device comes with a standard manufacturer's warranty. Please refer to the warranty card included in your package for specific terms and conditions, including warranty duration and coverage details. Keep your proof of purchase for warranty claims.

9.2 Customer Support

For any questions, technical assistance, or warranty claims, please contact Nulift customer support through the official website or the contact information provided on your warranty card. Please have your model number (NL04) and purchase details ready when contacting support.