

## Breo Breo Foot Massage & iPalm520e

# Breo Foot Massager Machine with Heat & iPalm520e Instruction Manual

Brand: Breo | Model: Breo Foot Massage & iPalm520e

### PRODUCT OVERVIEW

The Breo Foot Massager Machine with Heat and iPalm520e offers a great combination of hand and foot massage with pleasant heat. Designed to promote relaxation, it helps users achieve better sleep and recover from busy days or intense sports activities through powerful massage sessions.

### SAFETY INFORMATION

- Do not use while sleeping.
- Keep the cord away from heat sources and liquids to prevent electrical hazards.
- Always follow the manufacturer's instructions for safe use.

### PRODUCT FEATURES

#### Shiatsu Foot Massager with Heating

## Shiatsu Foot Massager With Heating



Rolling Relief



15 Minutes Timer



Air Compression



Soothing Heat



Deep Kneading




Adjustable Intensity



This image displays the key features of the Breo Shiatsu Foot Massager. It highlights Rolling Relief, a 15-minute timer, Air Compression, Soothing Heat, Deep Kneading, and Adjustable Intensity settings for a personalized massage experience.

### Full-Coverage Foot Relaxation

# Full-Coverage Relax

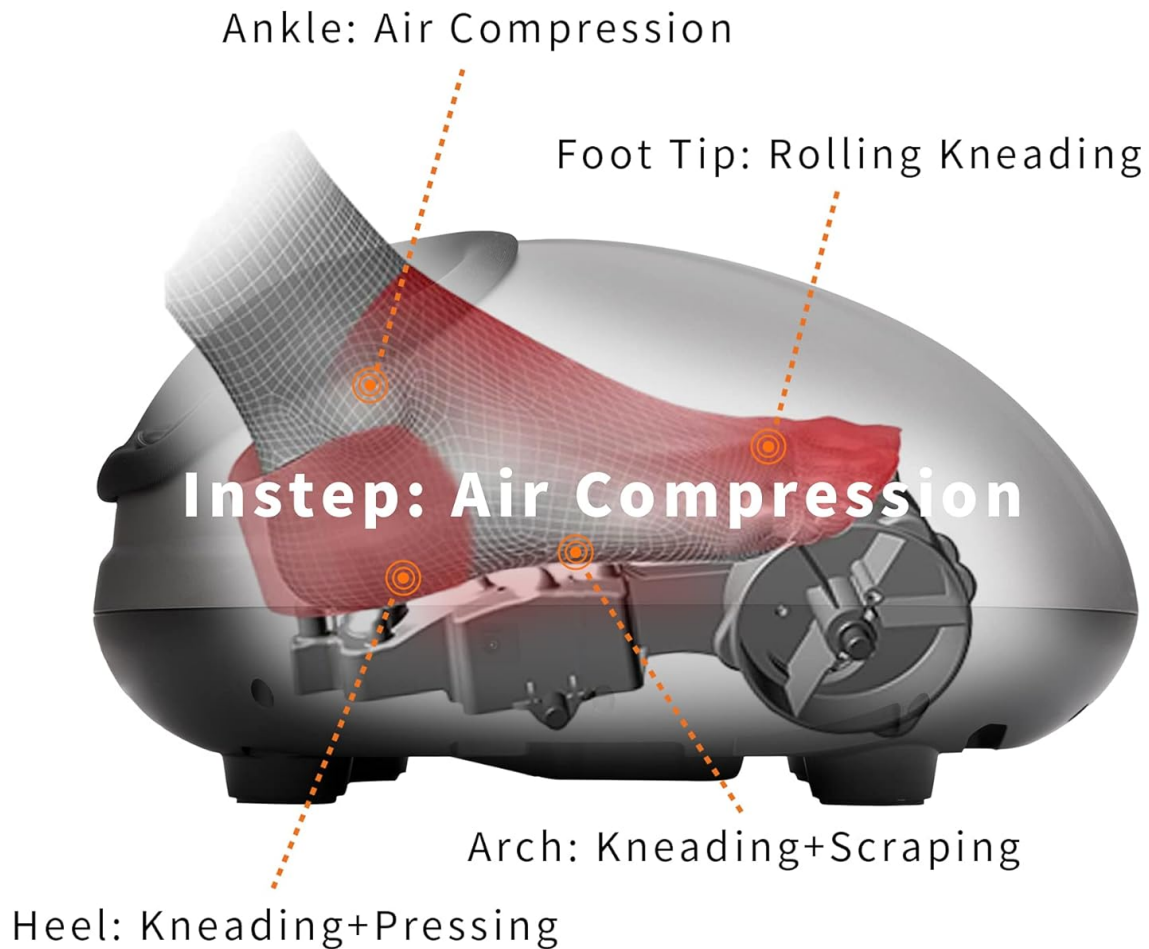

A detailed diagram illustrates how the massager provides full-coverage relaxation. It shows air compression targeting the ankle and instep, rolling kneading for the foot tip, and a combination of kneading and scraping for the arch and heel areas.

## Warm Soothing Heat Function

# Warm Soothing Heat


Improve Blood circulation & Revive Tired Muscles



This image shows the foot massager with its internal heating elements glowing red, indicating the warm soothing heat function. This feature is designed to improve blood circulation and revive tired muscles.

## Detachable & Washable Foot Sleeves



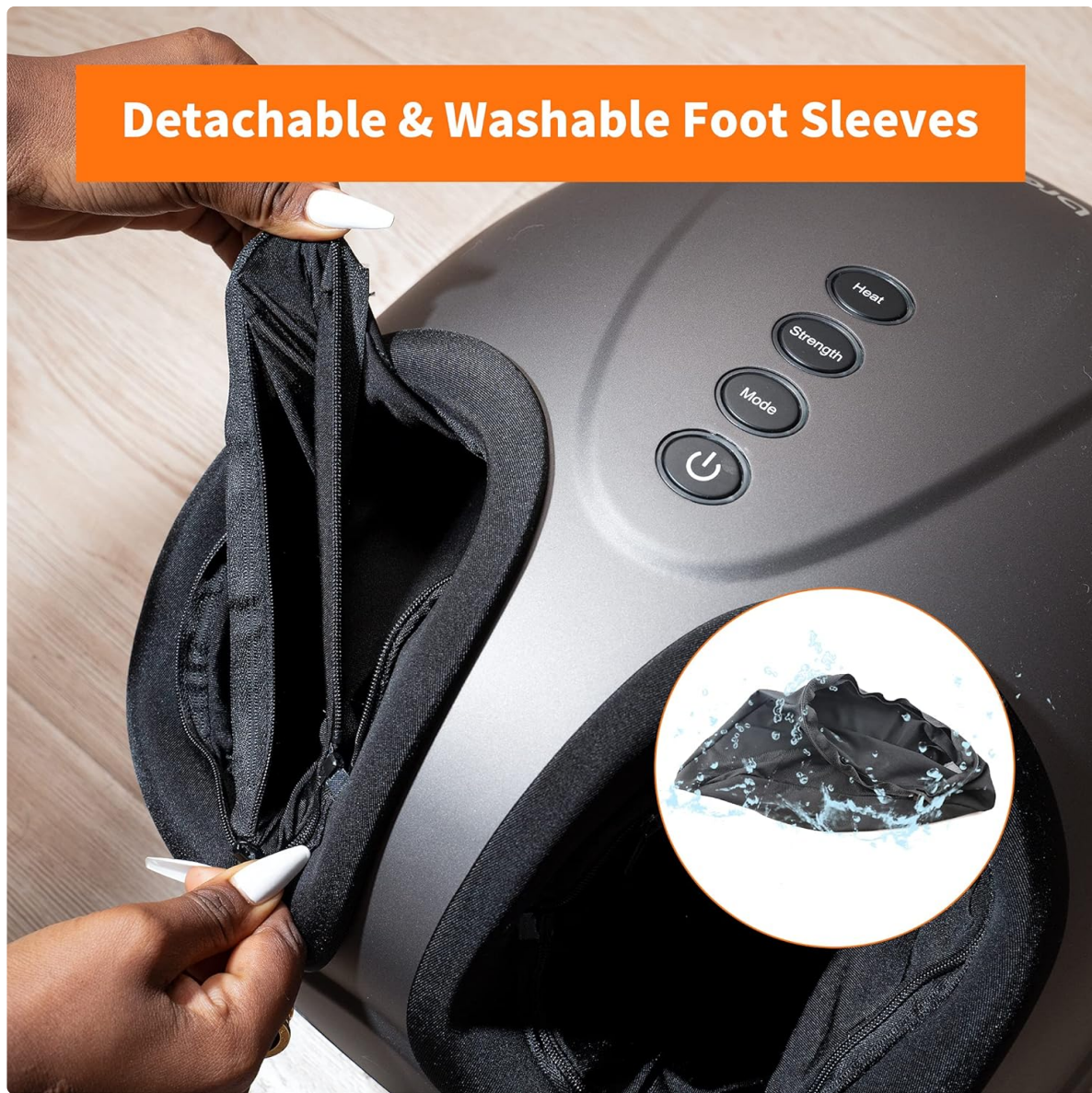

The image demonstrates the detachable and washable foot sleeves, highlighting the ease of maintenance and hygiene for the foot massager. A separate graphic shows a sleeve being washed.

### **iPalm520e Hand Massager with Air Pressure and Heat**

## Air Pressure and Heat Compression Function



### Air Pressure Massage

Squeeze the acupuncture points to promote blood circulation



### Massage Nail Design

Point to point acupuncture massage on your hand and finger



### Heat Compression

Comfortable temperature slowly makes the palm of your hand warm

This image showcases the iPalm520e hand massager, detailing its air pressure massage to squeeze acupuncture points, massage nail design for point-to-point massage on hand and finger, and heat compression for comfortable warmth.

## Portable and Wireless iPalm520e


# Portable And Wireless



The iPalm520e is shown as a portable and wireless device, featuring a Type-C charging port for convenient power. This design allows for use anywhere without being tethered to an outlet.

## Versatile Use of iPalm520e



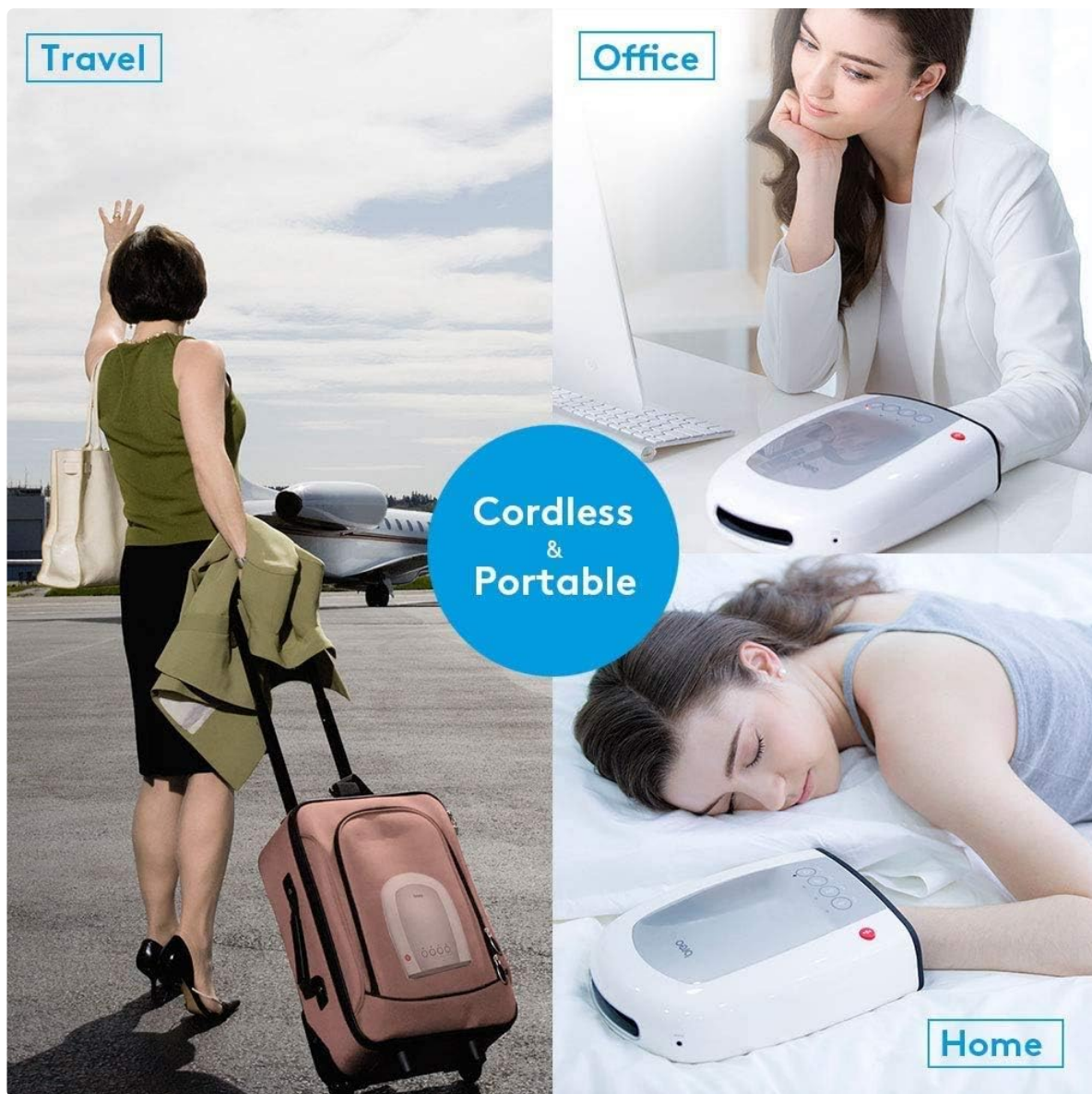

This image illustrates the iPalm520e's versatility, showing it being used in different environments such as during travel, in an office setting, and at home, emphasizing its cordless and portable nature.

## SETUP

1. Place the Breo Foot Massager on a flat, stable surface.
2. Ensure the power cord is securely plugged into the massager and a suitable power outlet.
3. For the iPalm520e, ensure it is charged via its Type-C port if using wirelessly.

## OPERATING INSTRUCTIONS

### Breo Foot Massager Operation

1. **Power On/Off:** Press the **Power** button to turn the massager on or off. The device will typically start with default settings.
2. **Heat Function:** Press the **Heat** button to activate or deactivate the soothing heat. The indicator light will turn red when active.
3. **Pressure Intensity:** Press the **Pressure** button to cycle through different air compression intensity levels (e.g., low, medium, high). The indicator light color may change to reflect the selected intensity.



- 4. **Kneading Intensity/Mode:** Press the **Kneading** button to adjust the intensity of the kneading massage or to cycle through different massage modes. The indicator light will change accordingly.
- 5. **Timer:** The massager typically has a built-in timer (e.g., 15 minutes) for automatic shut-off. Refer to the device's display for current time remaining.

Demonstration Video: Breo Foot Massager

Your browser does not support the video tag.

This video provides a visual guide on how to operate the Breo Foot Massager, demonstrating the use of the power, heat, pressure, and kneading buttons to customize your massage experience.

iPalm520e Hand Massager Operation

For the iPalm520e, place your hand into the device. Use the control panel or remote (if applicable) to select desired air pressure intensity, massage modes, and heat settings. The device is designed for intuitive use with clear button indicators.

MAINTENANCE

- **Cleaning Foot Sleeves:** The foot sleeves of the Breo Foot Massager are detachable via zippers and can be removed for washing. This ensures hygiene, especially with multiple users.
- **Exterior Cleaning:** Wipe the exterior of both massagers with a soft, damp cloth. Do not use abrasive cleaners or immerse the devices in water.
- **Storage:** Store the massagers in a cool, dry place away from direct sunlight when not in use.

SPECIFICATIONS

Feature	Detail
Brand	Breo
Model Number	Breo Foot Massage & iPalm520e
Use For	Feet (Foot Massager), Hands (iPalm520e)
Power Source	Corded Electric (Foot Massager), Wireless/Rechargeable (iPalm520e)
Special Feature	Heat, Kneading, Air Compression
Color	Black