

## MASDERM Conductive Gel Cream for RF & Microcurrent Devices

# MASDERM Face Conductive Gel Cream User Manual

Model: Conductive Gel Cream for RF & Microcurrent Devices

### PRODUCT OVERVIEW

The MASDERM Face Conductive Gel Cream is specifically formulated for use with radio frequency (RF), microcurrent, and cavitation facial machines and devices. This cream acts as a conductive medium, enhancing the effectiveness of these treatments while providing beneficial skincare properties. It is designed to leave the skin feeling softer, smoother, hydrated, and illuminated.

This product is recommended by aesthetics professionals and is suitable for both professional cabin use and at-home application.

### KEY INGREDIENTS

The MASDERM Conductive Gel Cream features a blend of high-quality natural ingredients to maximize skin hydration and effectiveness:

- **Hyaluronic Acid:** A powerful humectant that attracts and retains moisture, providing intense hydration to the skin.
- **Collagen:** A protein molecule that forms fibers essential for proper skin function, contributing to skin elasticity and firmness.
- **Bio-Retinol:** A vegan alternative to traditional Retinol, which helps to firm the skin.

## INGREDIENTES PRINCIPALES

### BIO-RETINOL:

Opción vegana frente al Retinol que ayuda a reafirmar la piel



### COLÁGENO:

Molécula proteica que forma fibras para el correcto funcionamiento de la piel



### ÁCIDO HIALURÓNICO:

Polisacárido que retiene agua para mantener la piel hidratada



*Image: Visual representation of the main ingredients: Bio-Retinol, Collagen, and Hyaluronic Acid, highlighting their molecular structures and benefits.*

## APPLICATION AND USAGE

This cream can be applied with manual massage or in conjunction with a radiofrequency device. For optimal results, ensure your skin is dry and clean before application.

### Manual Application:

Apply the cream to the desired area of the face and gently massage until absorbed. This can be done daily, morning and night.

### Application with Radiofrequency Device:

1. Ensure skin is clean and dry.
2. Apply a sufficient amount of MASDERM Conductive Gel Cream to the area to be treated.
3. Use your radiofrequency device over the cream for a maximum of 10 minutes per treated area.
4. It is recommended to use the device 3 days a week, alternating with manual massage on other days if desired.

## CÓMO APLICAR



**Aplica la crema cada día mañana y noche.  
Potencia sus resultados con el dispositivo\*  
de radiofrecuencia facial**

**(\*) 10min cada noche | Función RF - 3 días/semana**

*Image: Demonstrates the application of the cream and its use with a facial radiofrequency device.*

### **EXPECTED RESULTS AND SKIN SENSATION**

After applying the cream and performing the radiofrequency treatment, your skin will feel noticeably softer, smoother, and more hydrated. The cream also contributes to a more illuminated complexion. Users often describe the texture and scent as very pleasant.

MEJORAS VISIBLES

EN MENOS DE UN MES



**Flacidez**  
**Arrugas**  
**Sequedad**  
**Pómulos Caídos**  
**Papada**

• AVALADOS Y RECOMENDADOS POR DERMATÓLOGOS •

*Image: Before and after comparison illustrating visible improvements in skin texture and appearance within less than one month of use.*

## RESULTADOS EN MENOS DE 1 MES\*



\*Estudio realizado a 20 voluntarias durante 28 días (entre 40-60 años)

*Image: Clinical study results showing improvement in the perioral area (around the lips) after 28 days of consistent use.*

### QUALITY AND PROFESSIONAL ENDORSEMENT

MASDERM Conductive Gel Cream is widely used and trusted by aesthetics professionals. It is compatible with various devices including diathermy, radiofrequency, resistive, and capacitive currents.

- Tested under dermatological control.
- Not tested on animals.
- Paraben-free.
- Manufactured according to ISO 22716:2013 standards.


CON LA CONFIANZA DE LOS DERMATÓLOGOS



Dra. Isabel Pérez

Dermatóloga

“Si quieres mejorar la firmeza del rostro, una excelente opción es la crema de Masderm, yo me la aplico a diario y la recomiendo a mis pacientes”

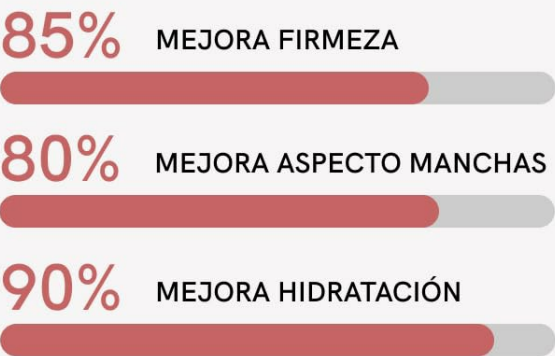

\*Estudio realizado a 20 voluntarias durante 28 días (entre 40-60 años)

Image: Dr. Isabel Pérez, a dermatologist, recommends the cream, highlighting improvements in firmness (85%), reduction of spots (80%), and hydration (90%) based on a study.

PRODUCT SPECIFICATIONS

Attribute	Detail
Brand	MASDERM
Item Volume	1000 Milliliters (16.9 fl. oz)
Product Dimensions (L x W x H)	3.94 x 3.94 x 7.09 inches
Item Weight	1.17 Pounds
Manufacturer	KERVLAN LAB S.L
Age Range (Description)	Adult
Active Ingredients	Hyaluronic Acid





*Image: Front view of the MASDERM RF Facial Cream bottle, showing the product name and volume.*



Image: Back view of the MASDERM RF Facial Cream bottle, displaying usage instructions in multiple languages and product identification codes.





Image: Close-up of the ingredients list and product information on the MASDERM RF Facial Cream bottle.

## STORAGE AND MAINTENANCE

To maintain the quality and effectiveness of your MASDERM Conductive Gel Cream, store it in a cool, dry place away from direct sunlight. Ensure the pump dispenser is securely closed after each use to prevent contamination and drying out.

Do not ingest. Keep out of reach of children.

## TROUBLESHOOTING

If you experience any issues with the product, please refer to the following common concerns:

- **Cream absorbs too quickly during device use:** While the cream is designed to be moisturizing, if it absorbs too quickly for your device, you may need to apply a slightly thicker layer or reapply as needed during the treatment session to maintain conductivity.
- **Skin irritation:** Discontinue use if irritation occurs. Consult a dermatologist if symptoms persist. Ensure your skin is clean before application.
- **Pump not dispensing:** Ensure the pump is unlocked and fully depressed. If the product has been stored for a long time, a small amount of air might be trapped; try priming the pump several times.


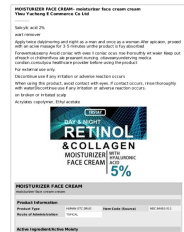


WARRANTY AND SUPPORT



For any questions, concerns, or support regarding your MASDERM Conductive Gel Cream, please contact the manufacturer, KERVLAN LAB S.L. Refer to the product packaging for specific contact details or visit the official MASDERM website.

Please retain your purchase receipt for any warranty claims, though specific warranty details for this consumable product are not provided in this manual.

© 2024 MASDERM. All rights reserved.

Related Documents - Conductive Gel Cream for RF & Microcurrent Devices

	<p><a href="#">JUNGLÜCK 2025 Advent Calendar: Your Guide to Skincare Discoveries</a></p> <p>Explore the JUNGLÜCK 2025 Advent Calendar with this comprehensive guide. Discover skincare products, learn about your skin type, and find the perfect routine for radiant, healthy skin.</p>
	<p><a href="#">TRSTAY Retinol &amp; Collagen Moisturizer Face Cream - Salicylic Acid &amp; Niacinamide</a></p> <p>Discover TRSTAY's Day &amp; Night Moisturizer Face Cream, formulated with Retinol, Collagen, Hyaluronic Acid, Salicylic Acid (2%), and Niacinamide. This topical treatment helps with wart removal and skin hydration. Product details, usage, and ingredients from Yiwu Yucheng E Commerce Co Ltd.</p>
	<p><a href="#">EvenSkyn Phoenix User Manual: Guide to Microcurrent Face Lifting</a></p> <p>Comprehensive user manual for the EvenSkyn Phoenix, a solar-powered microcurrent beauty device for lifting and toning the face, neck, and décolletage. Learn about its features, usage, and maintenance.</p>
	<p><a href="#">High Performance AquaTight™ Facial Treatment Protocol</a></p> <p>Detailed protocol for the Pevonia High Performance AquaTight™ Facial treatment, including product requirements, preparation steps, and a comprehensive application guide.</p>

	<p><a href="#">Prosana Lactic Enzyme Gel Peel: Professional Skincare Treatment Guide</a></p> <p>Comprehensive guide to the Prosana Lactic Enzyme Gel Peel, detailing its benefits for hyperpigmentation and acne scarring, application steps, skin type suitability, contraindications, and cost-profit analysis for wellness professionals.</p>
	<p><a href="#">Pevonia Aroma Aquatight™ Essentials Facial Treatment Protocol</a></p> <p>Detailed instructions and product list for the Pevonia Aroma Aquatight™ Essentials Facial treatment, including preparation and application steps for professional use.</p>