

ANLAN LM-2201

ANLAN Anti-Wrinkle Eye Device User Manual

Model: LM-2201 | Brand: ANLAN

1. INTRODUCTION

Thank you for choosing the ANLAN Anti-Wrinkle Eye Device. This device is designed to provide comprehensive eye care through its four distinct modes: EMS, Warm Compress, High-Frequency Vibration, and Light Therapy. Regular use can help address concerns such as dark circles, eye bags, and fine lines around the eyes. Please read this manual thoroughly before first use to ensure proper operation and maintenance.



Image: The ANLAN Anti-Wrinkle Eye Device, highlighting its primary function.

2. SAFETY INFORMATION

To ensure safe operation and prevent damage to the device or injury, please observe the following precautions:

- Do not use the device on broken, irritated, or infected skin.
- Avoid direct contact with eyeballs.
- Keep out of reach of children.
- Do not immerse the device in water or any other liquid.
- If you experience any discomfort, discontinue use immediately.
- Consult a medical professional before use if you have any pre-existing medical conditions, especially related to eyes or skin.
- Use only the provided charging cable.

3. PACKAGE CONTENTS

Verify that all items are present in the package:

- ANLAN Anti-Wrinkle Eye Device (Model: LM-2201)
- USB Charging Cable
- User Manual



Image: The ANLAN Anti-Wrinkle Eye Device shown with its retail packaging.

4. PRODUCT OVERVIEW

Device Components

The ANLAN Anti-Wrinkle Eye Device consists of the following main parts:

- **Massage Head:** The metallic tip that comes into contact with the skin, delivering various treatments.
- **LCD Display:** Shows current temperature, mode, and battery status.
- **Control Buttons:** For power, mode selection, and intensity adjustment.
- **Charging Port:** Located at the bottom of the device.

Control Buttons

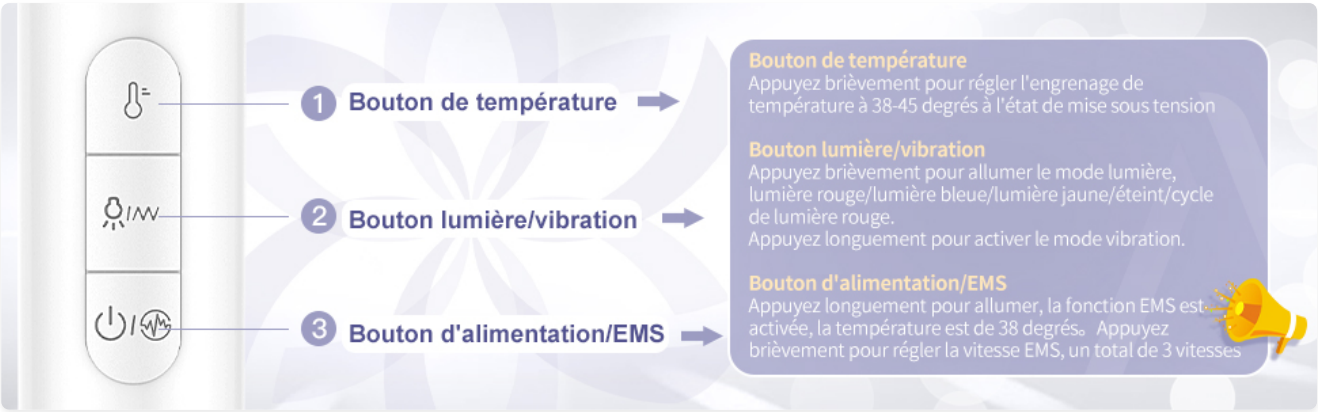


Image: Detailed diagram of the device's control buttons and their functions.

Button	Function
1. Temperature Button	Short press to adjust temperature (38-45°C) when powered on.
2. Light/Vibration Button	Short press to cycle through light modes (Red/Blue/Yellow/Off). Long press to activate vibration mode.
3. Power/EMS Button	Long press to power on/off. EMS function activates automatically at 38°C upon power on. Short press to adjust EMS intensity (3 levels).

5. SETUP

Charging the Device

Before first use, fully charge the device. The battery indicator on the LCD display will show charging status.

1. Connect the USB charging cable to the charging port at the bottom of the device.
2. Plug the other end of the USB cable into a compatible USB power adapter (not included) or a computer USB port.
3. The device will indicate charging on its LCD screen. Charging typically takes 2-3 hours.
4. Once fully charged, disconnect the device from the power source.

6. OPERATING INSTRUCTIONS

General Usage Steps



Image: Visual guide demonstrating the four steps for using the device.

1. Step 1: Preparation

After cleansing and removing makeup from your face, apply your preferred eye cream or serum to the area around your eyes. Select the desired mode on the device.

2. Step 2: Inner to Outer Eye

Gently move the massage head from the inner corner of the eye towards the outer corner in a semi-circular motion.

3. Step 3: Circular Motion

Move the device from the lower corner of the eye upwards in small circular motions.

4. Step 4: Multi-Area Use

This device can also be used on other parts of your face, including the lips, to address fine lines and improve skin condition.

Understanding the Modes

The device offers four primary modes that can be combined for a customized treatment:

- **EMS Mode:**

Stimulates muscle contraction, increases cellular activity, restores muscle elasticity, and helps tighten the eye contour. This mode can reduce the appearance of dark circles and eye bags.

EMS stimulates eye muscle movement, lifts the skin around the eyes, restores muscle elasticity, and reduces eye wrinkles

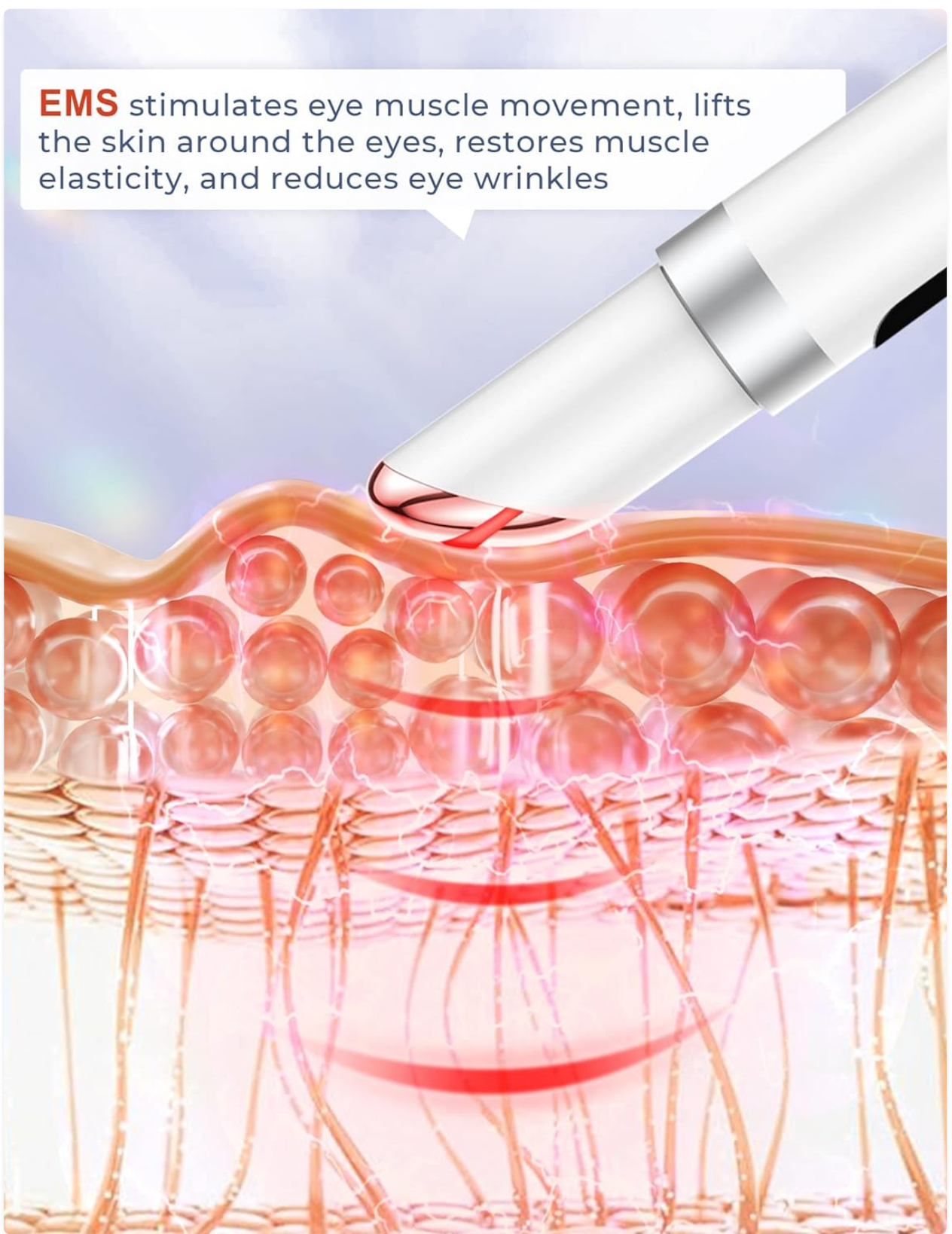


Image: Illustration of EMS stimulating eye muscle movement and reducing wrinkles.

- **Warm Compress Mode (38-45°C):**

Opens pores, accelerates blood circulation, and promotes the absorption of skincare nutrients, ensuring no product is wasted.

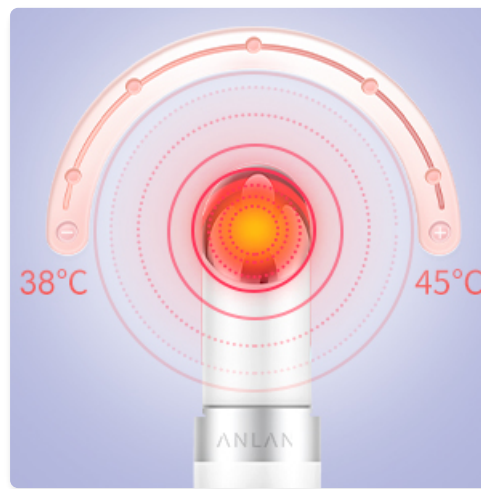


Image: Visual representation of the warm compress function and its temperature range.


- **High-Frequency Vibration ($11,000 \pm 2,000$ rpm):**

Relaxes muscles and alleviates fatigue around the eyes. This vibration helps in deeper penetration of skincare products.

HIGH FREQUENCY VIBRATION:

11000±2000rpm

Relax muscles and relieve fatigue

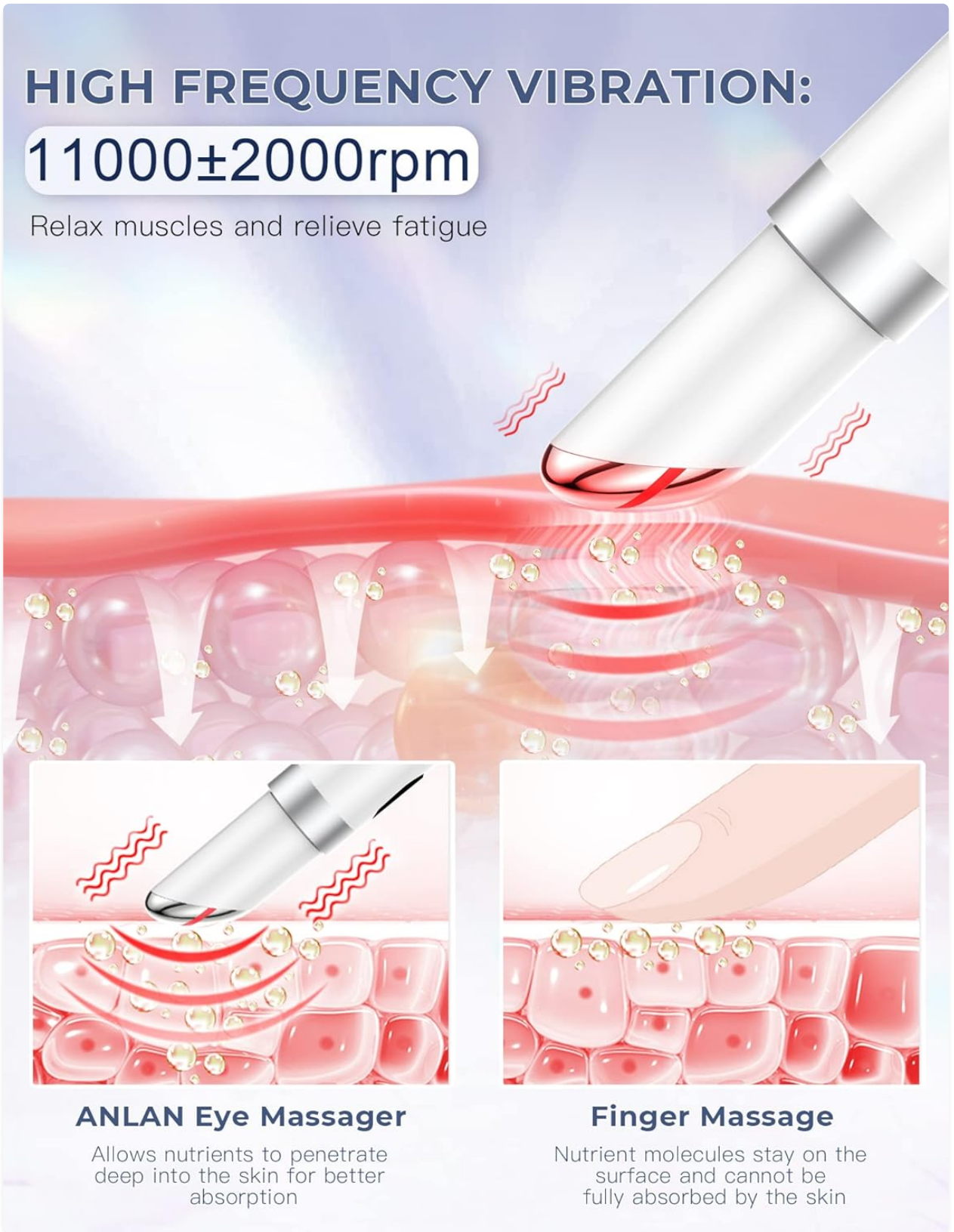

ANLAN Eye Massager

Allows nutrients to penetrate deep into the skin for better absorption

Finger Massage

Nutrient molecules stay on the surface and cannot be fully absorbed by the skin

Image: Illustration of high-frequency vibration relaxing muscles and aiding nutrient absorption.

- **Light Therapy Modes:**

The device incorporates three distinct light therapy options:



Image: Explanation of Red, Blue, and Yellow light therapy benefits.

- **Red Light (620 ± 10nm):** Promotes collagen regeneration.
- **Blue Light (470 ± 10nm):** Provides disinfection and bacteriostasis, inhibiting acne regeneration.
- **Yellow Light (590 ± 10nm):** Decomposes melanin to restore skin brightness and radiance.

7. MAINTENANCE

Cleaning

- After each use, wipe the massage head with a soft, damp cloth or a cotton pad.
- Do not use abrasive cleaners, alcohol, or harsh chemicals.
- Ensure the device is completely dry before storing.

Storage

- Store the device in a cool, dry place, away from direct sunlight and extreme temperatures.
- Keep the protective cap on the massage head when not in use.

8. TROUBLESHOOTING

If you encounter any issues with your device, please refer to the table below:

Problem	Possible Cause	Solution
Device does not turn on.	Low battery.	Charge the device fully.
No vibration.	Vibration mode not activated.	Long press the Light/Vibration button to activate.
Temperature not reaching desired level.	Device not in contact with skin or insufficient warm-up time.	Ensure good contact with skin. Allow a few moments for the device to heat up.
EMS sensation is too weak/strong.	EMS intensity setting.	Short press the Power/EMS button to adjust EMS intensity (3 levels).

If the problem persists after trying these solutions, please contact customer support.

9. SPECIFICATIONS

Feature	Detail
Model Number	LM-2201
Dimensions (L x W x H)	13 x 2.2 x 2.2 cm
Weight	170 grams
Power Source	Rechargeable Lithium-ion Battery (USB Cable)
Material	Plastic
Temperature Range	38-45°C
Vibration Frequency	11,000 ± 2,000 rpm
Light Wavelengths	Red: 620 ± 10nm, Blue: 470 ± 10nm, Yellow: 590 ± 10nm


10. WARRANTY AND SUPPORT

ANLAN products are manufactured to high-quality standards. This device comes with a standard manufacturer's warranty against defects in materials and workmanship. Please refer to your purchase documentation for specific warranty terms and duration.

For technical support, troubleshooting assistance, or warranty claims, please contact ANLAN customer service through the retailer where you purchased the product or visit the official ANLAN website for contact information.

© 2023 ANLAN. All rights reserved.

Related Documents

	<p>ANLAN Neck Beauty Device User Manual</p> <p>User manual for the ANLAN Neck Beauty Device, detailing its features, functions, usage instructions, and warranty information. The device utilizes positive and negative ions, vibration, warm massage, and red, blue, and green light therapy for skin cleansing, EMS care, and moisturizing.</p>
---	---

