

[Manuals.plus](#) /

> [Pursonic](#) /

> Pursonic Whitening USB Rechargeable Sonic Toothbrush User Manual

## Pursonic TBUSB200-12

# Pursonic Whitening USB Rechargeable Sonic Toothbrush User Manual

Model: TBUSB200-12

Brand: Pursonic

## INTRODUCTION

---

Thank you for choosing the Pursonic Whitening USB Rechargeable Sonic Toothbrush. This manual provides important information on the safe and effective use, maintenance, and care of your new sonic toothbrush. Please read it thoroughly before first use and keep it for future reference.



Image: The Pursonic Sonic Toothbrush, its packaging, and the included 12 brush heads.

## SAFETY INFORMATION

---

### Important Safety Precautions:

- Do not use for children.
- Be gentle with the soft bristles to avoid injury.
- Do not immerse the charging base in water.
- Consult your dentist before use if you have severe gum disease or oral surgery.

## PRODUCT FEATURES

---

- **38,000 Strokes Per Minute:** Powerful sonic vibrations for effective cleaning.
- **3 Brushing Modes:** Normal, Super, and Gumcare, catering to various oral care needs.
- **2-Minute Auto Timer:** Ensures dentist-recommended brushing time.
- **30-Second Interval Alert System:** Guides you to brush each quadrant of your mouth evenly.
- **Long-Lasting Battery:** One 6-hour charge provides up to 25 days of use (based on brushing twice a day for 2 minutes each time).

- **Smart LED Indicator:** Provides power shortage reminders and power-off protection.
- **12 Brush Heads Included:** A 3-year supply of Dupont nylon bristles for safe and comfortable brushing.



## Ultimate Technology Features

### Go Beyond Just Cleaning



2 minute timer with 30 second interval alert system



Smart LED indicator with power shortage reminder and power-off protection



38,000 strokes per minute



Intelligent memory function  
(Unit will restart on your last usage mode)



Super charged NI-MH battery 6000mAh, full charged 6 hours can last upto 25 days (twice a day based on 2 min use)

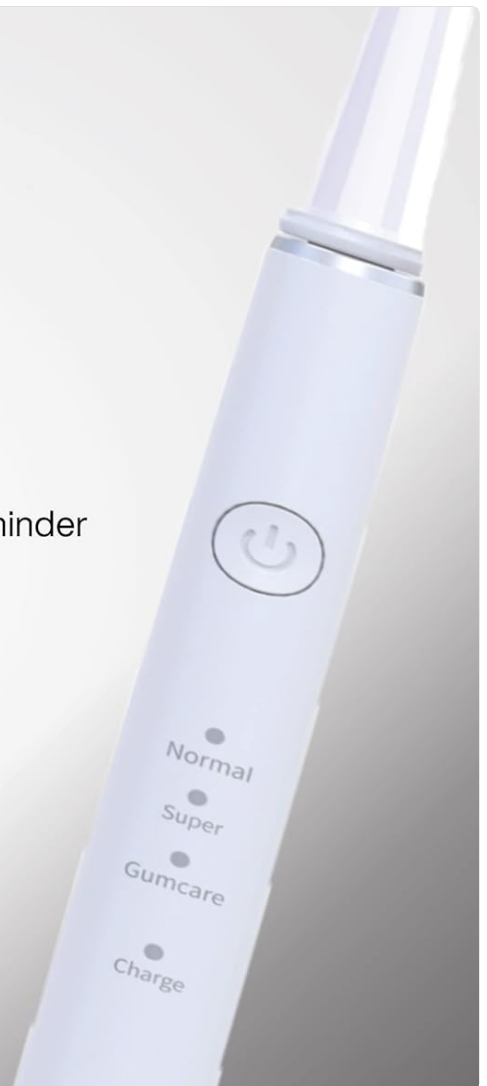


Image: Key technology features of the Pursonic Sonic Toothbrush, including timer, LED indicator, stroke count, and battery life.

## 3 BRUSHING MODES

USB RECHARGABLE SONIC TOOTH BRUSH



**NORMAL**  
Regular Cleaning



**SUPER**  
Eliminate stubborn stains and brighten your teeth.



**GUMCARE**  
Provide massage style brushing experience.



Comes with 12 brush heads and that it is 3 year supply

Image: Explanation of the three brushing modes: Normal for regular cleaning, Super for stain elimination, and Gumcare for massage.

## SETUP AND INITIAL USE

---

1. **Charge the Toothbrush:** Before first use, fully charge the toothbrush for approximately 6 hours. The LED indicator will show charging status.
2. **Attach Brush Head:** Align the brush head with the metal shaft on the toothbrush handle and push it down firmly until it clicks into place.
3. **Apply Toothpaste:** Apply a small amount of toothpaste to the bristles.



Image: The toothbrush handle alongside the included 12 brush heads, ready for attachment.

## OPERATING INSTRUCTIONS

---

1. **Power On/Off:** Press the power button to turn the toothbrush on or off.
2. **Select Brushing Mode:** Press the power button repeatedly to cycle through the three brushing modes: Normal, Super, and Gumcare. The indicator light next to the selected mode will illuminate.

3. **Brush Your Teeth:** Place the bristles against your teeth at a slight angle towards the gumline. Move the brush slowly across your teeth, allowing the sonic vibrations to do the work.
4. **Follow the Timer:** The toothbrush will pause briefly every 30 seconds to indicate when to move to the next quadrant of your mouth. After 2 minutes, the toothbrush will automatically turn off.



Image: A user demonstrating the proper way to hold and use the sonic toothbrush.

## MAINTENANCE AND CARE

---

- **After Each Use:** Rinse the brush head thoroughly under running water. Wipe the handle with a damp cloth.
- **Deep Cleaning:** Periodically remove the brush head and rinse both the brush head and the handle connection area.
- **Brush Head Replacement:** Replace the brush head every 3 months, or sooner if the bristles become frayed. The included 12 brush heads provide a 3-year supply.
- **Charging:** Recharge the toothbrush when the charge indicator light signals low battery.



## 12 Brush Heads is 3 Years Supply

No need for expensive brush head replacements



Image: A visual representation of the 12 brush heads included, emphasizing the long-term supply.

### TROUBLESHOOTING

Problem	Possible Cause	Solution
Toothbrush does not turn on.	Battery is depleted.	Charge the toothbrush for 6 hours.
Weak vibration.	Low battery or brush head not properly attached.	Recharge the toothbrush. Ensure brush head is firmly attached.
Brush head falls off during use.	Brush head not securely attached.	Push the brush head down firmly until it clicks into place.

### SPECIFICATIONS

**Model Number:** TBUSB200-12

**Strokes Per Minute:** 38,000

**Brushing Modes:** Normal, Super, Gumcare

**Timer:** 2-minute auto-timer with 30-second interval alerts

**Power Source:** Battery Powered (USB Rechargeable)

**Battery Life:** Up to 25 days on a single charge

**Bristle Firmness:** Soft

**Included:** 1 Pursonic Sonic Toothbrush, 12 Brush Heads, USB Charging Cable

## WARRANTY AND SUPPORT

---

For warranty information or technical support, please refer to the contact details provided with your product packaging or visit the official Pursonic website. Keep your purchase receipt for warranty claims.

## LEGAL DISCLAIMER

---

*Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.*