

[Manuals.plus](#) /

> [TOUCHBeauty](#) /

> TOUCHBeauty Facial Massager with Ion and LED Light Therapy, Model AG-1713 Instruction Manual

TOUCHBeauty AG-1713

TOUCHBeauty AG-1713 Facial Massager Instruction Manual

Model: AG-1713

1. PRODUCT OVERVIEW

The TOUCHBeauty AG-1713 is a versatile cosmetic device designed for comprehensive facial care, offering deep cleansing, enhanced absorption, anti-aging benefits, and relaxation. It integrates ion technology, microcurrent, LED light therapy, and massage functions to improve skin texture and appearance.

Key features include:

- **Positive Ion Cleansing:** Deeply cleanses pores, removing impurities.
- **Negative Ion Absorption:** Enhances the penetration and efficacy of skincare products.
- **Microvibration:** Gentle massage to stimulate skin.
- **LED Light Therapy:** Yellow LED light stimulates collagen production and brightens skin.
- **Interchangeable Heads:** Includes a facial beauty head, microcurrent scalp head, scalp massage head, and microcurrent facial head for various treatments.
- **Smart Timer:** Automatic shut-off after 2.5 minutes for safety and consistent treatment duration.
- **Rechargeable Battery:** USB rechargeable, providing up to 40 minutes of use on a single charge.

5 in 1 Beauty & Health Care Device



Image 1.1: The TOUCHBeauty AG-1713 device showcasing its five main functions: ION+/-, Microcurrent, 2 Colors Light Therapy, Scalp Massage, and Eye Massage.



Image 1.2: The device head highlighting its benefits for skin and relaxation, including deep cleansing, skin brightening, nourishment, head relaxation, and eye relief.

2. DEVICE COMPONENTS

Familiarize yourself with the main components of your TOUCHBeauty AG-1713 device:



Image 2.1: Labeled diagram of the device, showing the main body, interchangeable heads (Ion+/- Beauty Head, Microcurrent Scalp Head, Scalp Massage Head, Microcurrent Facial Head), control buttons (Power, ION+/-, LED), charging port, and accessories like the USB cable and cotton ring.

- **Device Body:** The main handheld unit.
- **Connection Point:** Where interchangeable heads attach.
- **Light Therapy Head:** Emits LED light for skin treatment.
- **Conductive Bars:** For microcurrent and ion functions.
- **Power Button:** Turns the device on/off.
- **ION+/- Microcurrent Button:** Selects ion and microcurrent modes.
- **LED Button:** Activates and cycles through LED light therapy modes.
- **USB Cable:** For charging the device.
- **Charging Port:** Located at the bottom of the device.
- **Cotton Ring:** Used to secure cotton pads for deep cleansing.
- **Interchangeable Heads:**
 - Ion+/- Beauty Head

- Microcurrent Scalp Head
- Scalp Massage Head
- Microcurrent Facial Head

3. SETUP AND CHARGING

3.1 Initial Charge

Before first use, fully charge the device. Connect the provided USB cable to the charging port at the bottom of the device and plug the other end into a compatible USB power adapter (not included) or computer USB port. The charging indicator light will show the charging status.

3.2 Attaching Heads

To attach an interchangeable head, align it with the connection point on the device body and gently twist or push until it clicks securely into place. To remove, gently twist or pull the head off.

4. OPERATING INSTRUCTIONS

The TOUCHBeauty AG-1713 offers multiple modes for various skin and scalp care needs. Each mode is designed for specific benefits.

4.1 General Usage Steps (Facial Treatment)

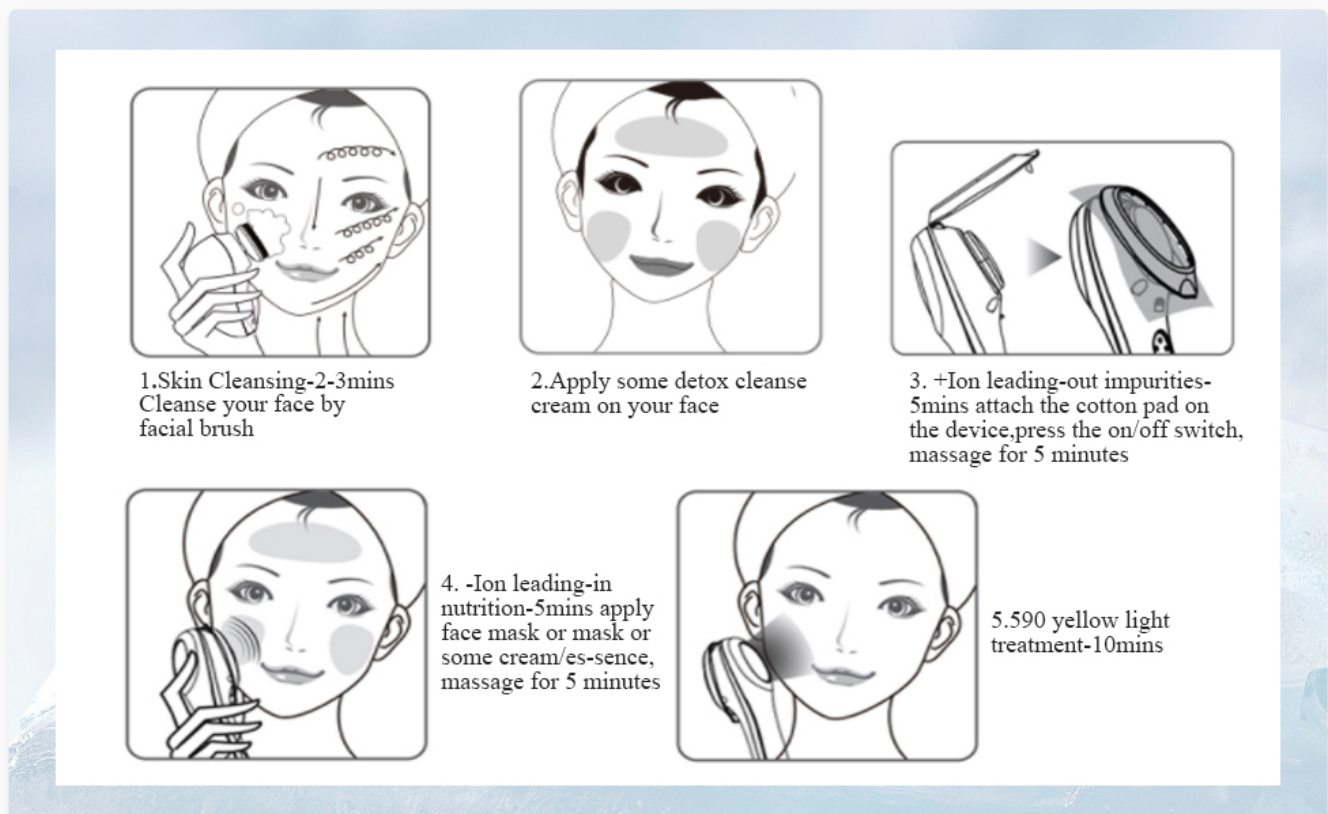


Image 4.1: A visual guide illustrating the recommended 5-step facial treatment routine using the device, from cleansing to light therapy.

- 1. Skin Cleansing (2-3 minutes):** Cleanse your face thoroughly using a facial brush or your preferred method.
- 2. Apply Detox Cleanse Cream:** Apply a suitable detox or cleansing cream to your face.
- 3. Positive Ion Leading-Out (5 minutes):** Attach a cotton pad with the cotton ring to the device. Press the ION+/- Microcurrent button to select the positive ion mode. Gently massage your face for 5 minutes to draw out

impurities.

4. **Negative Ion Leading-In (5 minutes):** Apply your favorite face mask, cream, or essence. Press the ION+/- Microcurrent button to select the negative ion mode. Massage for 5 minutes to enhance absorption.
5. **Yellow Light Treatment (10 minutes):** Press the LED button to activate the yellow light mode. Move the device gently over your face for 10 minutes.

4.2 Ion and Microcurrent Modes

The device utilizes positive and negative ions, along with microcurrent, to optimize skin health.

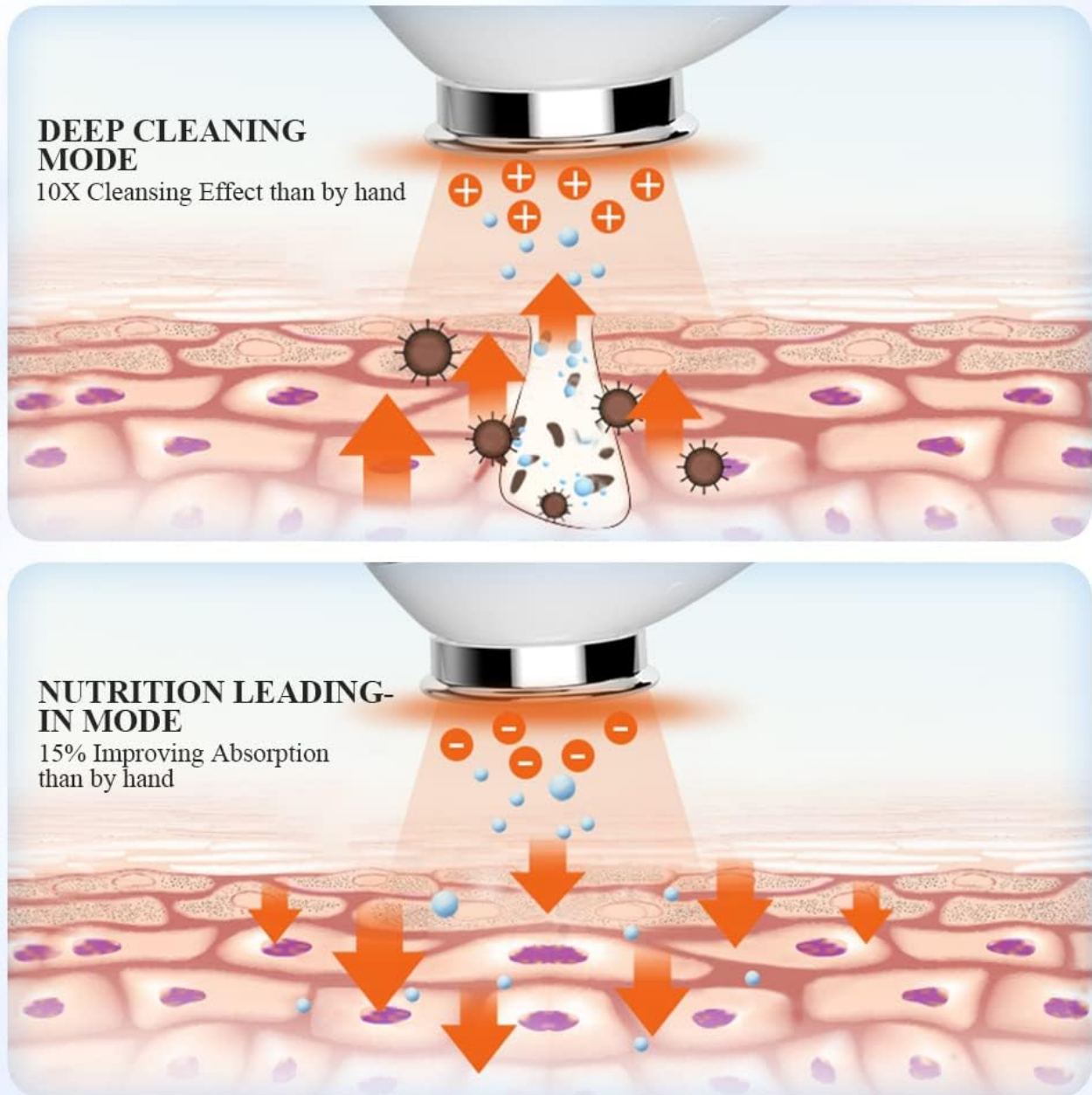


Image 4.2: Illustration of how positive ions facilitate deep cleaning by attracting impurities from pores, and how negative ions enhance nutrient absorption into the skin.

- **Deep Cleaning Mode (Positive Ions):** Use with a cotton pad and cleansing lotion. The positive ions attract negatively charged impurities and dirt from deep within the pores, providing a thorough cleanse.
- **Nutrition Leading-In Mode (Negative Ions):** Apply your preferred serum, cream, or mask. The negative ions help push active ingredients deeper into the skin, improving absorption and effectiveness.
- **Microcurrent:** This function works in conjunction with the ion modes to stimulate facial muscles and promote circulation.

4.3 LED Light Therapy

The device features two LED light modes for targeted skin benefits.

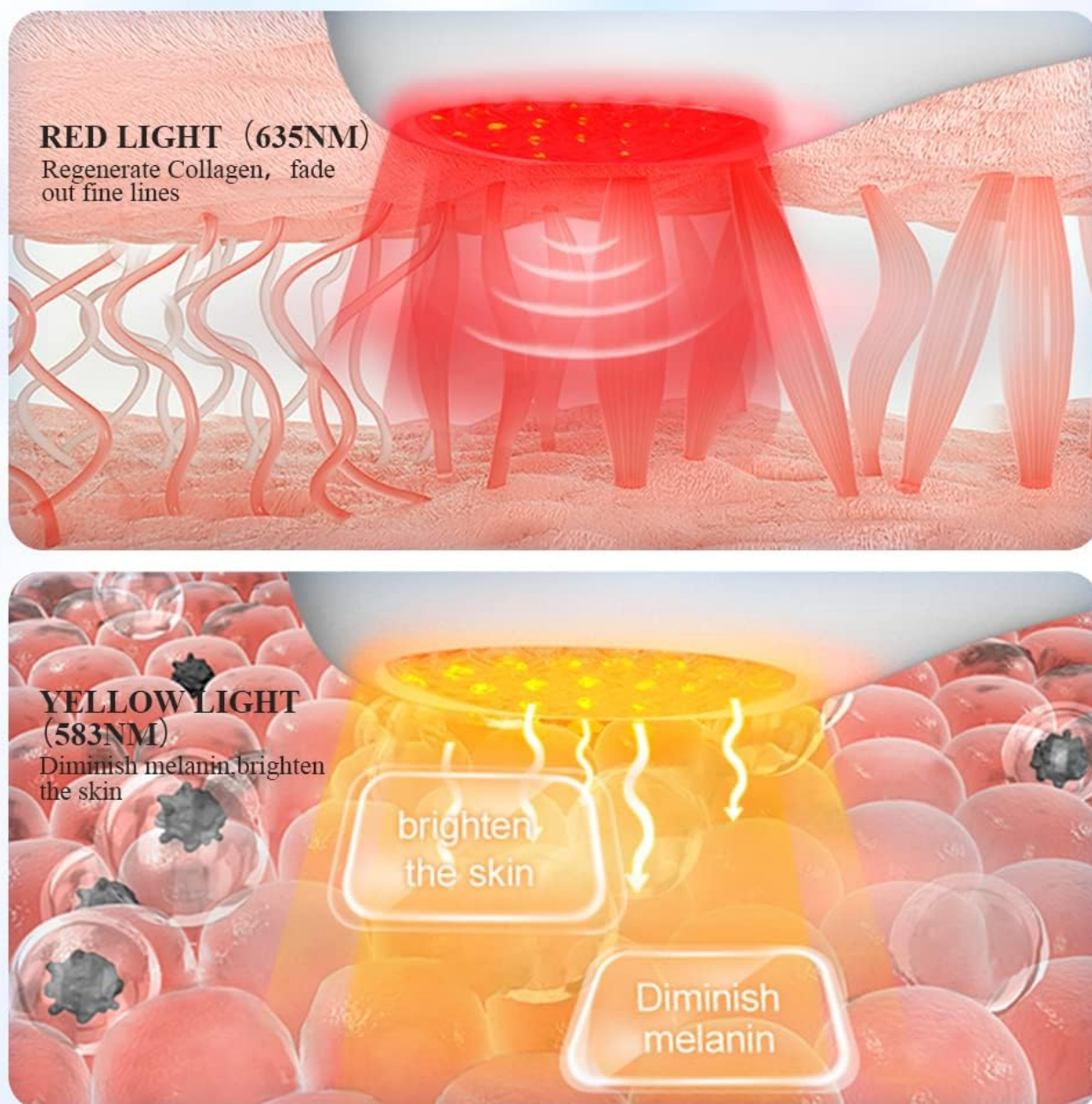


Image 4.3: Visual representation of the effects of Red Light (635nm) on collagen regeneration and fine line reduction, and Yellow Light (583nm) on melanin reduction and skin brightening.

- **Red Light (635nm):** Helps to regenerate collagen, reduce the appearance of fine lines, and improve skin elasticity.
- **Yellow Light (583nm):** Stimulates collagen production, strengthens cells, diminishes melanin, and brightens the skin tone.

4.4 Scalp and Eye Massage

Utilize the specialized heads for targeted massage.

Scalp Massager



Image 4.4: A user demonstrating the application of the scalp massager attachment for head relaxation and stimulation.

- **Scalp Massage Head:** Attach the scalp massage head and gently move it across your scalp to promote circulation and relaxation.
- **Eye Massage:** Use the appropriate head (e.g., Ion+/- Beauty Head or Microcurrent Facial Head) around the eye area with gentle movements to reduce puffiness and improve product absorption.

Smart Timer: The device is equipped with an automatic timer that turns off after 2.5 minutes of continuous operation in each mode. This ensures consistent treatment duration and prevents overuse. Simply restart the device if you wish to continue your session.

5. MAINTENANCE AND CARE

Proper maintenance ensures the longevity and hygiene of your device.

- **Cleaning the Device:** After each use, wipe the device head and body with a soft, damp cloth. For thorough cleaning, use a mild, non-abrasive cleanser. Ensure the device is completely dry before storing.
- **Do Not Submerge:** Do not immerse the device in water or any other liquid.

- **Storage:** Store the device and its attachments in a cool, dry place, away from direct sunlight and extreme temperatures. The provided storage box (Image 5.1) is ideal for this purpose.
- **Battery Care:** To maintain battery health, avoid fully discharging the device frequently. Charge it when the battery indicator is low.

Host size: 168*43*39mm
Size of storage box: 219*128*62mm



Image 5.1: The device and its dedicated storage box, designed for convenient and protective storage. Host size: 168*43*39mm, Size of storage box: 219*128*62mm.

6. TROUBLESHOOTING

This section addresses common issues you might encounter with your TOUCHBeauty AG-1713 device.

- **Device Not Turning On:**
 - Ensure the device is fully charged. Connect it to the charger and wait for the charging indicator to confirm.
 - Press and hold the Power button for a few seconds to ensure it registers.
- **No Vibration/Light in a Mode:**

- Confirm the correct head is securely attached.
- Ensure the device is powered on and the specific mode (ION+/-, LED) is selected using its respective button.
- The device has a smart timer and will turn off after 2.5 minutes. If it stops, simply restart the mode.

- **Reduced Effectiveness:**

- Ensure the device head is clean and free from residue.
- For ion modes, ensure proper contact with the skin and that a suitable skincare product (or cotton pad with lotion for positive ion mode) is used.
- Ensure the device is fully charged for optimal performance.

If you experience issues not covered here, please refer to the contact information for customer support.

7. SPECIFICATIONS

Attribute	Detail
Brand	TOUCHBeauty
Model Number	AG-1713
Product Dimensions (Host)	16.8 x 4.3 x 3.9 cm (6.61 x 1.69 x 1.54 inches)
Product Weight	790 g (1.74 lbs)
Storage Box Dimensions	21.9 x 12.8 x 6.2 cm (8.62 x 5.04 x 2.44 inches)
Power Source	USB Rechargeable Battery
Battery Life	Up to 40 minutes per charge
Charging Time	Less than 3 hours
LED Wavelengths	Red Light: 635nm, Yellow Light: 583nm
Automatic Timer	2.5 minutes per mode

8. WARRANTY AND SUPPORT

Specific warranty details for the TOUCHBeauty AG-1713 device are not provided in the product information. Please refer to the retailer or manufacturer's official website for the most up-to-date warranty policy and customer support contact information.

For technical assistance, product inquiries, or warranty claims, please contact TOUCHBeauty customer service through their official channels.

