

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [Oral-B](#) /
- › [Oral-B iO Series 3 Limited Electric Toothbrush Instruction Manual](#)

### Oral-B iO G3.2i6.1KD

# Oral-B iO Series 3 Limited Electric Toothbrush Instruction Manual

Model: iO G3.2i6.1KD

## 1. INTRODUCTION

---

The Oral-B iO Series 3 Limited Electric Toothbrush is designed to provide a superior cleaning experience. It features advanced iO technology with gentle micro-vibrations and a unique round brush head for effective plaque removal along the gumline. This manual provides essential information for the safe and optimal use of your new electric toothbrush.



Figure 1: Oral-B iO Series 3 Limited Electric Toothbrush with accessories.

## 2. PACKAGE CONTENTS

---

Your Oral-B iO Series 3 Limited package includes:

- 1 Oral-B iO Series 3 Limited Handle
- 2 Oral-B iO Brush Heads (1 Ultimate Clean, 1 Gentle Care)
- 1 Charging Base
- 1 Travel Case
- 1 Brush Head Refill Holder

For a visual guide on unboxing and included items, please refer to the video below:

Video 1: Unboxing and initial overview of the Oral-B iO Series 3 Limited Electric Toothbrush.

## 3. SETUP

---

### 3.1 Charging the Toothbrush

1. Plug the charging base into an electrical outlet.

2. Place the toothbrush handle on the charging base. The battery indicator will light up to show charging progress.
3. A full charge typically takes approximately 16 hours.

### 3.2 Attaching the Brush Head

Align the brush head with the shaft of the toothbrush handle and push it down until it clicks into place. Ensure it is securely attached before use.

## 4. OPERATING INSTRUCTIONS

---

### 4.1 Power On/Off

Press the power button to turn the toothbrush on or off.

### 4.2 Cleaning Modes

The Oral-B iO Series 3 Limited offers 3 cleaning modes:

- **Daily Clean:** For everyday cleaning.
- **Sensitive:** For gentle yet thorough cleaning of sensitive areas.
- **Whitening:** Designed to polish and whiten teeth.

Press the mode button (below the power button) to cycle through the available modes. The active mode will be indicated by an illuminated icon on the handle.

# Easily choose from 3 cleaning modes



SENSITIVE



DAILY  
CLEAN



WHITEN

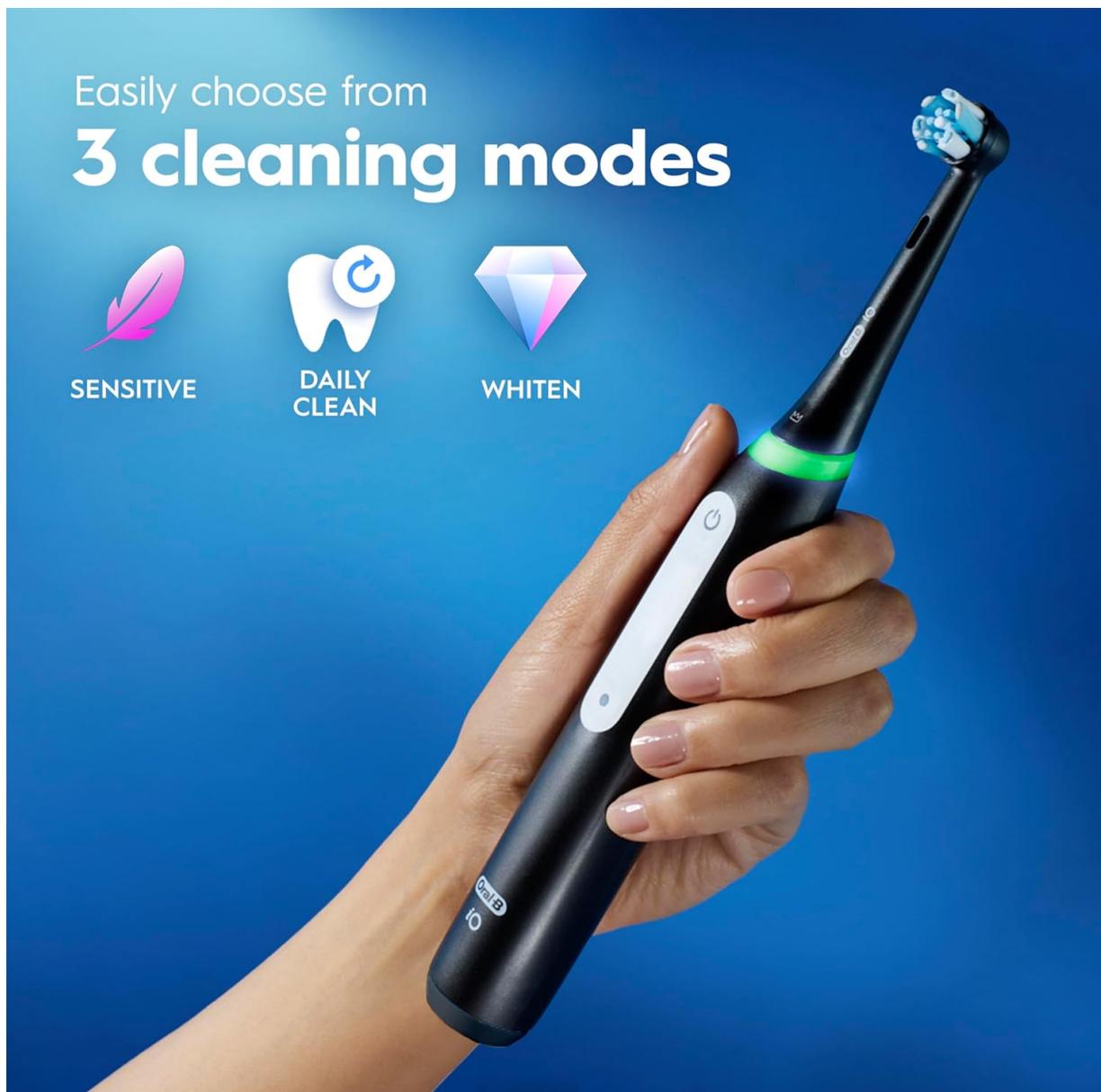


Figure 2: Icons indicating the three cleaning modes.

## 4.3 2-Minute Timer

The built-in 2-minute timer ensures you brush for the dentist-recommended duration. The toothbrush will briefly pause every 30 seconds to signal when to move to the next quadrant of your mouth, and will automatically turn off after 2 minutes.

## 4.4 Smart Pressure Sensor

The pressure sensor helps protect your gums and teeth from excessive brushing force. It provides real-time feedback:

- **Green Light:** Optimal brushing pressure.
- **Red Light:** Brushing too hard. Reduce pressure immediately.
- **Blue Light:** Brushing too softly. Apply slightly more pressure.



Figure 3: Pressure sensor light indicators.

For a detailed breakdown of the toothbrush's working modes and pressure sensor, watch this video:

Video 2: Demonstration of cleaning modes and pressure sensor functionality.

## 5. MAINTENANCE

### 5.1 Cleaning the Toothbrush

After each use, rinse the brush head thoroughly under running water. You can also detach the brush head and rinse both the head and the handle separately. Wipe the handle dry with a clean cloth.

### 5.2 Replacing Brush Heads

Dentists recommend replacing your brush head every 3 months, or sooner if the bristles become worn or frayed. The blue indicator bristles on your brush head will fade to white over time, signaling when it's time for a replacement.

## 6. TRAVEL CASE

The included travel case allows you to conveniently store and transport your toothbrush and up to two brush heads, keeping them hygienic and protected during travel.



Figure 4: Oral-B iO Series 3 Limited with travel case.

## 7. TROUBLESHOOTING

- **Toothbrush not turning on:** Ensure the toothbrush is fully charged.
- **Pressure sensor not working:** Check if the brush head is correctly attached. If issues persist, contact customer support.
- **Reduced battery life:** Ensure the charging base is properly connected and the toothbrush is placed correctly on it.

## 8. SPECIFICATIONS

<b>Model Number</b>	iO G3.2i6.1KD
<b>Power Source</b>	Battery Powered (1 Lithium Ion battery included)
<b>Special Features</b>	Brushing Timer, Pressure Sensor, Rechargeable, Travel Size
<b>Age Range</b>	Adult (Do not use on children under age 3)
<b>Manufacturer</b>	Procter & Gamble

## 9. WARRANTY AND SUPPORT

The Oral-B iO Series 3 Limited Electric Toothbrush comes with a 2-year warranty. Additionally, Oral-B offers a 30-day money-back guarantee. For further assistance or warranty claims, please refer to the contact information provided in your product packaging or visit the official Oral-B website.