



Manuals.plus /

- › Oral-B /
- › Oral-B Smart 4 4500 Electric Toothbrush User Manual

## Oral-B Smart 4 4500

# Oral-B Smart 4 4500 Electric Toothbrush User Manual

Model: Smart 4 4500 | Brand: Oral-B

## 1. INTRODUCTION AND OVERVIEW

The Oral-B Smart 4 4500 electric toothbrush is engineered to provide a superior cleaning experience, removing up to 100% more plaque than a manual toothbrush. It combines advanced cleaning technology with smart features to guide you towards healthier gums and a brighter smile. This manual provides essential information for the proper setup, operation, and maintenance of your new electric toothbrush.

### Key Features:

- **Superior Plaque Removal:** Round brush head designed for effective cleaning and healthier gums.
- **Smart Brushing Guidance:** Connects to the Oral-B app via Bluetooth for real-time feedback.
- **Visible Pressure Sensor:** Protects gums by alerting you when brushing too hard.
- **Gentle Whitening:** Removes surface stains for a whiter smile from day one.
- **Long-Lasting Battery:** Lithium-ion battery provides over 2 weeks of brushing on a single charge.
- **Professional Timer:** Ensures you brush for the dentist-recommended 2 minutes.



Figure 1: The Oral-B Smart 4 4500 electric toothbrush package, including the handle, charger, brush heads, travel case, and toothpaste.

## 2. SETUP

Follow these steps to prepare your Oral-B Smart 4 4500 for first use:

1. **Charge the Toothbrush:** Before initial use, fully charge the toothbrush. Place the handle on the charging unit. The battery indicator light will show charging status. A full charge can take up to 12 hours.
2. **Attach a Brush Head:** Firmly push a brush head onto the handle's shaft until it clicks into place.
3. **Download the Oral-B App (Optional):** For smart guidance and tracking, download the Oral-B app from your smartphone's app store. Enable Bluetooth on your phone and follow the in-app instructions to pair your toothbrush.



Figure 2: Components of the Oral-B Smart 4 4500, showing the handle, brush heads, charger, and travel case.



Figure 3: The Bluetooth symbol on the toothbrush handle, indicating connectivity for smart features.

### 3. OPERATING INSTRUCTIONS

To achieve optimal cleaning results, follow these operating guidelines:

1. **Prepare for Brushing:** Wet the brush head and apply a small amount of toothpaste.
2. **Power On and Select Mode:** Press the power button to turn on the toothbrush. Press the power button again to cycle through the available brushing modes (Daily Clean, Sensitive, Whitening). Select your preferred mode.
3. **Brush Your Teeth:** Guide the brush head slowly from tooth to tooth, spending a few seconds on each surface. Do not scrub hard. The round brush head is designed to clean each tooth individually.
4. **Utilize the Professional Timer:** The toothbrush will briefly stutter every 30 seconds to indicate when to move to the next quadrant of your mouth. A longer stutter will signal the completion of the dentist-recommended 2-minute brushing time.
5. **Monitor Pressure:** The visible pressure sensor will illuminate red if you are applying too much pressure. Reduce

pressure immediately to protect your gums and teeth.

6. **Power Off:** After completing your brushing session, press the power button to turn off the toothbrush.



Figure 4: Proper technique for brushing with the Oral-B electric toothbrush.

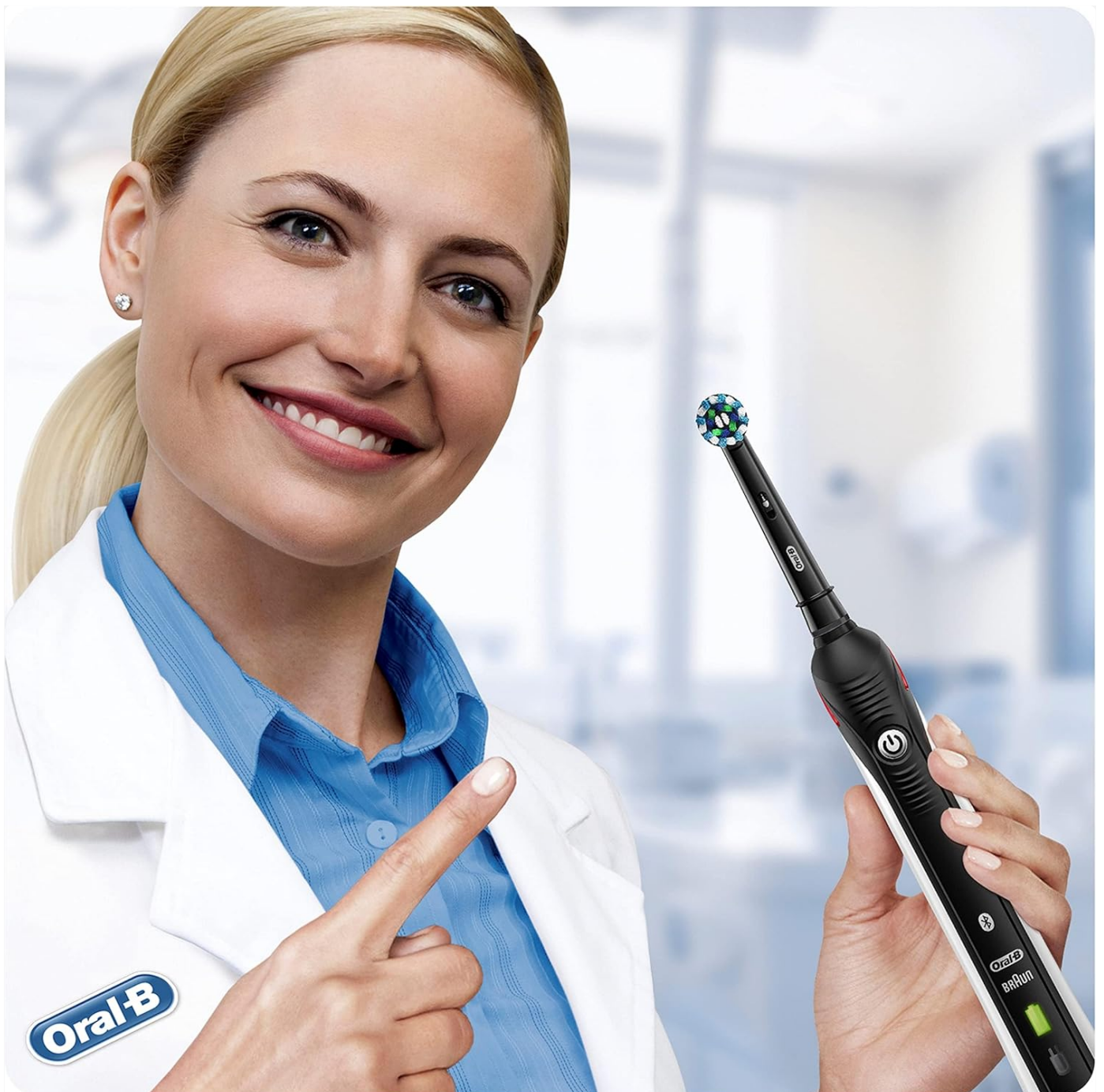


Figure 5: The red light on the handle indicates that too much pressure is being applied during brushing.

#### 4. MAINTENANCE

Proper maintenance ensures the longevity and effectiveness of your Oral-B Smart 4 4500 toothbrush:

- **Rinse Brush Head:** After each use, rinse the brush head thoroughly under running water with the handle switched on. Remove the brush head and clean the handle and brush head separately.
- **Clean Handle:** Wipe the handle with a damp cloth regularly to remove any toothpaste residue. Do not immerse the handle in water.
- **Replace Brush Head:** It is recommended to replace your brush head every three months, or sooner if the bristles are splayed or worn. Worn bristles can reduce cleaning effectiveness.
- **Recharge Battery:** Recharge the toothbrush when the battery indicator shows low power. The lithium-ion battery provides over two weeks of regular use on a full charge.

# SMART

## CARATTERISTICHE DI SPAZZOLAMENTO



Batteria al litio - fino a 2 settimane di ricarica



Timer Professionale



Sensore della pressione di spazzolamento visibile



App con guida intelligente via Bluetooth



3 Modalità di spazzolamento

Figure 6: A selection of Oral-B brush heads, highlighting the variety available for replacement.

## 5. TROUBLESHOOTING

If you encounter any issues with your Oral-B Smart 4 4500, refer to the following common troubleshooting tips:

- **Toothbrush Not Charging:**

Ensure the charging unit is securely plugged into a working power outlet. Verify that the toothbrush handle is correctly seated on the charging base. Check for any debris preventing proper contact.

- **Pressure Sensor Always Red:**

This indicates consistent excessive pressure. Try to brush more gently. The sensor is designed to protect your gums and teeth from damage.

- **Bluetooth Connectivity Issues:**

Ensure Bluetooth is enabled on your smartphone and the Oral-B app is open and updated. Try restarting both your smartphone and the toothbrush. If issues persist, try unpairing and re-pairing the device in the app settings.

- **Shorter Battery Life Than Expected:**

Battery life can vary based on usage frequency, brushing mode, and age of the battery. Ensure the toothbrush is fully charged before each use. If the issue persists after a full charge cycle, contact customer support.

- **Toothbrush Not Turning On/Off:**

Ensure the battery is charged. If the issue persists, try cleaning the power button area for any debris. If the problem continues, the unit may require service.

## 6. SPECIFICATIONS

Feature	Detail
Brand	Oral-B
Model Number	8006540747544
Power Source	Electric Cord
Special Features	Pressure Sensor, Rechargeable
Battery Type	1 Lithium-ion (included)
Brush Head Type	Round
Age Range	Adult
Number of Items	1 (toothbrush unit)
Package Dimensions	26.9 x 19 x 15 cm
Package Weight	770 grams

## 7. SAFETY WARNINGS

Please read and adhere to the following safety warnings to ensure safe operation and prevent injury:

- Do not use the electric toothbrush for any purpose other than oral hygiene.
- Supervise children closely when the toothbrush is in use to prevent choking hazards from small parts.
- It is recommended to replace the brush head every three months to maintain hygiene and effectiveness.
- Avoid applying excessive pressure during brushing. The visible pressure sensor will alert you if you are pressing too hard, which can damage tooth enamel and gums.
- Do not modify or attempt to repair the appliance. Refer all servicing to qualified personnel.
- Keep the charging unit away from water. Do not place or store it where it can fall or be pulled into a tub or sink.

## 8. WARRANTY AND SUPPORT

For detailed warranty information, product registration, or customer support, please refer to the documentation included in your product packaging. You may also visit the official Oral-B website for frequently asked questions, troubleshooting guides, and contact information for customer service.

Please retain your proof of purchase for any warranty claims.

