

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [Oral-B](#) /
- › [Oral-B Pro 100 CrossAction Battery Powered Toothbrush User Manual](#)

## Oral-B P\_1651

# Oral-B Pro 100 CrossAction Battery Powered Toothbrush User Manual

Model: P\_1651

## INTRODUCTION

The Oral-B Pro 100 CrossAction Battery Powered Toothbrush is designed to provide a superior clean compared to a regular manual toothbrush. Its perfectly angled bristles reach deep between teeth for effective plaque removal and healthier gums. This manual provides essential information for the setup, operation, maintenance, and specifications of your Oral-B Pro 100 toothbrush.

## SAFETY INFORMATION

- Do not use this product as a toy or without adult supervision.
- Small parts may pose a choking hazard.
- Do not chew or bite on the brush head.
- DO NOT USE ON CHILDREN UNDER AGE 3.
- Prior to use, ensure the brush head fits well.
- If the product is dropped, the brush head should be replaced before the next use, even if no damage is visible.
- Dispose of or recycle batteries in accordance with local regulations.

## PACKAGE CONTENTS

Your Oral-B Pro 100 CrossAction Battery Powered Toothbrush package includes:

- 1 Oral-B Pro 100 Toothbrush Handle (Black)
- 1 CrossAction Brush Head
- 2 AA Batteries





Figure 1: Oral-B Pro 100 Power Toothbrush package contents.

## SETUP

### 1. Inserting Batteries

1. Remove the battery cover from the bottom of the toothbrush handle by pressing the two small bumps on each side and pulling down.
2. Insert the two AA batteries, ensuring correct polarity (+/-).
3. Replace the battery cover, pushing it firmly until it clicks into place.



Figure 2: Inserting AA batteries into the toothbrush handle.

### 2. Attaching the Brush Head

1. Align the brush head with the metal shaft on the toothbrush handle.
2. Push the brush head down firmly until it snaps into place. Ensure the small plastic grommet is in place to secure the head.

## OPERATING INSTRUCTIONS

---

### 1. Turning On/Off

- Press the power button (marked with a universal power symbol) to turn the toothbrush on.
- Press the power button again to turn the toothbrush off.



Figure 3: Oral-B Pro 100 Toothbrush with power button.

## 2. Brushing Technique

- Apply a small amount of toothpaste to the brush head.
- Place the brush head against your teeth, ensuring it covers one tooth at a time.
- Move the brush head slowly from tooth to tooth, allowing the rotating action to clean each surface. The round brush head is designed for a tooth-by-tooth clean.
- The toothbrush features a 2-minute coaching timer, which helps ensure you brush for the dentist-recommended duration.

# Round cleans 100% better

than a regular  
manual toothbrush


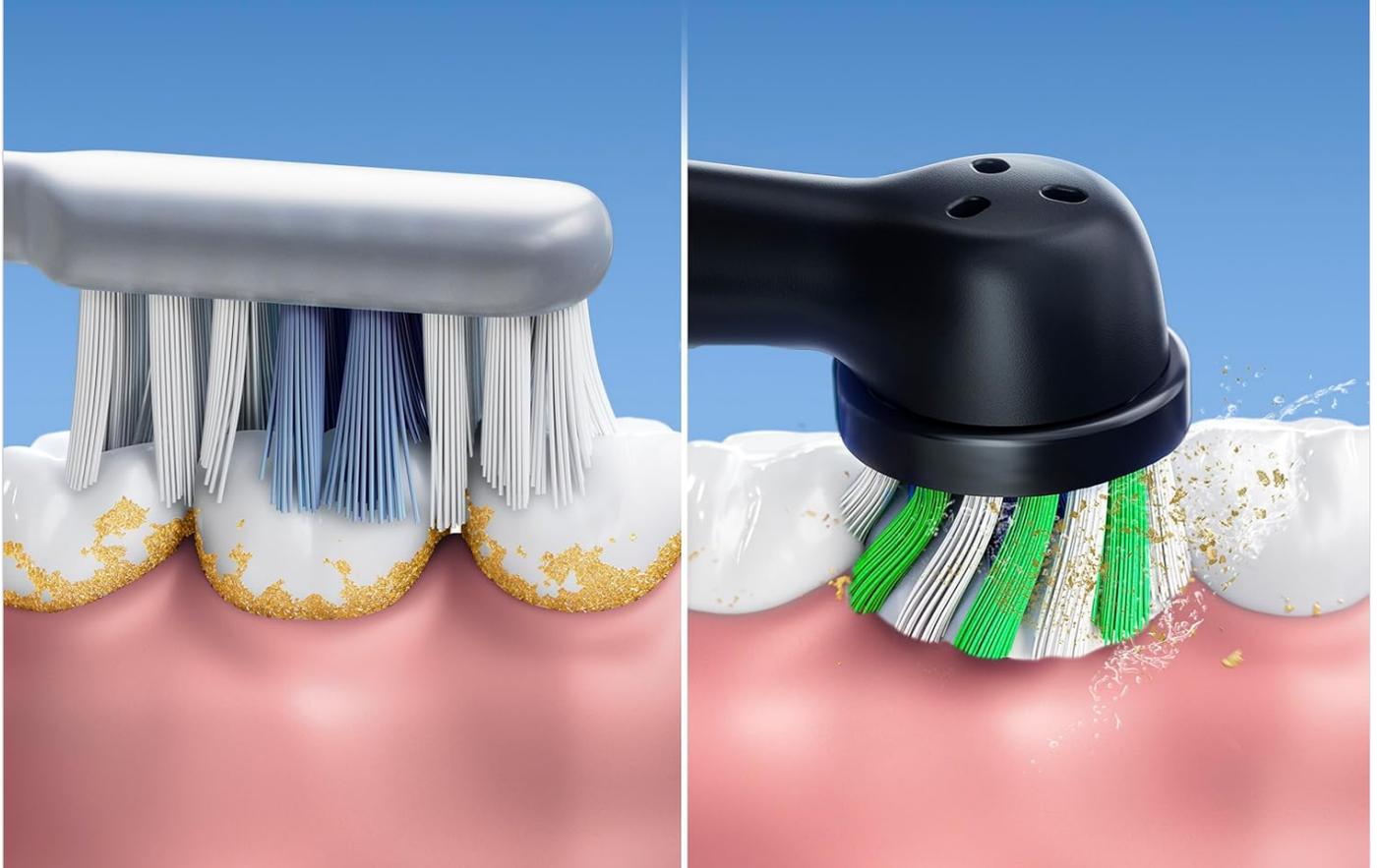


Figure 4: The round brush head provides superior cleaning.

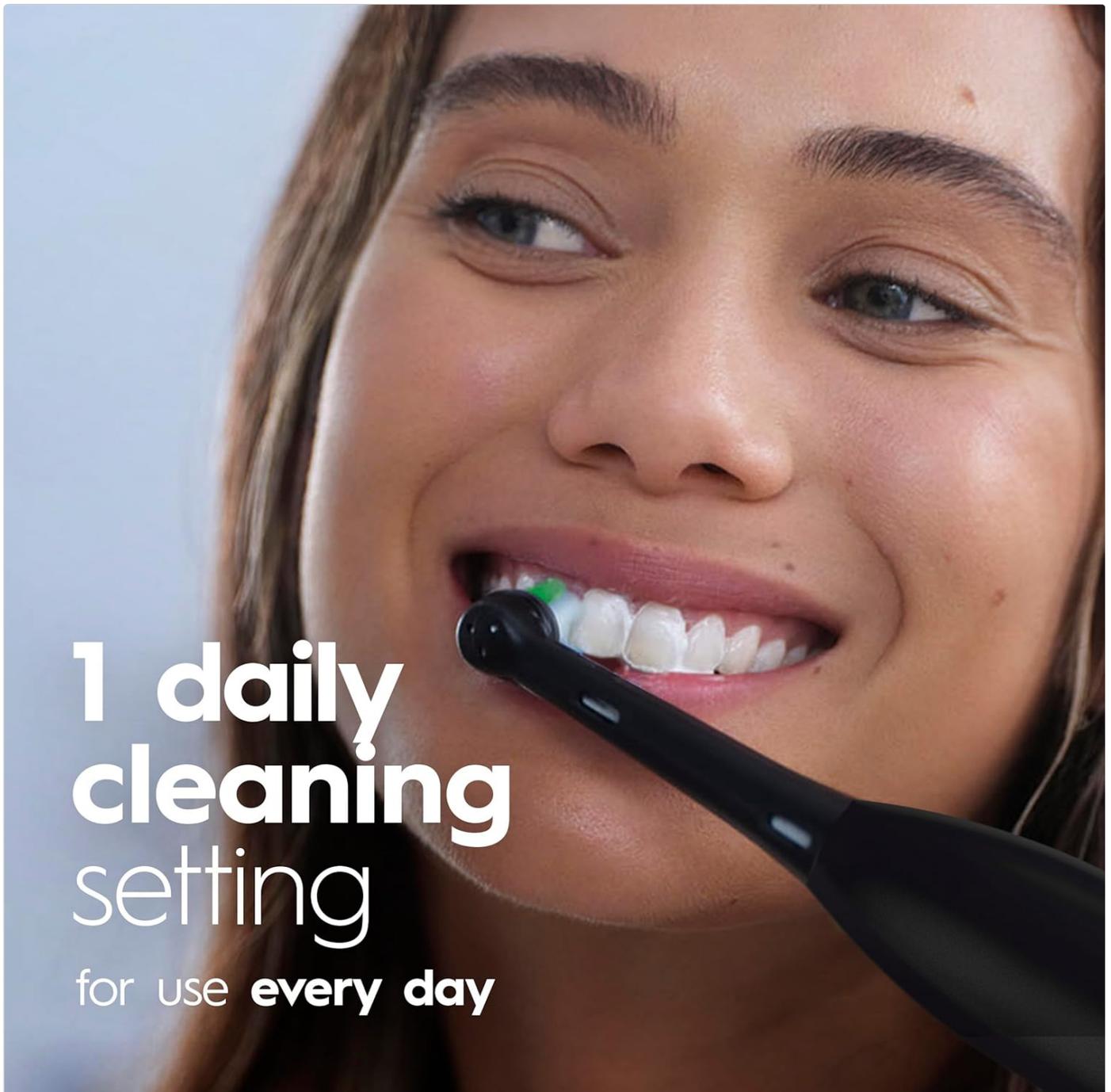

**1 daily  
cleaning  
setting**  
for use **every day**

*Figure 5: Daily cleaning with the Oral-B Pro 100.*

## **MAINTENANCE**

---

### **1. Cleaning**

- After each use, switch off the handle and remove the brush head.
- Rinse both parts separately with water.
- Wipe them dry before reassembling.

### **2. Replacing the Brush Head**

- Replace the brush head every 3 months, or sooner if the brush head becomes worn or no longer fits properly.
- The Oral-B Pro 100 is compatible with various Oral-B brush heads, including CrossAction, Sensitive Gum Care, FlossAction, and 3D White. Ensure you use genuine Oral-B replacement brush heads.



Only compatible with genuine  
**Oral-B** replacement  
brush heads

*Figure 6: Replace brush head every 3 months.*



Figure 7: Compatible with genuine Oral-B replacement brush heads.





Figure 8: Choose your perfect clean from various Oral-B brush heads.

## SPECIFICATIONS

| Feature | Detail |
|---------|--------|
|---------|--------|

|                        |   |
|------------------------|---|
| Brand                  | Oral-B  |
| Model Number           | P_1651  |
| Power Source           | Battery Powered (2 AA batteries included)       |
| Operation Mode         | Automatic                                       |
| Number of Settings     | 1   |
| Toothbrush Head Shape  | Round (CrossAction)                             |
| Revolutions per Minute | 9,000 (for superior plaque removal)             |
| Timer                  | 2-minute coaching timer                         |
| Included Components    | Toothbrush handle, 1 brush head, 2 AA batteries |
| Product Dimensions     | 2.08 x 1.77 x 9.37 inches                       |
| Item Weight            | 0.18 Kilograms (6.35 ounces)                    |
| Color                  | Black   |
| Handle Material        | Plastic   |
| Bristle Material       | Nylon   |

## TROUBLESHOOTING

### 1. Toothbrush not turning on or weak performance:

- Check if the batteries are correctly inserted with the proper polarity.
- Replace the two AA batteries with new ones. The included batteries are designed for consistent performance, but will eventually need replacement.

### 2. Brush head feels loose or ineffective:

- Ensure the brush head is firmly snapped onto the handle.
- If the brush head is worn (typically after 3 months of use), replace it with a new genuine Oral-B brush head.
- If the toothbrush was dropped, replace the brush head immediately, even if no visible damage is present.

## ADDITIONAL RESOURCES

For visual guidance on using and maintaining your Oral-B Pro 100 toothbrush, please refer to the following videos:

### Oral-B Battery Operated Toothbrush Review - Watch First

Your browser does not support the video tag.

*Description: A detailed review of the Oral-B Battery Operated Toothbrush, covering its features and performance.*

### Unboxing the EASY to Assemble OralB Pro100

Your browser does not support the video tag.

*Description: This video demonstrates the unboxing and easy assembly process of the Oral-B Pro 100 toothbrush.*

### **Why I bought the Oral-B Pro 100 Battery Powered Electric Toothbrush**

Your browser does not support the video tag.

*Description: A user shares their reasons for choosing the Oral-B Pro 100 Battery Powered Electric Toothbrush.*

### **Oral-B Pro 100 3D White, Battery Toothbrush Review**

Your browser does not support the video tag.

*Description: A review focusing on the Oral-B Pro 100 with a 3D White brush head.*

### **Watch for honest review**

Your browser does not support the video tag.

*Description: An honest review of the Oral-B Pro 100 Battery Powered Toothbrush.*

### **What you need to know about this Tooth Brush...**

Your browser does not support the video tag.

*Description: An overview of key information regarding the Oral-B Pro 100 Toothbrush.*

### **Is it good quality?**

Your browser does not support the video tag.

*Description: A video assessing the quality of the Oral-B Pro 100 Battery Powered Toothbrush.*

### **I like this tooth brush**

Your browser does not support the video tag.

*Description: A user's positive feedback and demonstration of the Oral-B Pro 100 toothbrush.*

### **Watch Before Buy Oral-B Oral B Pro 100 3D Battery Toothbrush**

Your browser does not support the video tag.

*Description: A quick overview of the Oral-B Pro 100 3D Battery Toothbrush for potential buyers.*

## WARRANTY AND SUPPORT

---

For warranty information or customer support, please refer to the official Oral-B website or contact their customer service directly. Details can typically be found on the product packaging or through the manufacturer's official channels.

