

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [RAEMAO](#) /

› [RAEMAO Massage Gun User Manual](#)

RAEMAO x6

RAEMAO Massage Gun User Manual

Model: x6 | Brand: RAEMAO

1. INTRODUCTION

The RAEMAO Massage Gun is designed to provide targeted deep tissue massage, helping to alleviate muscle stiffness, pain, and tension. Its powerful percussion action promotes blood circulation and aids in muscle recovery, making it ideal for athletes, individuals with office fatigue, or anyone seeking muscle relief. This manual provides essential information for safe and effective use of your device.

2. SAFETY INFORMATION

To reduce the risk of fire, electric shock, or personal injury, please follow these instructions carefully:

- For adult use only.
- Use only on a dry, clean surface of the body or clothing by lightly pressing and moving across the body for approximately 60 seconds per region.
- Use the massage gun only on the body's soft tissue as desired without producing pain or discomfort.
- Do not use on head or any hard or bone area of the body.
- Use only the attachments provided.
- Discontinue use if you experience pain, numbness, or tingling.
- Keep fingers, hair, or any other body part away from the shaft and backside of the applicator as pinching may occur.
- Do not place any objects into the ventilation holes of the massage gun.
- Do not immerse in water or allow water to enter the ventilation holes of the massage gun.
- Do not drop or misuse the device.
- Do not attempt to disassemble the device.
- Do not operate for longer than 10 minutes continuously to prevent overheating and preserve battery life.
- Consult a physician before using if you have any medical conditions.

3. FEATURES

- **Pain Relief:** Targets muscle stiffness, pain, and tension with precision.
- **Long Battery Life:** Up to 5 hours of operation (depending on speed), with a 10-minute auto-off function for safety.
- **8 Replaceable Massage Heads:** Designed to target specific areas of the body for various muscle massager needs.
- **10 Speed Settings:** Covers a wide range of muscle soreness problems.
- **Quiet Operation:** Equipped with an advanced brushless motor, operating at a quiet 40dB noise level.
- **Smart LED Screen:** Intuitive touch-sensitive LED screen to easily adjust speed and view battery percentage.
- **Ergonomic Design:** Reduces stress on wrists, hands, and arms.
- **Portable:** Lightweight design with a carrying case for easy storage and transportation.

4. PACKAGE CONTENTS

- RAEMAO Massage Gun (x1)
- Interchangeable Massage Heads (x8)
- USB Charging Cable (x1)
- Carrying Protective Case (x1)
- User Manual (x1)

5. PRODUCT OVERVIEW & FUNCTIONS



Figure 1: RAEMAO Massage Gun with included accessories.

LCD TOUCH SCREEN DISPLAY

Adjust the gear according to different degrees of muscle fatigue

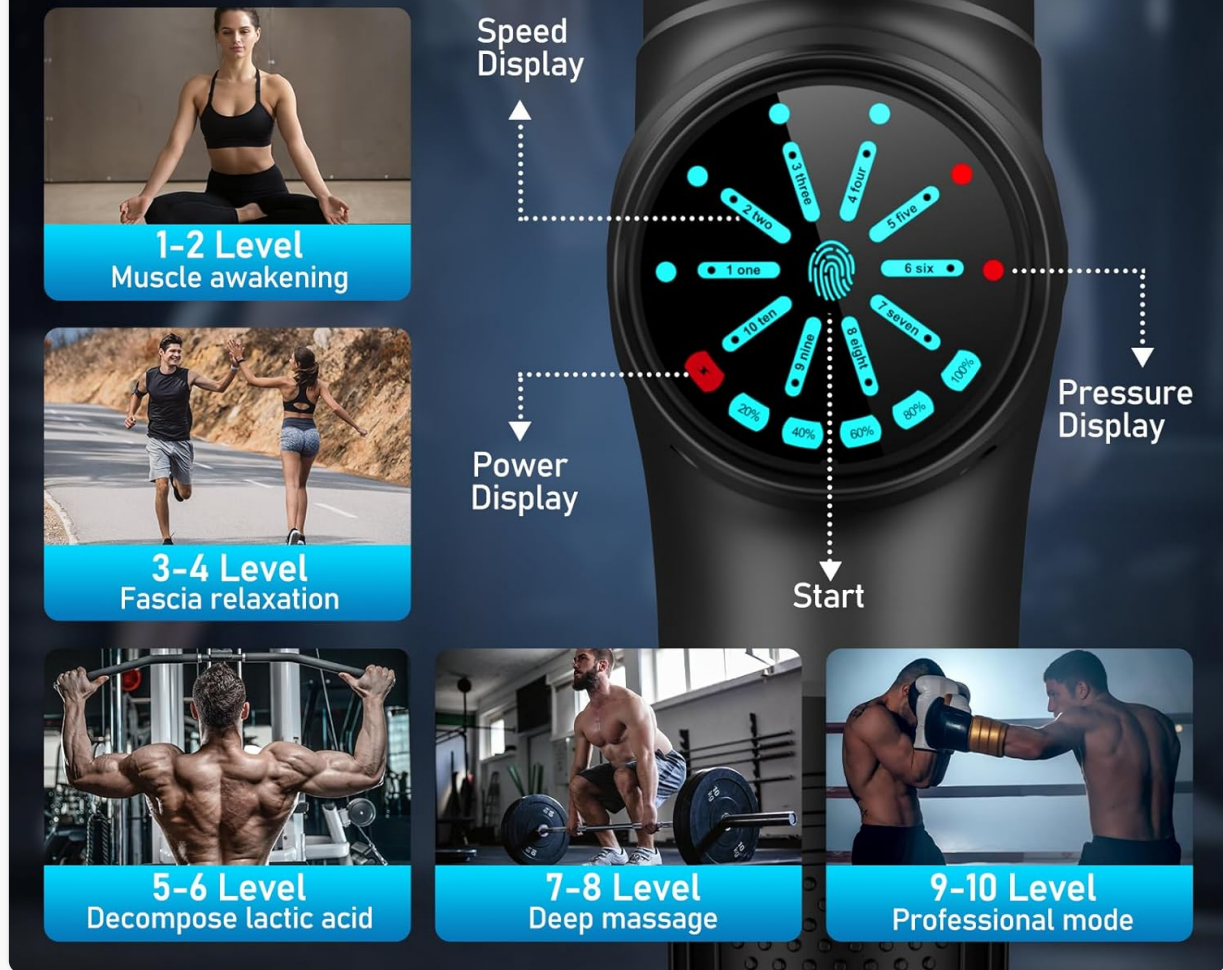


Figure 2: Intuitive LED touch screen for control.

The massage gun features an intuitive LED touch screen display on the top for speed and battery indication. The handle provides an ergonomic grip for comfortable use. The charging port is located at the bottom of the handle.

6. MASSAGE HEADS

The RAEMAO Massage Gun comes with 8 specialized massage heads, each designed for optimal use on different muscle groups and body parts. Refer to the diagram below for specific applications:

8 MASSAGE HEAD



Round Head
Suitable for thighs,
calves, back



Flat Head
Suitable for
all body parts



Wedge Head
Suitable for larger
muscles around Legs,
Chest, and Back



Fork Head
Suitable for neck,
spine and
achilles tendon



Convex Head
Suitable for quads,
gluteal



Y-Shape Head
Suitable for arm
and leg muscles



Thumb Head
Suitable for palm,
foot bottom



Bulge Head
Suitable for
abdominal muscles

Figure 3: Overview of the 8 interchangeable massage heads.

- **Round Head:** Suitable for thighs, calves, and back.
- **Flat Head:** Suitable for all body parts.
- **Wedge Head:** Suitable for larger muscles around legs, chest, and back.
- **Fork Head:** Suitable for neck, spine, and Achilles tendon.
- **Convex Head:** Suitable for quads and gluteal muscles.
- **Y-Shape Head:** Suitable for arm and leg muscles.
- **Thumb Head:** Suitable for palm and foot bottom.
- **Bulge Head:** Suitable for abdominal muscles.

To attach a massage head, align it with the port on the gun and firmly push it in. To detach, gently twist and pull the head out.

7. SETUP & OPERATING INSTRUCTIONS

Before first use, fully charge the massage gun. Connect the USB charging cable to the charging port at the bottom of the handle. The LED display will indicate charging status and battery level.

LONG BATTERY LIFE



Figure 4: Charging port and battery display.

To power on the device, slide the ON/OFF switch at the bottom of the handle to the 'ON' position. The LED screen will light up. Press the fingerprint icon on the screen to start the massage. To adjust the speed, tap the fingerprint icon to cycle through the 10 available speed levels. The LED display will show the current speed setting.

12MM AMPLITUDE

Activate the deep muscle groups, quickly relax full body



Figure 5: Using the massage gun for muscle relief.

The massage gun features a 10-minute auto-off function to prevent overheating and preserve battery life. After 10 minutes of continuous use, the device will automatically power off. To continue use, simply turn it back on.

8. CHARGING

Connect the USB charging cable (provided) to the charging port at the bottom of the massage gun's handle. The LED display on the screen will show the battery percentage, indicating when it's fully charged (100%). The battery can be recharged at any time, regardless of the current battery level. For safety, always ensure the machine is turned off when charging.

9. MAINTENANCE

- Use a slightly damp towel and a dry cloth to gently wipe the outer surface clean. Be sure to turn the massage gun off before cleaning.
- Do not attempt to clean internal surfaces; disassembling the device could result in device failure.
- Please turn off the power of the device and store it in the supplied case if you don't use it for a long time or for travel.

- Keep the device away from heat sources, liquids, and children.

10. TROUBLESHOOTING

Problem	Solution
Device does not turn on.	Ensure the ON/OFF switch at the bottom is in the 'ON' position. Check battery level on the LED screen and charge if necessary.
Massage head is not working properly.	Ensure the massage head is securely attached. Check for any obstructions around the shaft.
Device stops unexpectedly during use.	The device has a 10-minute auto-off safety feature. Turn it off and then back on to resume use. Check battery level.
Unusual noise or vibration.	Turn off the device immediately. Check for any foreign objects or damage to the massage head or shaft. Contact customer support if the issue persists.

11. SPECIFICATIONS

- **Model:** x6
- **Amplitude:** 12mm
- **Massage Heads:** 8 (interchangeable)
- **Working Noise:** 32-40dB
- **Charging Time:** About 3-4 hours
- **Motor Run Speed:** 1200-3200 RPM
- **Charging Voltage:** 100V-240V
- **Running Time:** Up to 5 hours (depending on speed level)
- **Battery:** 2500mAh Lithium Ion (included)
- **Adjustable Levels:** 10 speeds
- **Weight:** 14.11 ounces (400 Grams)
- **Dimensions:** 10.43 x 7.87 x 3.54 inches
- **Material:** Plastic

12. WARRANTY & SUPPORT

RAEMAO offers a 1-year warranty from the date of purchase against defects in material or workmanship. If within one year from the date of purchase, the product fails due to a defect, RAEMAO will repair or replace the product at no charge. This warranty does not cover damage caused by accident, abuse, improper handling, or unauthorized repairs.

For customer support, please visit the [RAEMAO Store on Amazon](#) or refer to the contact information provided in your product packaging.

13. OFFICIAL PRODUCT VIDEO

Your browser does not support the video tag.

This video provides a comprehensive overview of the RAEMAO Massage Gun's features and benefits, demonstrating its

quiet operation, adjustable speeds, and various massage heads in action.

© 2023 RAEMAO. All rights reserved.