

[Manuals.plus](#) /

› [Real Relax](#) /

› Real Relax Massage Chair Favor-06 Instruction Manual

Real Relax Favor-06

Real Relax Massage Chair Favor-06 Instruction Manual

Model: Favor-06 (BH-6201)

1. INTRODUCTION

Thank you for choosing the Real Relax Favor-06 Massage Chair. This manual provides essential information for the safe and effective operation, maintenance, and troubleshooting of your new massage chair. Please read this manual thoroughly before use and retain it for future reference.

2. SETUP AND ASSEMBLY

The Real Relax Favor-06 massage chair is designed for straightforward assembly. The main unit comes largely pre-assembled, requiring minimal steps to get started.

2.1 Unpacking and Placement

Ensure you have adequate space for the chair. The minimum doorway size required for entry is 29 inches. The product dimensions are 57.48"D x 28.34"W x 42.12"H.



Image 2.1.1: Minimum doorway size for the Real Relax Favor-06 Massage Chair. The chair requires a minimum doorway width of 29 inches for easy placement.

2026 Upgraded Favor 06



Phone Holder



Massage While Tiktok & Reading



Blood Oxygen/Heart Rate Measurement



Seat Back Pad

Image 2.1.2: Product dimensions of the Real Relax Favor-06 Massage Chair. The chair measures 57.48 inches deep, 28.34 inches wide, and 42.12 inches high.

2.2 Footrest Installation

The footrest is typically folded for shipping. To install:

1. Place the footrest on the floor.
2. Connect the power cord to reset the chair.
3. Connect the air tube and the wires from the footrest to the main unit.
4. Slide two white sleeves onto connector pins on both sides of the footrest.
5. Hang the footrest on the main unit.
6. Fold down the plastic latch plates and fasten them with screws.

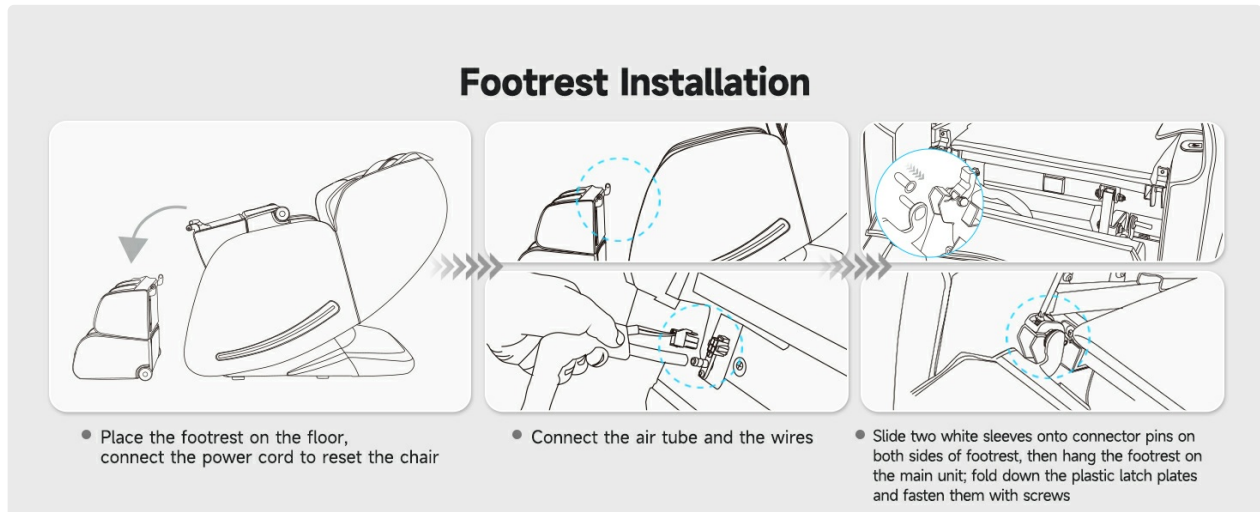


Image 2.2.1: Footrest installation guide. This image shows the step-by-step process of attaching the footrest to the main massage chair unit.

2.3 Power Connection

Once the footrest is securely attached, connect the main power cable to a suitable electrical outlet. Ensure the power switch on the chair is in the 'On' position.

3. OPERATING INSTRUCTIONS

Your Real Relax Favor-06 massage chair offers a variety of features and customizable settings for a personalized massage experience.

3.1 Control Options

The chair can be controlled via the upgraded white LCD touch screen, the dedicated phone controller, or the shortcut keys located on the armrest.

Product Dimension

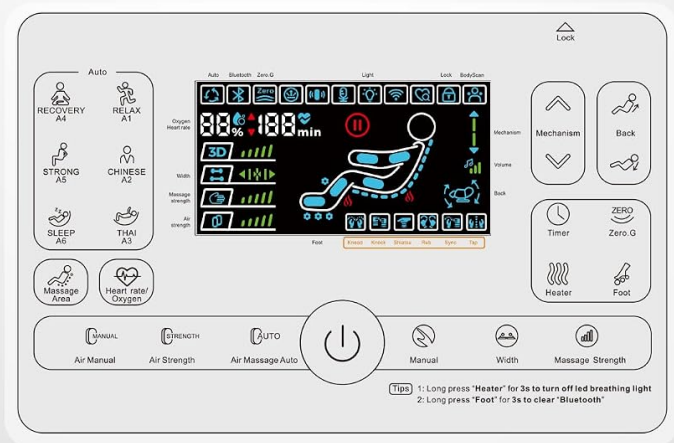


Image 3.1.1: Upgraded Phone Controller and Phone Holder. The new phone controller and integrated phone holder enhance user convenience during massage sessions.

APP Control

Control all functions of the massage chair through the app



Image 3.1.2: APP Control. Control all functions of the massage chair conveniently through the dedicated mobile application.


3.2 Massage Features

- **3D SL Track:** The advanced 3D SL track system massages your body from the neck down to the legs. You can choose fixed massage points, adjust the robot roller massage width to 3 levels, and select from 5 levels of adjustable intensity.

3D SL Track Massage Chair



5 levels of
adjustable intensity

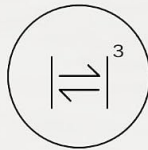

3 Levels of
width adjustment



Message Area
Customizable



3D Robot
Hand

Image 3.2.1: 3D SL Track Massage Chair. This diagram shows the extensive massage path covering neck, shoulders, back, waist, and hips, with adjustable intensity and width.

- **Zero Gravity:** Experience weightlessness with the one-button zero gravity design. This position elevates your legs above your heart, reducing pressure on your spine and enhancing the massage effect.

Zero Gravity Experience



“▶||” Want a break? Just **PAUSE**, use it as a normal recliner and chill

Image 3.2.2: Zero Gravity Experience. The chair reclines to a zero-gravity position, distributing body weight evenly and reducing spinal pressure for ultimate relaxation.

- **Full-Body Airbag Massage:** Enjoy a comprehensive massage with full-body airbags, offering 5 levels of adjustable pressure. Airbags target shoulders, arms, calves, and feet.

Full Body Massage

Full body airbag massage & foot roller massage



Image 3.2.3: Full Body Airbag Massage. Airbags provide compression massage to shoulders, arms, calves, and feet, with 3 adjustable intensity levels for foot rollers.

- **Foot Roller Massage:** Rollers located in the footrest provide a soothing massage to your soles, complementing the airbag compression.



Image 3.2.4: Foot Roller Massage. Rollers provide targeted massage to the soles of the feet, enhancing circulation and relieving tension.

- **Heating Function:** The chair includes a heating function for added comfort and muscle relaxation. A removable back pad is included to adjust massage intensity.



Image 3.2.5: Massage Chair in a bedroom setting. The chair features LED breath lights that can be turned off, adding to the ambiance.

- **Bluetooth & Speakers:** Connect your device via Bluetooth to enjoy music or podcasts through the integrated speakers during your massage.



Image 3.2.6: Bluetooth Speaker. Enjoy your favorite audio content through the chair's integrated Bluetooth speakers.

- **Heart & Blood Oxygen Measurement:** The chair features a sensor for heart rate and blood oxygen measurement. Note: This function will reset the massage chair while massaging. Test results are for reference only.

Blood Oxygen & Heart Rate Measurement



Image 3.2.7: Blood Oxygen & Heart Rate Measurement. The chair can measure your heart rate and blood oxygen levels, providing data for reference.

- **Child Lock:** For safety, the chair includes a child lock feature to prevent unintended operation.
- **Pause Button:** A convenient pause button allows you to temporarily stop the massage and use the chair as a normal recliner.

4. MAINTENANCE

To ensure the longevity and optimal performance of your massage chair, follow these maintenance guidelines:

- **Cleaning:** Spot clean the faux leather material as needed. Use a soft, damp cloth with mild soap for cleaning. Avoid harsh chemicals or abrasive cleaners.
- **General Care:** Keep the chair away from direct sunlight and heat sources to prevent material damage. Regularly check all connections and components for any signs of wear or looseness.

5. TROUBLESHOOTING

If you encounter any issues with your Real Relax Favor-06 massage chair, please refer to the following common troubleshooting tips:

- **Chair not powering on:** Ensure the power cord is securely plugged into both the chair and the wall outlet. Check that the power switch on the chair is in the 'On' position. Verify the wall outlet is functional.
- **Massage rollers not moving:** Ensure the chair is not in 'Pause' mode. Check if any objects are obstructing the rollers. Restart the chair.
- **Airbags not inflating:** Verify that all air tubes are securely connected, especially after footrest installation. Check the air strength settings on the controller.
- **Unusual noises:** Discontinue use and inspect the chair for any loose parts or obstructions. Contact customer support if the noise persists.

For further assistance, please contact Real Relax customer support.

6. SPECIFICATIONS

Feature	Detail
Brand	Real Relax
Model Name	Favor-06
Model Number	BH-6201
Color	Blue
Material	Faux Leather, Cotton (Fill Material)
Product Dimensions (D x W x H)	57.48" x 28.34" x 42.12"
Item Weight	174 Pounds
Maximum Weight Recommendation	400 Pounds
Height Limit	6.1 feet
Additional Features	3D SL Track, Child Lock, Phone Controller, Phone Holder, Seat Back Pad, Zero Gravity, App Control, Bluetooth, Foot Roller Massage, LED Breath Light, Heart & Blood Oxygen Measurement, Pause Button
Power	Electric

7. WARRANTY AND SUPPORT

Your Real Relax Favor-06 Massage Chair comes with a warranty card included in the accessory box. Please refer to the warranty card for detailed terms and conditions.

For any questions, technical support, or warranty claims, Real Relax offers 24-hour online chat support on workdays. Please feel free to contact us for assistance.