

Oclean Oclean Flow

Oclean Flow Sonic Electric Toothbrush User Manual

Model: Oclean Flow

Brand: Oclean

INTRODUCTION

This user manual provides comprehensive instructions for the Oclean Flow Sonic Electric Toothbrush. It covers everything from initial setup and daily operation to maintenance and troubleshooting, ensuring you get the most out of your device for superior oral hygiene. Please read this manual carefully before using the product.

WHAT'S INCLUDED

Your Oclean Flow Sonic Electric Toothbrush package includes the following items:

- 1 x Oclean Flow Electric Toothbrush Handle
- 7 x DuPont Replacement Brush Heads
- 1 x Type-C Charging Cable
- 1 x User Manual



Image: The Oclean Flow Electric Toothbrush, showing the main handle, a pack of 6 replacement brush heads, and an additional brush head attached to the toothbrush, along with the product packaging.

SETUP

Initial Charging

Before first use, fully charge your Oclean Flow toothbrush. The toothbrush boasts an impressive 180-day battery life on a full charge, eliminating the need for frequent recharging.

- Locate the Type-C charging port at the bottom of the toothbrush handle.
- Connect the provided Type-C charging cable to the toothbrush and a compatible USB power source (e.g., wall adapter, computer USB port).

- The indicator light on the handle will illuminate to show charging status. A full charge typically takes a few hours.
- **Note:** While the toothbrush uses a USB-C port, it is designed to charge optimally with the included USB-A to USB-C cable. Some USB-C to USB-C cables or chargers may not be compatible.

SUPER LONG BATTERY LIFE

180 days battery life

1 = **180**
full charge days uses



Image: A close-up view of the Oclean Flow toothbrush's USB-C charging port with the cable inserted, illustrating its long battery life.



Image: A visual representation emphasizing the 180-day battery life of the Oclean Flow toothbrush from a single full charge.

Attaching the Brush Head

Align the brush head with the metal shaft on the toothbrush handle and push it down firmly until it clicks into place.

OPERATING

Brushing Technique

For optimal cleaning, wet the brush head and apply a small amount of toothpaste. Place the brush head against your teeth at a 45-degree angle to the gum line. Move the brush in a gentle circular motion, allowing the sonic vibrations to clean your teeth effectively. Brush for two minutes, ensuring all surfaces of your teeth, gums, and tongue are cleaned.

POWERFUL SONIC CLEANING

UP TO
76,000

Movements Per Minute

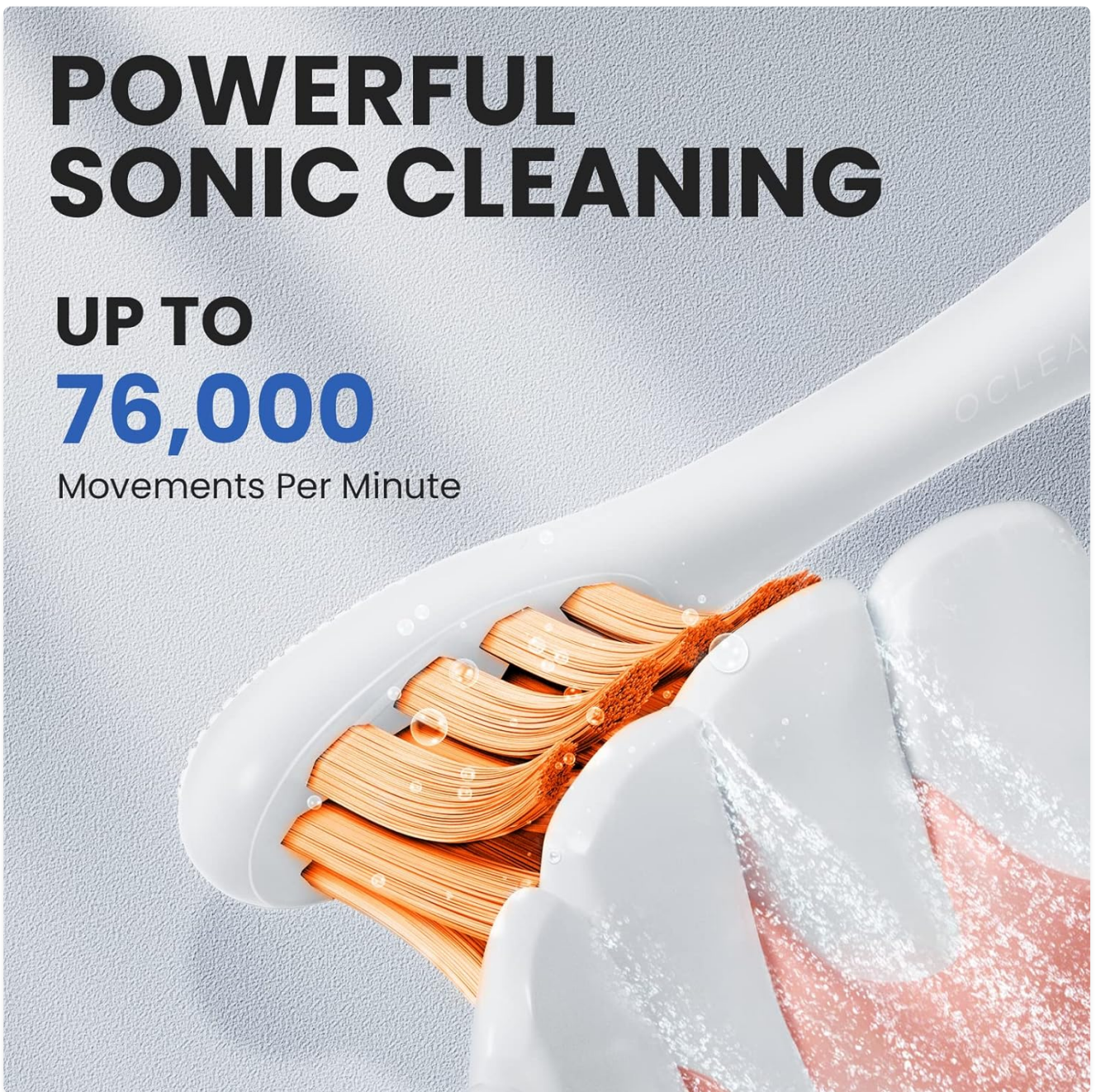


Image: A visual demonstrating the powerful sonic cleaning action of the Oclean Flow toothbrush, highlighting its 76,000 movements per minute.

Brushing Modes

The Oclean Flow offers 5 high-performance brushing modes to cater to various oral care needs. Press the power button to cycle through the modes:

- **Morning Mode:** For a superior cleaning experience.
- **Night Mode:** Designed for deep and thorough cleaning, effectively removing food residue.
- **Standard Cleaning Mode:** Ideal for daily general cleaning.
- **Whitening Mode:** Helps effectively remove stains and restore natural tooth whiteness.
- **Gentle Mode:** Recommended for individuals with sensitive gums or those new to electric toothbrushes.

5 BRUSHING MODES



Morning

For a superior cleaning

Brushing intensity: ★★★



Night

For a deep and thorough cleaning

Brushing intensity: ★★★★★



Standard Cleaning

For daily general cleaning

Brushing intensity: ★★★★★



Whitening

For a naturally whiter smile

Brushing intensity: ★★★★★



Gentle

For those who have sensitive gums

Brushing intensity: ★★



Image: An illustration of the five distinct brushing modes available on the Oclean Flow toothbrush, explaining their purpose and intensity.

Smart Timer and Zone Reminder

The Oclean Flow features a built-in 2-minute smart timer, which is the dentist-recommended brushing time.

Additionally, a 30-second zone switch reminder prompts you to move to a different quadrant of your mouth, ensuring a thorough and even clean.

2 MINS SMART TIMER



Image: A user demonstrating the Oclean Flow toothbrush, with visual aids explaining the 2-minute smart timer and 30-second zone change reminder for effective brushing.

MAINTENANCE

Cleaning the Toothbrush

The Oclean Flow Sonic Electric Toothbrush is IPX7-rated waterproof, allowing for easy cleaning under running water. After each use:

- Rinse the brush head thoroughly under running water.
- Rinse the toothbrush handle to remove any toothpaste residue.
- Wipe the handle dry with a soft cloth.

IPX7

WATERPROOF



Image: The Oclean Flow toothbrush fully submerged in water, illustrating its IPX7 waterproof capability for easy cleaning.

Replacing Brush Heads

It is recommended to replace your brush head every 3 months, or sooner if the bristles become frayed or worn. To replace:

- Pull the old brush head straight off the handle.
- Align the new brush head with the metal shaft and push it down firmly until it clicks into place.

TROUBLESHOOTING

Toothbrush Not Charging

If your toothbrush is not charging, please check the following:

- Ensure the charging cable is securely connected to both the toothbrush and the power source.
- Verify that the power source is active and functioning correctly.
- **Important:** The Oclean Flow toothbrush is designed to charge with the included USB-A to USB-C cable. While it has a USB-C port, some third-party USB-C to USB-C cables or power adapters may not be compatible. If you are using a different cable, try the one provided with the product.

Reduced Vibration or Power

If you notice a decrease in vibration intensity or power, it may indicate a low battery. Recharge the toothbrush fully as per the 'Initial Charging' instructions.

Brush Head Not Staying On

Ensure the brush head is pushed down firmly onto the metal shaft until you hear or feel a click. If the issue persists, the brush head or handle may be damaged. Contact customer support for assistance.

SPECIFICATIONS

Feature	Description
Brand	Oclean
Model	Flow
Power Source	Battery Powered
Battery Life	Up to 180 days on a full charge
Sonic Movements	Up to 76,000 movements per minute
Brushing Modes	5 (Morning, Night, Standard, Whitening, Gentle)
Waterproof Rating	IPX7
Timer	2-minute smart timer with 30-second zone reminder
Bristle Type	High-Density DuPont Bristles (Soft)
Charging Port	Type-C
Package Dimensions	11.14 x 6.57 x 3.74 inches
Item Weight	10.23 ounces

SAFETY INFORMATION

For adult use only. This toothbrush features soft bristles designed for gentle cleaning. Its IPX7 waterproof rating

ensures safe use in wet environments, such as the shower. Always ensure the charging port cover is securely closed when not charging to maintain water resistance.

WARRANTY AND SUPPORT

Oclean is committed to improving oral health through technology and innovation. Your Oclean Flow Sonic Electric Toothbrush comes with a **2-year service warranty**. This includes hassle-free replacement or return for eligible issues.

For any questions, concerns, or support needs, please contact Oclean customer service. We offer 7*24h customer support to assist you.



Image: A visual representation of Oclean's customer support, highlighting their commitment to customer care and warranty services.

You can also refer to the official Oclean store for more information and products:[Visit the Oclean Store](#)

For a downloadable version of the user manual, please click here:[Download User Manual \(PDF\)](#)