

Medivon AQUA_PRO

Medivon Aqua Pro Foldable Water Foot Massager Instruction Manual

Model: AQUA_PRO

1. INTRODUCTION

Thank you for choosing the Medivon Aqua Pro Foldable Water Foot Massager. This device is designed to provide a relaxing and therapeutic foot spa experience in the comfort of your home. It features bubble massage, acupressure rollers, infrared therapy, ozone sterilization, and temperature control for an optimal experience. Please read this manual carefully before use to ensure safe and proper operation, and retain it for future reference.

2. SAFETY INSTRUCTIONS

- Always unplug the device from the power outlet immediately after use and before cleaning.
- Do not use while bathing or showering.
- Do not place or store the device where it can fall or be pulled into a tub or sink.
- Do not immerse the device in water or any other liquid.
- Never reach for an appliance that has fallen into water. Unplug it immediately.
- Close supervision is necessary when this appliance is used by, on, or near children or individuals with disabilities.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Contact customer service for examination and repair.
- Keep the cord away from heated surfaces.
- Never drop or insert any object into any opening.
- Do not use outdoors.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.

- To disconnect, turn all controls to the 'OFF' position, then remove the plug from the outlet.
- This appliance is designed for household use only.
- Consult your doctor before using this massager if you have any medical conditions, including but not limited to diabetes, heart disease, pregnancy, or if you have a pacemaker.
- Do not use on open wounds, swollen areas, or inflamed skin.

3. PACKAGE CONTENTS

Please check the package for the following items:

- Medivon Aqua Pro Foldable Water Foot Massager Unit
- Power Cord
- Instruction Manual

4. SETUP

1. **Unfold the Unit:** Place the folded foot massager on a flat, stable surface. Gently pull up the sides of the unit until it fully expands and locks into place.
2. **Fill with Water:** Fill the foot massager with warm or cold water, depending on your preference, up to the indicated maximum fill line. Do not overfill.
3. **Connect Power:** Ensure your hands are dry. Plug the power cord into a standard electrical outlet.



Image 4.1: The Medivon Aqua Pro Foot Massager in its folded, compact form, ready for storage or initial setup.

5. OPERATING INSTRUCTIONS

Once the unit is set up and filled with water, you can begin your foot massage session.

1. **Power On:** Press the 'POWER' button on the control panel to turn on the device. The display will illuminate.
2. **Temperature Adjustment:** Use the 'TEMP +' and 'TEMP -' buttons to adjust the water temperature. The desired temperature range is typically between 35°C and 48°C. The device will actively heat and maintain the selected temperature.
3. **Mode Selection:** Press the 'MODE' button to cycle through different massage modes, which may include bubble massage, vibration, or a combination.
4. **Infrared Light:** Press the 'LIGHT' button to activate the infrared therapy function.
5. **Ozone Sterilization:** Press the 'OZONE' button to activate the ozone sterilization function. This helps maintain hygiene.
6. **Timer:** Some models may include a timer function. Refer to the control panel for specific timer options.

7. **Enjoy the Massage:** Place your feet into the water and relax. The integrated massage rollers will provide acupressure.
8. **Power Off:** After your session, press the 'POWER' button again to turn off the device. Unplug from the power outlet.



Image 5.1: An overhead view of the foot massager, showcasing the six massage heads and the intuitive control panel for various functions.

ROTLICHT-THERAPIE


RED-LIGHT THERAPY



Image 5.2: Feet immersed in the foot massager, demonstrating the red-light therapy feature in action.

6. FEATURES EXPLAINED

- **Bubble Massage:** Gentle bubbles are generated to stimulate the feet, promoting relaxation and improving circulation.
- **Acupressure Rollers:** Integrated rollers provide targeted pressure to key acupressure points on the soles of your feet, relieving tension and fatigue.
- **Infrared Therapy:** The red-light therapy function emits infrared light, which is believed to help soothe muscle aches, reduce inflammation, and promote healing.
- **Ozone Sterilization:** This feature uses ozone to sterilize the water, ensuring a hygienic environment for your foot spa and making it suitable for multiple users within a family.
- **Temperature Control (35-48°C):** The device actively heats and maintains the water at your desired temperature, providing consistent warmth throughout your massage session.

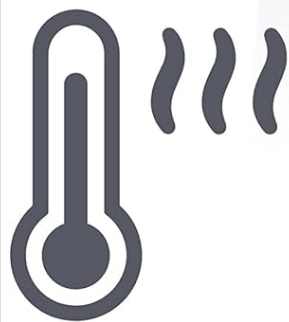

AKTIVE ERWÄRMUNG ACTIVE HEATING

Image 6.1: The side of the foot massager, highlighting its active heating capability to maintain water temperature.

7. MAINTENANCE AND CLEANING

1. **Drain Water:** After each use, unplug the device. Carefully drain the water using the integrated drainage spout.
2. **Clean the Unit:** Wipe the interior and exterior surfaces with a soft, damp cloth. Do not use abrasive cleaners or harsh chemicals.
3. **Dry Thoroughly:** Ensure the unit is completely dry before folding and storing.



Image 7.1: The foot massager in use and demonstrating the water drainage process.

8. STORAGE

The Medivon Aqua Pro is designed for convenient storage due to its foldable nature.

1. **Ensure Dryness:** Make sure the unit is completely dry before folding.
2. **Fold the Unit:** Gently press down on the sides to collapse the unit into its compact form.
3. **Store:** Store the folded unit in a cool, dry place, away from direct sunlight and out of reach of children.

EINFACHE LAGERUNG

EASY

STORAGE



Image 8.1: The foot massager in its folded state, easily stored in a drawer, highlighting its compact design.

9. TROUBLESHOOTING

- **Device does not power on:**
 - Ensure the power cord is securely plugged into both the device and a working electrical outlet.
 - Check if the power outlet is functional by plugging in another appliance.
- **No bubbles or massage function:**
 - Ensure the device is powered on and the specific massage mode is selected.
 - Check if the water level is within the recommended range.
- **Water not heating or temperature inconsistent:**
 - Verify that the temperature setting is correctly adjusted using the 'TEMP +' and 'TEMP -' buttons.
 - Allow a few minutes for the water to reach the desired temperature.

10. SPECIFICATIONS

Brand	Medivon
-------	---------

Model Number	AQUA_PRO
Dimensions (L x W x H)	20 x 15 x 30 cm (unfolded)
Weight	2 Kilograms
Power Source	Mains Powered (Not battery operated as per product details)
Temperature Range	35-48 °C
Special Features	Temperature Maintenance, Bubble Massage, Acupressure, Ozone Sterilization, Infrared Therapy

11. WARRANTY AND SUPPORT

For warranty information, technical support, or service inquiries, please refer to the warranty card included with your product or contact Medivon customer service directly. Keep your proof of purchase for any warranty claims.