

MLAY T14

MLAY T14 IPL Pulsed Light Hair Removal Device User Manual

Model: T14

1. INTRODUCTION

This manual provides essential information for the safe and effective use of your MLAY T14 IPL Pulsed Light Hair Removal Device. Please read it thoroughly before first use and keep it for future reference.

The MLAY T14 utilizes advanced IPL (Intense Pulsed Light) technology to target hair follicles, inhibiting hair growth. It features rapid flash speed, a high number of flashes, and an integrated ice cooling function for a comfortable experience.

2. SAFETY INFORMATION

Important Safety Precautions:

- Always wear the provided safety glasses during operation to protect your eyes from intense light.
- Do not use on tattooed areas, moles, or very dark skin tones (Fitzpatrick skin types V-VI).
- Avoid direct exposure to the eyes.
- Keep out of reach of children.
- Consult a doctor before use if you have any skin conditions, medical implants, or are pregnant/breastfeeding.
- Ensure the device is clean and dry before each use.
- Do not use on irritated, sunburned, or recently waxed skin.
- Perform a patch test on a small skin area 24 hours before full treatment to check for adverse reactions.

3. PACKAGE CONTENTS

Upon unpacking, please verify that all items listed below are present and in good condition:

- 1 x MLAY T14 IPL Main Unit (equipped with HR lamp head)
- 1 x 500,000 Flash Hair Removal Lamp Head (3.9 cm²)
- 1 x Safety Glasses
- 1 x Lens Cleaning Swab
- 1 x User Manual
- 1 x Power Adapter (100V-240V)
- 1 x Warranty Card



Image: The MLAY T14 IPL device, its packaging box, and a pair of protective glasses.

4. PRODUCT OVERVIEW

The MLAY T14 is an advanced IPL device designed for effective and comfortable hair removal. Key features include:

- **Fast Flash Speed:** 0.5 seconds per flash for quick treatments.
- **High Flash Capacity:** 500,000 flashes, sufficient for full body use over an extended period.
- **Ice Cooling Technology:** Maintains a skin contact temperature of 0-3°C to minimize discomfort and protect the skin.
- **5 Energy Levels:** Adjustable intensity from 8.97-17J to suit different skin sensitivities and hair types.
- **Replaceable Lamp Heads:** Designed for durability and reduced waste.
- **Multi-functional:** Can be used for hair removal (HR), skin rejuvenation (SR), and acne treatment (AC) with appropriate lamp heads.

MLAY T14

Épilateurs à Lumière Pulsée

3°C
Refroidissement



500000
flashes



5 niveaux
d'énergie



Technologie
IPL indolore



Zone de
traitement: 3,9cm²



Max. 17J
d'énergie



Image: A woman holding the MLAY T14 device, with icons highlighting its features: 3°C cooling, 500,000 flashes, 5 energy levels, painless IPL technology, 3.9 cm² treatment area, and max 17J energy.

3°C

=

Confort Glacial
Épilation sans Stress



Notre appareil T14

Refroidissement instantané

Sensation agréable

Idéal pour les peaux sensibles



Autres Appareils

Chaleur excessive

Irritation de la peau

Non adapté aux peaux sensibles

Image: A visual comparison showing the MLAY T14's 3°C glacial comfort and stress-free epilation versus other devices that may cause excessive heat and skin irritation.

IPL pulsed light

Comfortable and painless

Say goodbye to repetition. Delay growth. Deep hair removal



Image: A cross-section diagram illustrating how IPL pulsed light penetrates the epidermis and dermis to target hair follicles, delaying hair growth and providing deep hair removal.

5. SETUP

1. **Unpack:** Carefully remove all components from the packaging.
2. **Inspect:** Check the device and accessories for any signs of damage. Do not use if damaged.
3. **Clean Skin:** Ensure the treatment area is clean, dry, and free from any lotions, oils, or deodorants. Shave the hair in the treatment area before use; do not wax or pluck.
4. **Connect Power:** Plug the power adapter into the device and then into a suitable electrical outlet (100V-240V).
5. **Wear Safety Glasses:** Always put on the provided safety glasses before turning on the device.

6. OPERATING INSTRUCTIONS

6.1 Power On/Off

- Press and hold the power button for 2 seconds to turn the device on or off.
- The device will enter standby mode after 5 minutes of inactivity.

6.2 Adjusting Energy Levels

The MLAY T14 offers 5 energy levels. Start with the lowest level (Level 1) for your first treatment and gradually increase as your skin tolerates it. Higher energy levels generally yield better results.

- Press the power button briefly to cycle through energy levels (1 to 5). The selected level will be displayed on the screen.



Image: The MLAY T14 device illustrating its 5 energy levels, providing personalized epilation programs for different skin and hair types, with a maximum intensity of 17J.

6.3 Flash Modes

The device supports two flash modes: Manual and Auto.

- **Manual Mode:** Press the flash button for each individual flash. Ideal for small, sensitive areas like the face, underarms, and bikini line.
- **Auto Mode:** Press and hold the flash button for 2 seconds to activate. The device will automatically flash when the treatment window is in full contact with the skin. Suitable for larger areas like legs, arms, and back.



Image: The MLAY T14 device demonstrating its two flash modes (Auto and Manual) and recommended treatment times for various body

parts: Bikini (1 min), Underarms (30s), Arms (1.5 min), Legs (1.5 min), Face (45s), Back (1.5 min).

6.4 Treatment Procedure

1. Shave the treatment area.
2. Plug in the device and turn it on.
3. Select the appropriate energy level.
4. Place the treatment window firmly against your skin. Ensure full contact.
5. Press the flash button (Manual Mode) or glide the device over the skin (Auto Mode).
6. Move to the next area, ensuring slight overlap to avoid missing spots.
7. After treatment, apply a soothing lotion if desired.



Votre Corps Doux et Lisse en 6 Minute Chrono !



1 minute
Maillot



30 s
Aisselles



1,5 minute
Bras



1,5 minute
Jambes



45 s
Visage



1,5 minute
Dos

Image: A graphic illustrating estimated treatment times for various body parts using the MLAY T14: Bikini (1 minute), Underarms (30 seconds), Arms (1.5 minutes), Legs (1.5 minutes), Face (45 seconds), Back (1.5 minutes).

17J d'Intensité pour une Peau Lisse en 4 Semaines

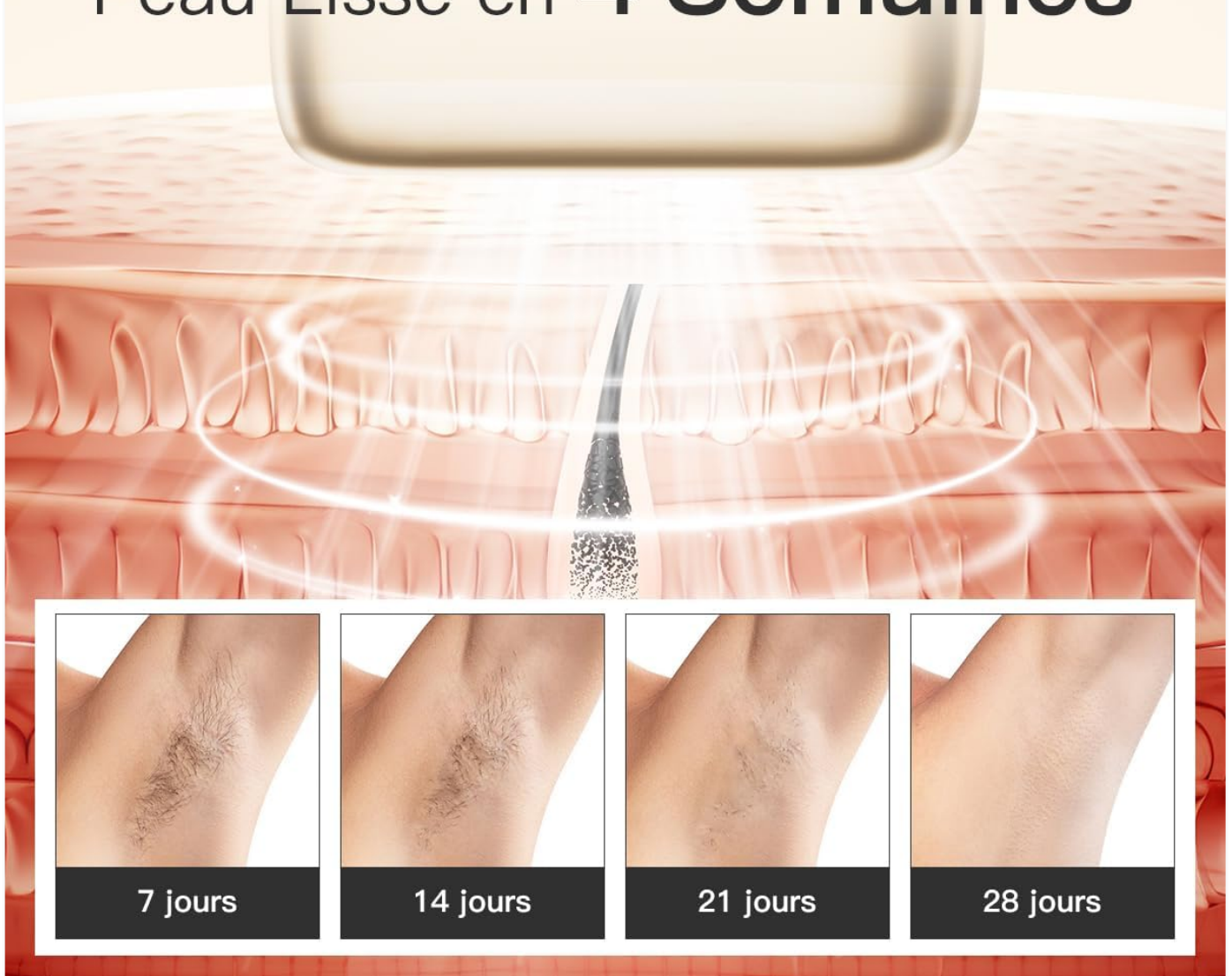


Image: A visual representation of hair reduction progress over 4 weeks, showing significant decrease in hair density after 7, 14, 21, and 28 days of consistent use.

6.5 Using Different Lamp Heads (If applicable)

The MLAY T14 supports different lamp heads for various functions:

- **HR (Hair Removal) Lamp Head:** For permanent hair reduction.
- **SR (Skin Rejuvenation) Lamp Head:** For improving skin texture and tone.
- **AC (Acne Treatment) Lamp Head:** For reducing acne and blemishes.

To change a lamp head, ensure the device is turned off and unplugged. Gently pull out the current lamp head and insert the desired one until it clicks into place.

Not Only A Hair Remove Machine But Also A Beauty Instrument



Remove Hair Without
Black Spot Left
500000 flash



Photo Rejuvenation
500000 flash



Acne Remove & Light Acne
500000 flash



Image: The MLAY T14 device demonstrating its versatility with different lamp heads for Hair Removal, Photo Rejuvenation, and Acne Removal, each with 500,000 flashes.

Tête de lampe
rajeunissante pour la peau

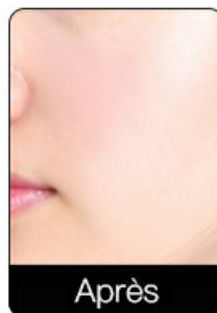
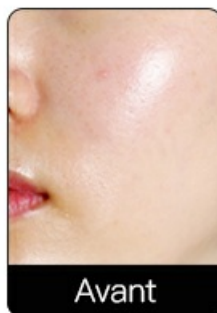


Image: Before and after comparisons for HR (Hair Removal), SR (Skin Rejuvenation), and AC (Acne Treatment) lamp heads, showing the potential results of each treatment.

7. MAINTENANCE

- **Cleaning:** After each use, ensure the device is turned off and unplugged. Use the provided lens cleaning swab or a soft, dry cloth to gently wipe the treatment window. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and moisture. Keep it in its original packaging or a protective case to prevent damage.

8. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Not properly plugged in; power button not pressed long enough.	Ensure power adapter is securely connected. Press and hold the power button for 2 seconds.
Device does not flash.	Treatment window not in full contact with skin; device in standby mode; lamp head expired.	Ensure full contact with skin. Press any button to wake from standby. Check remaining flashes on the display. Replace lamp head if necessary.
Skin feels too hot or uncomfortable.	Energy level too high for skin type; cooling function not active.	Reduce the energy level. Ensure the cooling function is active (indicated by snowflake icon). Move the device more quickly.
No visible hair reduction after several weeks.	Inconsistent use; energy level too low; hair color/skin type not suitable for IPL.	Ensure consistent treatment schedule. Increase energy level if comfortable. Consult the skin/hair type chart in the full manual (if available).

9. SPECIFICATIONS

Feature	Detail
Model	T14
Brand	MLAY
Technology	IPL (Intense Pulsed Light)
Flash Capacity	500,000 flashes
Flash Speed	0.5 seconds/flash
Energy Levels	5 levels (8.97-17J)
Cooling Technology	Ice Cooling (0-3°C)
Treatment Area	3.9 cm²

Feature	Detail
Power Input	100V-240V
Product Dimensions	11.7 x 4.9 x 12.5 cm
Weight	1.7 kilograms
Manufacturer	MLAY
Country of Origin	China

10. WARRANTY AND SUPPORT

The MLAY T14 IPL device comes with a **2-year warranty** from the date of purchase. This warranty covers manufacturing defects and malfunctions under normal use conditions.

For warranty claims, technical support, or any questions regarding your device, please contact MLAY customer service through the retailer where you purchased the product or visit the official MLAY website for contact information.

Please retain your proof of purchase for warranty validation.

© 2024 MLAY. All rights reserved.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question