

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [VRAIKO](#) /

› [VRAIKO Lily Neck Face Massager User Manual](#)

## VRAIKO LILY

# VRAIKO Lily Neck & Face Massager Instruction Manual

Model: LILY

Brand: VRAIKO

## PRODUCT OVERVIEW

---

The VRAIKO Lily Neck & Face Massager is a sophisticated skin rejuvenation device designed to enhance your skincare routine. It combines thermal therapy, triple-action LED light, and adjustable vibration to promote anti-aging, lifting, and tightening of sagging skin on both the neck and face, contributing to a radiant and glowing complexion.



Image: VRAIKO Lily Neck & Face Massager in Pearl White with its charging stand.

Key features include:

- **Triple Action Skin Rejuvenation:** Features three LED light modes (Blue, Green, Red) that offer distinct

dermatological benefits for skin beautification.

- **Enhanced Serum Absorption:** Provides gentle warmth (approx. 42°C to 45°C) to boost the absorption of serums, firming creams, or facial oils, maximizing their effectiveness.
- **Adjustable Vibrating Massage:** Offers three levels of vibration intensity (6000~8000, 8000~10000, and 10000~12000 times/min) for a customizable and comfortable massage experience.
- **Ergonomic Design:** The streamlined massage head is designed to comfortably glide along the contours of the neck, face, shoulders, and arms.

## SETUP

Before first use, ensure your VRAIKO Lily device is fully charged. The device comes with a USB-C charging cable and a custom stand for convenient charging and storage.

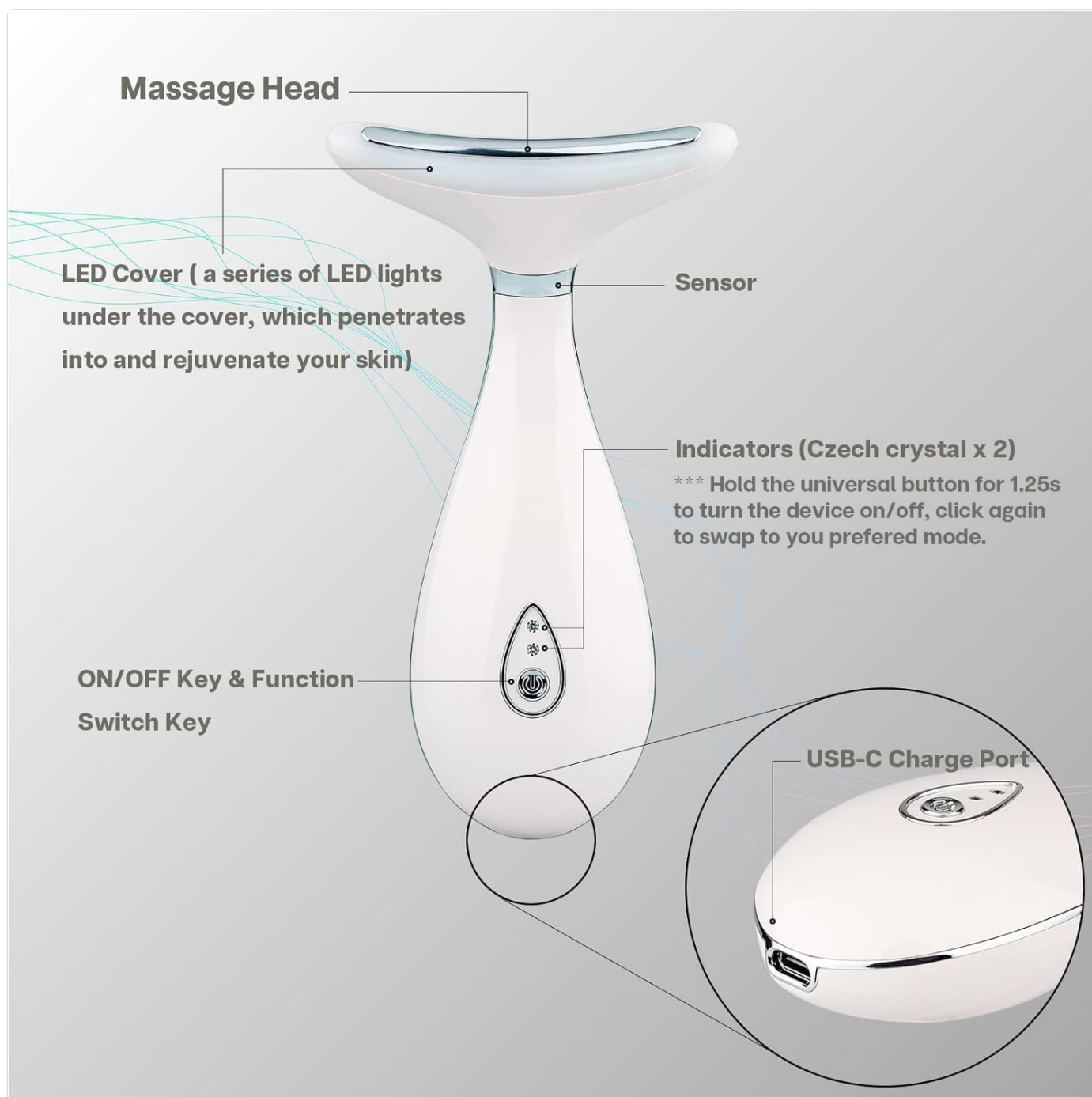


Image: Detailed diagram of the VRAIKO Lily device components.

- Connect the USB-C charging cable to the port at the base of the device.
- Plug the other end of the cable into a compatible USB power adapter (not included).

- Place the device on its custom stand while charging. The indicator lights will show charging status.
- A full charge typically takes a few hours. Once charged, the device is ready for use.

## OPERATING INSTRUCTIONS

---

Follow these steps for an effective skincare routine:

1. **Step 1: Prep Your Skin with Serums**

After facial cleansing, apply your preferred serums, firming creams, or facial oils to your neck and face. This helps the device glide smoothly and enhances product absorption.

2. **Step 2: Turn On and Select Mode**

Hold the universal button for 1.25 seconds to turn the LILY device on. Simply click the button again to switch between the different modes (Blue, Green, Red). Double-click the function key to adjust vibration intensity within each mode.

3. **Step 3: Apply and Glide**

Lightly press the top of the device to your skin. Glide in straight upward lines along your neck and jawline, or across each facial area you wish to treat. The ergonomic design allows for effortless application along contours.

# HOW TO USE

## Simple & Effortless



### Step1:

prep your skin with serums after facial cleansing.



### Step2:

Hold the universal button for 1.3 secs to turn the LILY on. Press the button again to pick the preferred mode for your skin.



### Step3:

Lightly press the top of the device to your skin. Glide straight upward lines with the heat and LED, along neck and jawline, or across each facial area you wish to treat.

Image: Visual guide for the 3-step usage process.

### LED Light Modes Explained:

# Multi-functional Beauty Device with Memory

starts with the mode at which you quit last time



Image: Explanation of the different LED light modes and their functions.

- **Blue Mode (6000~8000 times/min vibration, Normal temperature):** Ideal for tightening sagging skin and promoting collagen synthesis. Suitable for oily and sensitive skin.
- **Green Mode (8000~10000 times/min vibration, approx. 42°C):** Promotes skin surface microcirculation, fights edema and blackheads, and calms the skin.
- **Red Mode (10000~12000 times/min vibration, approx. 45°C):** Promotes blood circulation, reduces wrinkles and fine lines, and rejuvenates skin for a radiant complexion.

For optimal results, use the VRAIKO Lily 5-6 times a week as part of your daily routine, allowing for a rest day after each use.

## MAINTENANCE

To ensure the longevity and hygiene of your VRAIKO Lily device, proper maintenance is essential:

- After each use, gently wipe the massage head with a soft, damp cloth to remove any residue from skincare products.



- Do not submerge the device in water or use harsh chemicals for cleaning.
- Store the device on its custom stand in a cool, dry place away from direct sunlight.

## TROUBLESHOOTING

---

If you encounter any issues with your VRAIKO Lily device, please refer to the following common solutions:

- **Device not turning on:** Ensure the device is fully charged. Connect it to the USB-C charging cable and allow it to charge for at least 30 minutes before attempting to turn it on again.
- **Vibration not working:** Check if the device is powered on and in an active mode. Double-click the function key to cycle through vibration intensities.
- **LED lights not illuminating:** Verify the device is turned on and you have selected a mode. If the issue persists after charging, contact customer support.
- **Device not heating:** The thermal function is active in Green and Red modes. Ensure you are in one of these modes. The warmth is subtle and designed to enhance absorption, not to be intensely hot.

For persistent issues not covered here, please contact VRAIKO customer service.

## SPECIFICATIONS

---

Feature	Detail
Product Dimensions	6.3 x 3.54 x 1.18 inches
Item Weight	14.39 ounces
Item Model Number	LILY
Batteries	1 Lithium Ion battery (included)
Manufacturer	VRAIKO
Country of Origin	China
Color	Pearl White
Skin Type	All
Material	Acrylonitrile Butadiene Styrene (ABS)

## WARRANTY AND SUPPORT

---

The VRAIKO Lily Neck & Face Massager comes with a worry-free 12-month warranty from the date of purchase. This warranty covers manufacturing defects and issues arising from normal use.

For any questions, technical support, or warranty claims, please feel free to contact the VRAIKO customer service team. Contact information can typically be found on the product packaging or the official VRAIKO website.