

VRAIKO LILY

VRAIKO Lily Neck & Face Massager User Manual

Model: LILY

1. PRODUCT OVERVIEW

The VRAIKO Lily Neck & Face Massager is a sophisticated skin rejuvenation device designed to enhance your skincare routine. It combines thermal, triple action LED, and vibration technologies to provide comprehensive care for your neck and face, aiming to improve skin firmness, reduce signs of aging, and promote a radiant complexion.



Key Features:

- **Triple Action LED Technology:** Features blue, green, and red LED lights, each offering distinct benefits for skin rejuvenation based on dermatological theories.
- **Adjustable Vibrating Massage:** Provides gentle to strong vibrations (6000~12001 times/min) to stimulate blood flow and enhance product absorption.
- **Therapeutic Warmth:** Offers normal temperature, 42°C, and 45°C settings to boost the absorption of skincare products and provide a comfortable spa experience.
- **Ergonomic Design:** Streamlined massage head perfectly contours to the neck, face, shoulders, and arms for effortless and effective skincare.
- **Enhanced Absorption:** Designed to be used with serums, firming creams, or facial oils for improved results.

2. SAFETY INFORMATION

Warning:

- Avoid using the massager for prolonged periods beyond the advised duration.
- Keep out of reach of children. For external use only.
- Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

3. SETUP

Charging the Device:

1. Connect the provided USB-C Charging Cable to the device's charge port, located at the bottom.
2. Plug the other end of the USB-C cable into a compatible USB power adapter (not included).
3. The indicator lights will show charging status. Ensure the device is fully charged before first use.
4. The device comes with a custom stand for convenient storage and charging.

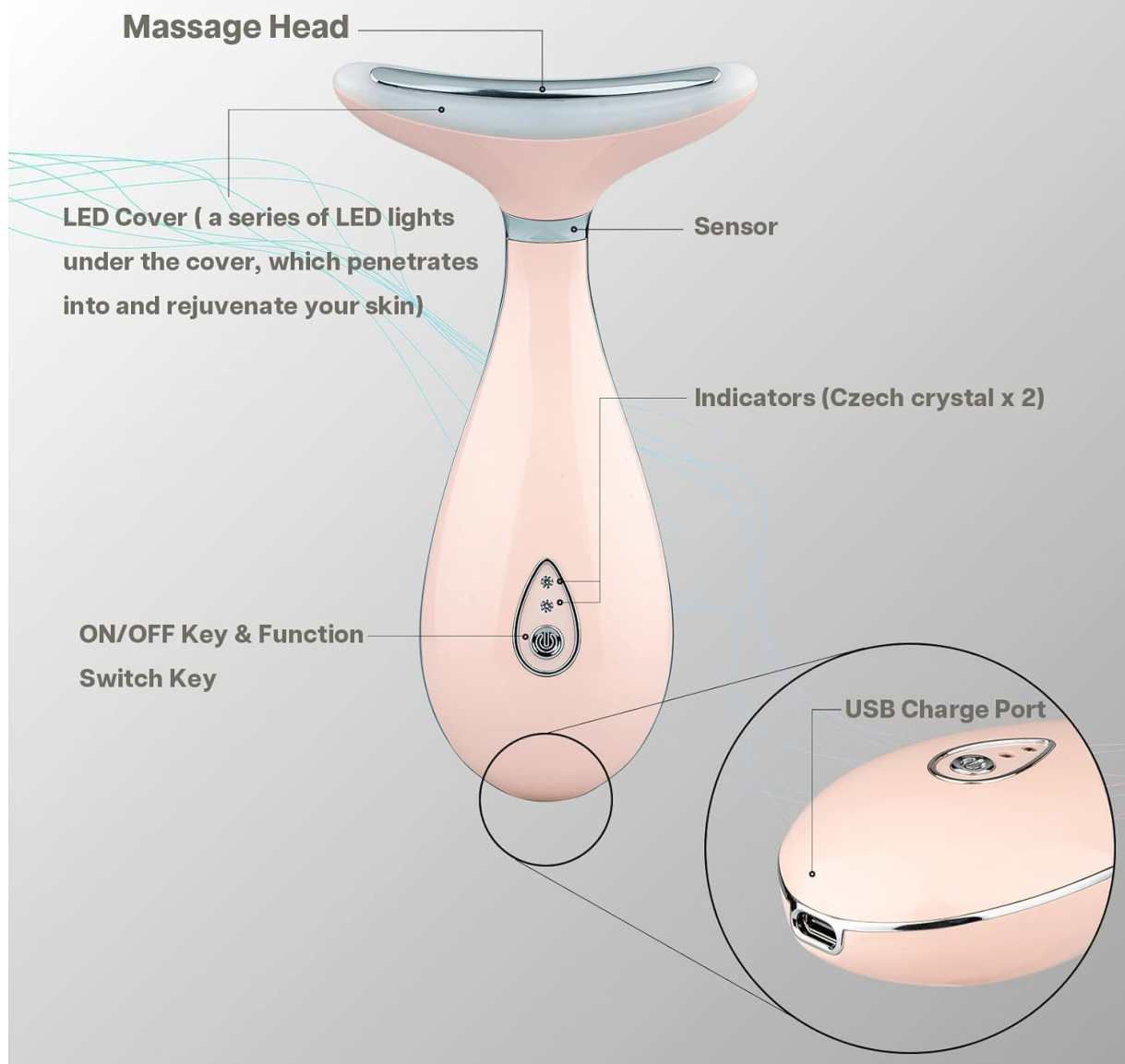


Image: Diagram illustrating the massager's components, including the USB charge port.

4. OPERATING INSTRUCTIONS

Usage Steps:

1. **Step 1: Prep Your Skin.** After facial cleansing, apply your preferred serums, firming cream, or facial oils to your skin.
2. **Step 2: Power On.** Hold the universal button for 1.25 seconds to turn the device On/Off. The device will remember the last mode used.
3. **Step 3: Massage.** Gently glide the massager along the contours of your neck, face, and jawline, moving from bottom to top. Allow the warmth to open pores and tone the skin.
4. **Step 4: Enhance Absorption.** After pores are opened, apply a facial mask to further boost the absorption of serums.

A daily routine of 2-6 minutes is recommended for effective skincare. Consistent use is key to achieving desired results. It is recommended to use the Lily 5-6 times a week, with a rest day after each use.

Understanding the Modes:

Multi-functional Beauty Device with Memory

starts with the mode at which you quit last time

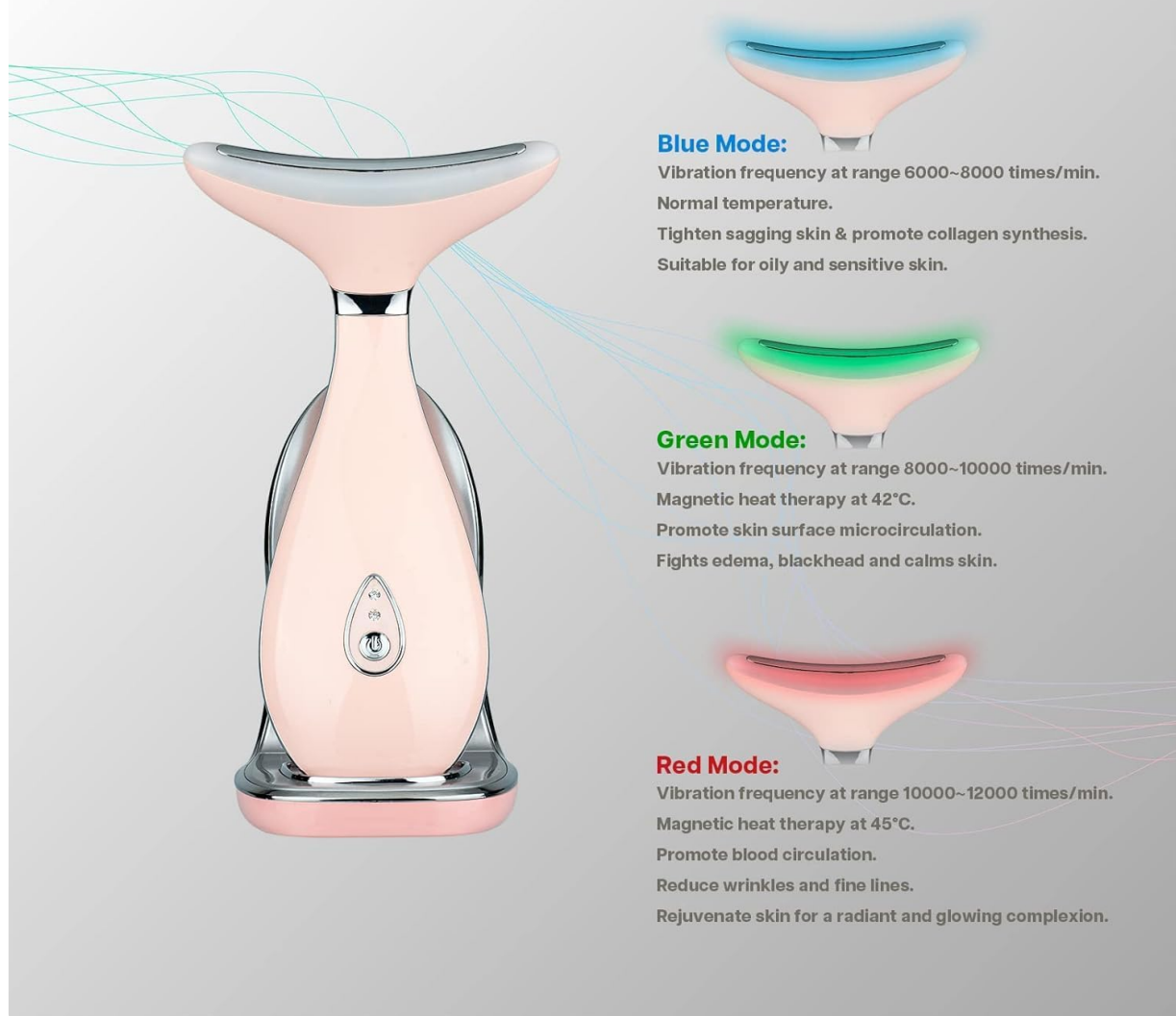


Image: Visual representation of the three distinct LED modes and their benefits.

• Blue Mode:

- Vibration frequency: 6000~8000 times/min.
- Temperature: Normal.
- Benefits: Tightens sagging skin, promotes collagen synthesis. Suitable for oily and sensitive skin.

• Green Mode:

- Vibration frequency: 8000~10000 times/min.
- Temperature: Approx. 42°C.
- Benefits: Promotes skin surface microcirculation, fights edema, blackheads, and calms skin.

• Red Mode:

- Vibration frequency: 10000~12001 times/min.
- Temperature: Approx. 45°C.
- Benefits: Promotes blood circulation, reduces wrinkles and fine lines, rejuvenates skin for a radiant and glowing complexion.

To adjust the vibrating intensity under each mode, simply double-click the function key to customize your most comfortable massage experience.

5. MAINTENANCE

Cleaning:

- After each use, gently wipe the massage head with a soft, damp cloth.
- Do not immerse the device in water or any other liquid.
- Avoid using abrasive cleaners or solvents that may damage the device's surface.

Storage:

- Store the device in a cool, dry place, away from direct sunlight and extreme temperatures.
- Utilize the custom stand provided for organized and safe storage.

6. TROUBLESHOOTING

If you encounter any issues with your VRAIKO Lily Neck & Face Massager, please refer to the following common troubleshooting tips:

- **Device not turning on:** Ensure the device is fully charged. Connect it to the USB-C charging cable and allow it to charge for at least 2 hours.
- **Vibration or LED not working:** Check if the device is powered on and if the correct mode is selected. Double-click the function key to cycle through modes.
- **Reduced performance:** Ensure the massage head is clean and free from any residue. Reapply serums or creams to ensure proper contact with the skin.

For persistent issues, please contact VRAIKO customer service.

7. SPECIFICATIONS

Feature	Detail
Product Dimensions	6.3 x 3.54 x 1.5 inches
Item Model Number	LILY
Batteries	1 Lithium Ion battery required (included)
Manufacturer	VRAIKO
ASIN	B0BMFS32FH
Country of Origin	China
Brand	VRAIKO
Color	Pink
Skin Type	All

Item Weight	0.42 Kilograms (14.82 ounces)
-------------	-------------------------------



Image: Dimensions of the VRAIKO Lily Neck & Face Massager.

8. WARRANTY AND SUPPORT

The VRAIKO Lily Neck & Face Massager comes with a worry-free 12-month warranty.

For any questions, concerns, or support needs, please feel free to contact the VRAIKO customer service team:

- **Email:** support@vraikocare.com
- **Website:** www.vraikobeaauty.com
- **Instagram:** [@vraiko_official](https://www.instagram.com/vraiko_official)

EN Warning! Keep out of reach of children. For external use only.
DE Warnung! Außer Reichweite von Kindern aufbewahren. Nur für äußere Anwendung.
FR Avertissement! Tenir hors de portée des enfants. Réservé à un usage externe.
ES Advertencia! Mantener fuera del alcance de los niños. Solo para uso externo.
IT Avvertenza! Tenere fuori dalla portata dei bambini. Solo per uso esterno.
NL Waarschuwing! Buiten bereik van kinderen houden. Alleen voor uitwendig gebruik.
PL Ostrzeżenie! Przechowywać poza zasięgiem dzieci. Tylko do użytku zewnętrznego.
SE Varning! Förvara utom räckhåll för barn. Endast för extern användning.

VRAIKO Beauty Technology Limited.

Website: www.vraikobeautey.com

Instagram: @vraiko_official

Support: support@vraikocare.com

Collaboration: brand@vraikocare.com


Manufacturer: Junlang Technology
(Shenzhen) Co., Ltd.

Address: Room 2007, Building 5,
Liangao Creative Park, Fuyuan 2nd Road,
Tangwei Community, Fu Hai Street,
Bao'an District, Shenzhen, China



APEX CE SPECIALISTS LIMITED

89 Princess Street, Manchester,
M1 4HT, UK
E-Mail : info@apex-ce.com

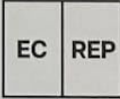


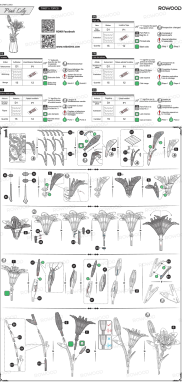




Apex CE Specialists GmbH

Ravolzhauser str.31a 63543
Neuberg Germany
E-Mail : Info@apex-ce.com



Image: Back of the product packaging with support details.

	<p>LILY 12W Microwave LED Sensor Lamp: Installation and Operation Guide</p> <p>Detailed guide for the LILY 12W Microwave LED Sensor Lamp, covering specifications, functions, installation, testing, and troubleshooting for optimal performance.</p>
	<p>Love in Lyrics: An Exploration of Supporting Textual Manifestation of Affection in Social Messaging</p> <p>This research paper explores a novel interactive system named Lily, designed to enhance affectionate communication in text messaging by suggesting relevant song lyrics. It details the system's design, implementation, and an empirical study with couples, offering insights into technology-mediated romantic interactions.</p>
	<p>Garmin Lily 2 Active Owner's Manual</p> <p>Comprehensive guide to the Garmin Lily 2 Active smartwatch, covering setup, activity tracking, health monitoring, connectivity, and troubleshooting.</p>
	<p>Rowood Pink Lily TW071 Assembly Instructions</p> <p>Step-by-step assembly guide for the Rowood Pink Lily TW071 3D wooden puzzle model. Includes a detailed part list and instructions for building the intricate floral model.</p>
	<p>Shape Club Programs - Lily of the Valley Wellness Luxury Hotel</p> <p>Explore the Shape Club programs at Lily of the Valley, a wellness luxury hotel. Discover tailored programs for weight loss, sport, detox, and better aging, designed by nutritionist Jacques Fricker.</p>
	<p>Garmin Lily 2 Användarhandbok</p> <p>Denna användarhandbok ger en detaljerad översikt över Garmin Lily 2 smartklocka. Lär dig hur du konfigurerar din enhet, utforskar dess olika widgets för hälso- och aktivitetsspårning, hanterar Bluetooth-anslutna funktioner och använder säkerhetsfunktioner. Handboken täcker även ämnen som träning, pulsfunktioner, klockinställningar och felsökning för att säkerställa en optimal användarupplevelse. För ytterligare information och support, besök garmin.com.</p>