

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [Sawgmore](#) /

› [Sawgmore Portable Oral Irrigator TH2D4 Instruction Manual](#)

Sawgmore TH2D4

Sawgmore Portable Oral Irrigator

Model: TH2D4 - User Manual

INTRODUCTION

Thank you for choosing the Sawgmore Portable Oral Irrigator. This device is designed to provide effective oral hygiene by using a high-pressure water pulse to clean between teeth and below the gumline, where traditional brushing cannot reach. Please read this manual thoroughly before use to ensure safe and optimal operation.

SAFETY INFORMATION

- Do not immerse the charging cable or adapter in water.
- Ensure the base cover is tightly closed before use to maintain IPX7 waterproof integrity.
- Keep out of reach of children.
- Do not direct the water jet at the eyes, nose, or other sensitive areas.
- Consult your dentist before use if you have had oral surgery or have severe gum disease.
- Use only with water or dentist-recommended mouthwash.

PRODUCT OVERVIEW

The Sawgmore Portable Oral Irrigator is a compact and powerful device for enhancing your daily oral care routine. It features multiple cleaning modes, a long-lasting battery, and a waterproof design for convenient use.



Figure 1: Contents of the Sawgmore Portable Oral Irrigator package, including the device, various jet tips, charging cable, and instruction manual.

Key Features:

- **Professional Dental Care:** High-pressure water pulse (1400 times/minute, 30-90 PSI) effectively removes food debris, massages gums, and cleans teeth.
- **Long Battery Life:** 1400 mAh lithium battery provides 30-40 days of use on a 4-hour charge.
- **IPX7 Waterproof Design:** Safe for use in the bathroom and easy to clean.
- **4 Cleaning Modes:** Normal, Soft, Pulse, and DIY customizable pressure.
- **Memory Function:** Remembers your last used mode.
- **4 Replaceable Jet Tips:** Includes standard tips, a tongue scraper, and a periodontal tip for comprehensive cleaning.
- **360-Degree Rotating Nozzle:** Reaches all areas of the mouth.
- **Portable Design:** Lightweight with a 220ml water tank, ideal for travel and home use.



Figure 2: The oral irrigator addressing common dental concerns such as gum bleeding, periodontal disease, tooth decay, bad breath, braces, periodontal pockets, dental implants, and dental calculus.

SETUP

Initial Charging:

Before first use, fully charge the device. Connect the USB charging cable to the irrigator's charging port and plug the other end into a compatible USB power source (e.g., computer, car charger, power bank, or power adapter).



Figure 3: The oral irrigator being charged. It can be powered by a computer, car charger, power bank, or a standard power adapter. A full charge takes approximately 4 hours and provides up to 30-40 days of usage.

The charging indicator light will show the charging status. A full charge typically takes 4 hours and provides 30-40 days of usage.

Attaching a Jet Tip:

1. Select the desired jet tip (standard, tongue scraper, or periodontal).
2. Insert the tip into the hole at the top of the irrigator until it clicks into place.
3. To remove, press the tip release button (usually located on the back or side of the device) and pull the tip out.

Filling the Water Tank:

1. Open the water tank cover located at the bottom or back of the device.
2. Fill the 220ml tank with lukewarm water. You may also use dentist-recommended mouthwash.
3. Close the water tank cover securely to prevent leakage.

OPERATING INSTRUCTIONS

Getting Started:

1. Lean over a sink and place the jet tip into your mouth, aiming it at your gumline.
2. Press the Power button to turn on the device.
3. Press the Mode button (M) to cycle through the available cleaning modes:

4 Pressure Settings

DIY more personal intensities choice




Normal

Deeply clean



Soft

For Beginners & Sensitive Teeth

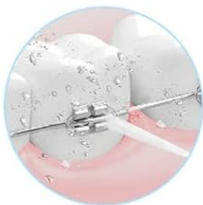

Pulse

Massage Gums



DIY

Set your personal pressure



Figure 4: The device's control panel highlighting the four pressure settings: Normal for deep cleaning, Soft for beginners and sensitive teeth, Pulse for gum massage, and DIY for custom pressure.

- **Normal:** Strong water pulses (1400 times/minute) for effective removal of food particles.
- **Soft:** Gentle pressure, ideal for beginners and sensitive teeth.
- **Pulse:** Rhythmic pulses for gum massage and improved circulation.
- **DIY:** Allows you to customize the water pressure to your preference.

The innovative memory function will save your last selected mode for convenience.

Flossing Technique:

- Start with the back teeth and work your way forward.
- Hold the handle at a 90-degree angle to your teeth.
- Guide the water jet along the gumline, pausing briefly between teeth.
- Allow the water to flow into the sink.
- Continue until you have cleaned all areas of your mouth.



Figure 5: A user demonstrating the proper way to hold and use the Sawgmore Portable Oral Irrigator for effective cleaning.

After use, press the Power button to turn off the device. Empty any remaining water from the tank.

CLEANING AND MAINTENANCE

Daily Cleaning:

1. After each use, empty the water tank.

2. Rinse the water tank with clean water.
3. Wipe the exterior of the device with a soft, damp cloth.
4. Remove the jet tip and rinse it under running water.

Deep Cleaning (Weekly/Bi-weekly):

- To prevent mineral buildup, fill the water tank with a mixture of 1 part white vinegar and 2 parts water.
- Run the device until the tank is empty.
- Rinse the tank and run the device with a full tank of clean water to flush out any remaining vinegar.

Storage:

Store the device in a cool, dry place. Ensure the water tank is empty and dry before storing. The compact design makes it easy to store or carry for travel.



Figure 6: The portable design of the Sawgmore Oral Irrigator allows it to be easily stored in a bag, making it convenient for travel.

Waterproof Feature:

IPX7 **Waterproof**

Note

Make sure the base cover is closed tightly!



Figure 7: The IPX7 waterproof rating ensures the device is protected against immersion in water, making it safe for use in the bathroom. Always ensure the base cover is closed tightly.

The device has an IPX7 waterproof rating, meaning it can withstand immersion in water up to 1 meter for 30 minutes. This makes it safe for use in the shower or bath. However, always ensure the base cover is tightly closed to maintain its waterproof integrity.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or no charge.	Charge the device fully (approx. 4 hours).
Weak water pressure.	Low battery, clogged tip, or mineral buildup.	Charge the device. Clean the jet tip. Perform a deep clean with vinegar solution.
Water leakage.	Water tank cover not closed properly, or tank/seal damaged.	Ensure the water tank cover is securely closed. If leakage persists, contact customer support.

Problem	Possible Cause	Solution
Device stops during use.	Low battery or water tank empty.	Recharge the device or refill the water tank.

SPECIFICATIONS

- **Brand:** Sawgmore
- **Model:** TH2D4
- **Power Source:** Battery Powered (1400 mAh Lithium Battery)
- **Charging Time:** Approximately 4 hours
- **Usage Time:** 30-40 days per full charge
- **Water Pressure:** 30-90 PSI
- **Water Pulse Frequency:** 1400 times/minute
- **Water Tank Volume:** 220 ml
- **Waterproof Rating:** IPX7
- **Material:** Plastic
- **Color:** Black
- **Included Tips:** 2 Standard Jet Tips, 1 Tongue Scraper, 1 Periodontal Tip

WARRANTY AND SUPPORT

This Sawgmore Portable Oral Irrigator comes with a standard manufacturer's warranty. Please refer to the warranty card included in your package for specific terms and conditions, including warranty duration and coverage details.

For technical support, troubleshooting assistance, or warranty claims, please contact Sawgmore customer service through the contact information provided on the product packaging or our official website. Please have your model number (TH2D4) and purchase date ready when contacting support.