

[manuals.plus](#) /

› [Beautifully](#) /

› [Beautifully B-Plasma PRO Facial Massager Instruction Manual](#)

## Beautifully B-Plasma PRO

# Beautifully B-Plasma PRO Facial Massager Instruction Manual

Model: B-Plasma PRO

## INTRODUCTION

---

The Beautifully B-Plasma PRO is an advanced skincare device designed for effective acne reduction, skin regeneration, and anti-aging care. It utilizes plasma technology and ultrasonic cavitation peeling to improve skin health and appearance. This manual provides essential information for safe and optimal use of your device.



Image: The Beautifully B-Plasma PRO device, showcasing its design and branding.

## IMPORTANT SAFETY INFORMATION

---

- Read all instructions before first use.
- Do not use on broken, irritated, or infected skin.
- Consult a physician before use if you have any skin conditions, allergies, or are pregnant.
- Keep the device dry and do not immerse it in water.
- Keep out of reach of children.
- Do not use if the device or charging cable is damaged.
- Use only the provided charging cable and power adapter.

## PACKAGE CONTENTS

---

---

The Beautifly B-Plasma PRO package includes the following items:

- Beautifly B-Plasma PRO Device
- Power adapter
- Charging cable
- User manual

## SETUP AND CHARGING

---

1. **Unpack:** Carefully remove all components from the packaging.
2. **Initial Charge:** Before first use, fully charge the device. Connect the charging cable to the device and the power adapter, then plug it into a suitable power outlet. The charging indicator light will show the charging status.
3. **Clean Device:** Wipe the device head with a soft, damp cloth before first use.



Image: Close-up view of the device, highlighting the power button and mode selection area.

## OPERATING INSTRUCTIONS

---

The Beautifly B-Plasma PRO offers two main modes: **Cavitation Peeling** and **Plasma Therapy**. Ensure your skin is clean before use.

### 1. Cavitation Peeling Mode

This mode uses ultrasonic vibrations to gently exfoliate the outer layers of the epidermis, removing dead skin cells and impurities. It helps to improve skin texture and prepare the skin for better absorption of skincare products.

1. **Preparation:** Dampen your face with water or a mild toner. The skin must be wet for effective peeling.
2. **Activate:** Press the power button to turn on the device. Select the "Peeling" mode.
3. **Application:** Gently glide the spatula-shaped head across your skin in an upward and outward motion, applying light pressure. Focus on areas with blackheads or rough texture.
4. **Duration:** Use for approximately 5-10 minutes, 1-2 times per week.



Image: A user demonstrating the application of the device in peeling mode on the cheek.

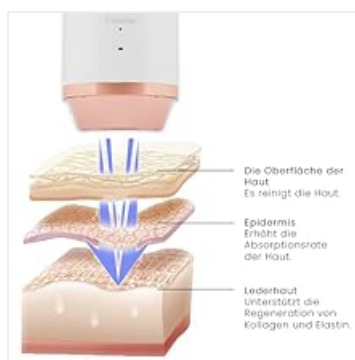


Image: Diagram illustrating how cavitation peeling cleanses the skin surface and enhances absorption of products into the epidermis.

## 2. Plasma Therapy Mode

Plasma therapy utilizes micro-discharges to generate plasma and ozone, which have disinfecting, bactericidal, and regenerative properties. This mode is effective for treating acne, stimulating collagen and elastin production, and

improving overall skin tension.

1. **Preparation:** Ensure your skin is clean and dry. You may apply a serum or cream suitable for plasma treatment.
2. **Activate:** Press the power button to turn on the device. Select the "Plasma" mode.
3. **Application:** Gently touch the tip of the device to the skin. You may feel a slight tingling sensation. Move the device slowly over the desired areas, such as acne spots, fine lines, or areas needing regeneration.
4. **Duration:** Use for approximately 5-10 minutes, 2-3 times per week.



Image: A user applying the device in plasma mode, highlighting its multi-functional capabilities including ozone sterilization.





Image: Diagram illustrating the effect of plasma therapy on acne-prone skin, leading to healthier skin.

## MAINTENANCE AND CARE

- **Cleaning:** After each use, wipe the device head with a soft, damp cloth or an alcohol wipe. Ensure no residue remains.
- **Drying:** Allow the device to air dry completely before storing.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and moisture.
- **Charging:** Recharge the device when the battery indicator is low to ensure optimal performance.

## TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery.	Charge the device fully.
Device not charging.	Charging cable or adapter issue; improper connection.	Ensure cable is securely connected. Try a different power outlet.
No sensation during use (Peeling mode).	Skin is not sufficiently wet.	Ensure skin is thoroughly dampened with water or toner.
No sensation during use (Plasma mode).	Improper contact with skin.	Ensure the device tip makes firm contact with the skin.

If you encounter issues not listed here, please contact customer support.

## PRODUCT SPECIFICATIONS

**Brand:** Beautifly

**Model Number:** 5903003701019

**Dimensions (L x W x H):** 8 x 7.5 x 7.5 cm

**Weight:** 235 Grams

**Color:** White

**Special Features:** Cleaning, anti-acne care, ozone sterilization

**Recommended Uses:** Acne reduction, skin cell reconstruction, skin regeneration, collagen structure reconstruction, ozone sterilization

**Skin Type:** Dry, All

**Scent:** Unscented

## WARRANTY AND SUPPORT

---

For warranty information and customer support, please refer to the documentation included with your purchase or contact Beautify customer service directly. Keep your proof of purchase for any warranty claims.

