

[Manuals.plus](#) /

› [HassoKon](#) /

› HassoKon Electric Teeth Cleaning Kit User Manual

HassoKon SOL026B

HassoKon Electric Teeth Cleaning Kit User Manual

Model: SOL026B

1. IMPORTANT SAFETY INFORMATION

Please read all safety information before using the HassoKon Electric Teeth Cleaning Kit. Failure to follow these instructions may result in injury or damage to the device.

- Keep away from children. This device is not a toy.
- Do not submerge the device in water beyond the recommended depth. The body is waterproof for rinsing, but not for full submersion.
- Consult a healthcare professional or dentist if you have pre-existing medical conditions related to gums or teeth, or if you experience discomfort during use.
- Ensure the device is properly charged before use to maintain optimal performance.
- This device is designed for cleaning teeth and removing stains/tartar. It is not a substitute for professional dental care.
- Use gentle pressure. Do not press too hard on teeth or gums. The intelligent identification technology is designed to prevent operation on soft tissues.

2. PRODUCT OVERVIEW AND COMPONENTS

The HassoKon Electric Teeth Cleaning Kit is designed for effective at-home oral hygiene. It includes various components to address different cleaning needs.



Figure 2.1: Kit Contents. This image displays the complete HassoKon Teeth Cleaning Kit, including the main electric cleaning device, a variety of interchangeable cleaning heads, a tongue scraper, an oral mirror, and a small wrench for attaching heads. The packaging box is

also visible.

2.1 Main Device Features

- **Intelligent Identification Technology:** The device is designed to only activate on hard surfaces (teeth) and stop vibrating when it touches soft tissues (gums), ensuring safe operation.
- **Built-in LED Light:** Provides illumination for better visibility during cleaning, especially in hard-to-reach areas.
- **4 Adjustable Modes:** Soft, Normal, Middle, and Strong modes cater to different sensitivity levels and cleaning requirements.
- **USB Fast Charging:** Convenient charging via USB with automatic power-off when fully charged.
- **IPX6 Waterproof:** The device body is highly water-resistant for easy rinsing and cleaning.

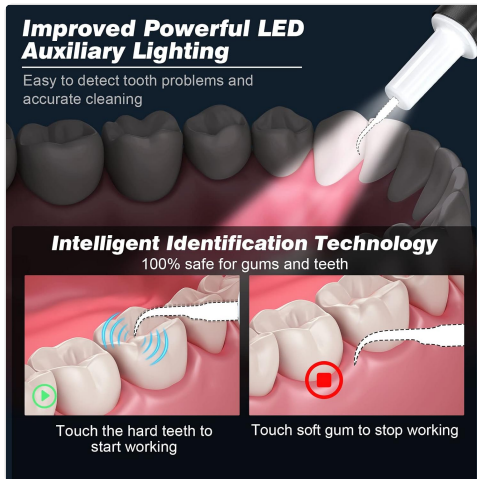


Figure 2.2: Intelligent Identification and LED Light. This image highlights the intelligent identification technology, showing how the device activates on teeth and pauses on gums. It also illustrates the integrated LED light for enhanced visibility during use.

2.2 Included Cleaning Heads and Accessories

- **Flat Head:** Ideal for cleaning deep and large calculus and stains on tooth surfaces.
- **Pointed Head:** Designed for precise cleaning of small calculus and plaque, particularly in the gums and between teeth.
- **Mouth Opener:** Helps keep the mouth open for easier access to all areas.
- **Oral Mirror:** Provides 360-degree observation of the oral cavity.
- **Tongue Scraper:** For effectively cleaning tongue coating.
- **Dental Tweezer:** Useful for clipping cotton or handling small items in the mouth.
- **Tartar Scraper:** For removing stubborn tartar on the surface of teeth.
- **Sickle Shape Scraper:** Designed to remove deep calculus between teeth and below the gumline.



Figure 2.3: Full Set Cleaning Kit. This image illustrates the different types of cleaning heads (Flat Head, Pointed Head) and additional dental tools (Mouth Opener, Oral Mirror, Tongue Scraper, Dental Tweezer, Tartar Scraper, Sickle Shape Scraper) that come with the HassoKon kit, along with their primary uses.

3. SETUP

3.1 Charging the Device

1. Before first use, fully charge the device.
2. Connect the USB charging cable to the device's charging port and a suitable USB power adapter (not included).
3. The charging indicator light will show the charging status. A full charge typically takes about 2 hours and provides up to 60 days of working time.
4. The device will automatically power off when fully charged to prevent overcharging.



Figure 3.1: USB Charging and Waterproofing. This image illustrates the USB charging capability of the HassoKon device and highlights its IPX6 waterproof rating, showing it being rinsed under water. It also displays the four cleaning modes and the intelligent identification feature.

3.2 Attaching a Cleaning Head

1. Select the appropriate cleaning head (Flat or Pointed) for your needs.
2. Align the cleaning head with the device's tip.
3. Use the provided small wrench to securely tighten the cleaning head onto the device. Ensure it is firmly attached to prevent detachment during use.



Figure 3.2: Attaching Cleaning Heads. This diagram demonstrates the process of attaching the professional oral-grade stainless steel cleaning heads (Pointed Head and Flat Head) to the HassoKon device using the provided spanner, ensuring a secure fit for effective cleaning.

4. OPERATING INSTRUCTIONS

4.1 Powering On and Mode Selection

1. Press the power button to turn on the device. The LED light will illuminate.
2. Press the power button repeatedly to cycle through the four intensity modes:
 - **SOFT (34KHz):** For sensitive teeth and normal cleaning.
 - **NORMAL (36KHz):** For daily oral cleaning.
 - **MIDDLE (38KHz):** For teeth polishing and stronger cleaning.
 - **STRONG (40KHz):** For huge plaque and stubborn stains.
3. Select the mode that best suits your comfort level and cleaning needs. Start with the SOFT mode if you are new to the device or have sensitive gums.



Figure 4.1: 4 Intensity Cleaning Modes. This image details the four adjustable intensity modes of the HassoKon device: Soft (34KHz), Normal (36KHz), Middle (38KHz), and Strong (40KHz), explaining their recommended uses for different oral care needs.

4.2 Cleaning Procedure

1. Ensure your mouth is clean. You may brush your teeth first.
2. Use the oral mirror and mouth opener for better visibility and access to all areas of your mouth.
3. Gently place the tip of the selected cleaning head onto the tooth surface where plaque or tartar is present. The device will only vibrate when it makes contact with a hard surface.
4. Move the cleaning head slowly and gently along the tooth surface, focusing on areas with stains or calculus. Do not apply excessive pressure. Let the device do the work.
5. Rinse your mouth with water periodically to remove dislodged debris.
6. After cleaning, turn off the device by pressing the power button.



Figure 4.2: Before and After Cleaning. This image provides a visual comparison of teeth before and after using the HassoKon Electric Teeth Cleaning Kit, demonstrating its effectiveness in removing stains and improving oral hygiene.

4.3 Using Other Accessories

- **Tongue Scraper:** Gently scrape the surface of your tongue from back to front to remove coating and improve breath freshness.
- **Dental Tweezer:** Use for handling small items or cotton during your oral care routine.
- **Tartar Scraper / Sickle Shape Scraper:** These manual tools can be used for more targeted removal of stubborn tartar, but exercise extreme caution to avoid damaging gums or enamel.



Figure 4.3: Accessory Usage. This image demonstrates the practical application of various accessories included in the HassoKon kit, such as the mouth opener, oral mirror, tongue scraper, dental tweezer, sickle shape scraper, and tartar scraper, illustrating their roles in a comprehensive oral hygiene routine.

5. MAINTENANCE

5.1 Cleaning the Device and Heads

1. After each use, remove the cleaning head from the device.
2. Rinse the cleaning head and the tip of the device under running water. The device body is IPX6 waterproof, allowing for easy rinsing.
3. Wipe the device body with a clean, dry cloth.
4. Ensure all parts are completely dry before storing.

5.2 Storage

- Store the HassoKon Teeth Cleaning Kit in a cool, dry place, away from direct sunlight and extreme temperatures.
- Keep all components, especially the cleaning heads, in the provided storage box to prevent damage and maintain hygiene.

6. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or no charge.	Charge the device fully using the provided USB cable.
Device does not vibrate on teeth.	Intelligent identification technology is active (touching soft tissue). Cleaning head not securely attached.	Ensure the cleaning head is in firm contact with a hard tooth surface. Re-attach the cleaning head securely with the wrench.
Reduced cleaning effectiveness.	Incorrect mode selected. Cleaning head worn out or dirty.	Switch to a stronger mode (Middle or Strong). Clean the cleaning head thoroughly or replace it if necessary.
Device makes unusual noise.	Cleaning head not properly attached.	Ensure the cleaning head is securely tightened with the wrench.

7. SPECIFICATIONS

- **Brand:** HassoKon
- **Model Number:** SOL026B
- **Power Source:** Corded Electric (USB Rechargeable)
- **Special Feature:** Rechargeable, LED Light, Intelligent Identification Technology
- **Product Benefits:** Removes Plaque, Tartar, Stains
- **Intensity Modes:** 4 (Soft, Normal, Middle, Strong)
- **Waterproof Rating:** IPX6

- **Charging Time:** Approximately 2 hours
- **Working Time (Full Charge):** Up to 60 days
- **Product Dimensions:** 24.5 x 9.7 x 4.5 cm
- **Item Weight:** 400 g

8. WARRANTY AND SUPPORT

For any questions, concerns, or issues regarding your HassoKon Electric Teeth Cleaning Kit, please contact our customer support. Our after-sales service experts are committed to assisting you.

Please refer to your purchase documentation for specific warranty details. If you are not satisfied with our product, please feel free to contact us. We will do our best to solve your problems.

Contact Information: Refer to the contact details provided with your product packaging or on the official HassoKon website.

© 2025 HassoKon. All rights reserved.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page’s manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question

Answered questions about this manual

Why is my HassoKon SOL026B teeth cleaning kit not lighting up or charging when plugged in?

If your HassoKon SOL026B is not lighting up or showing a green indicator when charging, it may be due to a faulty connection, a depleted battery, or a hardware failure. Please try the following...

[Read full answer](#)

