

Beautifly 5903003700302

Beautifly BLooky PRO Sonic Eye Massager User Manual

Model: 5903003700302

1. INTRODUCTION

The Beautifly BLooky PRO Sonic Eye Massager is an advanced anti-aging device designed to rejuvenate the delicate skin around the eyes. It integrates multiple technologies to promote collagen production, reduce signs of fatigue, and improve skin elasticity.

This manual provides essential information for the safe and effective use of your BLooky PRO device. Please read it thoroughly before first use and retain it for future reference.

2. SAFETY INFORMATION

Before using the BLooky PRO, please observe the following safety precautions:


- Do not use on broken, irritated, or infected skin.
- Consult a physician before use if you have any pre-existing medical conditions, especially related to skin, eyes, or circulation.
- Keep out of reach of children.
- Do not immerse the device in water. Clean with a damp cloth only.
- Use only the provided charging cable.
- Discontinue use if you experience any discomfort or adverse reactions.

3. PACKAGE CONTENTS

Verify that all items are present in your package:

- Beautifly BLooky PRO Sonic Eye Massager
- USB Charging Cable
- User Manual (this document)

4. PRODUCT OVERVIEW

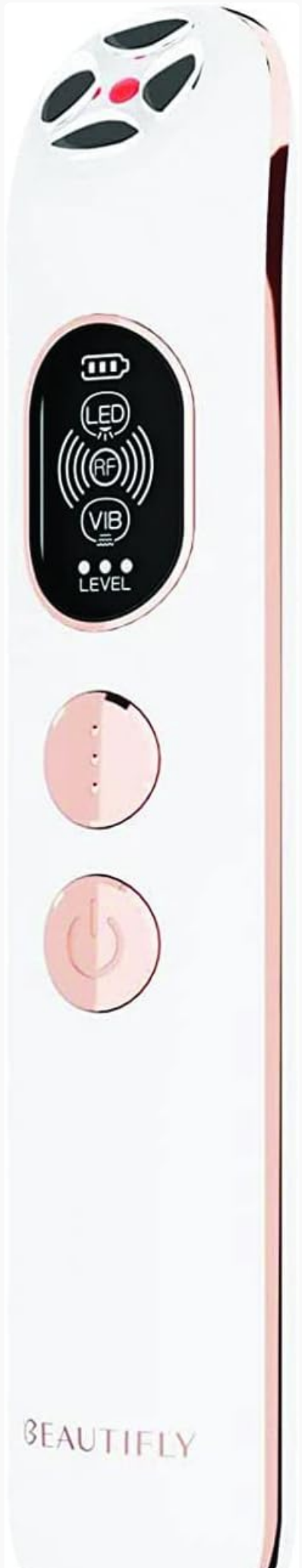

BEAUTIFLY



Image 1: Front view of the Beautify BLooky PRO Sonic Eye Massager, showing the massage head, display, and control buttons.

The BLooky PRO features an ergonomic design for comfortable handling and precise application. Key components include:

- **Message Head:** Hypoallergenic metal head for skin contact, incorporating RF, photon light, and vibration functions.
- **Display Screen:** Shows current mode, intensity level, and battery status.
- **Mode/Intensity Button:** Used to cycle through modes and adjust intensity.
- **Power Button:** Turns the device on/off.

5. SETUP

5.1 Charging the Device

1. Connect the USB charging cable to the charging port on the device.
2. Plug the other end of the USB cable into a compatible USB power adapter (not included) or a computer USB port.
3. The battery indicator on the display will show charging status. A full charge typically takes approximately 2-3 hours.
4. Unplug the device once fully charged.

5.2 Initial Use

Before the first use, ensure the device is fully charged. Clean the massage head with a soft, damp cloth.

6. OPERATING INSTRUCTIONS

6.1 Powering On/Off

- To turn on: Press and hold the Power Button for 2 seconds. The device will emit a beep and the display will light up.
- To turn off: Press and hold the Power Button for 2 seconds. The device will beep and the display will turn off.

6.2 Selecting Modes and Intensity

The BLooky PRO offers various modes combining RF, photon light, and vibration. It also has 3 intensity levels.

1. After powering on, press the Mode/Intensity Button to cycle through available modes. The display will indicate the active functions (LED, RF, VIB).
2. Once a mode is selected, press the Mode/Intensity Button again to adjust the intensity level (Level 1, 2, or 3).

6.3 Application



Image 2: A user applying the BLooky PRO to the under-eye area.



Image 3: A user applying the BLooky PRO to the cheek area, demonstrating broader application.

For optimal results, apply a thin layer of eye cream or serum to the desired treatment area before use. Gently glide the massage head over the skin around your eyes, moving in small, circular motions or upward strokes. Avoid direct contact with the eyeball. Recommended usage is 5-10 minutes per session, 2-3 times per week.

7. MAINTENANCE

7.1 Cleaning

After each use, wipe the massage head with a soft, slightly damp cloth. Do not use abrasive cleaners or immerse the device in water.

7.2 Storage

Store the device in a cool, dry place, away from direct sunlight and extreme temperatures. Keep it in its original packaging or a protective pouch to prevent damage.

8. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery.	Charge the device fully.
No vibration/light/RF.	Incorrect mode selected or device malfunction.	Ensure the correct mode is selected. If problem persists, contact customer support.
Skin irritation.	Sensitive skin or improper use.	Reduce intensity or frequency of use. Ensure skin is clean and moisturized. Discontinue use if irritation persists.

If you encounter issues not listed here, please contact Beautifly customer support.

9. SPECIFICATIONS

- **Model Number:** 5903003700302
- **Brand:** Beautifly
- **Dimensions (L x W x H):** 14.9 x 2.8 x 2.2 cm
- **Weight:** 54 Grams
- **Technology:** Multi-polar RF (400 KHz), Red Photon Light (650 nm), Sonic Vibration (12,000 rpm)
- **Intensity Levels:** 3
- **Battery:** 1 Lithium-ion (rechargeable)
- **Color:** White
- **Country of Origin:** Germany

10. WARRANTY AND SUPPORT

Beautifly products are manufactured to high-quality standards. For warranty information or technical support, please refer to the warranty card included with your purchase or visit the official Beautifly website. Please have your model number (5903003700302) and proof of purchase ready when contacting support.