

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

> [Quip](#) /

> [Quip Rechargeable Smart Sonic Toothbrush \(Model 920-00193\) Instruction Manual](#)

Quip 920-00193

Quip Rechargeable Smart Sonic Toothbrush (Model 920-00193) Instruction Manual

Your guide to optimal oral care with Quip.

1. PRODUCT OVERVIEW

The Quip Rechargeable Smart Sonic Toothbrush is designed to enhance your daily oral hygiene routine. It features timed sonic vibrations, a built-in 2-minute timer with 30-second pulses, and Bluetooth connectivity to track your brushing habits via the Quip app.



Image: The Quip Rechargeable Smart Sonic Toothbrush in its sleek black design.

2. SETUP AND FIRST USE

2.1 Charging Your Toothbrush

Before first use, fully charge your Quip toothbrush. The magnetic charging cable easily attaches to the base of the toothbrush. A single charge can last up to 3 months when used for 2 minutes, twice daily.

Magnetic cable, up to a 3-month charge

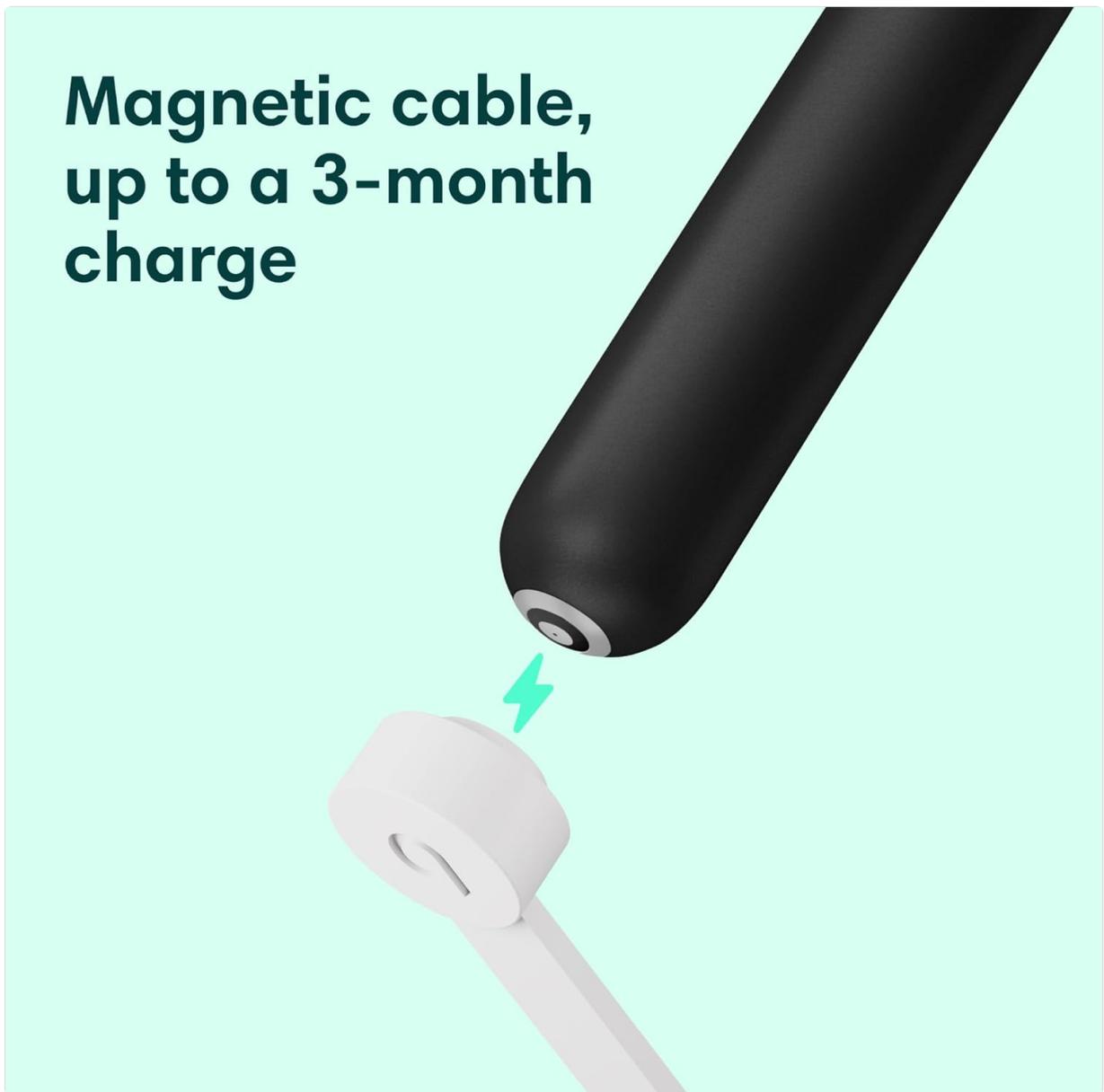


Image: The magnetic charging cable connecting to the base of the Quip toothbrush.

Video: This video demonstrates how to properly charge your Quip Smart Sonic Toothbrush using the magnetic cable.

2.2 Connecting to the Quip App

For advanced tracking and rewards, download the Quip app and connect your toothbrush via Bluetooth. The app monitors brushing time, strokes, coverage, and intensity, providing insights to improve your oral hygiene habits.

Earn rewards with the quip app

Auto-track brushing (time, strokes, coverage, intensity) to earn points and redeem rewards!

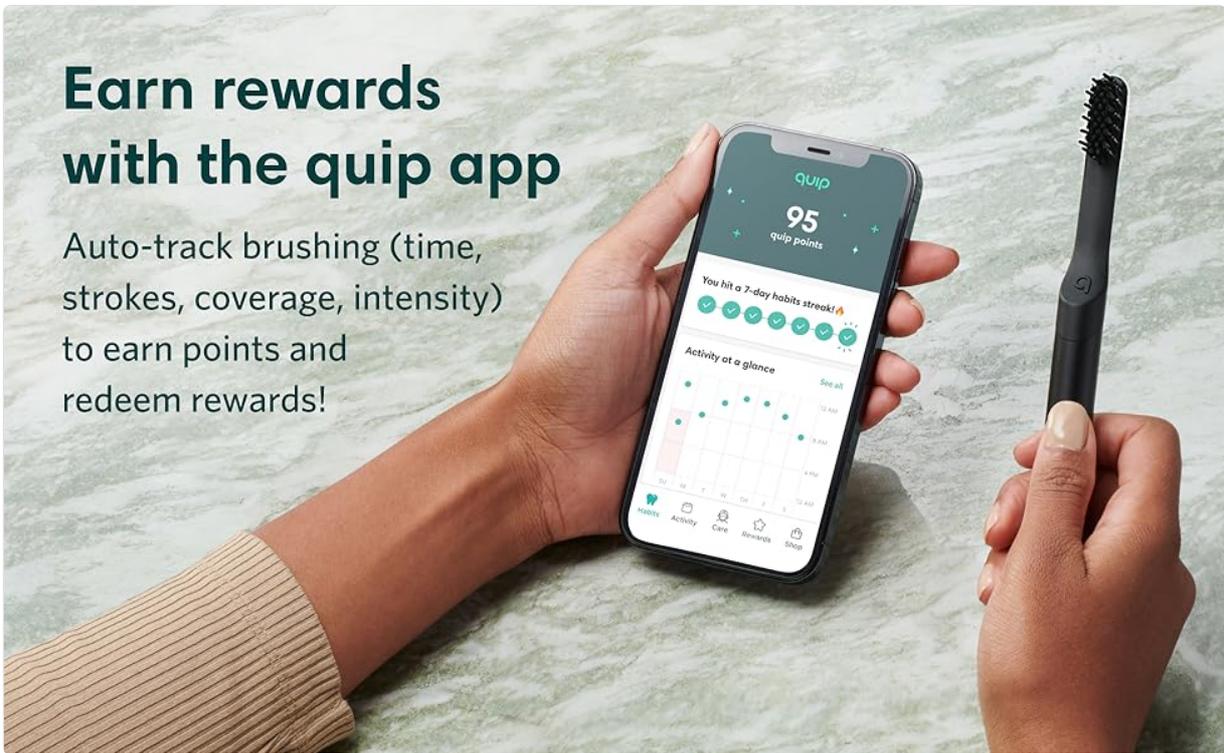


Image: A smartphone displaying the Quip app interface, showing brushing metrics and reward points.

3. OPERATING YOUR TOOTHBRUSH

The Quip toothbrush features timed sonic vibrations and soft bristles designed for effective cleaning of tight spaces and sensitive teeth and gums. Press the button on the handle to activate the brush.

Soft bristles, tongue scraper, safe for sensitive gums



Image: A close-up view of the Quip toothbrush head, highlighting its soft bristles and integrated tongue scraper.

3.1 Guided Brushing Technique

The toothbrush includes a 2-minute built-in timer for a complete routine, with 30-second pulses to guide you to switch brushing quadrants, ensuring an even clean across your mouth.

Video: This video demonstrates the proper use of the Quip Smart Sonic Toothbrush for effective cleaning.

4. MAINTENANCE

4.1 Replacing the Brush Head

Dentists recommend replacing brush heads every 3 months to maintain optimal cleaning performance. The Quip design allows you to keep the handle and replace only the brush head, promoting an eco-friendly approach.

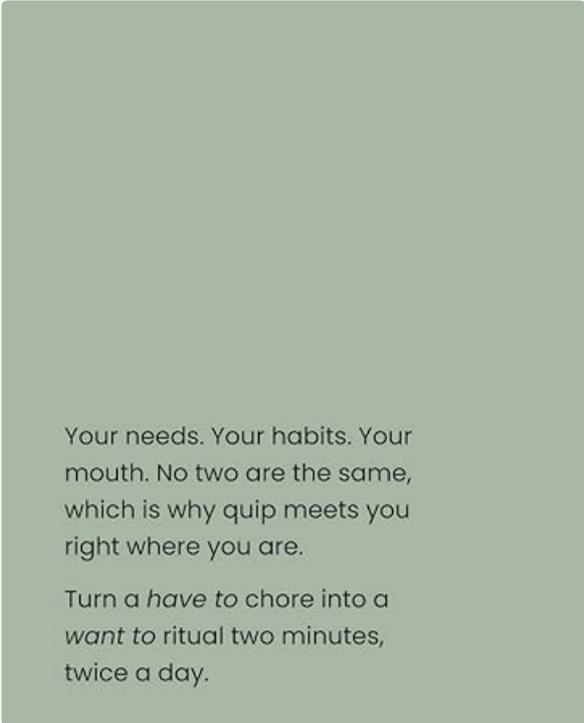


Image: An illustration highlighting the replaceable brush head feature of the Quip toothbrush, emphasizing its reusability.

4.2 Cleaning Your Toothbrush

After each use, rinse the brush head thoroughly under running water. Periodically wipe down the handle with a damp cloth to remove any residue. Ensure the toothbrush is dry before storing.

5. TROUBLESHOOTING

5.1 Toothbrush Not Vibrating

- **Check Charge:** Ensure the toothbrush is fully charged. Connect it to the magnetic charger and allow sufficient time for charging.
- **Re-seat Brush Head:** Remove the brush head and re-attach it firmly to ensure proper connection.
- **Contact Support:** If the issue persists after checking charge and re-seating the brush head, please contact Quip customer support for further assistance.

5.2 App Connectivity Issues

- **Bluetooth On:** Verify that Bluetooth is enabled on your smartphone.
- **App Permissions:** Ensure the Quip app has necessary Bluetooth permissions enabled in your phone settings.
- **Restart App/Toothbrush:** Close and reopen the Quip app, or turn the toothbrush off and on again, then attempt to reconnect.

6. SPECIFICATIONS

Feature	Detail
Product Dimensions	8.94 x 3.78 x 1.53 inches
Item Weight	6.24 ounces

Model Number	920-00193
Brand	Quip
Age Range	Adult
Special Feature	Brushing Timer, Bluetooth Connectivity
Power Source	Rechargeable
Bristle Firmness	Soft

7. SAFETY AND WARNINGS

Important: Use only the USB charging cable provided with an ETL or UL certified USB power supply, class III SELV output, 5V DC, 1A Max (power adapter not included).

Do not use for purposes other than oral hygiene. Replace brush heads every 3 months. Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

8. WARRANTY AND SUPPORT

For warranty information, product support, or to purchase replacement brush heads and accessories, please visit the official Quip website or contact Quip customer service. Keep your purchase receipt for warranty claims.