

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [Live Fine](#) /

› [Live Fine 28-Day Automatic Pill Dispenser User Manual](#)

### Live Fine LVFPDRF

# Live Fine 28-Day Automatic Pill Dispenser User Manual

Model: LVFPDRF

## INTRODUCTION

---

The Live Fine 28-Day Automatic Pill Dispenser is designed to simplify medication management by providing organized storage and timely reminders for prescriptions, vitamins, and supplements. This device features an upgraded LCD display, key lock, and customizable sound and light alerts to ensure accurate and consistent dosing for up to four weeks.



Figure 1: Live Fine 28-Day Automatic Pill Dispenser (Frosted Lid)

## KEY FEATURES

- **Easy Pill Organization & Dispensing:** Organizes, monitors, and dispenses medications for up to 4 weeks with audio and visual alerts (1 to 9 times daily).
- **9 Interchangeable Dosage Rings:** Provides various reminder options for improved accuracy and prevention of missed or double doses. 28 high-capacity slots minimize refills.
- **Upgraded LCD Display:** Large, bright screen shows current time, upcoming alarm, alert tone, volume control, battery indicator, and number of alarms set.
- **Adjustable Sound & Light-Up Alerts:** Offers 3 beeping patterns and 3 volume controls (Low, High, Silent). Interface and border lights illuminate, and tone sounds for 30 minutes until pills are dispensed.
- **Simple Setup, Programming & Security:** Battery-operated system with a secure lockable cover and manual key to prevent unauthorized access and accidental spills. Smooth sliding retrieval window for accessibility.

## SETUP GUIDE

### 1. Battery Installation

The dispenser requires 4 AA batteries (included). To install:

1. Locate the battery compartment on the underside of the dispenser.

2. Open the battery compartment cover.
3. Insert 4 AA batteries, ensuring correct polarity (+/-).
4. Close the battery compartment cover securely.

## **POWERED BY 4 AA** **BATTERIES INCLUDED**



Figure 2: Battery Compartment

## **2. Setting Time and Alarms**

The LCD display provides all necessary information for programming.

# EASY-TO-USE CONTROL PANEL WITH LARGE DISPLAY

SHOWING ALL THE INFORMATION YOU NEED!

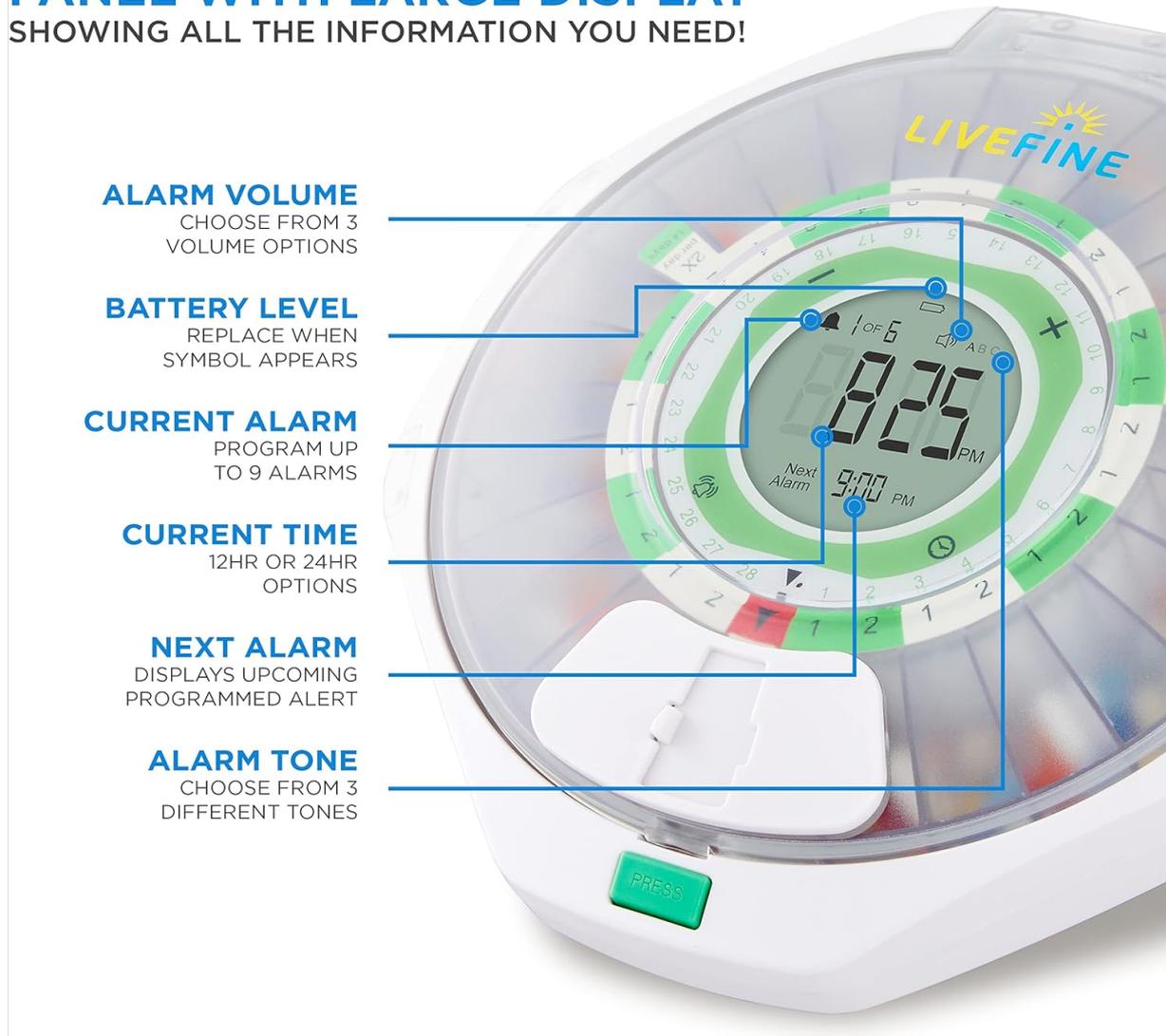



Figure 3: LCD Control Panel

- Use the **Menu** button to navigate through settings.
- Use the **+** and **-** buttons to adjust values.
- Set the current time (12HR or 24HR options available).
- Program up to 9 daily alarms.
- Select alarm tone from 3 different options.
- Choose alarm volume from 3 options (Low, High, Silent).

### 3. Loading Pills and Dosage Rings

The dispenser comes with 9 interchangeable dosage rings to accommodate various medication schedules.

**PILLS AT  
THE READY**  
EASILY ACCESSIBLE  
SLIDING WINDOW

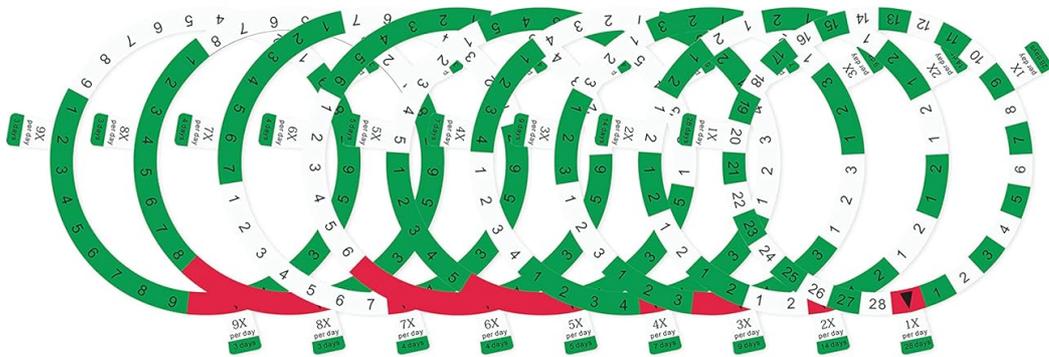

**NEVER LOSE TRACK OF YOUR SCHEDULE**  
9 DOSAGE RINGS INCLUDED

Figure 4: Interchangeable Dosage Rings

1. Unlock the lid using the provided key and open it.
2. Select the appropriate dosage ring based on your medication schedule (e.g., 1x per day for 28 days, 2x per day for 14 days).
3. Place the selected dosage ring onto the central mechanism.
4. Carefully fill each compartment with the prescribed pills according to the dosage ring and your schedule.
5. Ensure the pills are within the compartment boundaries to prevent jamming.



Figure 5: Loading the Pill Tray

## OPERATION

### 1. Receiving Alerts

When an alarm time is reached, the dispenser will activate both sound and light alerts. The lights around the display window will illuminate, and the chosen tone will sound for up to 30 minutes or until the dispenser is inverted to dispense pills.



Figure 6: Alarm and Light Alerts

## 2. Dispensing Pills

To dispense your medication:

1. When the alarm sounds, locate the sliding retrieval window.
2. Press the quick-press latch to open the unlocked cover.
3. Invert the dispenser to allow the pills from the current compartment to fall out into your hand or a container.
4. The alarm will stop once the dispenser is inverted.



Figure 7: Dispensing Pills

### 3. Security Lock

The dispenser features a key lock to secure the lid, preventing unauthorized access and accidental spills. Always use the provided key to lock the lid after loading pills or adjusting settings.



**INCLUDED LOCKING  
KEY AND SECURE LID**  
KEEP PILLS SAFE AND  
PREVENT SPILLS

Figure 8: Key Lock Mechanism

## MAINTENANCE

- **Cleaning:** Wipe the exterior with a soft, damp cloth. Do not submerge the dispenser in water. Ensure the pill compartments are dry before refilling.
- **Battery Replacement:** Replace batteries when the low battery indicator appears on the LCD screen.
- **Storage:** Store the dispenser in a cool, dry place away from direct sunlight and extreme temperatures.

## TROUBLESHOOTING

Problem	Possible Cause	Solution
Dispenser not turning on.	Dead or incorrectly installed batteries.	Check battery polarity or replace with new AA batteries.
Alarm not sounding.	Alarm volume set to silent or alarm not programmed.	Check volume settings on the LCD display. Ensure alarms are correctly programmed.

Problem	Possible Cause	Solution
Pills not dispensing or getting stuck.	Overfilled compartments or pills too large/small for slot.	Ensure compartments are not overfilled. If pills are consistently sticking, try a different dosage ring or consult with a healthcare professional about pill size. Ensure the dispenser is fully inverted.
Lid not locking.	Lid not properly aligned or key mechanism obstructed.	Ensure the lid is closed flush with the base. Check for any debris obstructing the lock mechanism.

## SPECIFICATIONS

- **Model Number:** LVFPDRF
- **Dimensions:** 10.16 x 10.16 x 3.07 inches
- **Weight:** 2.36 Pounds
- **Power Source:** 4 AA batteries (included)
- **Lid Type:** Frosted Lid
- **Features:** Automatic, Portable



## WARRANTY & SUPPORT

For detailed information regarding product warranty, please refer to the documentation included with your purchase or contact Live Fine customer support.

A comprehensive User Guide (PDF) is available for download: [Live Fine Pill Dispenser User Guide](#)

## Related Documents - LVFPDRF

	<p><a href="#">LiveFine Smart Pill Dispenser User Guide: Setup, Operation &amp; Warranty</a></p> <p>Comprehensive user guide for the LiveFine Smart Pill Dispenser (LVFWIFIDC, LVFWIFIDF). Learn how to set up, operate, maintain, and understand the warranty for your automatic medication organizer.</p>
	<p><a href="#">LiveFine Automatic Pill Dispenser User Guide</a></p> <p>Comprehensive user guide for the LiveFine Automatic Pill Dispenser, covering setup, operation, safety precautions, and warranty information. Learn how to program alarms, use dosage rings, and maintain your device for optimal medication management.</p>
	<p><a href="#">LiveFine Automatic Pill Dispenser LVFRAPD5 User Guide</a></p> <p>User guide for the LiveFine Automatic Pill Dispenser (Model LVFRAPD5), detailing features, setup, operation, and maintenance for organizing and dispensing medications.</p>

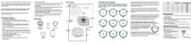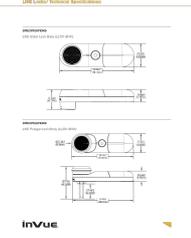	<p><a href="#">LiveFine Automatic Pill Dispenser User Guide: Setup, Operation, and Maintenance</a></p> <p>Comprehensive user guide for the LiveFine Automatic Pill Dispenser. Learn how to set up, install batteries, program alarms, use dosage rings, and maintain your device for safe and effective medication management.</p>
	<p><a href="#">LiveFine Automatic Pill Dispenser Limited Warranty Information</a></p> <p>Details the limited warranty provided by C&amp;A Marketing, Inc. for the LiveFine Automatic Pill Dispenser, including coverage, exclusions, limitations, and contact information.</p>
	<p><a href="#">InVue LIVE Locks Technical Specifications</a></p> <p>Comprehensive technical specifications for InVue's LIVE Lock product line, detailing various lock bodies, accessories, and replacement adhesives with dimensions and order codes.</p>